

CFX Qualification Plate — 96-Well Plate Format

Catalog #	Description
1845098	CFX Qualification Plate , predispensed 96-well plate for performance validation of CFX Duet, CFX Opus 96, CFX Opus 96 Dx, CFX Opus Deepwell, CFX Opus Deepwell Dx, CFX96, CFX96 Touch, CFX96 Touch Deep Well, and CFX Connect Real-Time PCR Systems

For research purposes only.

Description

The CFX Qualification Plate is an optimized assay predispensed in a 96-well plate format for performance validation of the CFX Duet, CFX Opus 96, CFX Opus 96 Dx, CFX Opus Deepwell, CFX Opus Deepwell Dx, CFX96, CFX96 Touch, CFX96 Touch Deep Well, and CFX Connect Real-Time PCR Systems. The plate includes template DNA, SYBR® Green Supermix, primer mix, and nuclease-free water. The CFX Qualification Plate can be incorporated into an instrument qualification procedure.

- Demonstrates a two-fold discrimination with a 99.7% confidence level
- Provides a streamlined workflow with predispensed plate, predefined thermal cycling protocol, and plate templates
- Provides one-click pass/fail performance result with CFX Manager Software, version 3.1, or CFX Maestro Software
- Includes criteria for pass/fail performance result

Contents

One 96-well plate.

Storage

Store at -20°C.

Running the Plate

Note: Always wear gloves when handling the CFX Qualification Plate.

1. Thaw the CFX Qualification Plate in the aluminum bag at room temperature for 5–10 min or until contents are thawed.

Important: The plate contents are light sensitive. Do not thaw the plate outside of the aluminum bag.

2. Once the contents are thawed, centrifuge the plate at 2,000 x g for 2 min.
3. Place the plate into your CFX instrument.
4. Launch CFX Manager Software, version 3.1, or CFX Maestro Software.
5. Click the **Run** tab in the menu bar and select **Qualification Plate Run** (Figure 1).
6. In the Run Setup window, under Start Run on Selected Block(s), make sure the checkbox corresponding to the CFX instrument being used for the Qualification Plate run is selected (Figure 2).
7. Enter the plate's barcode number and click **Start Run** (Figure 2).

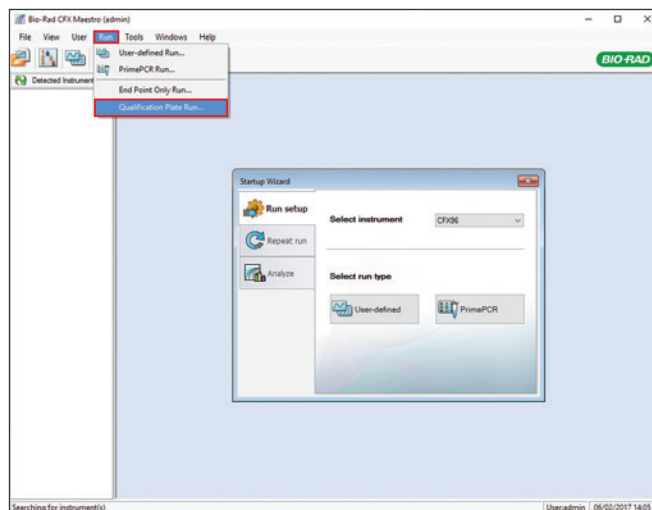


Fig. 1. Selecting Qualification Plate Run in the Run dropdown menu.

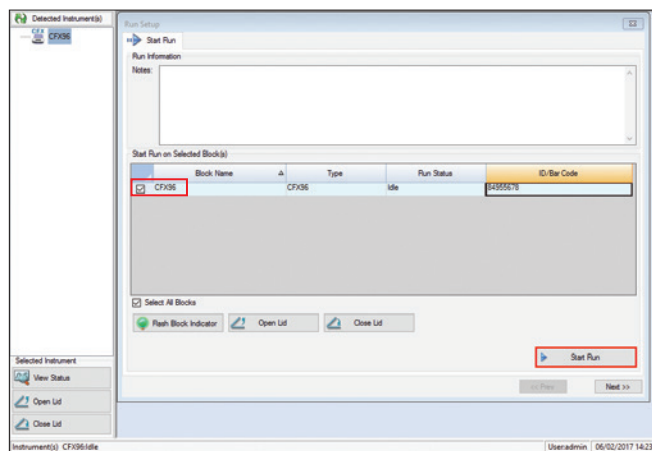


Fig. 2. Selecting Start Run in the Run Setup window.



Analyzing the Data

Click **Tools** in the menu bar. Select **Qualification Plate Report** for a pass/fail performance result (Figure 3).

Note: Up to six wells per unknown group and one well per no template control (NTC) and standard (STD) group can be excluded from analysis to meet passing criteria. Any NTC wells with a nonspecific product melt peak may be excluded. After excluding wells, repeat the step above for a new performance result. If your plate fails after well exclusion, contact Technical Support at 1-800-424-6723, option 2.

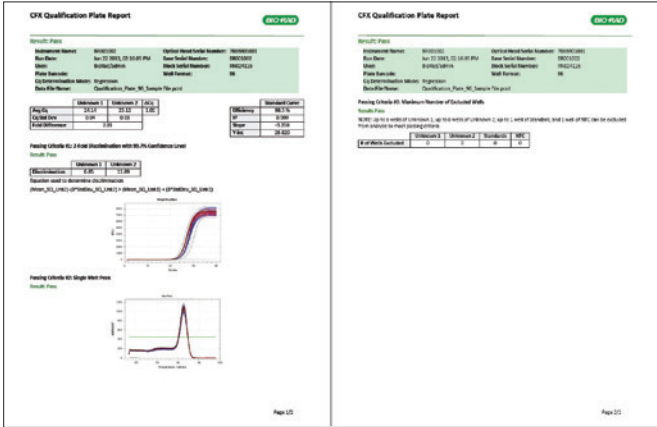


Fig. 3. CFX Qualification Plate Report window.

Passing Criteria

1. Two-Fold Discrimination with 99.7% Confidence Level

$$\begin{aligned} &(\text{Mean_SQ_Unk2}) - (3 \times \text{Stdev_SQ_Unk2}) > \\ &(\text{Mean_SQ_Unk1}) + (3 \times \text{Stdev_SQ_Unk1}) \end{aligned}$$

2. Single Product Melt Peak

Definitions

- Mean_SQ_Unk2 = the average starting quantity of unknown 2
- Stdev_SQ_Unk2 = the standard deviation of the average starting quantity of unknown 2
- Mean_SQ_Unk1 = the average starting quantity of unknown 1
- Stdev_SQ_Unk1 = the standard deviation of the average starting quantity of unknown 1

Note: Quantification cycle (Cq) determination mode is automatically set to regression. CFX96 Qualification Plate data must be analyzed in regression mode to determine pass/fail results.

BIO-RAD is a trademark of Bio-Rad Laboratories, Inc. All trademarks used herein are the property of their respective owner. © 2024 Bio-Rad Laboratories, Inc.

SYBR is a trademark of Thermo Fisher Scientific Inc. Bio-Rad Laboratories, Inc. is licensed by Thermo Fisher Scientific Inc. to sell reagents containing SYBR Green I for use in real-time PCR, for research purposes only.