

Catalog #	Description
161-4013	SureBeads Protein A Magnetic Beads , 3 ml
161-4023	SureBeads Protein G Magnetic Beads , 3 ml
161-4916	SureBeads Magnetic Rack , 16 tube holder
161-4011	SureBeads Protein A Magnetic Beads , 1 ml
161-4021	SureBeads Protein G Magnetic Beads , 1 ml
161-4833	Immunoprecipitation Starter Kit Protein A/G , 1 x 3 ml SureBeads protein A, 1 x 3 ml SureBeads protein G, and 1 x magnetic rack 16 tube holder
161-4813	Immunoprecipitation Starter Kit Protein A , 2 x 3 ml SureBeads protein A and 1 x magnetic rack 16 tube holder
161-4823	Immunoprecipitation Starter Kit Protein G , 2 x 3 ml SureBeads protein G and 1 x magnetic rack 16 tube holder

SureBeads are manufactured with magnetic bead technology from JSR Life Sciences.

SureBeads™ Magnetic Beads for Immunoprecipitation




Catalog #	Description
161-4013	SureBeads Protein A Magnetic Beads , 3 ml
161-4023	SureBeads Protein G Magnetic Beads , 3 ml
161-4916	SureBeads Magnetic Rack , 16 tube holder

BIO-RAD

Protocol


1

Add SureBeads protein A or G.

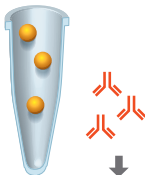

2

Add capture antibody, incubate, and magnetize beads to remove unbound antibody.

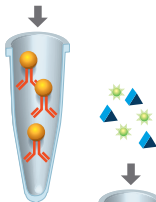


3

Add sample containing target protein and incubate.




4

Magnetize beads, remove supernatant, and wash unbound protein fractions.

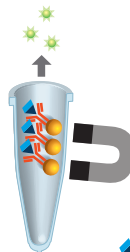

5

Add elution buffer, magnetize beads, and collect purified target protein.

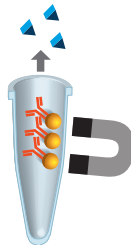

For detailed protocol and step-by-step directions, go to www.bio-rad.com/immunoprecipitation