



iQ-Check™ *Aspergillus*

Quick Guide



37 ± 1°C



48 ± 3 hr



- Enrich the sample in room-temperature buffered peptone water (BPW) supplemented with 0.3 g/L chloramphenicol for 48 ± 3 hr at 37 ± 1°C

- Cannabis and cannabis-infused products: 1:10 dilution of the sample in BPW + chloramphenicol
- Concentrates and cannabis-infused fats and oils: 1:10 dilution of the sample in BPW + chloramphenicol with 1% v/v of Tergitol 7 (or equivalent surfactant)

- After enrichment, transfer 1 ml of the enriched sample to a 1.5 ml screwcap tube

Note: If homogenization of samples using a 1:10 enrichment is difficult, a 1:30 enrichment may be used instead.




5 min



Supernatant

- Centrifuge at 10,000–12,000 x g for 5 min
- Discard all the supernatant

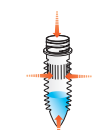

3 ± 1 min



Be sure the lysis reagent is constantly stirring to keep it in suspension.

- Add 200 µl of the lysis reagent (reagent A + reagent F) to the pellet
- Resuspend by pipetting up and down 5 times
- Grind at high speed for 3 ± 1 min in a vortex disruptor



1,300 rpm


95–100°C



10–15 min




- Incubate the tubes at 95–100°C for 10–15 min in a thermoshaker at 1,300 rpm
- If a thermoshaker with a tube adaptor is not available, transfer the entire sample volume to a deep well plate and incubate at 95–100°C for 10–15 min at 1,300–1,600 rpm

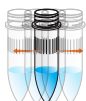


5 min



2 min

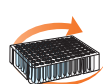

- For tubes, vortex at high speed, then centrifuge at 10,000–12,000 x g for 5 min
- For deep well plates, centrifuge at 2,250 x g for 2 min or allow to settle for 30 min undisturbed

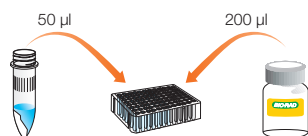


- If needed, to purify the DNA, combine 50 µl of the DNA extract with 200 µl of iQ-Check Purification Reagent* in a new deep well plate


* Sample treatment with iQ-Check Purification Reagent is recommended for inhibitory matrices such as cannabis flower/plant material and chocolate.

Note: Alternatively, dilute the DNA extract 1:10 in sterile water.

- Allow the samples to settle for at least 5 min but not more than 1 hr at room temperature
Avoid pipetting the beads. Do not vortex before collecting the samples.

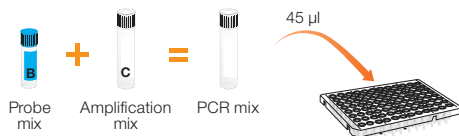

- Prepare the PCR mix (see PCR Mix Calculation Guide)
- Distribute 45 µl/well in a PCR microplate

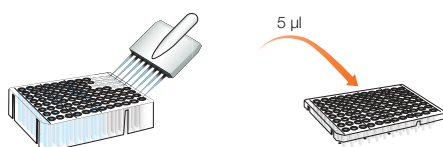

- Add 5 µl each of positive and negative controls
- Transfer 5 µl of the extracted DNA from the deep well plate to the microplate
Do not vortex before collecting the sample.
Eliminate any bubbles.
- Seal the microplate



- Launch CFX Manager IDE Software
- Create the plate setup
- Start the amplification by clicking **Run**

PCR Mix Calculation Guide

To find the correct volumes to use when preparing the PCR mix, add the total number of samples and controls to be analyzed and find the corresponding volumes of reagent B and reagent C in the table.

Total Number of Samples and Controls	Probes Reagent B, μ l	Amplification Mix Reagent C, μ l
1	5	40
2	11	86
3	16	130
4	22	173
5	27	216
6	32	259
7	38	302
8	43	346
9	49	389
10	54	432
11	59	475
12	65	518
13	70	562
14	76	605
15	81	648
16	86	691
17	92	734
18	97	778
19	103	821
20	108	864
21	113	907
22	119	950
23	124	994
24	130	1,000
25	135	1,100
26	140	1,100
27	146	1,200
28	151	1,200
29	157	1,300
30	162	1,300
31	167	1,300
32	173	1,400

Total Number of Samples and Controls	Probes Reagent B, μ l	Amplification Mix Reagent C, μ l
33	178	1,400
34	184	1,500
35	189	1,500
36	194	1,600
37	200	1,600
38	205	1,600
39	211	1,700
40	216	1,700
41	221	1,800
42	227	1,800
43	232	1,900
44	238	1,900
45	243	1,900
46	248	2,000
47	254	2,000
48	259	2,100
49	265	2,100
50	270	2,200
51	275	2,200
52	281	2,200
53	286	2,300
54	292	2,300
55	297	2,400
56	302	2,400
57	308	2,500
58	313	2,500
59	319	2,500
60	324	2,600
61	329	2,600
62	335	2,700
63	340	2,700
64	346	2,800

Total Number of Samples and Controls	Probes Reagent B, μ l	Amplification Mix Reagent C, μ l
65	351	2,800
66	356	2,900
67	362	2,900
68	367	2,900
69	373	3,000
70	378	3,000
71	383	3,100
72	389	3,100
73	394	3,200
74	400	3,200
75	405	3,200
76	410	3,300
77	416	3,300
78	421	3,400
79	427	3,400
80	432	3,500
81	437	3,500
82	443	3,500
83	448	3,600
84	454	3,600
85	459	3,700
86	464	3,700
87	470	3,800
88	475	3,800
89	481	3,800
90	486	3,900
91	491	3,900
92	497	4,000
93	502	4,000
94	508	4,100
95	513	4,100
96	518	4,100

Visit [bio-rad.com/Cannabis](https://www.bio-rad.com/Cannabis) for more information.

Visit bio-rad.com/FoodScience for more information.

BIO-RAD is a trademark of Bio-Rad Laboratories, Inc. IQ-CHECK is a trademark of Bio-Rad Europe GmbH in certain jurisdictions. All trademarks used herein are the property of their respective owner. © 2025 Bio-Rad Laboratories, Inc.



**Bio-Rad
Laboratories, Inc.**

Life Science
Group

Website bio-rad.com **USA** 1 800 424 6723 **Australia** 61 2 9914 2800 **Austria** 00 800 00 24 67 23 **Belgium** 00 800 00 24 67 23 **Brazil** 4003 0399
Canada 1 800 361 1808 **China** 86 21 6169 8500 **Czech Republic** 00 800 00 24 67 23 **Denmark** 00 800 00 24 67 23 **Finland** 00 800 00 24 67 23
France 00 800 00 24 67 23 **Germany** 00 800 00 24 67 23 **Hong Kong** 852 2789 3300 **Hungary** 00 800 00 24 67 23 **India** 91 124 4029300 **Israel** 0 3 9636050
Italy 00 800 00 24 67 23 **Japan** 81 3 6361 7000 **Korea** 82 080 007 7373 **Luxembourg** 00 800 00 24 67 23 **Mexico** 52 555 488 7670
The Netherlands 00 800 00 24 67 23 **New Zealand** 64 9 415 2280 **Norway** 00 800 00 24 67 23 **Poland** 00 800 00 24 67 23 **Portugal** 00 800 00 24 67 23
Russian Federation 00 800 00 24 67 23 **Singapore** 65 6415 3188 **South Africa** 00 800 00 24 67 23 **Spain** 00 800 00 24 67 23 **Sweden** 00 800 00 24 67 23
Switzerland 00 800 00 24 67 23 **Taiwan** 886 2 2578 7189 **Thailand** 66 2 651 8311 **United Arab Emirates** 36 1 459 6150 **United Kingdom** 00 800 00 24 67 23