

Opticon™ 2 Software Tutorial


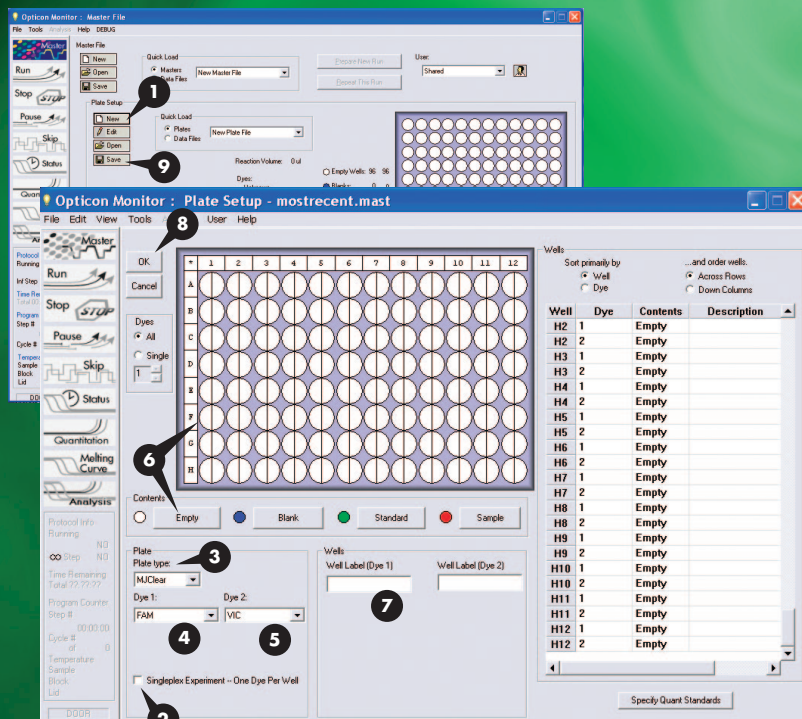


Plate Setup


1. Select **New** from the Plate Setup option bar.
2. If performing a singleplex experiment, check this option.
3. Select plate type from the options listed.
4. Specify Dye 1.
5. Specify Dye 2.
6. Select wells from the Plate diagram and categorize them as **Empty**, **Blank**, **Standard** or **Sample**.
7. Fill in the well label for dye(s) used.
8. When the contents for the entire plate are defined, select **OK**.
9. Select **Save**.

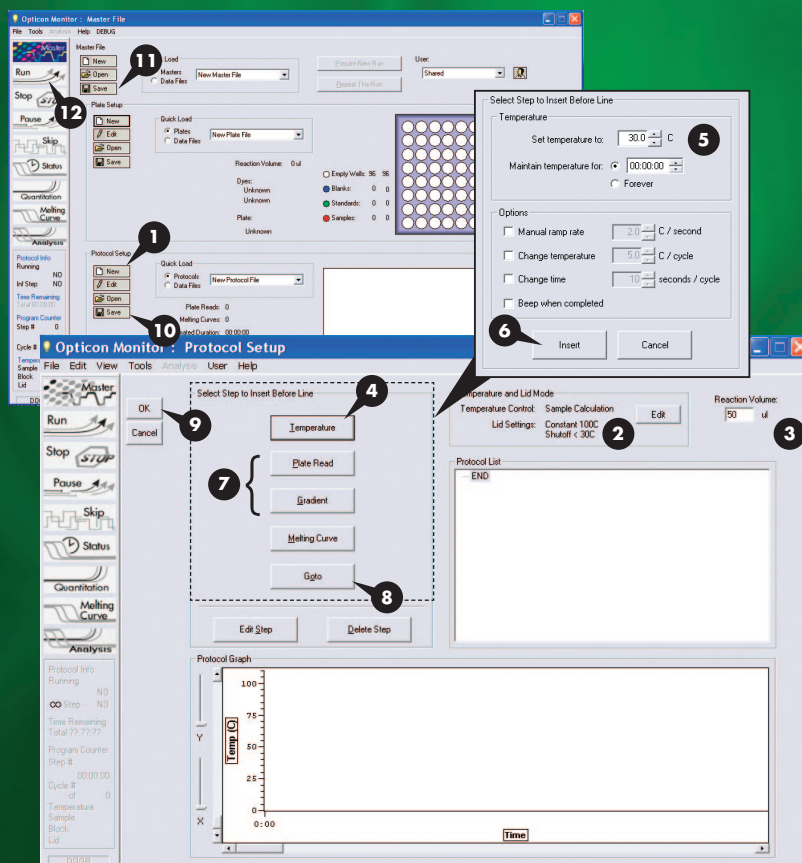

Protocol Setup


1. Select **New** from the Protocol Setup option bar.
2. Select temperature and lid control modes.
3. Enter reaction volume.
4. Select **Temperature**.
5. Fill in the temperature (°C) and hold time.
6. Press **Insert**.
7. Select **Plate Read**, **Gradient** or **Melting Curve** when appropriate. Press **Insert** to save function.
8. Use **Go To** to cycle steps.
9. When the protocol is completely defined, select **OK**.
10. Select **Save**.
11. Select **Save** from the Master File options bar to save both your plate setup and protocol setup under one name.
12. Select **Run**.

Opticon 2 Data Analysis Tutorial

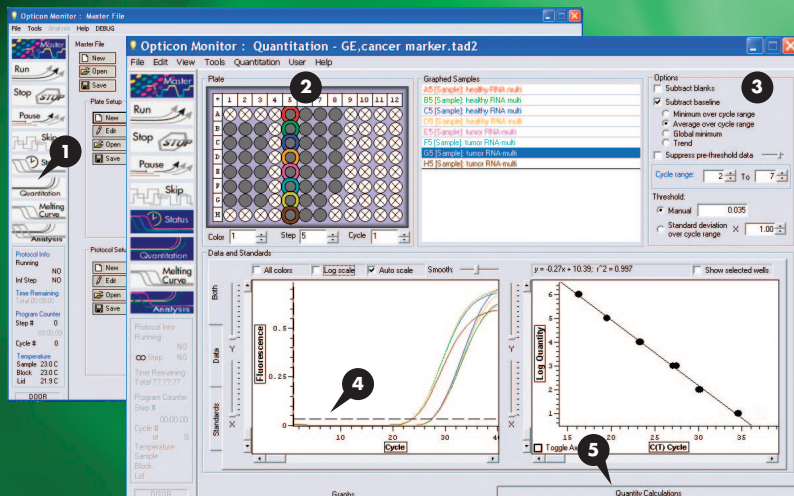


Quantitation

1. Select the **Quantitation** option.
2. Highlight wells for inclusion.
3. Select baseline option.
4. Click and drag the C(t) line as needed.
5. Click **Quantity Calculations** to display these values. Select the **Copy to Clipboard** option under the **Quantitation** drop-down menu to export data.

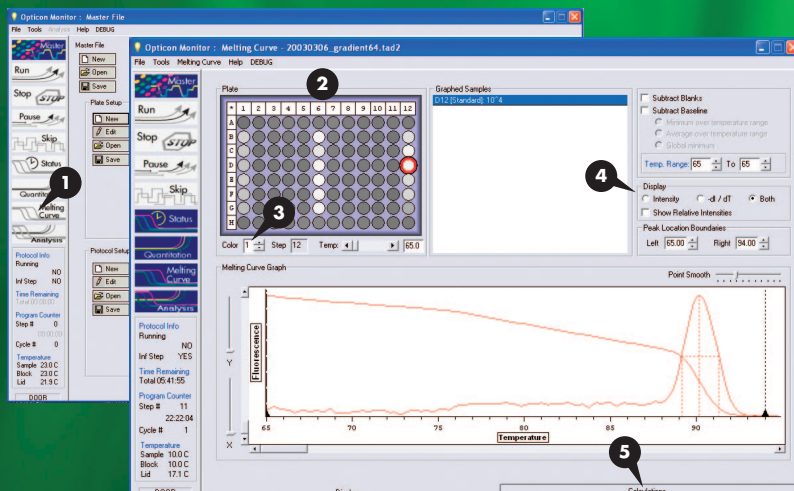

Melting Curves

1. Select the **Melting Curve** option.
2. Highlight wells for inclusion.
3. Select **Color 1** or **Color 2** to display melting curve data.
4. Use the **Display** panel to specify what to graph.
5. Click the **Calculations** tab to display the calculation screen. Select the **Copy to Clipboard** option under the **Melting Curve** drop-down menu to export data.

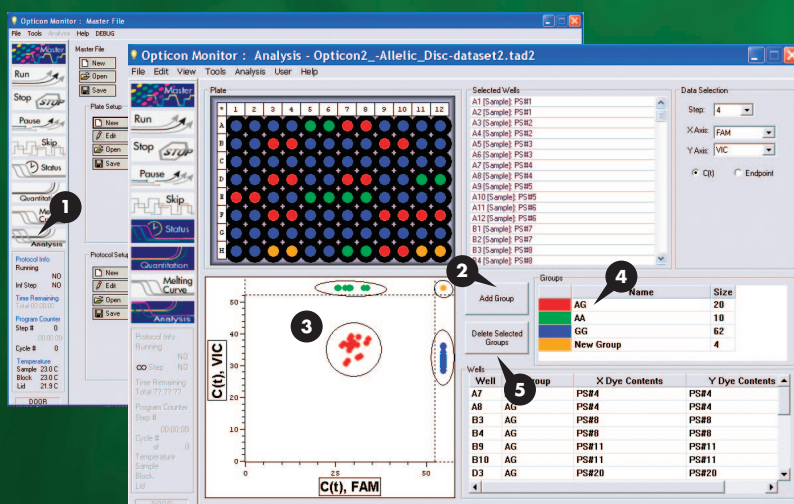

Analysis

1. Select the **Analysis** option.
2. Click **Add Group**.
3. Left-click the mouse on one side of a cluster of points. Move the cross-hairs to the other side of the cluster and left-click again to form an ellipse. Move the cross-hairs perpendicular to the ellipse to expand it and enclose all points in the cluster. Left-click again to complete the set.
4. Rename the group by highlighting the group name and left-clicking. Type the new name and press **Enter**.
5. To delete a group, highlight the group name and click **Delete Selected Groups**.

