

PDQUEST™

User Guide for Version 6.1
Windows and Macintosh

PDQUEST User Guide

Bio-Rad Technical Service Department

Phone: (800) 424-6723, option 2, option 3

(510) 741-6576

Fax: (510) 741-5802

E-mail: LSG.TechServ.US@Bio-Rad.com (U.S.)

LSG.TechServ.Intl@Bio-Rad.com (International)

Notice:

No part of this publication may be reproduced or transmitted in any form or by any means, electronic or mechanical, including photocopy, recording, or any information storage or retrieval system, without permission in writing from Bio-Rad.

PDQUEST and **The Discovery Series** are trademarks of Bio-Rad Laboratories. All other trademarks and registered trademarks are of their respective companies.

WASTE text engine © 1993–1997 Marco Piovanelli

Limitations of Liability:

Bio-Rad is not responsible for the misinterpretation of results obtained by following the instructions in this guide. Whenever possible, you should contact the Technical Services Department at Bio-Rad to discuss your results. As with all scientific studies, we recommend that you repeat your experiment at least once before making any significant conclusions for presentation or publication.

Copyright © 1999 by Bio-Rad Laboratories. All rights reserved.

Table of Contents

1. Introduction	1-1
1.1. Overview of PDQUEST	1-1
1.2. Overview of 2-D Gel Electrophoresis	1-2
1.3. Digital Data and Signal Intensity	1-4
1.4. Steps Involved in PDQUEST	1-5
1.5. Computer Requirements	1-10
1.6. Installing PDQUEST	1-12
1.7. Starting PDQUEST	1-14
1.8. Contacting Bio-Rad	1-22
2. General Operation	2-1
2.1. The Graphical Interface	2-1
2.2. Windows and Subwindows in PDQUEST	2-5
2.3. Mouse-assignable Tools	2-9
2.4. Keyboard Shortcuts	2-10
2.5. File Commands and Functions	2-10
2.6. Preferences	2-22
2.7. File Paths	2-26
2.8. Mouse and Mouse Macros	2-26
2.9. User Settings	2-27
3. Viewing and Editing Images	3-1
3.1. Magnifying and Positioning Tools	3-1
3.2. Hiding Overlays	3-4
3.3. Advanced Viewing Tools	3-4
3.4. Assigning and Interchanging Images	3-6
3.5. Density Tools	3-8

PDQUEST User Guide

3.6. Windows and Subwindows Commands	3-10
3.7. Colors	3-14
3.8. Transform	3-17
3.9. Cropping Images	3-23
3.10. Flipping and Rotating Images	3-26
3.11. Filtering Images	3-29
3.12. Invert Data	3-35
4. Calibrating Gels	4-1
4.1. Calibrating Gels	4-1
4.2. Merging Multiple Exposures	4-11
5. Detecting and Editing Spots	5-1
5.1. Selecting the Detection Parameters	5-1
5.2. Detecting Spots on All Gels	5-10
5.3. Gel Images and Gel Spots Files	5-14
5.4. Displaying Spots	5-14
5.5. Spot Information	5-17
5.6. Editing the Results of Spot Detection	5-21
6. Creating and Editing Matchsets	6-1
6.1. Creating a Matchset	6-2
6.2. Matchset Standard	6-6
6.3. Editing a Matchset	6-7
6.4. Landmarking	6-9
6.5. Auto-Matching	6-10
6.6. Checking and Editing Matches	6-11
6.7. Removing All Matches from an Image	6-21
6.8. Colony Matching	6-21
7. MrpI Values and Spot Quantitation	7-1
7.1. MrpI Standards	7-1

Contents

7.2. Replicate Groups	7-4
7.3. Normalization	7-7
8. Spot Numbering and Annotations	8-1
8.1. Standard Spot Numbers	8-1
8.2. Database Spot Numbers	8-3
8.3. Finding Spots	8-4
8.4. Annotations	8-5
9. Databasing Spots	9-1
9.1. Analysis Sets	9-1
9.2. Higher Level Matchsets	9-15
10. Graphs and Reports	10-1
10.1. Histogram Graphs	10-1
10.2. A-B Comparison	10-7
10.3. Scatter Plots	10-8
10.4. Quantity Reports	10-11
10.5. Quantity Graphs	10-12
10.6. Standard Image Reports	10-14
11. Printing and Exporting	11-1
11.1. Printing	11-1
11.2. Exporting	11-4
Appendix A. Gel Doc 1000/2000	A-1
Appendix B. Chemi Doc	B-1
Appendix C. GS-700 Imaging Densitometer	C-1
Appendix D. GS-710 Imaging Densitometer	D-1

PDQUEST User Guide

Appendix E. Fluor-S Multimager	E-1
Appendix F. Fluor-S MAX Multimager	F-1
Appendix G. Personal Molecular Imager FX	G-1
Appendix H. Molecular Imager FX	H-1
Appendix I. Cross-Platform File Exchange	I-1
Appendix J. Techniques	J-1

Preface

1. About This Document

This user guide is designed to be used as a reference in your everyday use of PDQUEST™. It provides detailed information about the tools and commands of PDQUEST for the Windows and Macintosh platforms. Any platform differences in procedures and commands are noted in the text.

This guide assumes you have a working knowledge of your computer operating system and its conventions, including how to use a mouse and standard menus and commands, and how to open, save, and close files. For help with any of these techniques, see the documentation that came with your computer.

This guide uses certain text conventions to describe specific commands and functions.

Example	Indicates
File > Open	Choosing the Open command under the File menu.
Dragging	Positioning the cursor on an object and holding down the left mouse button while you move the mouse.
CTRL+S	Holding down the Control key while typing the letter s.
Right-click / Left-click / Double-click	Clicking the right mouse button / Clicking the left mouse button / Clicking the left mouse button twice.

Some of the illustrations of menus and dialog boxes found in this manual are taken from the Windows version of the software, and some are taken from the

PDQUEST User Guide

Macintosh version. Both versions of a menu or dialog box will be shown only when there is a significant difference between the two.

2. Bio-Rad Listens

The staff at Bio-Rad are receptive to your suggestions. Many of the new features and enhancements in this version of PDQUEST are a direct result of conversations with our customers. Please let us know what you would like to see in the next version of PDQUEST by faxing, calling, or e-mailing our Technical Services staff. You can also use Solobug (installed with PDQUEST) to make software feature requests.

1. Introduction

1.1 Overview of PDQUEST

PDQUEST is a software package for imaging, analyzing, and databasing 2-D electrophoresis gels.

The software runs in a Windows or Macintosh environment and has an easy-to-use graphical interface with standard pull-down menus, toolbars, and keyboard commands.

PDQUEST collects data from gels scanned with one of several Bio-Rad imaging systems. An image of a sample is captured using the controls in the imaging device window and displayed on your computer monitor. Image processing and analysis operations are performed using commands from the menus and toolbars.

Once a gel has been scanned, advanced algorithms are available to remove background noise, gel artifacts, and horizontal or vertical streaking from the image. PDQUEST then uses a spot segmentation facility to detect and quantify protein spots. Protein spots that change over time can be traced, quantified, displayed on-screen, and exported to other applications for statistical analysis.

The software can be used to simultaneously analyze up to 10,000 spots on a given 2-D gel. PDQUEST can track and report all of the protein pattern and concentration changes in the sample, making feasible both qualitative and quantitative analysis of the data generated from 2-D gels.

The spot matching and databasing facilities make it possible to objectively compare hundreds of different gels. Histograms allow you to quickly compare the amount of the same spot in each gel within an experimental series. Spots can be grouped into user-defined sets. Higher level comparisons and analysis can be performed on master images that you generate from your experimental gels.

PDQUEST User Guide

Scan files can be transferred into or exported out of Discovery Series applications. Scans can be converted into TIFF format for easy compatibility with other applications. Match data can be exported into spreadsheets for easy analysis.

1.2 Overview of 2-D Gel Electrophoresis

2-D gel electrophoresis combines the techniques of isoelectric focusing with SDS polyacrylamide gel electrophoresis. Since these two separation methods rely on independent properties of proteins, chemical (IEF) and physical (SDS), this procedure can resolve complex biological samples and polypeptide mixes. The extraordinary resolution capabilities of SDS-PAGE make it useful for studying complex systems in great detail.

2-D gels generate high resolution separation of proteins. The isoelectric focusing (IEF) dimension is most often run in cylindrical tube gels. After focusing, the cylindrical gel is placed on top of an SDS-PAGE slab gel which is then run as usual. Proteins can be visualized by metabolic radiolabeling or a variety of staining methods.

Gene Expression Studies

2-D gel electrophoresis is an excellent technique for the study of differential gene expression under various growth conditions. Since the expression and regulation of individual proteins can be detected, 2-D gels are an indirect way to monitor gene activity. They allow for the investigation of quantitative as well as qualitative changes in cellular protein expression.

The environmental conditions of a cell can be changed in order to determine optimal growth conditions as well as monitor the cell's response to different stresses. Environmental conditions that can stress the cell include changes in temperature, pH, and nutrient availability. Some examples of the chemical stresses that can be placed on a cell include drug and hormone treatments. Since protein structure and function are the direct result of gene expression, the loss or change of a protein as detected by a 2-D gel can be extrapolated back to events occurring at the DNA level.

Introduction

Many questions encountered when genes are inserted and expressed in bacteria, yeast, and other cell types can be answered with 2-D analysis: Is the cell making the protein? Is the cell's progeny making the protein? Is the protein being made but not secreted? Have mutations occurred?

Medical Applications

2-D gels can also have important applications in medical research. For example, this technique can be used to verify the presence of specific protein markers that are linked to genetic diseases and disease states. Used in conjunction with other tests, 2-D gels can be part of medical screening procedures associated with mutations and teratology linked to genetic damage.

Growth factors are being studied for their role in the regulation of cell growth. 2-D gels can be used to evaluate quantitative and qualitative changes in cellular proteins in response to growth factor stimulation.

2-D gels allow the visualization of proteins whose expression is altered as the result of cell transformation, introducing oncogenes into the host genome. Assessment of phosphorylation, sulfation, or other secondary modifications could reveal functional protein pathways affected by oncogene expression. This information could contribute to a better understanding of cell growth and regulation.

Sample Experiments

Valuable information can be learned by exposing cells to a set of specific experimental conditions and subsequently examining their biological response.

A preliminary in vitro experiment that is useful when beginning 2-D gel work is to radiolabel a cell's proteins to steady-state with ^{35}S -methionine. Running gels of these proteins will help to determine the commonly expressed proteins in the cell under normal conditions.

Subsequently, other amino acids such as ^3H -leucine, ^3H -proline, ^3H -lysine, etc., can be used to radiolabel the proteins to steady-state. This will label any proteins that do not contain methionine and were, therefore, not detected in

PDQUEST User Guide

the first experiment, and will provide preliminary information on the amino acid composition.

Other experiments complementary to 2-D gels include: cell fractionation procedures, and post-translational modifications (phosphorylation, methylation, etc.). Data from these experiments can be added to a database, accumulating information about these proteins.

1.3 Digital Data and Signal Intensity

The Bio-Rad imaging devices supported by PDQUEST are light and / or radiation detectors that convert signals from biological samples into digital data. PDQUEST then displays the digital data on your computer screen, in the form of gray scale or color images.

A data object as displayed on the computer is composed of tiny individual screen pixels. Each pixel has an X and Y coordinate, and a value Z. The X and Y coordinates are the pixel's horizontal and vertical positions on the image, and the Z value is the signal intensity of the pixel.

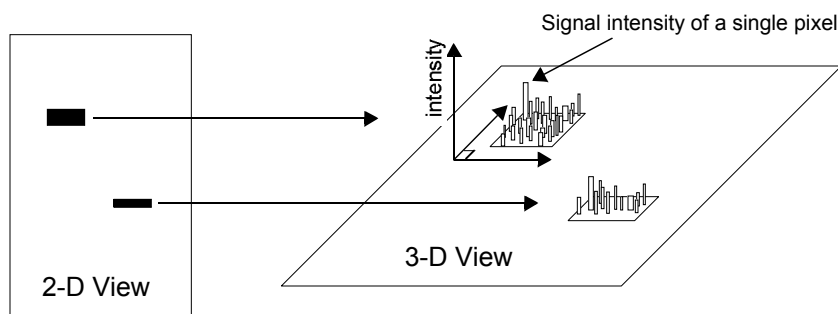


Fig.1-1. Representation of the pixels in two digitally imaged bands in a gel.

For a data object to be visible and quantifiable, the intensity of its clustered pixels must be higher than the intensity of the pixels that make up the background of the image. The total intensity of a data object is the sum of the

Introduction

intensities of all the pixels that make up the object. The mean intensity of a data object is the total intensity divided by the number of pixels in the object.

The units of signal intensity are Optical Density (O.D.) in the case of the GS-700 and GS-710 densitometers, the Gel Doc and Chemi Doc with a white light source, or the Fluor-S and Fluor-S MAX MultiImagers with white light illumination. Signal intensity is expressed in counts when using the Personal FX or FX, or in the case of the Gel Doc, Chemi Doc, Fluor-S, or Fluor-S MAX when using the UV light source.

1.4 Steps Involved in PDQUEST

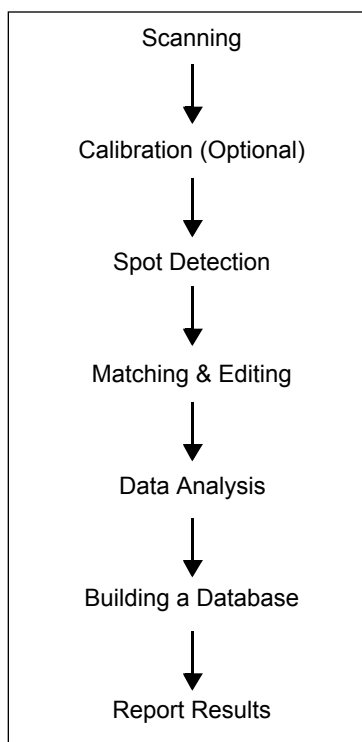


Fig.1-2. Steps in the Analysis of 2-D Gels.

PDQUEST User Guide

Scanning

Before you can use PDQUEST to analyze a gel, you need to capture its image into your computer. This may be done with one of the several Bio-Rad instruments supported by this software: the Molecular Imager FX and Personal Molecular Imager FX; the GS-700 and GS-710 Imaging Densitometers; the Gel Doc 1000/2000 and Chemi Doc gel documentation systems; and the Fluor-S and Fluor-S MAX MultiImagers.

The resulting images can be stored on your hard disk, network file server, or removable storage media.

PDQUEST supports the imaging of the following types of gels and films:

Radiolabeled (^{35}S , ^3H , ^{14}C , ^{32}P , ^{33}P , ^{125}I) gel films, wet, silver, and Coomassie-stained gels, and duplicate contact films from stained gels.

Calibration (Optional)

Calibration strips can be used by PDQUEST to create a calibration curve that establishes the relationship between film density and radioactivity per unit area in the gel.¹ Calibration strips consisting of segments of known radioactivity, are processed simultaneously with ^{35}S , ^3H , or ^{14}C radiolabeled gels. The data from each calibration strip is used by PDQUEST to compare the average optical density of the individual segments of the calibration strip with the average disintegrations per minute (DPM) of each sampled segment, and assemble a calibration curve based on that data. If the curve is acceptable, you save the experiment and thereby save the associated calibration information.

Spot Detection

Once your gel has been scanned, the next step is detecting and quantifying the protein spots. If you have multiple exposures of the same gel image, and those exposures are calibrated, the spot detection routine in PDQUEST “merges” the available information in each exposure to produce a single **Gel**

1. See Garrels, J.I., Farrar, J.T., and Burwell, C.B., *J. Biol. Chem.*, **254**, 7961-7977 [1979]

Introduction

Image. If you have uncalibrated gels or single exposures, a Gel Image is still made, but the information comes from a single 2-D scan.

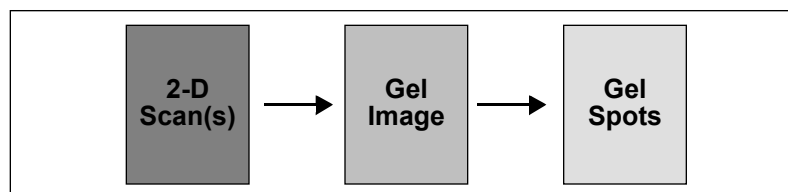


Fig.1-3. Images created and processed during spot detection.

The advantage of having multiple gel exposures is that you can obtain more quantitative information from the proteins on the gel. On a dark exposure, faint spots are quantifiable but dark spots are saturated and cannot be quantitated. On a light exposure, faint spots may not be detectable. Merging lets you quantitate faint spots from the dark exposures and uses light exposures to quantitate dark spots.

The merged Gel Image is processed to remove streaks and background noise. Spots are then detected and quantified from the processed Gel Image, and a new, synthetic, image called a **Gel Spots file** is created. Spot quantitation is stored as part of the Gel Spots file.

Spot detection and quantitation are performed using 2-dimensional Gaussian modeling. This method of detection will accurately detect and quantify overlapped spots as well as isolated spots. The end result is a Gel Spots file that represents each detected spot as an ideal Gaussian with associated DPM or PDU (protein data unit) values.

Spot detection reports image data in units of optical density (O.D.). O.D. is calculated by dividing the incident light intensity by the transmitted light intensity, and taking the log (base 10) of this quotient:

Optical Density (O.D.) = $\log_{10}(\text{Intensity } i / \text{Intensity } t)$, where:

Intensity ***i*** = incident light intensity, and

Intensity ***t*** = transmitted light intensity.

PDQUEST User Guide

Thus, an O.D. of 0.0 means that all the light is transmitted through the film, whereas an O.D. of 2.0 means that only 1 percent ($\log_{10} [100/1]$) of the light is transmitted through the film.

Spot detection requires several steps—detecting spot centers and quantifying spots are the last two steps in the overall process. The Auto-Detect Spots function (Quantify > Auto-Detect Spots) performs all the necessary steps automatically. These steps are outlined below:

<i>Create Gel Image</i>	Using calibration strip data, you can merge multiple scans derived from different film exposures of the same gel into a single gel image. This increases the effective dynamic range of the film such that both very faint and very dark spots can be accurately quantified. Each pixel of a 2-D scan is originally assigned an O.D. value based upon a step tablet calibration. For radiolabeled gels that include calibration strips, O.D. values are converted to DPMs; for uncalibrated or non-radiolabeled gels, linear interpolations of O.D.s are used to express spot quantitation. These images are called O.D. Gels.
<i>Initial Smoothing</i>	Gel Images are passed through an initial smoothing algorithm to remove densitometer noise.
<i>Background Subtraction</i>	After the initial smoothing pass, background density is subtracted from the image. This includes vertical and horizontal streaks as well as film fog and overall background density caused by the carrier medium (film, gel matrix, or blot matrix).
<i>Final Smoothing</i>	One or more additional smoothing passes are necessary to remove extraneous spots at or near the background level. The exact number of passes is empirically determined (e.g., 2 passes for radiolabeled gels).

Introduction

Locate Spots in Gel Image

This step locates the center and position of each recognizable spot in the Gel Image. This is the spot detection step of the process.

Fit and Quantify Spots

This step fits ideal Gaussian distributions to spot centers. Quantitation for each spot is determined by calculating the volume under the Gaussian.

The final fitting step of spot detection produces a gel spots file for each scan. A gel spots file is a list of all the detected spots, their X and Y coordinates, and their associated Gaussian parameters (X-sigma, Y-sigma, and height). PDQUEST uses this information to store gel data and to display a synthetic representation of the gel image on the computer screen.

Matching and Editing

After you have spot-detected a set of gels, you will probably want to compare the protein spots from one gel with those from another to examine changes in protein levels under experimental conditions. Groups of gels can be edited and matched to one another in a “matchset.” Matchsets consist of gel spots files and can include gel images. In a matchset, protein spots are matched to each other, enabling you to compare their quantities.

One member of your matchset can be used as a template for creating a matchset Standard. The Standard combines all the unique, matchable spots in a matchset into a single file.

Data Analysis

PDQUEST provides a variety of analytical tools to help you determine which sets of protein spots are statistically and scientifically meaningful.

Since it is virtually impossible to analyze all your experimental data at once, PDQUEST offers methods for breaking down the information into smaller, more manageable pieces. This is done by creating sets of spots called “Analysis Sets.” All the spots in an Analysis Set satisfy qualitative, quantitative, or statistical criteria that you specify.

PDQUEST User Guide

Building a Database

As the amount of protein data that you collect increases, it will become important to store that information in an organized and accessible format. PDQUEST provides several functions for building and analyzing a protein spot database, enabling you to store and compare data from various experiments using higher level matchsets.

Report Results

When your analysis is complete, you can print your experimental data or export it to another system for further analysis.

1.5 Computer Requirements

PDQUEST will run on a PC as a Windows 95/98 or Windows NT 4.0 32-bit application or on a Macintosh PowerPC.

The amount of computer memory required for using PDQUEST is mainly determined by the file size of the images you will scan and analyze. Images scanned at high resolution can be quite large. For this reason, we recommend that you archive images on a network file server or removable storage media.

PC

The following are the **minimum** system requirements for installing and running PDQUEST on a PC:

Operating system:	Windows 95/98. Recommended: Windows NT 4.0.
Processor:	Pentium 166.
RAM:	64 MB for Gel Doc, Chemi Doc, and Fluor-S imaging systems. 128 MB for FX, Personal FX, and densitometers.
Hard disk space:	3 GB. Recommended: Removable storage media (such as an Iomega Jaz drive) or a network file server.

Introduction

Monitor:	17" monitor, 1024 x 768 resolution, 256 colors.
SCSI:	Required for all Bio-Rad imaging devices except the Gel Doc and Chemi Doc. Adaptec SCSI card with EZ-SCSI software.
Printer:	Optional.

Macintosh

Note: The default amount of memory assigned to PDQUEST on the Macintosh is 64 MB. If the total RAM in your Macintosh is 64 MB or less (the minimum recommended amount is 128 MB), you should reduce the amount of memory assigned to PDQUEST to 10 MB less than your total RAM. With the PDQUEST application icon selected, go to File > Get Info in your Finder to reduce the memory requirements for the application. See your Macintosh computer documentation for details.

The following are the **minimum** system requirements for installing and running PDQUEST on a Macintosh:

Operating system:	System 8.0 or higher.
Processor / Model:	PowerPC Macintosh 9500.
RAM:	64 MB for Gel Doc, Chemi Doc, and Fluor-S imaging systems. 128 MB for FX, Personal FX, and densitometers.
Hard disk space:	3 GB. Recommended: Removable storage media (such as an Iomega Jaz drive) or network file server.
Monitor:	17" monitor, 1024 x 768 resolution, 256 colors.
SCSI:	Required to run all Bio-Rad imaging devices except the Gel Doc and Chemi Doc. Macintosh has built-in SCSI.
Printer:	Optional.

PDQUEST User Guide

1.6 Installing PDQUEST

The procedures for installing PDQUEST on a PC and a Macintosh are described in separate sections below.

1.6.a Installing on a PC

PDQUEST can be installed on a PC from a CD-ROM or you can download the installation program from the Internet. Make sure that Windows is up and running on your computer before attempting to install PDQUEST.

Installing from a CD-ROM

Insert the Discovery Series™ CD-ROM into the CD-ROM drive on your computer. The installer will start automatically.

Downloading from the Internet

You can download PDQUEST from the Internet using Bio-Rad's Web site. Go to www.bio-rad.com/984248.html and select Download PDQUEST from the list of Windows applications. Follow the instructions to download the installer onto your computer. Then you can double-click on the Setup.exe icon on your desktop to begin running the installer.



Fig.1-4. Setup.exe icon (Windows).

The Installation Program

The PDQUEST installation program will guide you through a series of screens. The installer will create a default directory tree under Program Files on your hard disk called Bio-Rad / The Discovery Series / Bin (you can select a different directory if you wish). The main PDQUEST program will be placed

Introduction

in the Bin directory. Additional directories for storing user profiles and sample images will also be created

The installer will place a PDQUEST icon on your desktop and create a Discovery Series folder in the Programs folder on your Windows Start menu.

After installation is complete, the installer will ask you if you want to open PDQUEST.

1.6.b Installing on a Macintosh

PDQUEST can be installed on a Macintosh from a CD-ROM or you can download the installation program from the Internet.

Installing from a CD-ROM

Insert the Discovery Series CD-ROM into the CD-ROM drive on your Macintosh. Double-click on the CD icon on your desktop, then double-click on the Install PDQUEST icon in the CD window.



Fig.1-5. PDQUEST installation program icon (Macintosh).

Downloading from the Internet

You can download PDQUEST from the Internet using Bio-Rad's Web site. Go to www.bio-rad.com/984248.html and select Download PDQUEST from the list of Macintosh applications. Follow the instructions to download the installer and place it on your desktop. Then you can double-click on the installer icon to begin running the installer.

Installed Files and Folders

The installer will guide you through a series of screens. The installer will create a folder on your hard drive called PDQUEST, which contains the main

PDQUEST User Guide

application and associated sample images (you can select a different folder if you wish). The installation will also create a folder called The Discovery Series in the Preferences folder in your System Folder; this contains the on-line Help and various system files.

Once installation is complete, the folder containing the PDQUEST application will appear open on your desktop.

1.7 Starting PDQUEST

PC

The installation program creates a PDQUEST application icon on your desktop. To start the program, double-click on this icon.

The installer also creates a Discovery Series folder in the Programs folder on your Windows Start menu. You can start PDQUEST by selecting the icon in this folder.



Fig.1-6. PDQUEST application icon.

Macintosh

After installation, the PDQUEST folder will be open on your desktop. To start the program, double-click on the application icon inside the folder.

1.7.a Software License

When PDQUEST first opens, you should see a dialog box titled “Welcome to PDQUEST,” which shows the current status of your software license.

Introduction

PDQUEST has a software licensing system. You can have full use of the software for 30 days free of charge, after which you must purchase the software and obtain a password for continued use.

Most users will be able to start PDQUEST and begin the free trial period immediately after installation. If you open the software and see a “Free Trial” message, you can skip to section 1.7.c, Trial Period.

Unable to Obtain Authorization

If you attempt to start PDQUEST and receive an “Unable to obtain authorization” message, you need a hardware security key (HSK) to run PDQUEST. HSKs for both PC and Macintosh are included with the full CD-ROM package, or are available from Bio-Rad for downloaded software. If you receive this message, turn off your PC or Macintosh, attach the appropriate HSK as described in the following section, and restart the computer and PDQUEST.

Note: Network License holders: Your network administrator is responsible for the software license. If you have difficulty starting and running PDQUEST, contact your administrator.

1.7.b Attaching the Hardware Security Key (If Necessary)

Before attaching the HSK, you should first attempt to start PDQUEST without the HSK. If PDQUEST opens successfully, you can skip this section.

Note: Some parallel port devices such as zip drives may be incompatible with HSKs. Please check with your peripherals vendor.

PDQUEST User Guide

PC



Fig.1-7. PC Hardware Security Key

Before proceeding, please turn off your computer. If you have a printer attached to your computer's parallel port, please turn that off as well.

The HSK attaches to the parallel port on the back of your PC. If a printer cable is attached to this port, disconnect it. After you have attached the HSK, you can attach the printer cable to the key itself and restart your computer and printer.

The code for the PC hardware security key is EYYCY. This is printed on the key itself.

You will also need to install the system driver that allows your PDQUEST software to "read" the HSK.

Note: Windows NT users must be in the local administrator group to install the HSK driver.

To install the driver, click on your Windows Start menu and select Programs > The Discovery Series. (This directory will not exist until you install PDQUEST.) Select Install HASP Hardware Security Key driver to begin installation.

Note: Windows 95/98 users must reboot their computer after installing the HSK. Windows NT users do not have to reboot.

Introduction

Macintosh

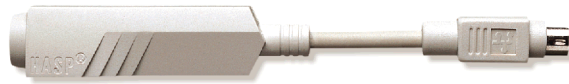


Fig.1-8. Macintosh Hardware Security Key

Before proceeding, please turn off your Macintosh.

The Macintosh HSK must be inserted in the Apple Desktop Bus (ADB) path. The ADB port is located on the back of your Macintosh.



Fig.1-9. Apple Desktop Bus icon on back of Macintosh.

The HSK can be inserted at any point in the ADB path—between the computer and the keyboard, between the keyboard and the mouse, between the keyboard and the monitor, etc. After you have attached the HSK, you can restart your computer.

The code for the Macintosh HSK is QCDIY. This code is printed on the key itself.

1.7.c Trial Period

When you start PDQUEST, the Software License screen will open. The type of screen you see will depend on whether you have an HSK, Network License, or neither.

Free Trial (No HSK or Network License)

If you do not have an HSK attached and do not have a Network License, the “Free Trial” screen will be displayed.

PDQUEST User Guide

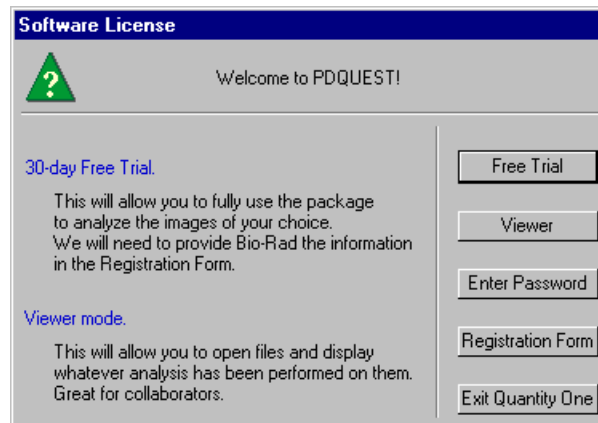


Fig.1-10. Main Registration Screen without HPK

Click on the Free Trial button. This will open the Software License Registration Form (see next section).

Temporary License (With HSK or Network License)

If you have an HSK attached or are using a Network License, the Software License screen will reflect the fact that you receive a 30-day temporary license ("Your license will expire on _____").

Introduction

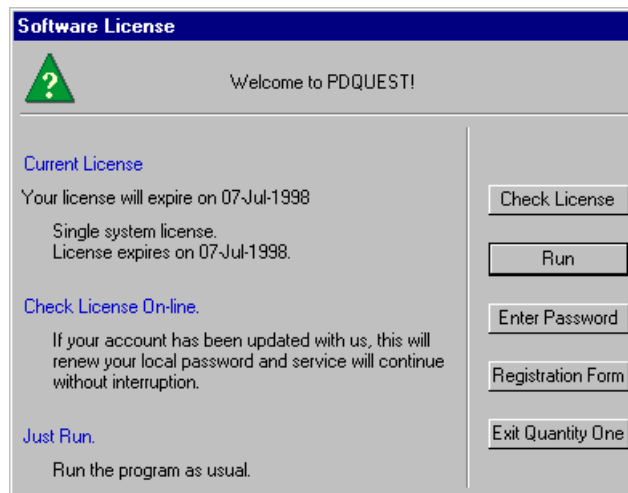


Fig.1-11. Main Registration Screen with HPK.

Click on the Run button to begin using the software.

If you are using an HSK, some time during the 30-day license period, click on the Registration Form button to register your software. If your 30-day period has expired, a Free Trial button will appear when you open the software. Click on this button to open the Software License Registration Form and register your software.

If you are using a Network License, any time during the initial 30 days click on the Check License button. Your full Network License will be activated.


Note: Network License holders: If your Network License is not activated when you click on Check License, notify your network administrator.

1.7.d Software Registration

Note: Network License holders do not need to register their software, and can skip the following section.

To obtain a full individual software license, you will need to fill out in the information in the Software License Registration Form.

PDQUEST User Guide

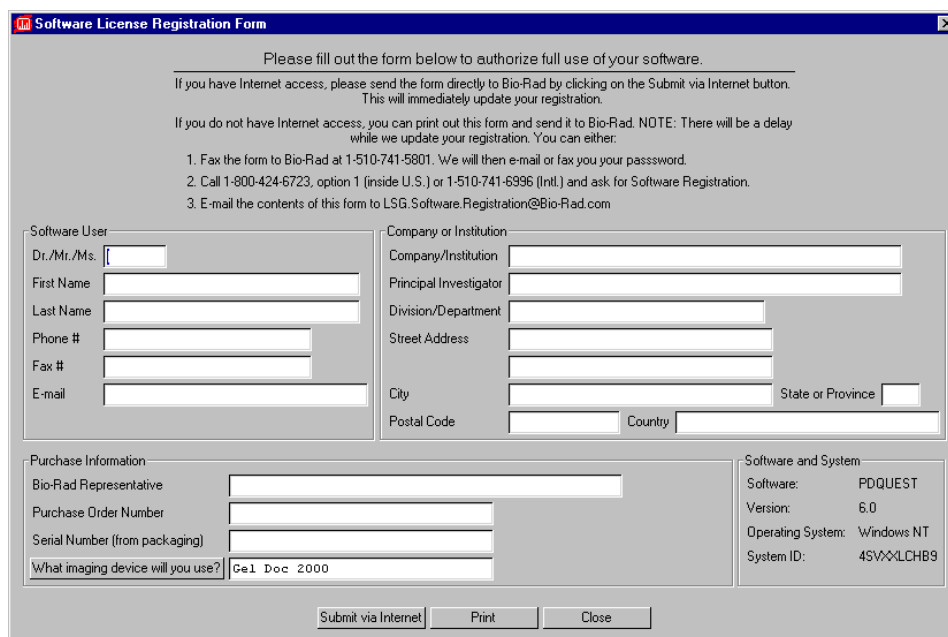

The image shows a screenshot of a software license registration form titled "Software License Registration Form". The form is designed for users to provide their contact and company information to authorize software use. It includes instructions on how to submit the form (via Internet, fax, or mail) and a list of three submission options. The form is divided into several sections: "Software User" (with fields for name, phone, fax, and email), "Company or Institution" (with fields for company name, principal investigator, division, street address, city, state/province, postal code, and country), "Purchase Information" (with fields for Bio-Rad representative, purchase order number, serial number, and imaging device), and "Software and System" (with fields for software name, version, operating system, and system ID). At the bottom, there are buttons for "Submit via Internet", "Print", and "Close".

Please fill out the form below to authorize full use of your software.

If you have Internet access, please send the form directly to Bio-Rad by clicking on the Submit via Internet button. This will immediately update your registration.

If you do not have Internet access, you can print out this form and send it to Bio-Rad. NOTE: There will be a delay while we update your registration. You can either:

1. Fax the form to Bio-Rad at 1-510-741-5801. We will then e-mail or fax you your password.
2. Call 1-800-424-6723, option 1 (inside U.S.) or 1-510-741-6996 (Intl.) and ask for Software Registration.
3. E-mail the contents of this form to LSG.Software.Registration@Bio-Rad.com

Software User

Dr./Mr./Ms.

First Name

Last Name

Phone #

Fax #

E-mail

Company or Institution

Company/Institution

Principal Investigator

Division/Department

Street Address

City State or Province

Postal Code Country

Purchase Information

Bio-Rad Representative

Purchase Order Number

Serial Number (from packaging)

What imaging device will you use?

Software and System

Software: PDQUEST

Version: 6.0

Operating System: Windows NT

System ID: 4SVXQLCHB9

Fig.1-12. Software License Registration Form.

Note: If you do not yet have a Purchase Order Number or Software Serial Number, you may leave these fields blank to receive a trial license.

Registering by Internet

If you have Internet access on your computer (the same computer on which you loaded PDQUEST), you can register quickly and easily.

In the Software License Registration Form, click on the Submit via Internet button.

Your information will be submitted automatically over the Internet, and a password will be generated automatically and sent back to your computer. Simply continue to run PDQUEST as before.

Introduction

This password will be good for 30 days. During this period, if you have already submitted your Software Serial Number, click on Check License in the Software License screen to update your license. (To access the Software License screen from within PDQUEST, select Help > Register.)

If you have not yet submitted your Software Serial Number, open the form again, enter the serial number, and resubmit it over the Internet. In 1–2 days, click on the Check License button to update your license.

Registering by Fax or E-mail

If you do not have Internet access, click on the Print button in the Software License Registration Form and fax the form to Bio-Rad at the number listed on the form. Alternatively, you can enter the contents of the form into an e-mail and send it to Bio-Rad at the address listed in the Registration Form.

If you include your Software Serial Number in the form, Bio-Rad will contact you by fax or e-mail in 2–3 days with a full license.

If you do not include your Software Serial Number in the form, Bio-Rad will contact you in 2–3 days with a temporary 30-day license. Then, some time during the 30-day period, enter the serial number in the form and fax the revised form to Bio-Rad. Bio-Rad will contact you by fax or e-mail in 2–3 days with a full license.

Entering a Password

If you fax or e-mail your registration information, you will receive a password from Bio-Rad. You must enter this password manually.

To enter your password, click on Enter Password in the Software License screen. If you are not currently in the Software License screen, select Register from the Help menu in PDQUEST.

PDQUEST User Guide

The 'Enter Password' dialog box contains the following information:

- Registration Information:**
 - To register, please contact Bio-Rad during US/Pacific business hours:
 - Phone: (800) 424-6723 opt 2 opt 3 (in the U.S.)
(510) 741-6576
 - Fax: (510) 741-5802
 - e-mail: LSG.TechServ.US@Bio-Rad.com (in the U.S.)
LSG.TechServ.Intl@Bio-Rad.com (International)
 - Please provide the information in the completed registration form.
- Current license:**
 - Single system license.
 - No time limit.
 - Version limit = 6.9.
- System ID:**
 - 4SVXXLCHB9
- Password:**
 - PDQUEST []
 - Buttons: OK (yellow), Not OK (yellow)
- Navigation:**
 - Buttons: Enter, Cancel

Fig.1-13. Enter Password screen.

In the Enter Password screen, type your password next to the PDQUEST prompt.

Once you have entered the correct password, the OK light next to the password field will change to green and the Enter button will activate. Click on Enter to run PDQUEST.

1.8 Contacting Bio-Rad

Bio-Rad technical service hours are from 8:00 a.m. to 4:00 p.m., Pacific Standard Time in the U.S.

Phone: (800) 424-6723, option 2, option 3
(510) 741-6576

Fax: (510) 741-5802

E-mail: LSG.TechServ.US@Bio-Rad.com (in the U.S.)
LSG.TechServ.Intl@Bio-Rad.com (International)

For software registration, phone:

(800) 424-6723, option 1, or (510) 741-6996

2. General Operation

2.1 The Graphical Interface

2.1.a Menu Bar

PDQUEST has a standard menu bar with pulldown menus that contain all the major features and functions available in the software.

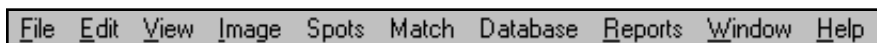


Fig.2-1. Menu bar.

The pulldown menus contain groups of related functions:

- File—General file commands (Load, Save, Print), imaging device acquisition windows (Gel Doc, Chemi Doc, GS-700, GS-710, Fluor-S, Fluor-S MAX, Personal FX, FX).
- Edit—Calibration, Normalization, Preferences, File Paths.
- View—General image viewing tools (Zoom Box, Grab), quick analysis tools (Density in Box, Plot Density), image assign functions.
- Image—Transform function, image processing tools (Crop, Flip, Filter).
- Spots—Spot detection, editing, and displaying tools.
- Match—Spot matching and matchset functions.
- Database—Analysis set and annotation functions, standards and replicate groups functions, spot analysis tools.
- Reports—Graphs and reports.
- Window—Tile windows, Imitate Zoom, Configure Subwindows.
- Help—On-line help, keyboard layout, software registration.

PDQUEST User Guide

Below the menu bar is the main toolbar, containing some of the most commonly used commands. Next to the main toolbar are the status boxes, which provide information about cursor selection and toolbar buttons.

2.1.b Main Toolbar

The main toolbar in PDQUEST appears below the menu bar. It includes buttons for the main file commands (Open, Save, Print) and essential viewing tools (Zoom Box, Grab, etc.), and buttons to open the secondary toolbars.



Fig.2-2. Main toolbar.

Tool Help

If you hold the cursor over a toolbar icon, the name of the command will pop up below the icon. This utility is called Tool Help. Tool Help appears on a time delay basis that can be specified in the Preferences dialog box under Edit > Preferences. You can also specify how long the Tool Help will remain displayed.

2.1.c Status Boxes

There are two status boxes in PDQUEST. These appear to the right of the main toolbar.

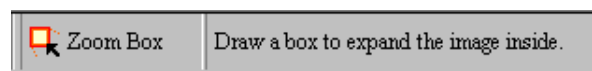


Fig.2-3. Status boxes.

The first box displays any function that is assigned to the mouse (see section 2.3, Mouse-assignable Tools). If you select a command such as Zoom Box, the

General Operation

name and icon of that command will appear in this status box and remain there until another mouse function is selected or the mouse is deassigned.

The second status box is designed to supplement Tool Help (see above). It provides additional information about the toolbar buttons. If you hold your cursor over a button, a short explanation about that command will be displayed in this second status box.

2.1.d Secondary Toolbars

PDQUEST has secondary toolbars that contain icons for groups of associated functions. You can open these toolbars from the main toolbar or from the View > Toolbars pulldown menu.

The secondary toolbars can be toggled between vertical, horizontal, and expanded formats by clicking on the resize button on the toolbar itself.

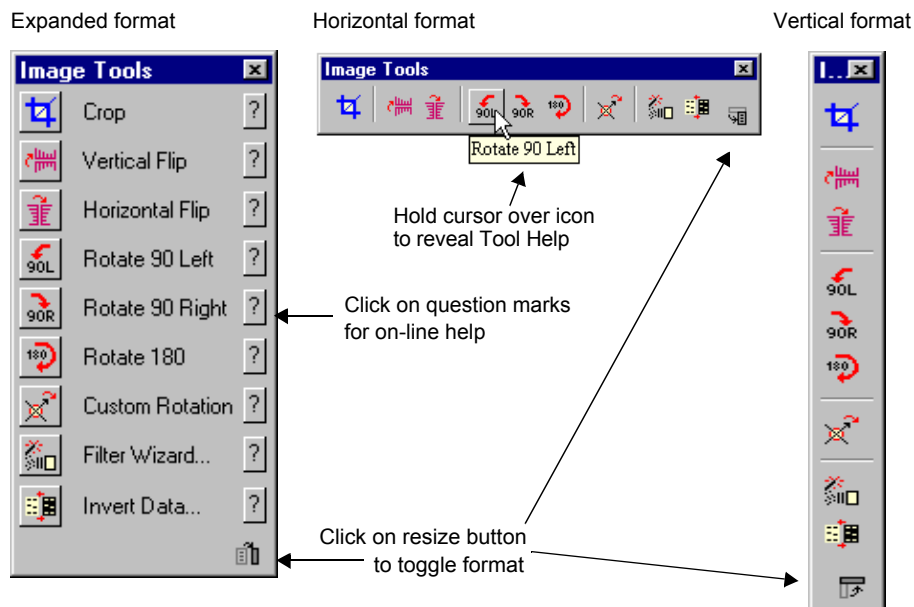


Fig.2-4. Secondary toolbar formats and features.

PDQUEST User Guide

The expanded toolbar format shows the name of each of the commands. Click on the “?” button next to the name to display on-line Help for that command.

2.1.e Quick Guide

The PDQUEST Quick Guide (accessed from the Help menu or main toolbar) is a tool palette that is organized to guide you through the major functions of the software. Like the secondary toolbars, it can be toggled between vertical, horizontal, and expanded formats by clicking on the resize button.

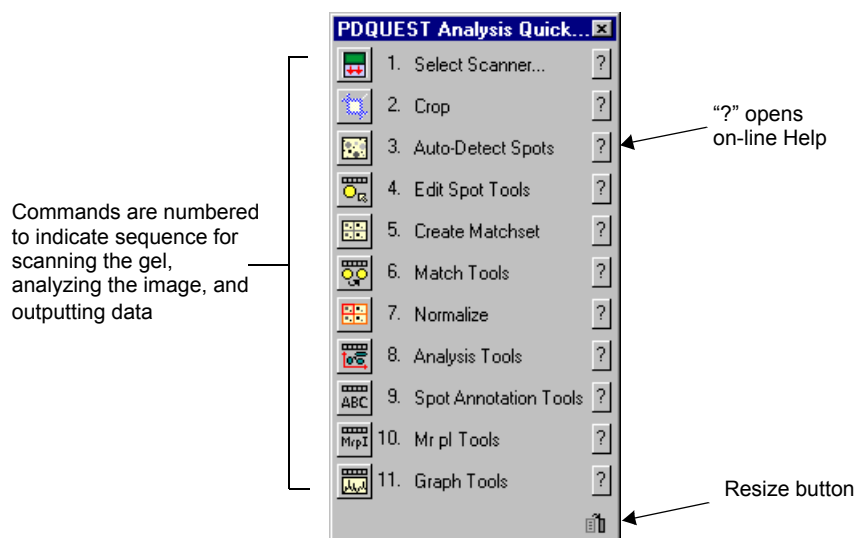


Fig.2-5. PDQUEST Quick Guide

In the expanded format, the Quick Guide commands are numbered as well as named, so that the order of operation is clear. Simply follow the steps and the Quick Guide will lead you through the analysis. Note that many of the buttons open toolbars of related commands (e.g., Match Tools).

As with the secondary toolbars, you can click on the “?” next to the name of a function to display the Help text.

2.2 Windows and Subwindows in PDQUEST

PDQUEST handles image windows differently than other, standard Windows applications.

In a standard Windows application, you open a single file and it appears in a single window. But in PDQUEST, you can “load” more than one image into the same window. You can then display your images in subwindows of the main window, as well as switch between different images in the same window / subwindow.

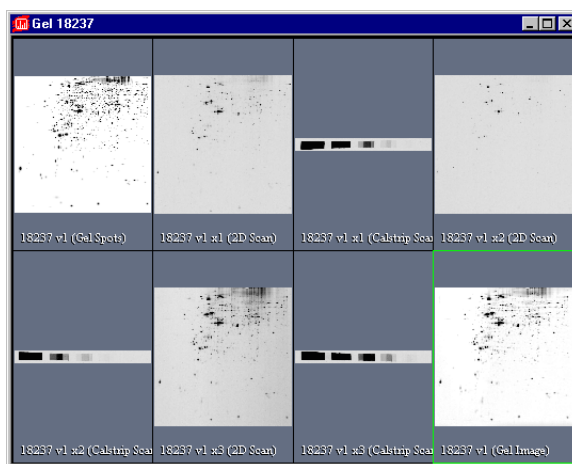


Fig.2-6. A single image window with seven subwindows.

This grows out of the fact that, in PDQUEST, a single gel can have multiple images, all of which have the same root name. Subwindows are also necessary when working with matchsets, which are composed of multiple images.

PDQUEST User Guide

	*18237 v1	Gel Image
	*18237 v1	Gel Spots
	*18237 v1 x1	2D Scan
	18237 v1 x1	Calstrip Scan
	18237 v1 x2	2D Scan
	*18237 v1 x2	Calstrip Scan
	18237 v1 x3	2D Scan
	18237 v1 x3	Calstrip Scan

Fig.2-7. Gel 18237 has eight image files, all of which can be loaded into the same window.

You start with a 2D Scan, which you convert into a Gel Image, which you convert into a Gel Spots image. You can also have multiple exposures of a 2D Scan, as well as Calstrip Scans associated with those exposures. All of these files can be loaded into a single window.

In addition, related images can be grouped into matchsets, which are displayed in one window. And the individual images in those matchsets can consist of both Gel Spots images and Gel Images.

You start by loading each image individually into a single window. Then PDQUEST has various functions that allow you to view the different images that are loaded in that window.

The Assign command under the View menu allows you view and select from the list of images that can be displayed in a given window.

General Operation

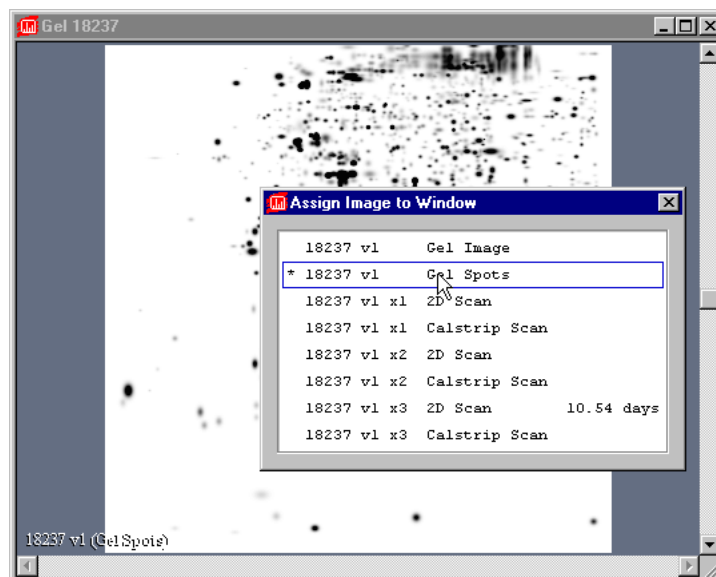


Fig.2-8. Assign Image lists the images loaded in an image window and allows you to select the one to display.

The Configure Subwindows command under the Window menu allows you to display the different images assigned to a given window in different combinations of subwindows.

PDQUEST User Guide

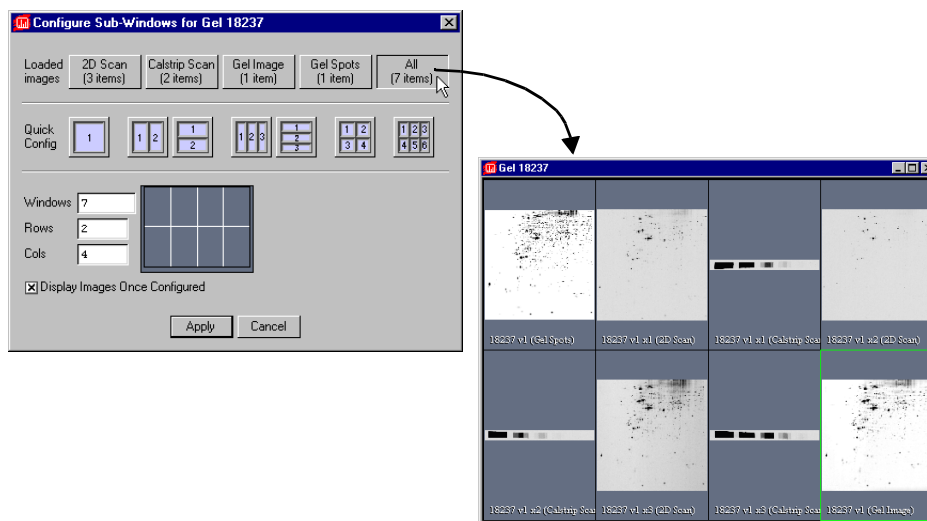


Fig.2-9. Configuring subwindows.

These windows and subwindows commands are described in greater detail in Chapter 3.

Note: Some commands in PDQUEST become mouse-assignable if you are displaying multiple subwindows. See the next section.

Navigating Subwindows

When multiple subwindows are displayed, you activate a subwindow by simply moving your cursor over it. A green border appears around the subwindow, and the border shifts as you move your cursor over a different subwindow. To execute a command in the active subwindow, you can click in the subwindow (for mouse-assignable functions) or use a keyboard key, toolbar button, or menu command.

Note: When selecting toolbar and menu commands, be careful that you don't inadvertently activate a different subwindow as you are moving your cursor to select the command. For this reason, it is advisable to use the

General Operation

mouse-assignable commands and keyboard functions when multiple subwindows are displayed.

If you prefer to activate subwindows by clicking on them, select the Click to Focus Tile option under Edit > Preferences. With this option selected, you must first click on the subwindow to activate it, then perform the command.

2.3 Mouse-assignable Tools

PDQUEST has a number of commands that don't perform actions right away, but instead "assign" a function to your mouse (e.g., Zoom Box, Density at Cursor, Make Spot Crosshair). You first select one of these tools from a menu or toolbar, then execute the command by clicking or dragging on the image.

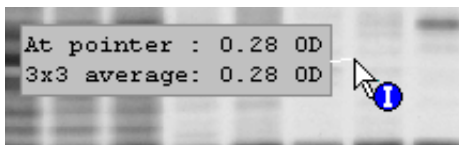


Fig.2-10. Example of a mouse-assignable tool: Density at Cursor

When a "mouse-assignable" function is selected, the cursor appearance changes. The name and icon of the function appear in the status box next to the main toolbar (see section 2.1.c).

To *deassign* a function from the mouse, click on the toolbar button of the assigned function or click in the status box displaying the assigned function. You can also double-click on the Hide Overlays button.

Note: Some commands in PDQUEST *become* mouse-assignable if you are displaying multiple subwindows. For example, if you select Spots > Spot Crosshairs with a single image displayed in a window, crosshairs will simply be displayed. However, if you select Spots > Spot Crosshairs with multiple subwindows displayed, the command will be assigned to your mouse and you must click on a subwindow to display the crosshairs in that subwindow.

2.4 Keyboard Shortcuts

Many of the functions in PDQUEST may be executed using keyboard commands (e.g., pressing the F1 key will activate the View Entire Image command). A list of key combinations and their associated commands will be displayed if you select Keyboard Layout from the Help menu.

The pulldown menus also list the shortcut keys for the menu commands.

Note: Mouse-assignable commands (see above) behave differently if you assign them using the keyboard versus selecting them from the menus or toolbars. For example, to use the Zoom Box command as a keyboard command, position your cursor on the image where you want to begin to create the magnifying box, then press the F2 key. The command is assigned to your mouse and immediately activated, so that when you move your cursor over the image, the zoom box is created. When you click the mouse button once, the defined region is magnified and the tool is automatically deassigned from your mouse.

2.5 File Commands and Functions

This section describes the basic file commands and functions of PDQUEST. These are selected from the File menu.

Note: Print and Export commands are described in Chapter 11

General Operation



Fig.2-11. File menu.

Load

You use the File > Load command to load a previously saved image or a matchset into a window. It is similar to “Open” in a standard Windows application, except that it doesn’t always open a new window for each image. Images with the same root name are loaded into the same window (see section 2.2, above).

You can also click on the Load button on the main toolbar.

Note: This version of PDQUEST will load any images and matchsets created with earlier versions of PDQUEST.

PDQUEST User Guide

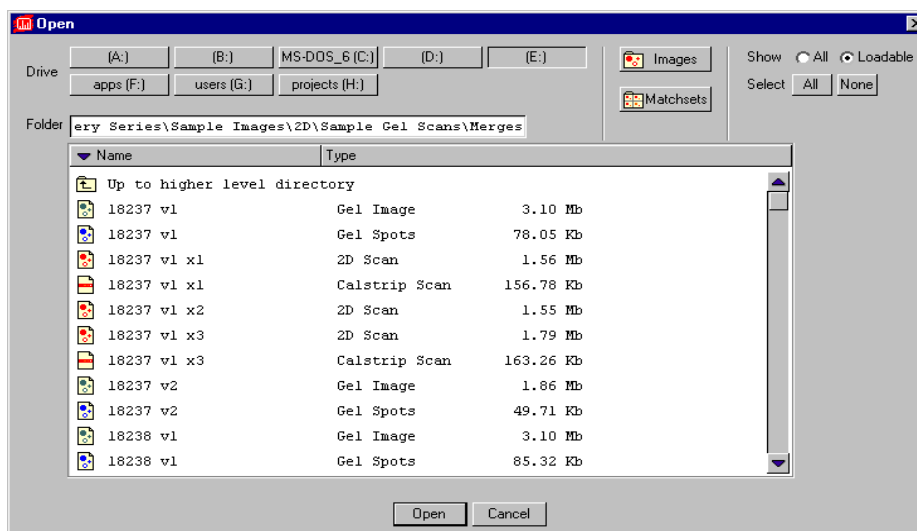


Fig.2-12. Load dialog box.

The Load dialog box differs from the standard Windows Open dialog in several key ways:

- Loaded images are marked with an asterisk.
- File types are listed next to each image name.
- Matchsets are listed as single files, even though they are really collections of images.

In the dialog box, load an image by double-clicking on its name or selecting the file name and clicking on the Open button.

To open multiple files, hold down the CTRL key and click on each file name. To open all the files in a directory, click on the Open button with no file names selected. You will be prompted to open all files.

Note: Remember that files with the same root name will be loaded into the same window.

General Operation

To navigate through other directories, click on the appropriate folder name or click on the “Up to higher level directory” icon. You can also type the full file path into the Folder field.

To look for files on other drives, click on the appropriate drive letter button at the top of the dialog.

Note: The drive buttons at the top of the dialog box are only for mapped drives. If you need to access an unmapped (network) drive, you can type the drive name directly into the Folder field using UNC notation (e.g., \\remote server name\shared directory name\...\file name).

To go directly to the folder where your images are stored, click on the Images button. To go to the folder containing your matchsets, click on Matchsets. The locations of these folders are specified in the File Paths dialog box, accessed under the Edit menu.

Click on the Desktop button to access your Windows or Macintosh desktop.

You can show all the files or only loadable files by selecting the appropriate button next to the Show prompt. You can also Sort the files by name or type.

An image created in the Windows version of PDQUEST can be opened in the Macintosh version, and visa versa. However, you must add the appropriate extension to your Macintosh files to open them in Windows (e.g., .gsc for gel scans, .gim for gel images, .csc for calstrip scans, etc.).

PDQUEST comes with a selection of sample images. In Windows, these files may be found in the Discovery Series / Sample Images / 2D directory. On the Macintosh, they are stored in the Sample Images folder in the PDQUEST folder. These sample images are related, and can be used to create a matchset.

Loading TIFF Images

The Load command can also be used to import TIFF images from other software applications.

TIFF images do not contain all of the tagged information that would normally be included in a Discovery Series image file (e.g., scanner type, date, color, etc.). For this reason, the File Info and Scan Report commands will not display all such tagged information from TIFF files.

PDQUEST User Guide

There are many types of TIFF formats that exist on the market. Not all are supported by the Discovery Series. There are two broad categories of TIFF files that are supported:

1. 8-bit Grayscale. Most scanners have an option between line art, full color, and grayscale formats. Select grayscale for use with the Discovery Series software. In a grayscale format, each pixel is assigned a value from 0 to 255, with each value corresponding to a particular shade of gray. Hewlett Packard™, Microtek™, and Sharp™ each make scanners that produce compatible 8-bit grayscale images.
2. 16-bit Grayscale. Bio-Rad's Molecular Imager (storage phosphor) systems use 16-bit pixel values to describe intensity of scale. Molecular Dynamics™ imagers also use 16-bit pixel values. The Discovery Series understands these formats and can interpret images from both Bio-Rad and Molecular Dynamics storage phosphor systems.

Note: PDQUEST can import 8- and 16-bit TIFF images from both Macintosh and PC platforms.

TIFF files that are *not* supported include:

1. 1-bit Line Art. This format is generally used for scanning text for optical character recognition or line drawings. Each pixel in an image is read as either black or white. Because the software needs to read continuous gradations to perform gel analysis, this on-off pixel format is not used.
2. 24-bit Full Color or 256 Indexed Color. These formats are frequently used for retouching photographs and are currently unsupported in the Discovery Series, although most scanners that are capable of producing 24-bit and indexed color images will be able to produce grayscale scans as well.
3. Compressed Files. The software does not read compressed TIFF images. Since most programs offer compression as a selectable option, files intended for compatibility with the Discovery Series should be formatted with the compression option turned off.

General Operation

Unload

To unload a single image from a window or subwindow, select the window or subwindow, then select File > Unload. If the image is in a subwindow, the image will close but the larger window will remain open; if the image is in its own window, the entire window will close. If the file has changed since it was last saved, you will be prompted to save the changes before closing.

Close

File > Close will close the entire active window, including any subwindows. If any of the images in the window has changed, you will be prompted to save the changes before closing.

Close All

File > Close All closes all open images in all windows. If you have made changes to any of the images, a message box will open for each unsaved image giving you the option of saving those changes.

Save

Save will save a new or old scan, scan set, or matchset to your hard disk, network drive, or other storage media.

In Windows, new PDQUEST images will be given a .gsc extension when they are first saved.

Save As

Save As can be used to save scans, scan sets, and matchsets under different names and different directories.

Note that in the case of scan sets (i.e., multiple exposures, versions, and/or calstrips of the same gel, all with the same root name), all the different files in the set must be loaded when you select Save As if you want the root name to be changed for all of them.

PDQUEST User Guide

Save All

Save All on the main toolbar or File menu saves all images, scan sets, and /or matchsets that are currently open.

Change Version/Exposure

Change Version /Exposure allows you to change the version and /or exposure number of a scan file.

Note: The root name of the file will remain unchanged.

You can also change the type of the file (to 1-D, 2-D, Calstrip, or DNA scan).

Select Change Version /Exposure from the File menu.

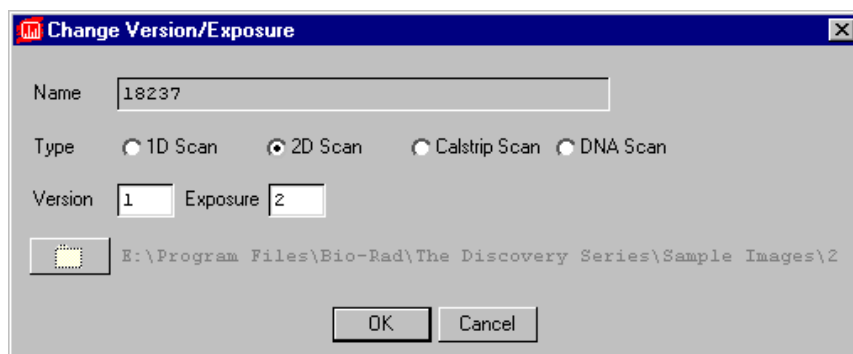


Fig.2-13. Change Version/Exposure dialog box.

In the dialog box, enter a new version and /or exposure number in the appropriate field(s), and click on OK.

Revert to Saved

Revert to Saved will reload the last saved version of the image you are working on. This is a quick way to undo any alterations you may have made to an image since you last saved it.

General Operation

Any changes you have made since last saving the file will be lost. (Also, all open dialog boxes will be closed.) A message box will warn you before completing the command.

File Info

File Info displays general information about your image, including the scan date, scan area, number of pixels in the image, data range, and the size of the file. Gel Spots files include the number of spots, total data in the spots, and other information. There is also a field where you can type in a file description or comments.

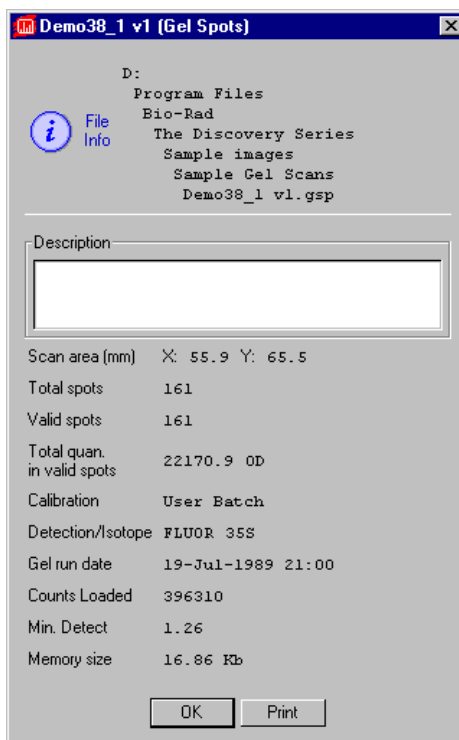


Fig.2-14. File Info box.

PDQUEST User Guide

To print the file information, click on the Print button.

Change the Image Dimensions

You can change the dimensions of certain images using the File Info dialog box.

Note: This feature is only available for images captured by a camera (such as the Gel Doc or Fluor-S) or imported TIF images in which the dimensions are not already specified.

When you select File > File Info for these types of images, the File Info dialog box will include fields for changing the image dimensions.

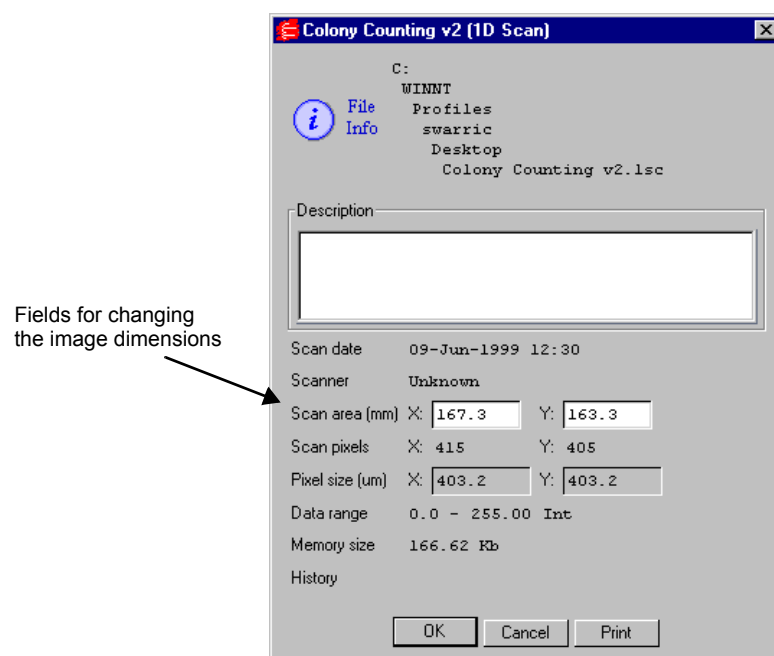


Fig.2-15. File Info dialog box with fields for changing the image dimensions.

General Operation

Enter the new image dimensions (in millimeters) in the appropriate fields. Note that the pixel size in the image (in micrometers) will change to retain the same number of pixels in the image.

Reduce File Size

Scan files can be quite large, and computer systems do not have unlimited memory or storage space. If you are having difficulty loading or storing a particular scan, you may want to reduce the size of the file by reducing the number of pixels in the image. (You can also trim unneeded parts of an image to reduce its memory size. See section 3.9, Cropping Images.)

Note: You cannot reduce the file size of Gel Images or Gel Spots files

This function is comparable to scanning at a lower resolution, in that you are increasing the size of the pixels in the image, thereby reducing the total number of pixels and thus memory size.

Note: Reducing the file size of an image will result in some loss of resolution. In most cases this will not affect quantitation. In general, as long as the pixel size remains less than 10 percent of the size of the objects in your image, changing the pixel size will not affect quantitation.

Select File > Reduce File Size to open the Reduce File Size dialog box.

The dialog box shows you the size of the pixels in the image (Pixel Size: X by Y microns), the number of pixels in the image (Pixel Count: X by Y pixels), and the Memory Size of the image.

As you increase the size of the pixels, the pixel count will decrease, as will the memory size. You can increase the pixel size in either dimension (see the following figure for an example).

PDQUEST User Guide

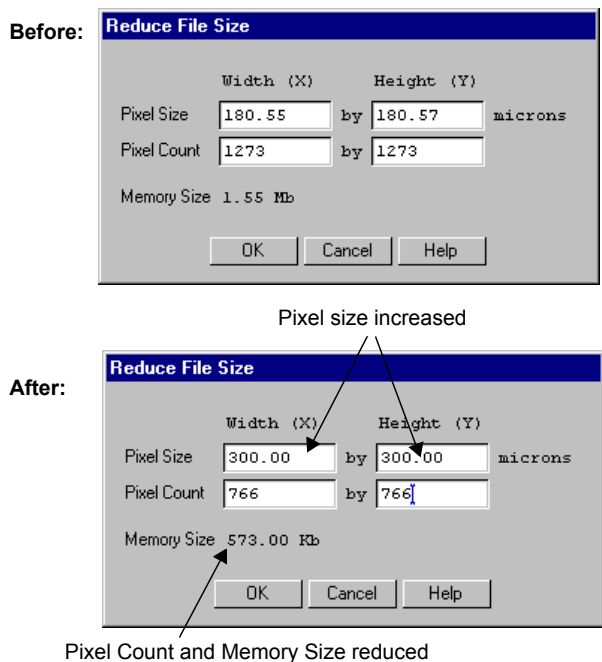


Fig.2-16. Reduce File Size dialog box, before and after pixel size increase.

Note: Asymmetric pixel reduction (i.e., making pixels smaller in one dimension than the other) is not recommended for 2-D gels, since you want to resolve spots in both dimensions.

When you are finished, click on the OK button.

A pop-up box will give you the option of reducing the file size of the displayed image or making a copy of the image and then reducing the copy's size.

Reducing the file size is an irreversible process. For that reason, we suggest you make a copy of the image and reduce its file size. That way, if you lose too much resolution, you can simply delete the copy and try again. Once you are happy with the reduced image, don't forget to delete the original. The goal is to save space!

General Operation

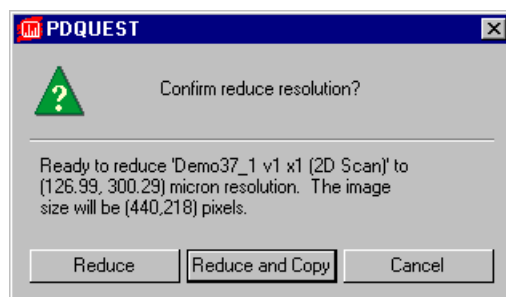


Fig.2-17. Confirm Reduce File Size pop-up box.

If you choose to make a copy of the image, you will be asked to enter a name for the new copy before the operation is performed.

Imaging Device Acquisition Windows

The File menu contains a list of Bio-Rad imaging devices supported by the Discovery Series software. These are:

1. Gel Doc 1000 / 2000
2. Chemi Doc
3. GS-700 Imaging Densitometer
4. GS-710 Imaging Densitometer
5. Fluor-S MultiImager
6. Fluor-S MAX MultiImager
7. Personal Molecular Imager FX
8. Molecular Imager FX

Selecting a name in the File menu will open the acquisition window that allows you to scan using that instrument.

See the individual appendices on the imaging devices for more details.

PDQUEST User Guide

Exit

File > Exit quits the application. You will be prompted to save any unsaved changes to images or matchsets.

2.6 Preferences

The Preferences dialog box (accessed under the Edit menu) allows you to customize basic features of your system, such as displays and toolbars.

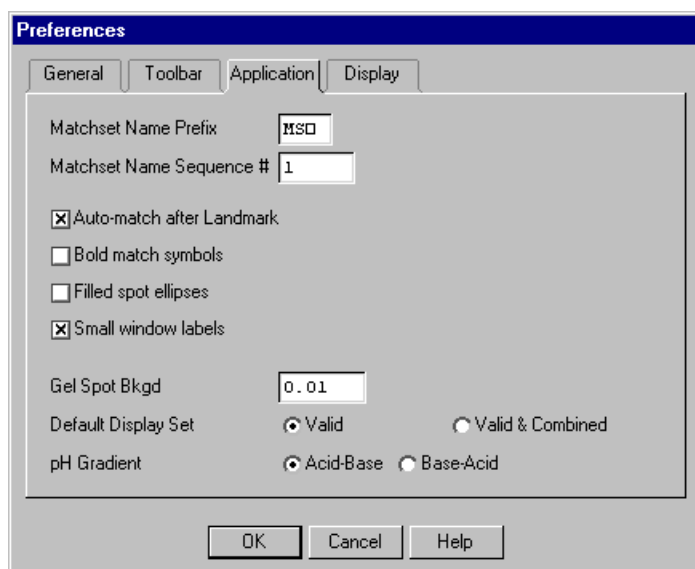


Fig.2-18.

Click on the appropriate tab to access groups of related preferences. After you have selected your preferences, click on OK to implement them. Some of your preferences will only be implemented by exiting and restarting the program. A message will notify you if this is necessary.

General Operation

2.6.a General

Click on the General tab to access the following preferences.

Start Maximized

In the Windows version, the Start Maximized checkbox determines whether PDQUEST occupies the entire computer screen when first opened. If this is unchecked, PDQUEST's menu and status bars will appear across the top of the screen and any toolbars will appear "floating" on the screen.

Multiple Mouse Assignment Mode

Multiple Mouse Assignment Mode is designed for backward compatibility with previous versions of PDQUEST. It is not documented in this manual.

2.6.b Toolbar

These settings determine the behavior and positioning of the secondary toolbars and Quick Guide. Click on the Toolbar tab to access these preferences.

If the Show Quick Guide checkbox is selected, the Quick Guide will open automatically when you open PDQUEST.

If Align Quick Guide with Document is selected, the Quick Guide will pop up flush with the edge of your documents. Otherwise, it will appear flush with the edge of the screen.

Quick Guide Placement and Toolbar Placement determine which side of the screen the Quick Guide and toolbars will first open.

The Placement Behavior setting determines whether a Quick Guide or toolbar will always pop up in the same place and format (Always Auto), or whether they will pop up in the last location they were moved to and the last format selected (Save Prior).

The Toolbar Orientation option buttons specify whether toolbars will first appear in a vertical, horizontal, or expanded format when you open PDQUEST.

PDQUEST User Guide

The Guides Always on Top checkbox determines whether Quick Guides will always be displayed on top of images, or whether they can be hidden behind images.

Tool Help Delay allows you to specify the amount of time (in seconds) the cursor must remain over a toolbar icon before the Tool Help appears. First-time users may want to specify a short delay to learn the names of the toolbar functions, while experienced users can specify a longer delay once they are familiar with the icons.

Persistence determines how long the Tool Help will linger on the screen after you move the cursor off a toolbar icon.

2.6.c Application

Click on the Application tab to access these preferences.

Enter up to three letters for the Matchset Name Prefix and a number for the Matchset Name Sequence #.

Each matchset that you create will be identified by a prefix followed by a unique sequence number. The prefix remains the same for all the matchsets while the sequence number increases by 1 each time a new matchset is created.

Once the prefix and the sequence number have been initially set, you probably will not need to change them.

If Auto-match after Landmark is selected, PDQUEST will automatically attempt to match the spots in your images when you create a landmark.

If Bold Match Symbols is selected, landmarks, manual matches, added spots, etc. will appear larger and bolder on the screen.

If Filled Spot Ellipses is checked, marked spots (e.g., spots marked using the Saturated Spots command) will be displayed as filled ellipses. If unchecked, they will be displayed as unfilled ellipses.

If Small Window Labels is checked, any text in the windows/subwindows will be displayed in a smaller font.

General Operation

If you want to apply a uniform background to Gel Spots images so they more closely resemble the original gels, enter a background optical density (O.D.) value in the text box next to the Gel Spot Bkgd prompt.

Specify whether you want to see only valid spots or valid spots along with combined spots by clicking on the appropriate button next to the Default Display Set prompt.

Finally, indicate the orientation in which your gels are run by clicking on either the Acid-Base or Base-Acid button next to the pH Gradient prompt.

2.6.d Display

Click on the Display tab to access these preferences.

Zoom %

The Zoom % field allows you to specify the percentage by which an image moves closer or farther away when you use the Zoom functions. This percentage is based on the size of the image.

Pan %

Pan % determines the percentage by which the image moves side to side or up and down when you use the arrow keys. This percentage is based on the size of the image.

Jump Cursor on Alert (Windows only)

When an alert box pops up in PDQUEST for Windows, you can set your cursor to automatically go to the OK button in the box by selecting Jump Cursor on Alert.

Click to Focus Tile

When multiple subwindows are displayed, you activate a subwindow by simply moving your cursor over it. To activate subwindows by clicking on

PDQUEST User Guide

them, select the Click to Focus Tile option. With this option selected, you first click on the subwindow to activate it, then perform the command.

Auto 'All Same'

Auto 'All Same' duplicates the zoom and positioning commands made in one subwindow in all displayed subwindows automatically.

2.7 File Paths

Edit > File Paths allows you to map the shortcut buttons in the Load dialog. You can set the shortcut paths to the locations on your hard disk or file server for the Matchsets, Images, Temporary, and Log directories.

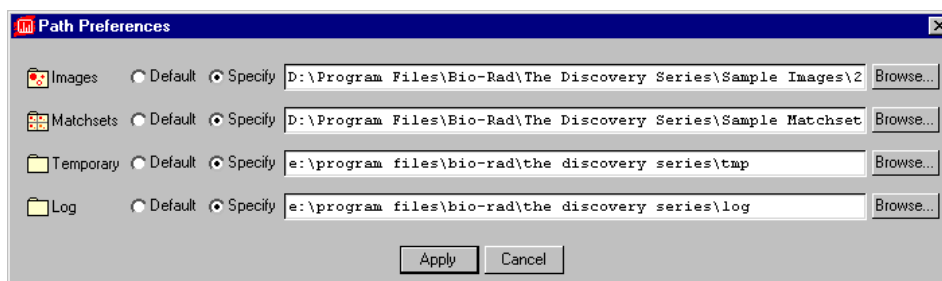


Fig.2-19. Path Preferences dialog box.

To change a shortcut, click on the Specify button and type the desired path in the text box, or click on the Browse button to search for the directory you want to map to.

2.8 Mouse and Mouse Macros

These functions are designed for backward compatibility with previous versions of PDQUEST. They are not recommended for new users and are not documented in this manual.

2.9 User Settings

If PDQUEST is on a workstation that has multiple users, each user can have his or her own set of PDQUEST preferences and settings.

In multiple-user situations, PDQUEST's preferences and settings are associated with individual user names. On a PC with Windows 95, 98, or NT 4.0, your user name is the name you use to log onto the computer. On a Macintosh, your user name is the Owner Name under the File Sharing control panel.

Note: If you do not log onto your PC under Windows 95/98 or do not have a Owner Name on your Macintosh, then you do not have a user name and your preferences and settings will be saved in a generic file.

PDQUEST User Guide

3. Viewing and Editing Images

This chapter describes the tools for magnifying and optimizing images and commands for configuring windows and subwindows. This chapter also describes the tools for cropping, flipping, and rotating images and filtering image noise.

These tools may be found on the View, Edit, Image, and Window menus.

Note: See section 2.2, Windows and Subwindows in PDQUEST, for background information on how windows and subwindows work.

3.1 Magnifying and Positioning Tools

The magnifying and positioning tools in PDQUEST are located on the View menu; some of these functions are also found on the main toolbar.

These commands will only change how the image is displayed on your computer screen. *They will not change the underlying data.*

PDQUEST User Guide

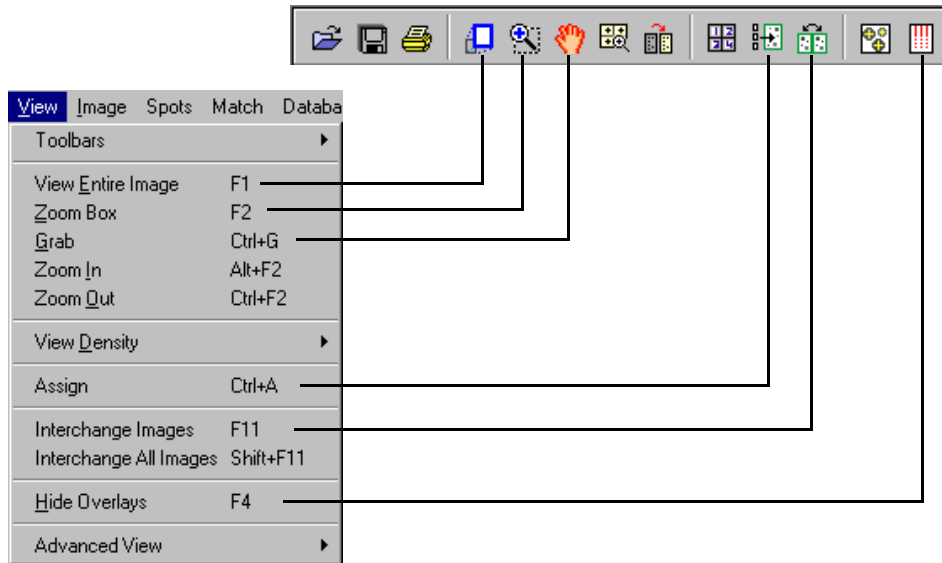


Fig.3-1. Viewing functions on View menu and main toolbar.

Zoom Box

Zoom Box is a mouse-assignable tool (see section 2.3) that allows you to select a small area of the image to magnify so that it fills the entire image window.

First, click on the Zoom Box button on your main toolbar or select View > Zoom Box. Your cursor arrow will change to a cross. Then drag the cursor on the image to enclose the area you want to magnify, and release the mouse button. The area of the image you selected will be magnified to fill the entire window.

Viewing and Editing Images

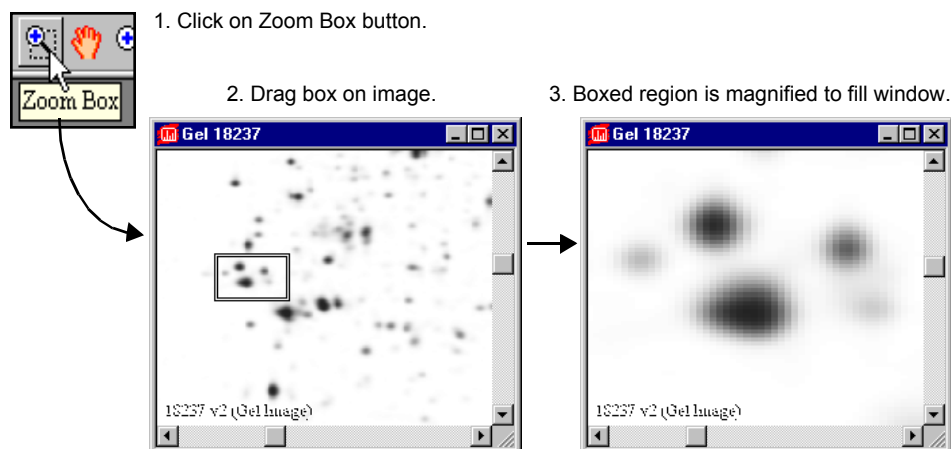


Fig.3-2. Zoom Box tool.

Zoom In/Zoom Out

These tools work like standard magnifying tools in other applications.

Click on the Zoom In or Zoom Out button on your main toolbar (or select from the View menu). Your cursor will change to a magnifying glass. Click on an area of the image to zoom in or zoom out a defined amount, determined by the setting in the Edit > Preferences dialog box.

Grab

This tool allows you to change the position of your image in the image window. Select View > Grab or click on the Grab button. Your cursor will change to a “hand” symbol. Dragging the Grab symbol on the image will move the image in any direction.

Arrow Keys

You can also move your image inside the image window by using the ARROW keys on your keyboard. Click on an arrow button to shift the image

PDQUEST User Guide

incrementally within the window. The amount the image shifts is determined by the Pan % setting in the Preferences dialog box.

View Entire Image

If you have magnified part of an image or moved part of an image out of view, View Entire Image returns to the original, full view of the image.

3.2 Hiding Overlays

To hide all plots, graphs, info boxes, and spot overlays and symbols on an image, select Hide Overlays from the main toolbar or View menu.

Note: Clicking once on Hide Overlays will conceal the overlays. Clicking twice will deassign any function that has been assigned to the mouse.

3.3 Advanced Viewing Tools

The Advanced View submenus contain a number of advanced positioning and viewing tools for use on windows and subwindows.

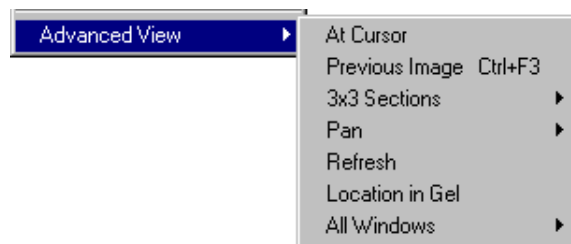


Fig.3-3. Advanced View submenu.

If you select At Cursor then click on the image, the image will be centered in the window at the point where you click.

Viewing and Editing Images

Previous Image (CTRL+F3) will “undo” the last zooming or positioning command.

3x3 Sections will magnify and position one of nine different quadrants of the image within the window. Select 3x3 Sections to display a list of available quadrants (Upper Left, Upper Middle, Middle Right, etc.).

Pan will move the image up, down, left, or right by increments within the window. You can also pan by using the arrow keys on the keyboard (see above). The amount the image shifts is determined by the Pan % setting in the Preferences dialog box.

Refresh refreshes the overlays and redisplay the image and overlays on the screen.

If you have magnified or repositioned a portion of the image, Location in Gel allows you to see which region of the image you are currently examining. Selecting Location in Gel displays two overlays in the window. The yellow overlay represents in the entire image, while the white overlay represents the region currently displayed within the window.

All Windows Viewing Commands

The All Windows submenu accesses advanced viewing features that can be applied to all displayed subwindows within a window. (See section 2.2 for background information on windows and subwindows.) Many of these commands are similar to those described above, except that they will change the display of all the subwindows at once. Unique commands are described below.

PDQUEST User Guide

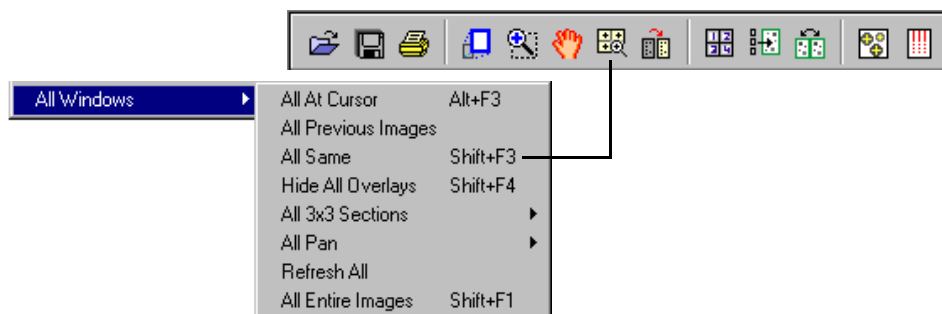


Fig.3-4. All Windows submenu.

Hide All Overlays conceals all the overlays in the displayed subwindows.

All Same applies the magnification and positioning of a single subwindow to all the subwindows. First select All Same from the submenu or main toolbar, then click on the subwindow whose view you want to duplicate.

Alternatively, select the subwindow you want to imitate and use the SHIFT+F3 keyboard command. See also Imitate Zoom, below.

Note: All Same will only work on comparable images with similar dimensions.

All Entire Images (SHIFT+F1) returns all displayed subwindows to their unmagnified, unrepositioned display.

3.4 Assigning and Interchanging Images

Because in PDQUEST you can load more than one image into a window or subwindow, you can use the Assign and Interchange commands to select the which loaded image you want to display. (See section 2.2 for background information on windows and subwindows.)

Assigning Images

To assign an image to a window, select View > Assign (CTRL+A) or click on the button on the main toolbar and click on the desired window or subwindow. This will open a pop-up box listing all the images currently loaded.

Viewing and Editing Images

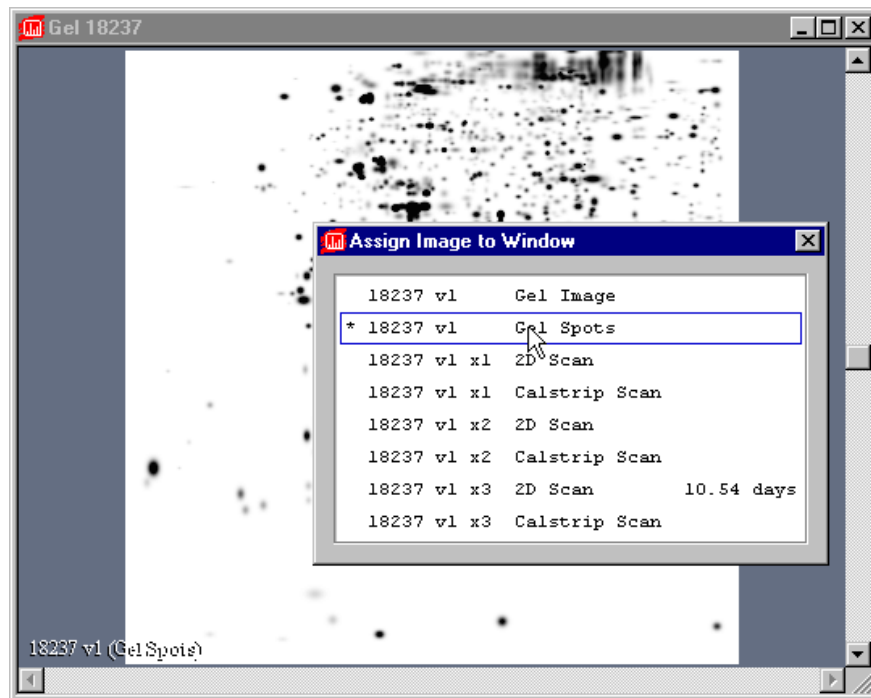


Fig.3-5. Example of an Assign pop-up list.

Select the image that you want to display. The image will appear in the window that you chose.

Auto-Assign opens a dialog box in which you can select the type of file to assign to a particular subwindow.

PDQUEST User Guide

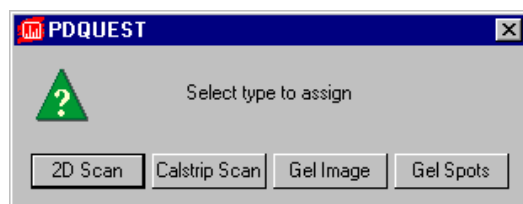


Fig.3-6. Auto-assign dialog box.

Interchanging Images

View > Advanced Assign > Interchange Images (F11) switches between Gel Images, Gel Spots images, and 2D scans within a particular subwindow.

View > Advanced Assign > Interchange All Images (SHIFT+F11) switches between Gel Images and Gel Spots images in all the displayed subwindows of a given matchset.

3.5 Density Tools

The Density tools (under the View > View Density submenu and on the Density Tools toolbar) are designed to provide a quick measure of the signal intensity of the data in your scan or Gel Image files.

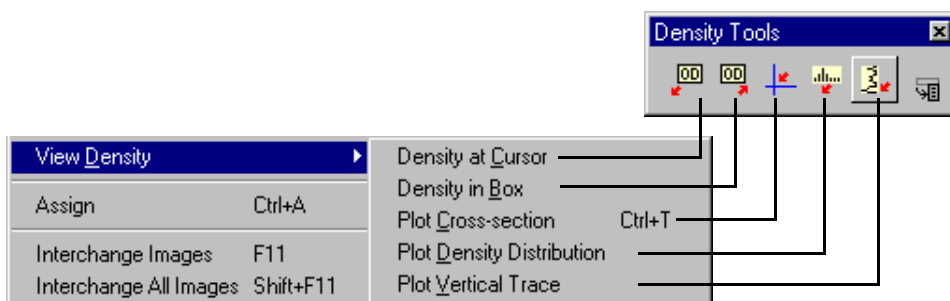


Fig.3-7. Density submenu and toolbar.

Viewing and Editing Images

Density at Cursor

Density at Cursor displays the signal intensity of the pixel on the image where you click your mouse. It also shows the average intensity for a 3 x 3 pixel box centered on that point.

To use this tool, click on the Density at Cursor button on the Density toolbar or select View > View Density > Density at Cursor. Then click on the point of interest on the image.

Density in Box

Density in Box displays the average and total intensity within a boxed region on the image. Select Density in Box from the Density toolbar or View menu, then define the region you want to measure by dragging your cursor across the image.

Plot Density Distribution

Plot Density Distribution displays a histogram of the signal intensity distribution for the part of the image displayed in the image window. The average intensity is marked in yellow on the histogram.

The histogram will appear along the right side of the image. Use the Zoom functions to display the histogram for magnified regions of your image.

Plot Cross-section

Plot Cross-section displays an intensity trace of a cross-section of the image centered on the point where you click your mouse. The horizontal trace is displayed along the top of the image, and the vertical trace is displayed along the side of the image.

The intensity at the point you clicked on is displayed, as is the maximum intensity along the lines of the cross-section. Dragging the mouse with this function selected will continuously update the display.

Note: The sampling width for the density traces is one pixel.

PDQUEST User Guide

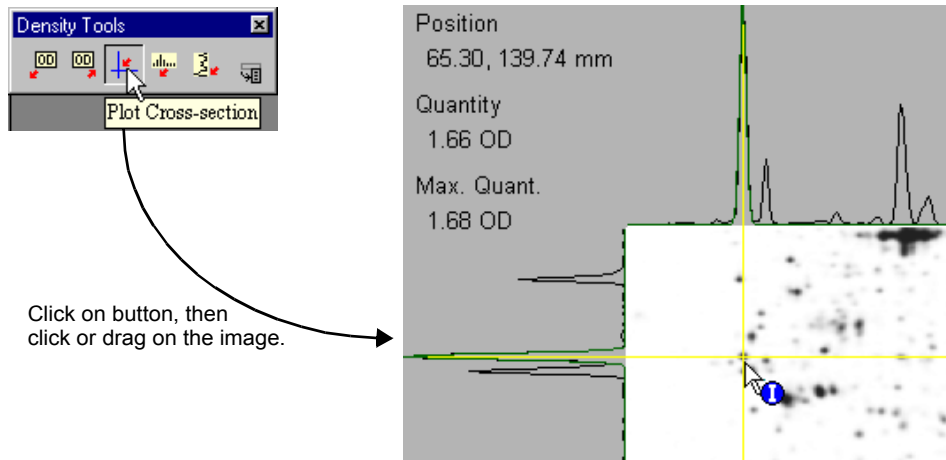


Fig.3-8. Plot Cross-section tool.

Plot Vertical Trace

Plot Vertical Trace plots an intensity trace of a vertical cross-section of the image centered on the point where you click your mouse. Dragging the mouse with this function selected will continuously update the display.

Note: The sampling width for the density traces is one pixel.

3.6 Windows and Subwindows Commands

The commands in the Window menu control how your windows and subwindows are configured.

Viewing and Editing Images

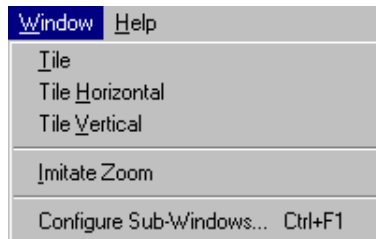


Fig.3-9. Window menu.

3.6.a Tiling Windows

If you have more than one image open, the Tile commands on the Window menu allow you to arrange your images neatly on the screen.

Tile will resize your windows and arrange them on the screen left to right and top to bottom.

Tile Vertical will resize your windows and arrange them side-by-side on the screen.

Tile Horizontal will resize your windows and stack them top-to-bottom on the screen.

3.6.b Imitate Zoom

Imitate Zoom applies the magnification and positioning of the active window or subwindow to all the displayed windows or subwindows.

First make sure that the window you want to imitate is active, then select Imitate Zoom from the Window menu. The other windows will change their view to match that of the active window

Note: Imitate Zoom will only work on comparable images with similar dimensions.

PDQUEST User Guide

3.6.c Configure Subwindows

If you have more than one image loaded into a window (i.e., in the case of scan sets or matchsets), you can control how the window is partitioned into subwindows using the Configure Subwindows command.

Select Configure Subwindows from the Window menu to open the Configure Subwindows dialog box.

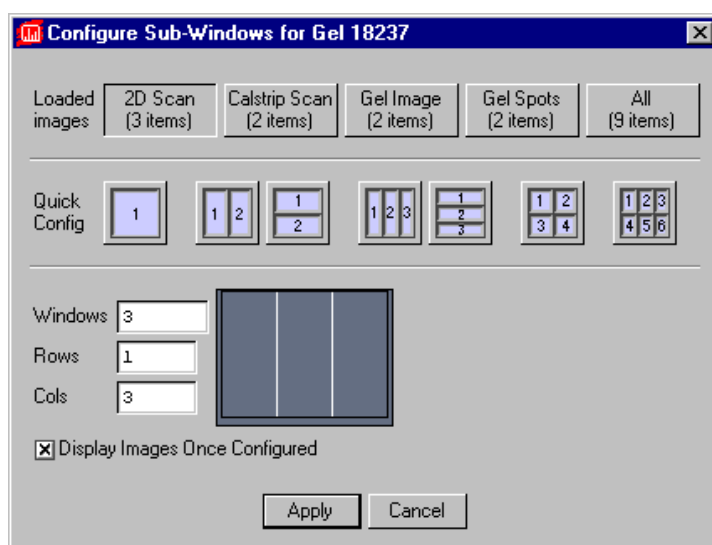


Fig.3-10. Configure Subwindows dialog box.

Loaded Images

The buttons next to the Loaded images prompt indicate the types of images that are loaded (2D scans, Gel Spots images, Gel Images, etc.) and indicate how many of each type are loaded. Click on a button to restrict your display to all the images of that type. (Click on All to not restrict by type.)

Viewing and Editing Images

Quick Config

The buttons next to Quick Config will immediately partition the window into the subwindow configuration depicted on the button. You do not need to click on the Apply button. (Note that if you have elected to display images of only one type using the Loaded Images button, only those types of images will be displayed in the subwindows.)

Custom Configuration

Alternatively, you can create the configuration of your choice by typing values in the Windows, Rows, and Cols fields. Windows specifies how many subwindows to create. Rows and Cols specify how to layer the subwindows.

Display Images Once Configured

Note: For most users, this checkbox should remain selected. Deselect this checkbox only if you have a slow computer and are displaying a lot of images.

If you have a matchset with many members and do not want to take time to refresh your subwindows after reconfiguring, deselect the Display Images Once Configured checkbox. Then, when you reconfigure your subwindows, the images will not immediately appear in the windows.

Click in a subwindow to make the image in that subwindow appear. Now perform a command in that window such as Zoom In. Select View > Advanced View > All Windows > All Same to reproduce that command in all the other subwindows and refresh the images in them with the appropriate zoom factor.

Apply

Click on the Apply button to implement your configuration selections.

PDQUEST User Guide

3.7 Colors

Edit > Colors opens a dialog box in which you can adjust the colors of the image, windows, buttons, etc.

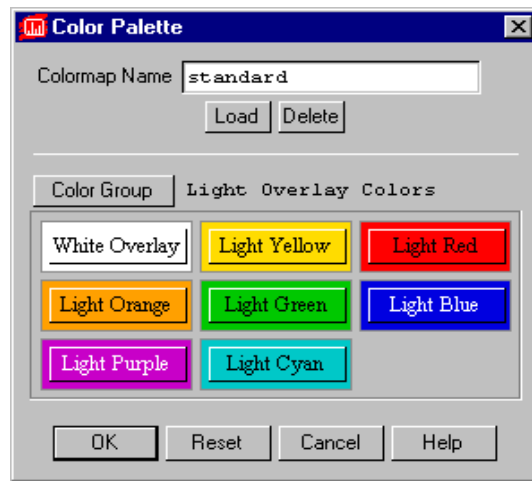


Fig.3-11. Colors dialog box.

Selecting a Color Group

Within the dialog box, the Color Group button allows you to select the colors of a particular group of objects (e.g., pop-up boxes, image colors, etc.). Click on the button to open the list of objects you can change.

Viewing and Editing Images

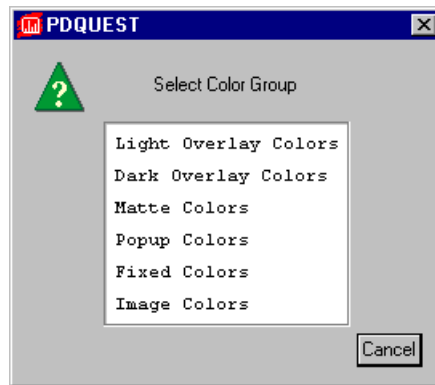


Fig.3-12. List of Color Groups.

Click on a color group in the list to select it.

Changing a Color

After you have selected the color group to change, click on the specific color button to change. The Color Edit dialog box will open, allowing you to adjust the red-green-blue (RGB) values of the color you selected.

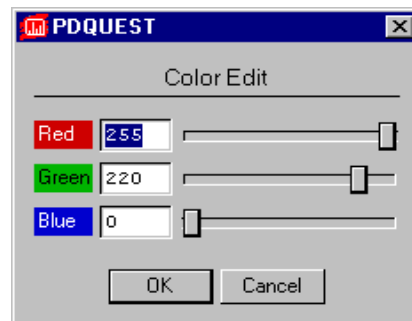


Fig.3-13. Color Edit dialog box.

Drag the sliders or enter a value in the fields. The color of the button will change with your adjustments.

PDQUEST User Guide

Saving/Selecting a Defined Set of Colors

After you have changed the colors within color groups, you can save these settings for future use on other images. The Colormap Name field displays the name of a defined set of colors and color groups. There are several standard colormaps included with PDQUEST, or you can create your own.

To select a predefined colormap, click on the Load button.

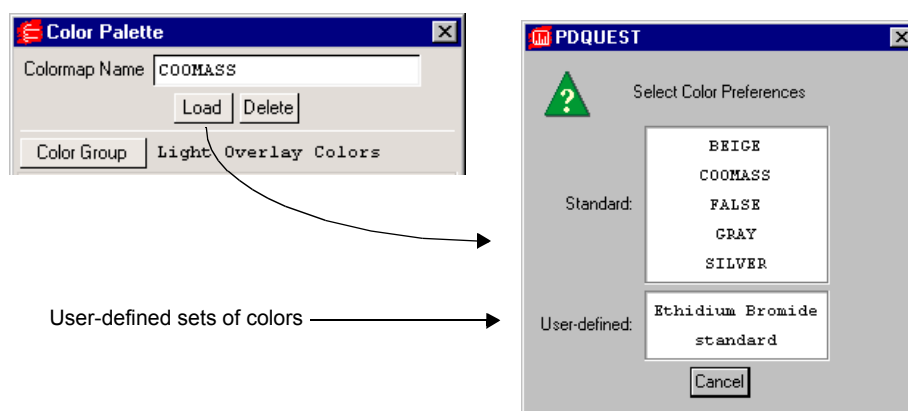


Fig.3-14. Selecting a Colormap.

From the list displayed, click on the set of colors you want to apply.

To create your own colormap, adjust the colors within the color groups as described above and type in a new colormap name. Click on OK to apply your changes.

To remove a colormap, click on the Delete button. Select the colormap to be deleted from the displayed list. A pop-up box will ask you to confirm the deletion.

If you change your mind about applying any changes you make, click on the Cancel button. If you want to return to the Standard colormap, click on the Reset button. All colors will be returned to their default values.

3.8 Transform

If features in your image are indistinct or fine details appear to be lost in background noise, you can use the Transform functions to optimize your displayed image. To open the Transform dialog box, select Image > Transform, or click on the Transform button on the main toolbar.

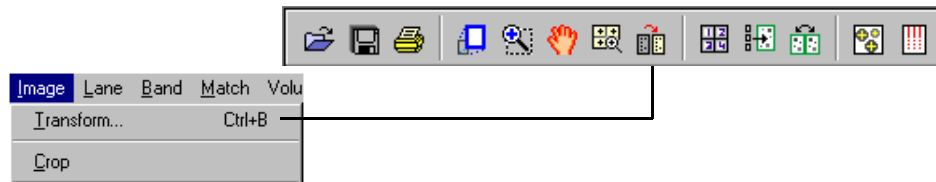


Fig.3-15. Transform command.

The Transform dialog box contains a Preview Window, a Frequency Distribution histogram, a Transform Plot, and three main methods of optimizing your image: Auto-scale, High and Low sliders, and a Gamma slider. You can use these controls to adjust the way the software transforms raw data into visual data.

PDQUEST User Guide

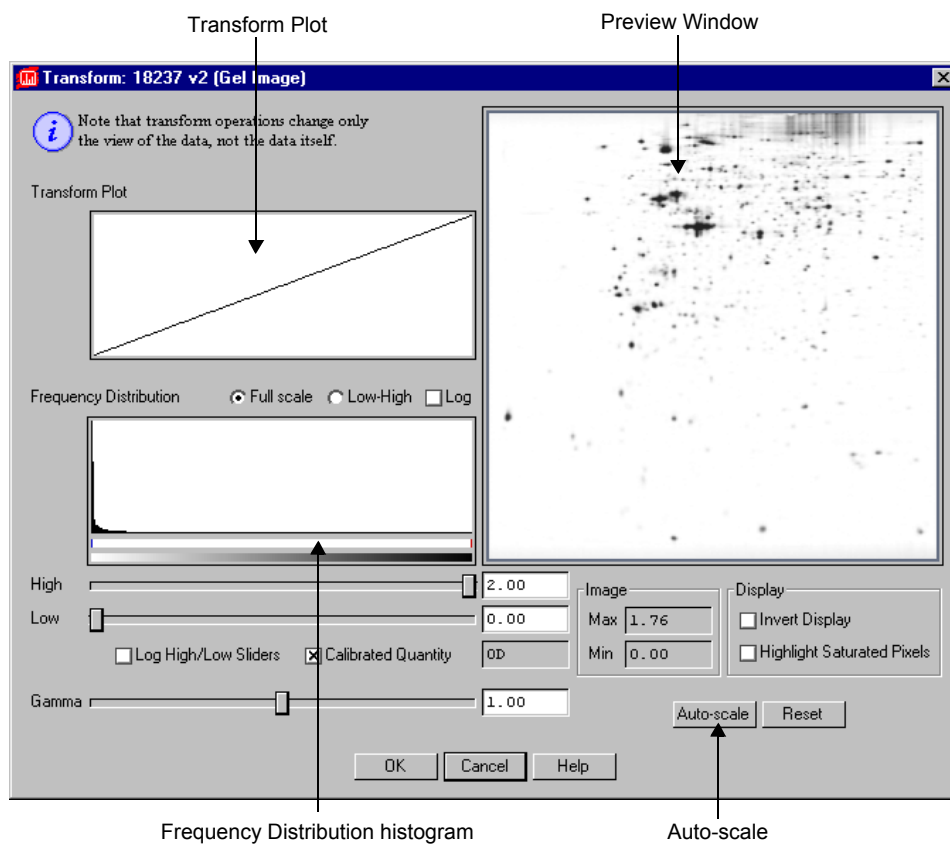


Fig.3-16. Transform dialog box.

Note: Any changes made with the Transform controls will only affect how the image is displayed on your screen. *They will not affect the underlying data.*

Viewing and Editing Images

3.8.a Transform Subwindows

Preview Window

The Preview Window in the Transform dialog shows a smaller view of the same image that is displayed in the main image window. Changes in the Transform controls are automatically reflected in the Preview Window. They are only applied to the main image when you click on OK.

You can use view tools like Zoom Box and View Entire Image in the Preview Window just as you can in the main image window, to focus on particular regions of interest. You can also use Grab and the ARROW keys to move the image within the Preview Window.

Frequency Distribution Histogram

The Frequency Distribution histogram shows the total data range in the image and the amount of data at each point in the range. In a typical scan, there is a signal spike at the left (“gray”) end of the histogram due to background noise.

Transform Plot

The Transform Plot is a logarithmic representation of how the raw pixel data are mapped to the pixels of your computer screen.

3.8.b Transform Controls

Auto-scale

Clicking on the Auto-scale button will optimize your image automatically. The lightest part of the image will be set to the minimum intensity (e.g., white), and the darkest will be set to the maximum intensity (e.g., black). This enhances minor variations in the image, making fine details easier to see. You can then “fine-tune” the display using the High, Low, and Gamma sliders described below.

PDQUEST User Guide

High/Low Sliders

If Auto-scale doesn't give you the appearance you want, you can use the High and Low sliders to redraw the image yourself. Dragging the High slider handle to the left will make weak signals appear darker. Dragging the Low slider handle to the right will reduce background noise.

As you drag the sliders, the slider markers on the Frequency Distribution histogram will move. Everything to the left of the Low marker will be remapped to minimum intensity, while everything to the right of the High marker will be remapped to maximum intensity. Using the histogram, you can position the markers at either end of the data range in your image, and use the low slider to cut off the "spike" of background noise.

You can also type specific High and Low values in the text boxes next to the sliders. Clicking anywhere on the slider bars will move the sliders incrementally.

Log High /Low Sliders changes the feedback from the slider handles, so that when you drag them, the slider markers move a shorter distance in the histogram. This allows for finer adjustments when your data is in a narrow range.

Viewing and Editing Images

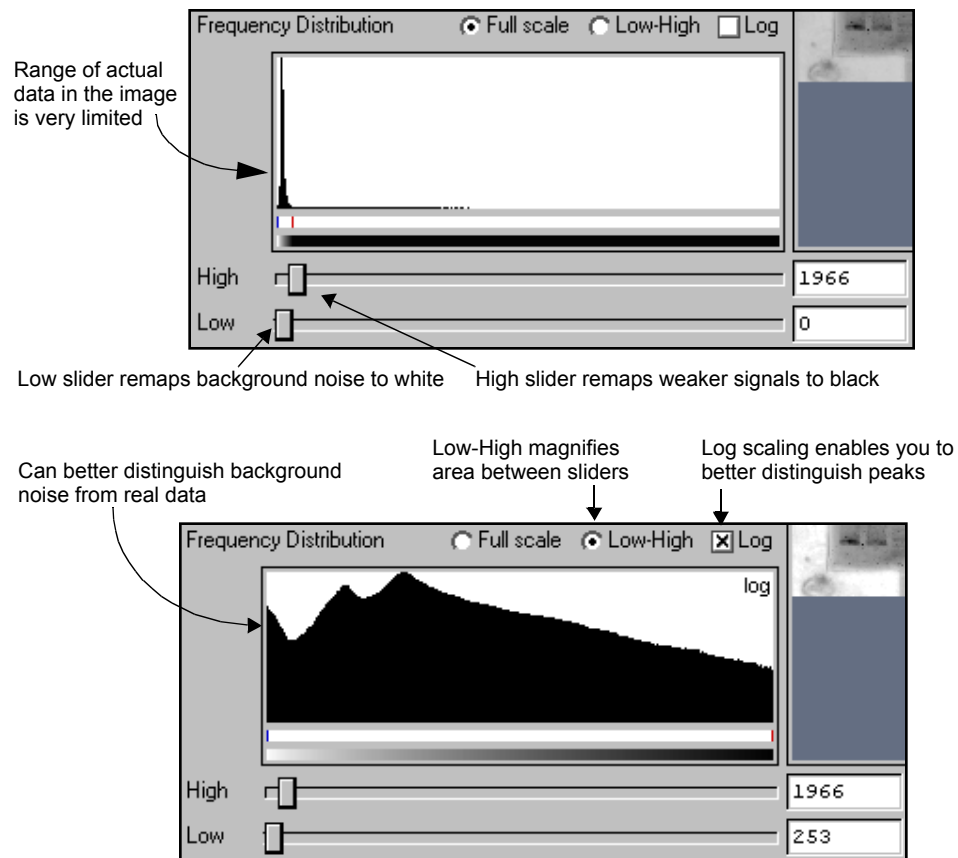


Fig.3-17. Two views of the Frequency Distribution histogram.

Gamma Slider

Some images may be more effectively visualized if their data are mapped to the computer screen in a nonlinear fashion. Adjusting the Gamma slider handle expands or compresses the contrast range at the dark or light end of the range, and this is reflected in the Transform Plot and Preview Window.

PDQUEST User Guide

3.8.c Other Features

Full Scale/Low-High

The Full Scale/Low-High option buttons adjust how the range of data in the image is displayed in the Frequency Distribution histogram and Transform Plot. They do not change how the data is displayed in the image window.

Selecting Full Scale adjusts the Frequency Distribution and Transform Plot displays so they show the full intensity range of the image.

Selecting From Low to High magnifies the range between the Low and High sliders. This makes it easier to view your data if it does not occupy the full intensity range of the image.

Log

The Log checkbox changes the way the data is displayed in the histogram so you can better discern subtle changes in signal intensity.

Image Max/Min and Units

Image Max and Min display the range of signal intensity in the image.

The image units are determined by the type of scanner used to create the image. In the case of densitometers, you can select Calibrated Quantity to display your image units in O.D.s.

Invert Display

The Invert Display checkbox flips light bands on a dark background to dark bands on a light background, and visa versa. Once, again, the actual data will not change—only the image.

Highlight Saturated Pixels

When the Highlight Saturated Pixels checkbox is selected, areas of the image with saturated signal intensity are highlighted in red.

Viewing and Editing Images

Reset

If at any time you want to return to an unmodified view of the scan data, click on Reset.

Transform All

If you have multiple subwindows open, you can transform all of them at once by clicking on the All button at the bottom of the dialog.

3.9 Cropping Images

If you are merging multiple exposures and/or matching spots between images, you will need to trim your scans to make them all the same size and shape. The tools for doing this are located on the Image menu.

The crop tools can also be used to reduce the memory size of a scan.

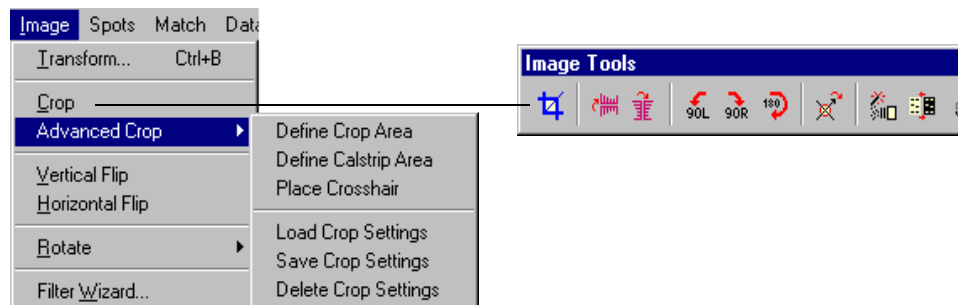


Fig.3-18. Crop commands.

The Crop command is automatically assigned after you complete a scan. Alternatively, you can click on the Crop button on the Image toolbar or select Image > Crop. Your cursor appearance will change to a Crop symbol.

PDQUEST User Guide

3.9.a Drawing a Crop Box

Define the region to be cropped by dragging the cursor across the image, creating a box. Everything outside the box will be deleted.

Information about the physical dimensions of the crop area (given in millimeters and number of pixels) is listed at the bottom of the crop box, as is information about the memory size of the image inside the crop area.

1. To *reposition* the crop box, position your cursor at the center of the box. The cursor appearance will change to a multidirectional arrow symbol. You can then drag the box to a new position.
2. To *resize* the box, position your cursor on a box border or corner. The cursor appearance will change to a bidirectional arrow. You can then drag that border or corner in or out, resizing the box.
3. To *redraw* the box, position your cursor outside the box. The cursor appearance will change back to the Crop tool, and you can draw another box, replacing the one you just drew.

3.9.b Advanced Crop Commands

If you are merging multiple exposures of the same gel (see Chapter 4) and / or matching spots between different gels, the images must all have the same dimensions, which means that they must be cropped to exactly the same size. The advanced crop commands allow you to save the crop box size and position from one scan and apply it to the others. They also allow you to crop images of gels and images of calstrips out of the same scan.

These commands are located on the Image > Advanced Crop submenu.

Define Crop Area/Define Calstrip Area

Define Crop Area allows you to draw a crop box as described in the previous section.

Define Calstrip Area allows you to draw a second crop box on the same scan. This is useful for scans the contain both an image of the gel and an image of the calstrip.

Viewing and Editing Images

First draw your crop box on the image, then select Define Calstrip Area and draw a second crop box around the calstrip portion of the image. When you crop the calstrip crop box (see next section), the resulting image will be labeled and saved as a calstrip scan.

Note: Be careful to define your regular crop box before using Define Calstrip Area. Also, use the Define Calstrip Area command to crop only images of calstrips, as the resulting cropped image will be defined as a calstrip image (i.e., with a .csc extension in Windows).

Place Crosshairs

Place Crosshair allows you to mark a reference point within the crop box. Select Place Crosshair, then click inside the box on an easily recognizable spot that is present on each exposure you scan. The crosshair enables you to correctly align the crop box to enclose the same region on all your images.

Crop Settings

With a crop box and crosshair displayed, select Save Crop Settings. Enter a name for your crop settings in the pop-up dialog box. Click on the Apply button to save your settings.

You can recall saved crop settings in other scan windows using the Load Crop Settings function. To delete crop settings, select Delete Crop Settings and choose the settings to be deleted from the list.

3.9.c Completing the Crop

Once you are satisfied with your crop box and/or calstrip box, position your cursor within the box slightly off-center. The cursor appearance will change to a scissors symbol. You can then click on the mouse to perform the crop.

A pop-up box will ask you whether you want to: (1) crop the original image, (2) save a copy of the area inside the crop box as a separate image, keeping the original image intact, or (3) cancel out of the cropping operation.

PDQUEST User Guide

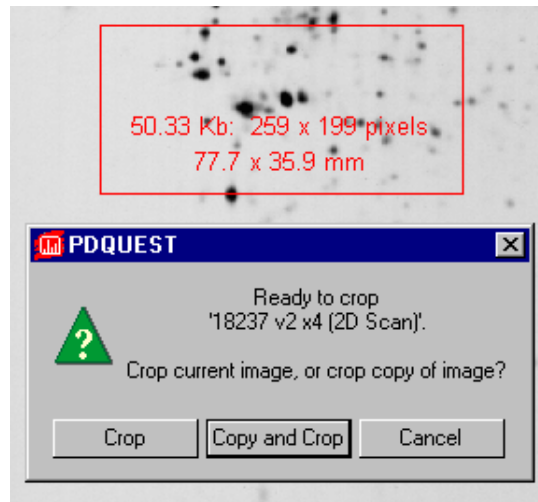


Fig.3-19. Crop box and pop-up Crop dialog.

If you select Crop, your cropped image will be displayed in the image window. If you select Crop and you have defined both a regular crop box and a calstrip crop box, both will be displayed as subwindows within the image window.

If you select the Copy and Crop button, a dialog box will open in which you can enter the name of the cropped image, its version number, exposure number, and /or file type. Click on OK, and the image window will display both the original image and the cropped image in subwindows.

3.10 Flipping and Rotating Images

The matching functions in PDQUEST will tolerate minor distortions in the shape and orientation of images. However, if your images are grossly misaligned, you may have to flip and /or rotate them. These commands are located on the Image menu and toolbar.

Viewing and Editing Images

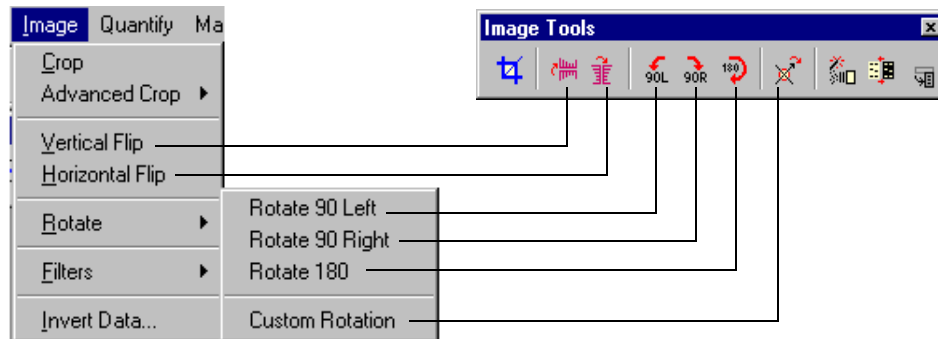


Fig.3-20. Flipping and rotating commands.

Note: These actions will erase any analysis you have performed on an image. You will be asked to confirm the flip / rotation before the command is executed.

Flipping

To flip the image right-to-left, select Horizontal Flip from the Image menu or toolbar. To flip the image top-to-bottom, select Vertical Flip.

90° Rotations

Select Rotate 90 Left, Rotate 90 Right, or Rotate 180 from the Image > Rotate menu or Image toolbar to perform the specified rotation. You will be asked to confirm your choice before the command is executed.

Custom Rotation

If you need to rotate your image in increments other than 90°, you can use the Custom Rotation command.

Select Custom Rotation from the Image > Rotate menu or Image toolbar. A green “plus” sign will appear next to your cursor. Click on the image you want to rotate and a circular overlay with an orange arrow will appear. A

PDQUEST User Guide

small dialog box also will open, indicating the angle of rotation in degrees and radians.

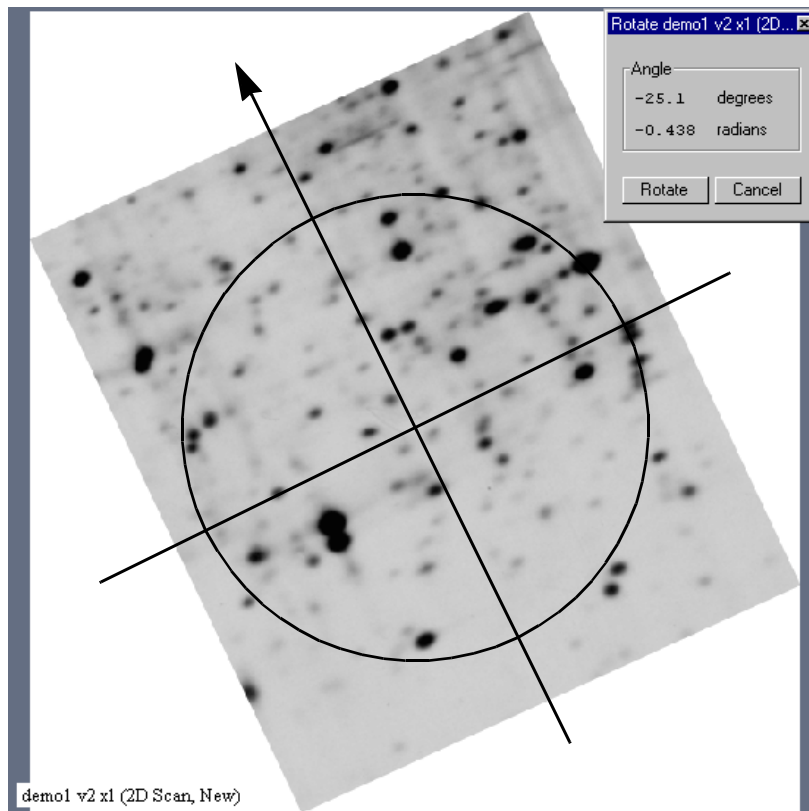


Fig.3-21. Custom rotation; the arrow points in the direction of the new top of the image.

To perform the rotation, position your cursor on the arrowhead and drag. As you drag, the arrow will rotate and the angle in the dialog box will change. Position the arrow so that it points in the direction of the new top of the image. You can “fine-tune” your rotation as much as you like.

Viewing and Editing Images

Note: If you want to center your arrow on a particular point on the image (e.g., to align along a particular lane), you can use the F3 key command. Position your cursor on the point on the image you want to center on, and press the F3 key. The image will move so that the center of the arrow and your cursor point are aligned.

To complete the rotation, click on the Rotate button in the small dialog box. A subwindow will open containing the rotated image. You will then have the option of renaming your new image and changing the version number.

If you are not satisfied with your rotated image, simply delete it and start over.

Note: Because an image is composed of square or rectangular pixels, Custom Rotation has to perform some minor smoothing on the image to turn it at a non-90° angle. In addition, any analysis performed on the image cannot be rotated and will be lost.

3.11 Filtering Images

Filtering is a process that removes small noise features on an image while leaving larger features (like data) relatively unaffected. PDQUEST has a wide range of filters for removing different types of noise from images. Depending on the nature of your data, you will probably need to use only one or two of the available filters. However, you should experiment with several different filters before selecting the ones that work best for your images.

The filtering commands are located on the Image menu and toolbar.

PDQUEST User Guide

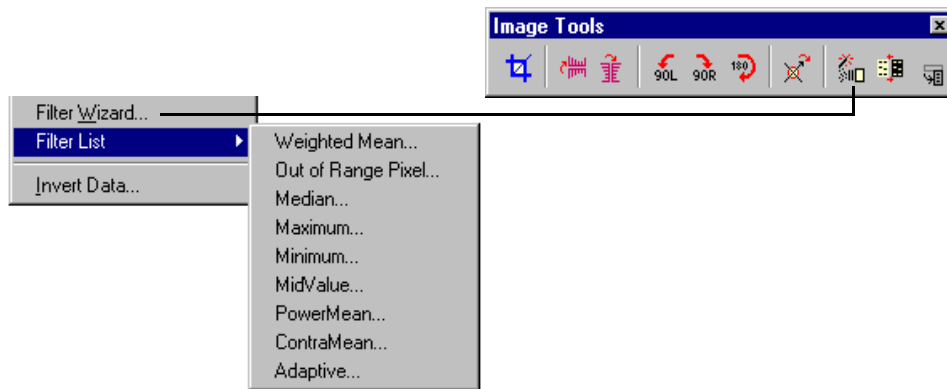


Fig.3-22. Filtering commands.

Note: Since filtering is an irreversible process, you will be asked if you want to create a copy of the original image before you filter. If you are experimenting with various filters, you should create copies of your image and compare them side-by-side. If you filter the original image and save it, *you cannot return to the original, unfiltered state.*

3.11.a Filter Wizard

The Filter Wizard is designed to guide you through the filter selection process. First, you will identify the type of noise in your image. Next, you will select the size of the filter to use on that noise. Finally, you will filter the image.

To open the Wizard, select Image > Filter Wizard or click on the Filter Wizard button on the Image Tools toolbar.

Viewing and Editing Images

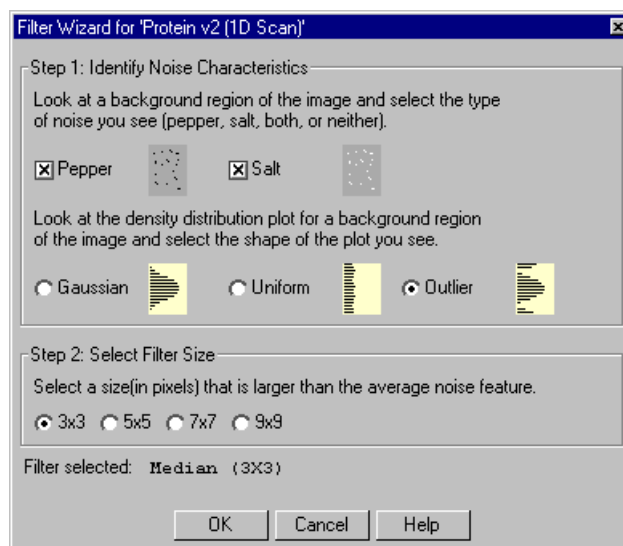


Fig.3-23. Filter Wizard dialog box.

Note: When using the Filter Wizard, it is useful to display a vertical plot of the density distribution in the image. The shape of this plot will help you determine the type of noise in the image. Select View > Plot Density > Plot Density Distribution; the plot will be displayed alongside the image.

Step I: Identify Noise Characteristics

The first step in the Wizard is to identify the type of noise in your image. The Wizard allows you to select from a variety of different noise types and combinations.

To determine the type of noise in your image, first magnify a background region of the image (i.e., a region where there is no data) using the Zoom Box command on the main toolbar. Note that as you magnify a region, the density distribution histogram changes to reflect only the signal within the region you have magnified.

After examining the magnified region, select one, both, or neither of the check boxes described below.

PDQUEST User Guide

- **Salt.** This type of noise appears as specks that are lighter than the surrounding background. The density distribution histogram of this type of noise displays noise peaks at the high end of the range. This type of noise is common in electronic cameras with malfunctioning pixels. It can also be caused by dust or lint in the imaging optics or scratches on photographic film. Salt is a type of outlier noise (see below).
- **Pepper.** This type of noise appears as specks that are darker than the surrounding background. The distribution histogram of this type of noise displays noise peaks at the low end of the range. Its causes are similar to those of salt noise. Pepper is a type of outlier noise (see below).

Next, select one of the following option buttons to describe additional features of your noise.

- **Gaussian.** The distribution histogram of this type of noise has a Gaussian profile, usually at the bottom of the data range. This type of noise is usually an electronic artifact created by cameras and sensors, or by a combination of independent unknown noise sources.
- **Uniform noise.** This type of noise appears in the histogram as a uniform layer of noise across the data range of the image.
- **Outlier noise.** This category of noise includes salt and pepper noise (see above). The distribution histogram of this type of noise displays noise peaks at the high and low ends of the range.

After you have identified the type of noise, go to Step 2.

Step 2: Select Filter Size

Image noise is filtered by means of a filtering window (or kernel), which is measured in pixels. This filtering window slides across the image, processing the pixels within it.

The available filter dimensions range from 3 x 3 pixels to 9 x 9 pixels. To select an appropriate size, magnify a background region of your image so that you can see the individual pixels. The filter size you select should be larger than the average noise feature but smaller than your data features.

Viewing and Editing Images

Note: A smaller filter will alter your image less than a larger filter. Large filters can result in better suppression of noise, but can also blur desirable features in the image.

Step 3: Begin Filtering

After you have completed your selections, the filter name and size will be displayed at the bottom of the Filter Wizard dialog box.

To begin filtering, click on the OK button. Because filtering is an irreversible process, a pop-up box will give you the option of filtering the original image, creating a copy of the image to filter, or cancelling out of this operation.

If you choose to Copy and Filter, you will be asked to enter a new name and / or version number for the new copy before the operation is performed. Once this information has been entered, the filtering operation will be performed.

3.11.b Selecting a Filter Directly

If you already know the type and size of filter you want, you can select it directly by selecting Image > Filters > Filter List. The pulldown menu includes all the filters available in the Discovery Series.

The types of filters in the Discovery Series are:

- **Weighted Mean.** This filter is useful for reducing Gaussian noise; it has little effect on salt and pepper. It calculates the weighted mean of the pixels within the filtering window and uses it to replace the value of the pixel being processed.
- **Out of Range Pixel.** This filter is useful for suppressing salt and pepper noise; its effect on Gaussian noise is minimal. This filter calculates the mean of the pixel values in the filtering window, including the pixel being processed. If the difference between the mean and the individual pixel value is above a certain threshold, then the individual value is replaced by the mean.
- **Median.** Also useful for suppressing salt and pepper noise, this filter has little effect on Gaussian noise. It calculates the median value of the pixels within the filtering window and uses it to replace the value of the pixel

PDQUEST User Guide

being processed. The median filter produces very little blurring if a small-sized window is selected.

- **Maximum.** This filter is useful for eliminating pepper noise in an image (it worsens the effect of salt noise). It replaces the value of the pixel being processed with the maximum value of the pixels within the filtering window.
- **Minimum.** This filter replaces the value of the pixel being processed with the minimum pixel value within the filtering window. This filter is useful for eliminating salt noise in an image (it worsens the effect of pepper).
- **MidValue.** This filter is useful for suppressing uniform noise within an image; however, it worsens the effect of pepper and salt. This filter replaces the value of the pixel being processed with the mean of the maximum and minimum pixel values within the filtering window.
- **Power Mean.** This filter is useful for suppressing salt and Gaussian noise within an image (it worsens the effect of pepper noise). It replaces the value of the pixel being processed with the power mean of the pixel values within the filtering window.
- **ContraMean.** This filter is useful for suppressing pepper and Gaussian noise within an image (it worsens the effect of salt). It replaces the value of the pixel being processed with the contra-harmonic mean of the pixel values within the filtering window.
- **Adaptive.** This filter is useful for suppressing Gaussian noise and salt and/or pepper within an image. If your image contains a mix of salt and pepper, select this filter.

To begin filtering, select a filter type from the pull-down list. A pop-up box will ask you to select a filter size.

Viewing and Editing Images

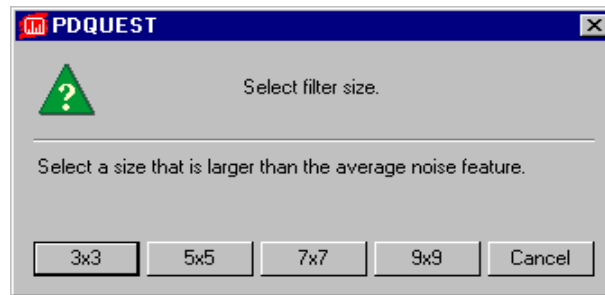


Fig.3-24. Selecting a filter size.

The available filter dimensions range from 3 x 3 pixels to 9 x 9 pixels. (See the previous section for guidance on selecting a size.)

Because filtering is an irreversible process, a pop-up box will give you the option of filtering the original image, creating a copy of the image to filter, or cancelling out of this operation.

If you choose to copy and filter, you will be asked to enter a new name and / or version number for the new copy before the operation is performed. Once this information has been entered, the filtering operation will be performed.

3.12 Invert Data

The Invert checkbox in the Transform dialog box (section 3.8.c) inverts only the appearance of the image. In some cases, however, you may want to invert the actual image data.

If your image has light spots on a dark background (i.e., the signal intensity of the background is greater than the signal intensity of the sample), you will need to use the Invert Data function before you can analyze the image.

This function is reversible, so if you change your mind, you can always switch it back.

To invert your image data, select Image > Invert Data. You may need to use the Transform function to adjust the appearance of your inverted image.

PDQUEST User Guide

4. Calibrating Gels

You can use calibration strips (calstrips) to calibrate your scans for accurate quantitation. See also Appendix J for a discussion of calibration.

Calstrips are required for both gel calibration and merging multiple exposures. If you aren't using calstrips, please proceed to the next chapter.

4.1 Calibrating Gels

The commands for calibrating gels are located on the Edit > Calibrate Gels submenu.

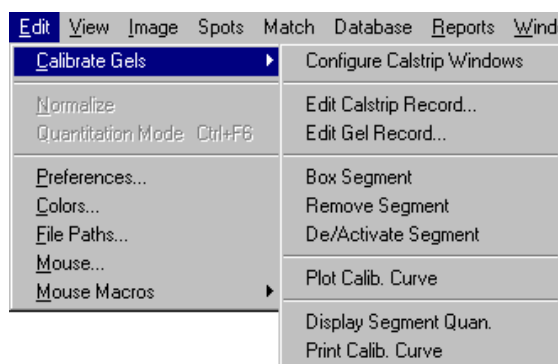


Fig.4-1. Edit > Calibrate Gels submenu.

Creating a Calstrip File

If calstrip images are currently part of your gel images, you will have to crop them and save them as separate files for use in PDQUEST. To do so, select Advanced Crop > Define Calstrip Area from the Image menu, draw a box around your calstrip, then crop it and save it as a calstrip file. See section 3.9.b for instructions.

PDQUEST User Guide

Configuring Calstrip Subwindows

If you are displaying multiple calstrips linked to multiple exposures (see section 4.2, Merging Multiple Exposures), Configure Calstrip Windows on the Edit > Calibrate Gels submenu will automatically arrange the calstrip images in subwindows.

4.1.a Calibration Strip Data

Values and units for each segment of your calstrip must be entered into a file called a Calstrip Record prior to calibration.

To do this, select Edit Calstrip Record from the Edit > Calibrate Gels submenu. The Calstrip Record dialog box will open.

Calibrating Gels

Calstrip Record Editor

Name:

Description:

Detection: ☒ Silver Stain ☐ Fluorography ☐ Syber Fluorescent Stain
☐ Coomassie Stain ☐ Autoradiography
☐ Gelcode Stain

Isotope: ☒ 3H ☒ 14C ☒ 32P ☐ 33P ☐ 35S ☐ 125I

Counting date:

Units: ☐ DPM ☒ OD ☐ PSL ☐ Other
☐ CPM ☐ HLU ☐ RD ☐ INT
☐ ODu ☐ RDu ☐ CNT

Number of segments: 0

Segment values in quantity/cm²:

1:		7:		13:		19:	
2:		8:		14:		20:	
3:		9:		15:		21:	
4:		10:		16:		22:	
5:		11:		17:		23:	
6:		12:		18:		24:	

Fig.4-2. Calstrip Record Editor dialog box.

Enter a title for the data record next to the Name prompt. You can enter descriptive information about the calstrip in the Description box.

Next, select the Detection method /stain that was used to make the calstrip.

For Fluorography and Autoradiography detection methods, select the Isotope used. Also enter the date the calstrip was made next to the Counting date prompt. The format in which the date should be entered is dd-mmm-yyyy (e.g., 01-Nov-1993 or 01-11-1993). The first three letters of the name of the month or the number of the month may be entered.

Specify the units for the quantitative values by clicking on one of the buttons next to the Units prompt or typing the units in the text box next to the Other button.

PDQUEST User Guide

Finally, enter the Segment values for each segment of the calstrip. These values should correspond to the calstrip segments as they appear from left to right (i.e., if the darkest segments are on the left, the values will decrease).

Range of Segment Values and Spot Quantities

The range of segment values should be approximately 2.0 to 2×10^4 . The minimum value the Calstrip Record will accept is 2.0. A typical maximum is 2×10^4 , but this value may vary. The range should span four orders of magnitude.

Typical spot quantities will range from 10^{-2} to 5×10^3 , depending on the size of the spot.

An example of segment values (DPM/cm²) for an ³⁵S calstrip are:

Segment #:	1. 278579.0	7. 14892.0	13. 827.0
	2. 159376.0	8. 9207.0	14. 516.0
	3. 100299.0	9. 5791.0	15. 315.0
	4. 67078.0	10. 3670.0	16. 134.0
	5. 40895.0	11. 2220.0	17. 93.0
	6. 24686.0	12. 1396.0	18. 44.0

Such a calstrip will produce spot quantitation values that typically range from 2.4 DPM to 502.6 DPM.

If you are working with Coomassie gels and you want to make a calstrip, you should begin with a stock solution and make approximately 10 serial dilutions, each of which will be used to make one segment of the calstrip.

For example, you might begin with a 50 µg / ml stock solution of BSA and make segments of the calstrip starting at 50 µg and progressing down to 5.0 ng in constant ratio steps.

In order to enter numerical values in the Calstrip Record that will produce interpretable spot quantitation values, convert all the numbers to ng (i.e., multiply all µg values by 10^3).

For the example mentioned above, the values would range from 50,000 ng to 5.0 ng.

Calibrating Gels

Saving/loading Calstrip Records

When you have completed the information in the dialog box, click on the Save button at the bottom of the box.

To access the information in an existing Calstrip Record, click on the Load button. A list of saved Calstrip Records will be displayed. Select the one whose data you wish to see by clicking on it. The data associated with that record will be displayed on the form.

To delete a saved Calstrip Record, click on the Delete button.

Click on the Cancel button to close the Calstrip Record Editor form.


4.1.b Scan Information

After the calibration data has been entered, you will also need to enter some information about the gel that goes with it.

Start by loading the 2-D and calstrip scans. Display the calstrip scans by using the Configure Calstrip Windows function (Edit > Calibrate Gels > Configure Calstrip Windows).

Click on Gel Record in the Calibrate Gels submenu and click on the calstrip scan. The Gel Record dialog box will open, in which you can enter information about the calstrip.

PDQUEST User Guide

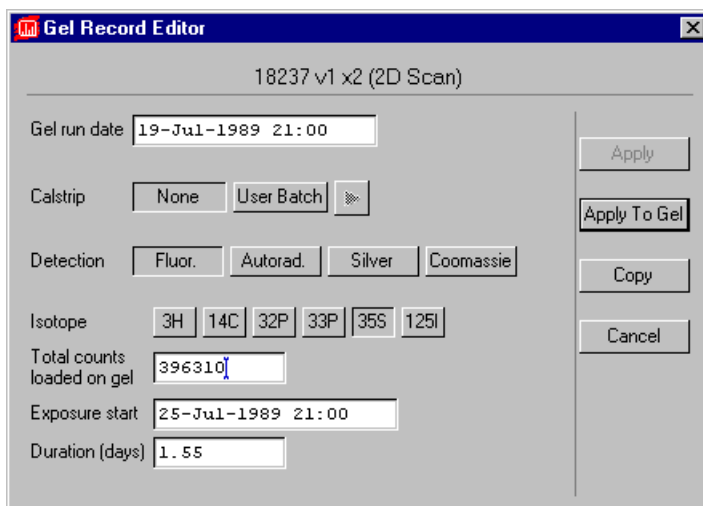

The screenshot shows a window titled "Gel Record Editor" with a close button (X) in the top right corner. The window contains several input fields and buttons. At the top, it says "18237 v1 x2 (2D Scan)". Below this, there are fields for "Gel run date" (19-Jul-1989 21:00), "Calstrip" (None, User Batch, and a right-pointing arrow button), "Detection" (Fluor., Autorad., Silver, Coomassie), "Isotope" (3H, 14C, 32P, 33P, 35S, 125I), "Total counts loaded on gel" (396310), "Exposure start" (25-Jul-1989 21:00), and "Duration (days)" (1.55). On the right side of the window, there are four buttons: "Apply", "Apply To Gel", "Copy", and "Cancel".

Fig.4-3. Gel Record Editor.

Enter the date on which the gel was run in the text box next to the Gel run date prompt.

Next to the Calstrip prompt, click on the button that corresponds to the kind of calibration associated with this gel. In most cases, you or someone in your lab has made the calstrips, so the button labeled User Batch is the appropriate choice.

Next to the User Batch prompt, click on the arrow. From the list displayed, select the calibration record name for the calstrip. The Detection and Isotope fields are automatically obtained from the specified Calstrip Record.

Click on the appropriate button to indicate the type of detection used for your gel, if not already defined from the batch file.

Click on the Isotope button to indicate the isotope used in the gel. This is usually, though not necessarily, the same isotope as the one used in the calstrip.

Next, enter the number of counts loaded on the gel.

Calibrating Gels

Enter the day on which you began exposing the gel to film in the text box next to the Exposure start prompt.

Note: Most calibration problems arise from incorrectly entering the dates used to calculate radioactive decay factors for the isotopes. Wrong values can make the calstrip appear to have decayed to nothing.

Next to Duration (days), enter the number of days that the film was exposed.

Once this data has been entered, you will want to apply it to the 2-D scan and the calstrip scan. To do so, click on the Apply button.

If you want to apply different calibration data to your 2-D scan than to your calstrip scan (if, for example, your 2-D scan and calstrip scans use different isotopes), use the Apply Here button to apply the gel record to the currently displayed image only (2-D or calstrip scan).

Copying Gel Record Information

If you have scanned several exposures of the same gel, most of the calibration data will be the same for all of them (except for the exposure start date and length). Instead of entering all the information for each exposure, you can enter it once, and then copy it between exposures using the Copy button.

Clicking on Copy will open a list of available gel records. Select the record that you would like to apply to the current image and click the mouse.

Click on the Cancel button to exit the operation.

4.1.c Sampling the Calibration Strip Segments

The next step in calibrating your gels involves sampling the densities in the segments of the calstrip to establish a correlation between the sampled intensities recorded by the scanner and the counts entered into the Calstrip Record.

If you have not already done so, display the calstrips. If you have gels with multiple exposures, use the Configure Calstrip Windows function to automatically divide the screen into enough windows to display all the calstrips for a given gel.

PDQUEST User Guide

Flip Calstrip (if necessary)

If the orientation of the scanned calstrip does not correspond to the order in which the segment activities are entered in the Calstrip Record, a message will be displayed, indicating that the calibration density decreases while the counts increase (or vice versa).

To change the orientation of the calstrip so it corresponds to the Calstrip Record segment activities, use the Horizontal Flip function (under the Image menu). If the dark segments were on the left, they will now be on the right, and vice versa. Again, we recommend that you enter your data from highest to lowest activity and display the calibration strip with the darkest segment on the left.

Box Segments in a Calstrip

Once the calstrip is in the correct orientation, go back to the Edit > Calibrate Gels submenu and select Box Segment.

Go to the first (leftmost) segment of your calstrip. Sampling boxes must be drawn beginning with the leftmost segment of the calstrip and moving right.

With Box Segment selected, drag a box on the calstrip. When you release the mouse button, the area will be sampled.

Note: When defining your box, be sure to avoid sampling regions that are very close to the borders of the segment or any irregular areas (i.e., areas that are overly dark or light).

The number of each segment will be displayed in the upper left corner of each sampling box.

Above the box will be a number between 0.00 and 2.00. This represents the O.D. value for the segment.

Beneath each box is the number of counts associated with that segment from the Calstrip Record.

If you make a mistake, simply redraw a new box over the old one. Repeat this procedure for all the segments of the calstrip.

Calibrating Gels

Note: All segments entered in the calstrip record must be boxed starting from the left. The leftmost segment corresponds to Segment 1 on the Calstrip Record. Your calstrips and calstrip data should be designed so that either (1) the darkest segments are on the left and the count values decrease, or (2) the lightest values are on the left and the count values increase.

You can draw more boxes than the number of segments defined in the Calstrip Record. However, they are automatically and permanently deactivated.

4.1.d Plot the Calibration Curve

Once you have created the sampling boxes, select Plot Calib. Curve to display the calibration curve.

A graph will be displayed that correlates the optical density (usually between 0.00 and 2.00 O.D.) with DPMs or other specified units.

The Y axis represents the optical density (0.00–2.00) in each sampled segment, while the X axis represents the corresponding DPM or other unit per image unit area calculated from the calibration record data.

PDQUEST User Guide

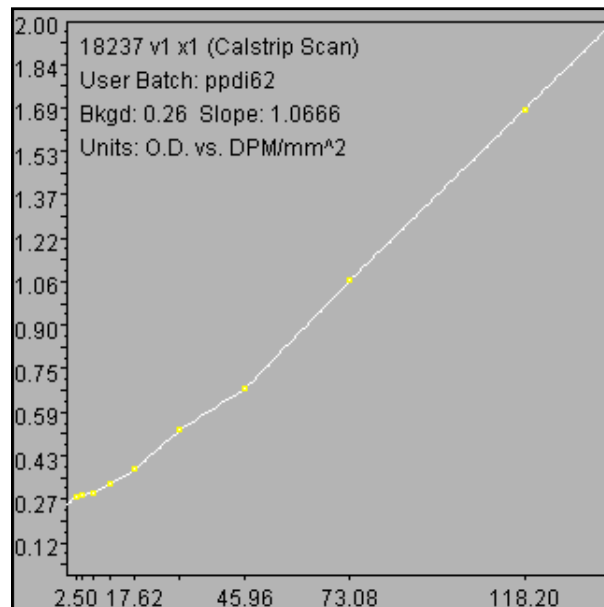


Fig.4-4. Example of a calibration curve.

Calibration Curve Problems

Find errant data points and their corresponding segments using the function Display Segment Quan.

Above each sampling box will be displayed the average O.D. value within the box, and just below the box will be displayed the corresponding DPM (or other specified unit). This should tell you which sampling boxes correspond to errant data points on the graph.

Usually, problems arise when a box straddles the boundary between two adjacent segments or the density is very uneven within the segment due to bleeding from an adjacent segment. These sorts of problems should be remedied by redoing the sampling boxes.

Use the function Remove Segment to remove an errant sampling box.

Calibrating Gels

Draw a new sampling box on the segment using Box Segment, as described above.

Place the new box carefully to avoid repeating the problem. For example, keep the box away from segment boundaries. If the segment shows bleeding, place the smallest box possible in an area of the segment with representative O.D. Try not to sample blotchy areas.

Replot the calibration curve using the function Plot Calib. Curve.

If there is still a problem, you can opt not to include that particular segment in the calculation of the calibration curve. Select De / Activate Segment and click on the problem segment. The segment will still exist, so its counts in the corresponding Calstrip Record do not need to be changed or deleted. However, the segment will not be included in the calculation of the calibration curve.

If you change your mind and want to reinclude the segment, select De / Activate Segment again and click on the segment. The segment will be reactivated and included in calibration curve calculations.

4.1.e Repeat for Each Calstrip

If you have multiple film exposures of the same gel, repeat the steps outlined above for each exposure. The calibration strip batch number for each film exposure of a single gel should be the same since the same physical calibration strip was used to make each exposure.

4.2 Merging Multiple Exposures

If you are working with 8-bit images and are using calstrips on your scans, you can merge multiple exposures of the same gel into a single image with a greater dynamic range.

Note: You can only merge 8-bit images with calstrips. If you are not working with 8-bit images or are not using calstrips, you cannot merge exposures.

PDQUEST User Guide

To merge multiple exposures, make sure all your exposures and calstrips are loaded into the same image window, then simply auto-detect spots as you would normally (see the following chapter).

Creating Different Exposure Files

To identify different gel scans as different exposures of the same gel, open the scan files, select Change Version / Exposure from the File menu, and enter consecutive exposure numbers for the different scans. Make sure that the root names of the scan files are the same.

Note: You can give different exposures the same root name when you are scanning them, or you can change them later using the Save As command. After you change the exposure numbers, you will end up with files named RootName x1, RootName x2, RootName x3, etc.

Why Merge Exposures?

Because 8-bit scanners have a limited dynamic range, dark spots may saturate on long exposures, while faint spots may be missed on short exposures, and either extreme can be missed on medium duration exposures. Therefore, it is frequently necessary to combine the quantitative data from a number of 8-bit exposures of the same gel exposed for different lengths of time into a “merged” synthetic image. Merging creates such an image automatically as you spot detect.

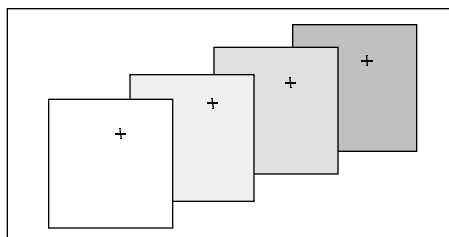


Fig.4-5. Diagrammatic representation of merging.

Calibrating Gels

Cropping and Aligning Exposures of the Same Gel for Merging

Before you can merge multiple exposures, you must crop and align them so that they are exactly the same size and position. To do this, you will use the Advanced Crop commands on the Image menu. See section 3.9.b for instructions.

PDQUEST User Guide

5. Detecting and Editing Spots

If you are satisfied with your calibration data, or if you are working with uncalibrated gels, you are now ready to detect and quantitate the protein spots on your gel. This chapter describes how to use PDQUEST to detect and quantify spots in 2-D protein gels that have been scanned and prepared according to the instructions in the previous chapters.

5.1 Selecting the Detection Parameters

The Spot Detection Wizard simplifies the process of creating and saving different spot detection parameters for different types of gel scans. In general, you can use a specific set of parameters to detect spots on all gels that use the same stain and have been run and scanned in the same manner.

To select the detection parameters for a particular type of gel, make sure that a representative image is open and active, then select Spot Detection Wizard from the Spots menu or PDQUEST quick guide.

PDQUEST User Guide

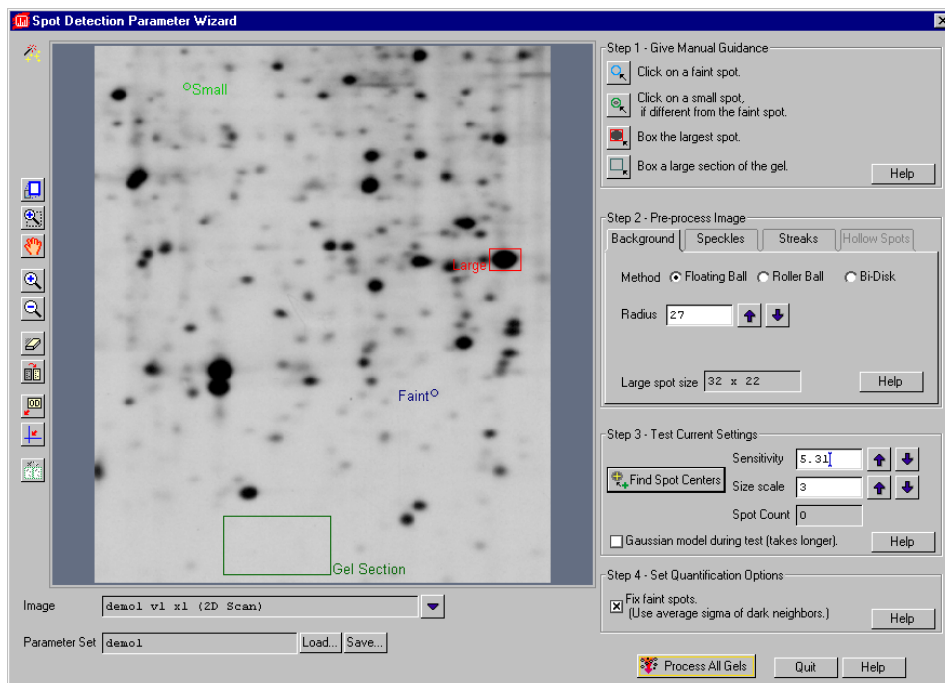


Fig.5-1. Spot Detection Wizard.

The Wizard has a preview window showing the currently loaded image and a control panel that is designed to guide you through the parameter selection process.

Loading Other Images

You can select parameters for the currently displayed image, or you can load a different open image by clicking on the arrow button next to the Image field under the preview window.

Detecting and Editing Spots

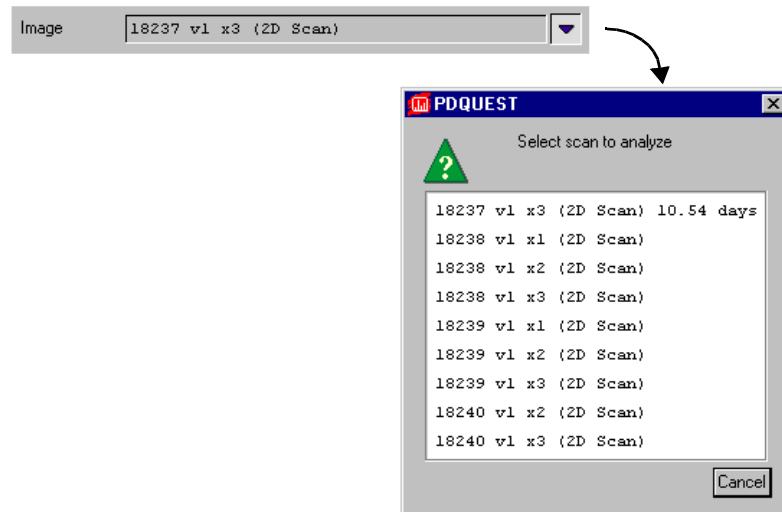


Fig.5-2. Selecting a new image in the Spot Detection Wizard.

A pop-up box will display a list of open images. Click on the image you want to load into the Wizard.

Note: If you load a new image into the Spot Detection Wizard without finding spot centers or saving the parameters in the current image, these parameters will be lost. Be sure to save any parameters you want to keep (see section 5.1.e) before loading new images.

Loading Saved Spot Detection Parameters

If you have saved spot detection parameters (see section 5.1.e) that you want to test on the current image, you can load them into the Wizard by clicking on the Load button next to the Parameter Set field below the preview window.

A pop-up box will display a list of saved parameters. Click on the parameter set you want to load. When the parameters are loaded, you can skip to Step 3, Test Current Settings.

PDQUEST User Guide

5.1.a Step 1. Give Manual Guidance

In this step, you mark features on your gel scan that automatically guide the selection of the proper spot detection parameters.

1. Click on the Faintest Spot button, then click on the faintest spot on the scan image that you want to detect. This will set the sensitivity of the spot detection.
2. If the smallest spot you want to detect is smaller than the faintest spot, click on the Smallest Spot button, then click on the smallest spot. This will set the size scale of the spot detection. In most cases, the smallest and faintest spots are the same, and you can skip this step.
3. Click on the Largest Spot button, then draw a box around the largest spot on the image that you want to detect. This sets the radius of the background subtraction rolling ball and the streak removal rolling disk.
4. Finally, click on the Box Large Section button, then draw a box around a background region of the image. You should box the largest region you can find that contains no spots or streaks. This sets the correct type of noise filter.

Detecting and Editing Spots

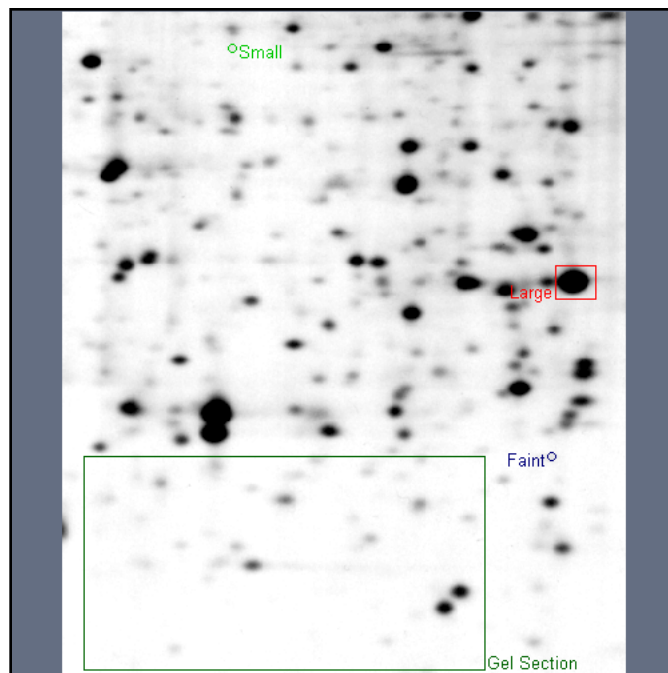


Fig.5-3. Image with large, small, and faint spots and background region defined.

5.1.b Step 2. Preprocess Image

Most of the settings in Step 2 are set automatically based on your selections in Step 1. They specify the levels of background subtraction, noise smoothing, and streak removal.

Note: Initially, you should accept the default settings under the Background and Speckles tabs and proceed to Step 3. If your image has streaks, select streak removal under the Streaks tab. If Find Spot Centers in Step 3 results in too much or too little background subtraction, noise smoothing, or streak removal, return to Step 2 and adjust these settings.

PDQUEST User Guide

Background

Click on the Background tab to access the settings for subtracting background intensity from the image.

Background subtraction is required for spot detection. The default method, Floating Ball, and default ball radius size are highly recommended. Radius size is calculated automatically based on the largest spot you want to detect (selected in Step 1 and displayed in the field at the bottom of the tab).

If Floating Ball results in too much background subtraction or undesirable image artifacts, you can select the Roller Ball method. Typically, the roller ball method results in more image artifacts than the floating ball.

The Bi-Disk method is slightly faster than the floating or roller balls, and will perform adequately for images with very smooth, even backgrounds. However, bi-disk will perform less well for images with uneven backgrounds.

Speckles

Click on the Speckles tab to access the settings for removing random image noise.

Image smoothing in your scans is optional but highly recommended. The default state of the Perform Smoothing checkbox is checked. To deactivate this function, deselect the checkbox by clicking in it.

The Filter Type is initially set based on the section of the image you boxed in Step 1. The Kernel Size is set based on the faintest/smallest spot you want to detect (also selected in Step 1).

Note: You should first attempt to find spot centers (Step 3) using these default settings before attempting to adjust them.

If you know the type of noise in your images, you can select a different filter type and kernel size. Section 3.11, Filtering Images, provides additional guidance on filter selection.

- **Median.** This filter is useful for suppressing salt and pepper noise. It has little effect on Gaussian noise. This filter produces very little blurring if a small-sized window is selected.

Detecting and Editing Spots

- **Weighted Mean.** This filter is useful for reducing Gaussian noise. It has little effect on salt and pepper noise.
- **Power Mean.** This filter is useful for suppressing salt and Gaussian noise within an image (it worsens the effect of pepper noise).
- **ContraMean.** This filter is useful for suppressing pepper and Gaussian noise within an image (it worsens the effect of salt).

To adjust the Kernel Size, magnify a background region of your image using the Zoom Box tool so that you can see the individual pixels. The filter size you select should be larger than the average noise feature but smaller than the smallest spots.

Streaks

Protein gels frequently display streaking in areas of high spot concentration. Click on the Streaks tab to access the streak removal settings.

To remove vertical and / or horizontal streaks from your image, select the appropriate checkbox(es). The disk radius size is determined automatically based on the size of the largest spot you want to detect (selected in Step 1 and displayed in the field at the bottom of the tab).

The default disk radius size is conservative. If, after you find spot centers (Step 3), you would like to remove more streaking from the image, reduce the size of each disk radius by about half and find spot centers again.

5.1.c Step 3. Test Current Settings

Click on the Find Spot Centers button to identify the spot centers in your image based on the parameters you selected in Steps 1 and 2.

Note: Find Spot Centers provides a quick view of how the current settings will work on your image. It is designed so that you can find spot centers, adjust the parameters, and repeat the operation. After you have determined the optimal settings using this procedure, you are ready to save the parameters and fully process your image for analysis (see below).

PDQUEST User Guide

When you click on Find Spot Centers, noise smoothing, background subtraction, and streak removal (if selected) will be performed, and the spot centers will appear marked with crosshairs.

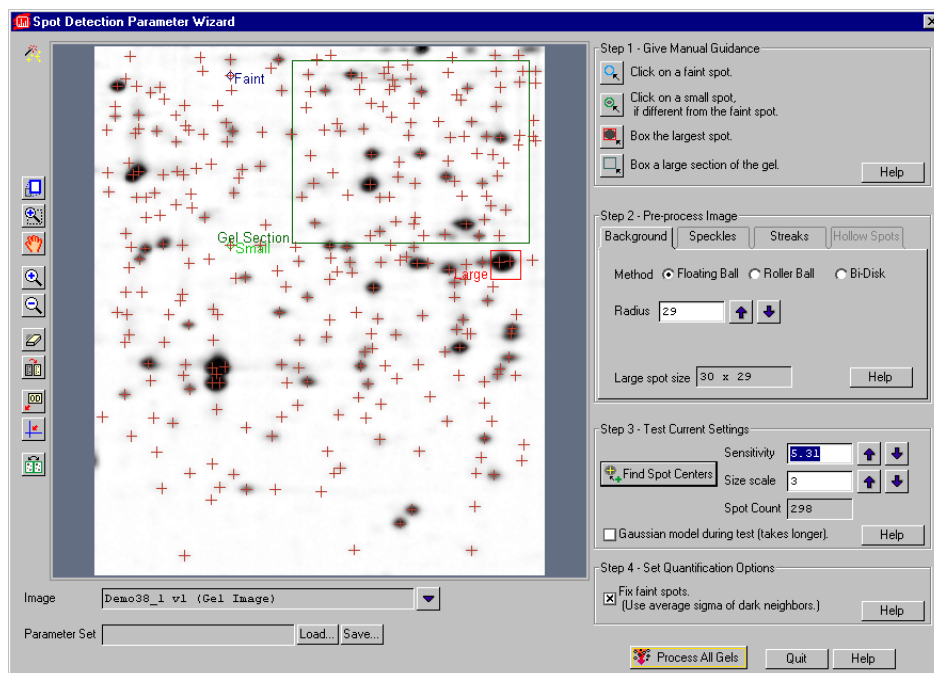


Fig.5-4. Spot crosshairs marked in the Spot Detection Wizard.

Review the results of spot detection by studying the crosshairs on the image. The number of spots detected is displayed in the Spot Count field.

If too much or too little background subtraction, noise smoothing, or streak removal has been performed, return to Step 2 and adjust the settings as described above.

Gaussian Model During Test

If this checkbox is selected, two-dimensional Gaussian fitting will be performed on all spots when you find spot centers. It will take longer to

Detecting and Editing Spots

detect spots; however, spot detection will be more accurate for overlapping spots.

Note: Fewer spots may be found if Gaussian Model During Test is selected, because large spots that would be identified as overlapping spot clusters can be correctly resolved with Gaussian modeling.

Sensitivity and Size Scale Adjustments

Adjustments to these settings will automatically redetect spots.

The Sensitivity setting determines the minimum signal intensity in the image that will be detected as a spot. The higher the sensitivity value, the more spots will be detected. The initial sensitivity value is set based on the faint spot selected in Step 1.

If too few spots are detected, increase the sensitivity by typing a number directly in the field or clicking on the adjustment arrows. If faint noise features or background artifacts have been detected as spots, you can reduce the sensitivity setting or perform additional speckle removal (Step 2).

The Size Scale field helps distinguish trends in signal intensity (spots) from random intensity fluctuations (noise). It uses the size of objects in the image (in pixels) to determine the nature of those objects. The initial size scale setting is based on the small spot you selected in Step 1.

If a gel image has noticeable noise features, a larger size scale may be preferable. If smoothing has adequately removed noise, small values are preferred so that small spots will be detected. You can enter odd values from 3 to 9 in the field or adjust the value using the arrow buttons.

5.1.d Step 4. Set Quantification Options

Fix Faint Spots is recommended for more accurate quantitation of faint spots. This corrects spots below a certain intensity over background (e.g., 0.2 O.D.) by making their X and Y sigma values the same as those of the spots in the same Y-coordinate region of the film. This overcomes the difficulty of obtaining an accurate fit when very few pixels are available.

PDQUEST User Guide

The Fix Faint Spots checkbox is selected by default. If you do not want to perform this operation, deselect the box by clicking in it.

5.1.e Saving the Parameters

After you have determined the best detection settings for your image, you are ready to save the parameters for use on this and other images.

Next to the Parameter Set field, click on the Save button, and enter an appropriate name for the parameter set (e.g., the type of stain used in the gel) in the pop-up box. Click on OK to save the set.

The parameter set will now be available for processing all gel scans.

Processing Gels

Click on the Process All Gels button to close the Spot Detection Wizard and open the Auto-Detect Spots control panel.

If you have not saved the current detection parameters, you will be prompted to do so before you exit the Wizard.

The Quit button will close the Wizard without opening the Auto-Detect Spots form. If you have not found spot centers, any unsaved parameters will be lost.

5.2 Detecting Spots on All Gels

After you have created and saved one or more sets of detection parameters using the Spot Detection Wizard, you are ready to detect spots on your images.

Note: Only gels that have been stained and scanned in the same manner should use the same detection parameters.

If you click on Process All Gels in the Wizard, the Auto-Detect Spots control panel will be displayed. You can also open the control panel directly by selecting Detect Spots on All Gels from the Spots menu.

Detecting and Editing Spots

Note: If you are detecting spots on calibrated gels, the gel scans and corresponding calstrip scans must both be loaded for proper calibration. Also, if you have multiple exposures of the same gel, make sure they are loaded for proper merging.

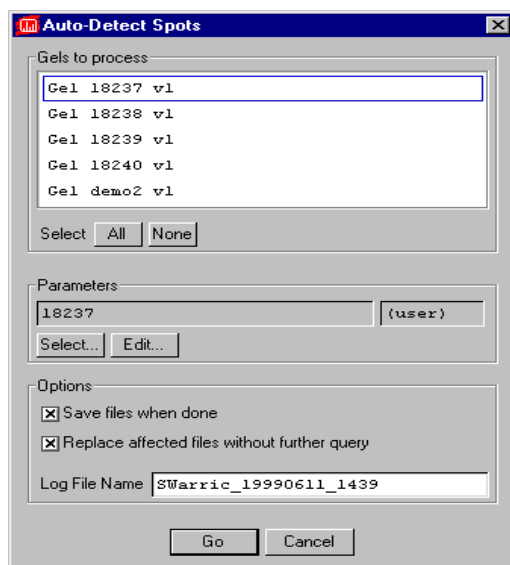


Fig.5-5. Auto-Detect Spots control panel.

All currently loaded gel scans will be displayed in the panel.

Note: Only the root name and version number of each gel scan will be displayed. Calstrip scans and multiple exposures will not be listed separately.

Select the scan or scans you want to process by clicking on them individually, or clicking on the All or None buttons.

Select the appropriate set of detection parameters by clicking on the Select button and selecting from the pop-up list of saved parameters.

PDQUEST User Guide

Note: If you clicked on Process All Gels in the Spot Detection Wizard, the detection parameters selected in Wizard will be displayed in the Parameters field.

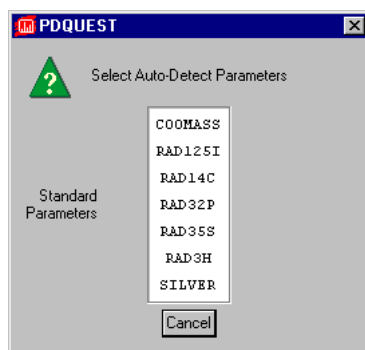


Fig.5-6. List of parameters for spot detection.

The parameters you select will be displayed in the Parameters field.

Click on Edit to return to the Spot Detection Wizard.

To automatically save the results of image processing, select the “Save files when done” checkbox.

If these gel scans have been processed before, specify whether you want to replace the existing files by selecting or deselecting the “Replace affected files without further query” checkbox.

You can review the results of spot detection in a log file that is created when you process your gel(s). The Log File Name field displays a default name for the log file. Accept this name or enter a new name.

Click on the Go button to proceed with spot detection.

A status screen will display the progress of spot detection.

Detecting and Editing Spots

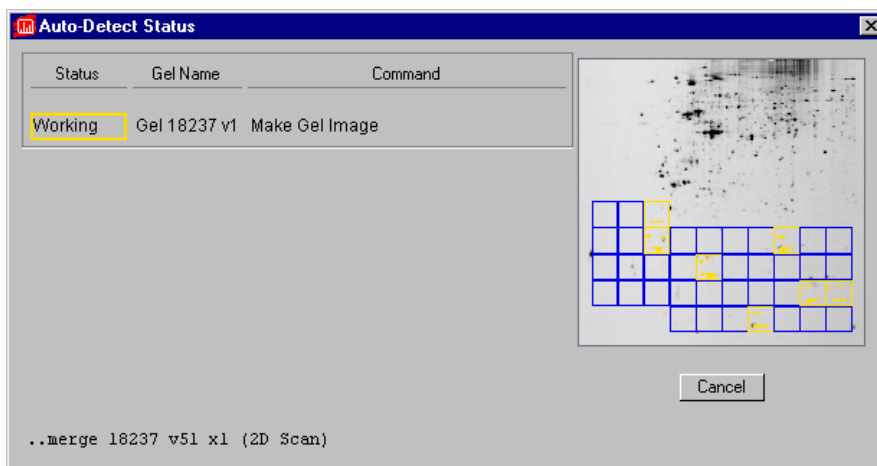


Fig.5-7. Displaying the progress of spot detection.

Continuously updated messages display information about the individual steps as they are performed. A preview image shows a graphical depiction of each step as it occurs.

Click on Cancel to interrupt spot detection before it is finished. Any steps already performed will be lost.

When detection is complete, a pop-up box will ask if you want to proceed with creating a matchset (see Chapter 6). You can create a matchset immediately, or review the results of spot detection and create a matchset later.

If you did not select "Save files when done" in the Auto-Detect control panel, save the results of your autodetection by clicking on the Save button on the main toolbar.

Spot Detection Report

To review the log file created during spot detection, select Spot Detection Report from the Spots menu. This report documents when spot detection began, which gels were detected (or failed to detect), which steps were performed, and when spot detection was completed.

5.3 Gel Images and Gel Spots Files

When you autodetect spots in a gel scan, you create two new files with the same root file name as the scan: a Gel Image file and a Gel Spots file.

The Gel Image file contains the merged and processed image. The Gel Spots file contains a synthetic image with Gaussian representations of the spots in your Gel Image file. The Gaussian fitting is essential for the purposes of higher-level spot analysis. Therefore, it is essential that the “synthetic” spots in the Gel Spots file accurately match the “real” spots in the Gel Image file.

When autodetection is complete, both the Gel Spots and Gel Image files will be loaded into the same image window with the scan that was used to create them.

Note: Use the Assign tool and Configure Subwindows command on the main toolbar to arrange your images in the window. Use Interchange Images (F11) on the main toolbar to toggle between the Gel Image and Gel Spots file. See Chapter 3 for more information about viewing commands.

5.4 Displaying Spots

The tools for displaying and editing spots are located on the Spots menu and Edit Spot toolbar

Detecting and Editing Spots

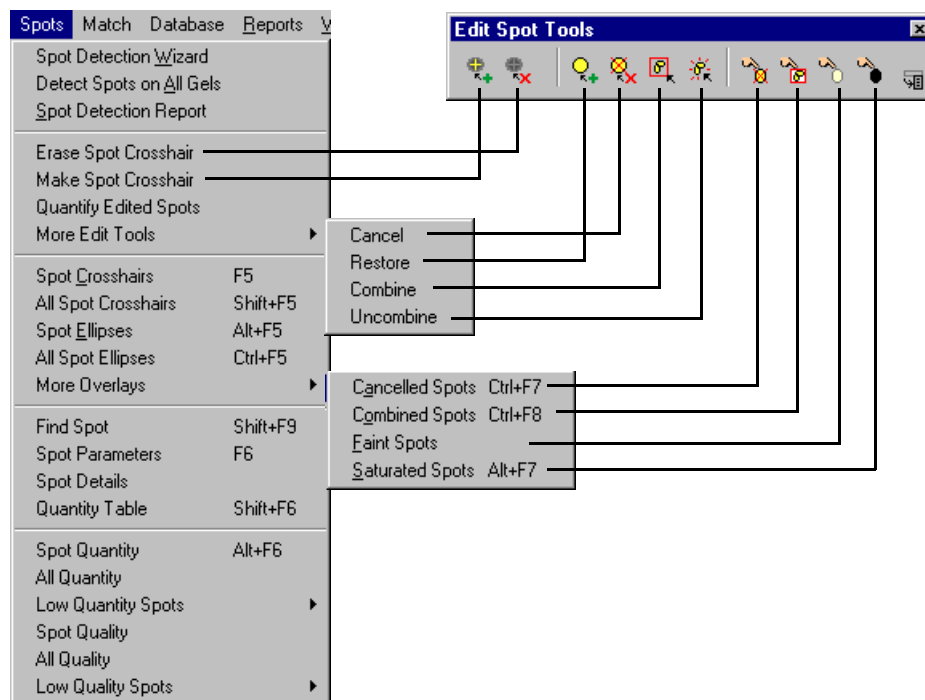


Fig.5-8. Spots menu and Edit Spot toolbar.

Spot crosshairs mark the exact centers of spots in the Gel Spots file. The crosshairs in the Gel Image are exact duplicates of those in the Gel Spots file.

To display spot crosshairs in a subwindow you click in, select Spot Crosshairs from the Spots menu (F5). To display crosshairs in all subwindows, select All Spot Crosshairs (SHIFT+F5).

PDQUEST User Guide

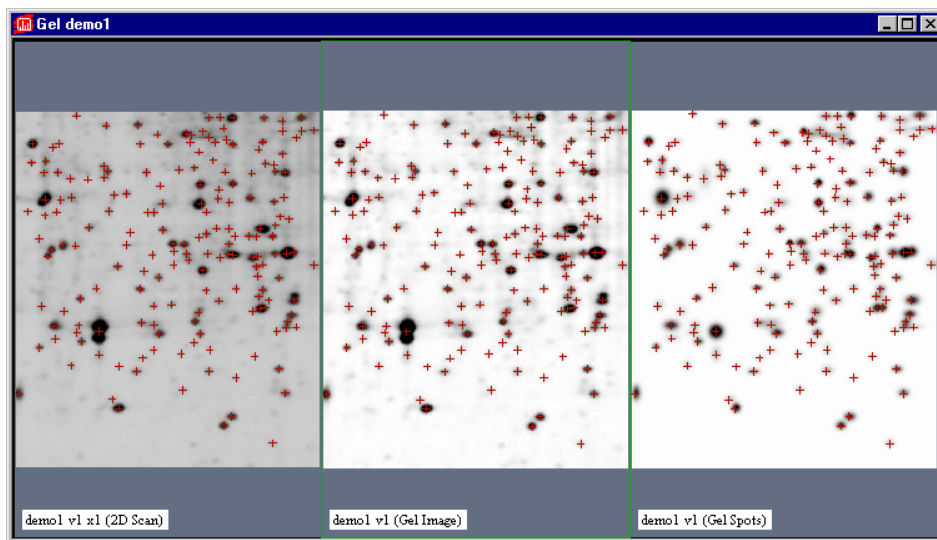


Fig.5-9. Spot crosshairs displayed in a Gel Scan, Gel Image, and Gel Spots file.

Each spot in the Gel Image should have a crosshair on top of it to indicate that this intensity has been detected as a spot.

Spot Ellipses on the Spots menu (ALT+F5) shows the shape of each spot generated by Gaussian fitting in the subwindow you click in. All Spot Ellipses (CTRL+F5) shows the spot ellipses in all the subwindows.

Automatic spot detection should accurately detect a high percentage of the spots you can see on the gel. Due to gel artifacts, etc., some spots may not be properly detected and fitted by Auto-Detect. If this is the case, you can manually add and delete spots using the commands described in the following section.

Viewing Tips

When studying spots on an image, you should magnify selected regions of the image using the Zoom Box tool on the main toolbar. You can navigate between zoomed sections using the ARROW KEYS on the keyboard.

Detecting and Editing Spots

Also, you can magnify and review the image section by section by using the View > Advanced View > 3x3 Sections function.

When checking for spot detection errors, use the Image > Transform command to change the display of the Gel Image. This can reveal spots that have gone undetected.

Use Hide Overlays on the main toolbar and View menu to conceal any overlays you have created. It can also be used to erase plots, traces, and information boxes that are created by other functions.

These and other viewing tools are described in detail in Chapter 3.

5.5 Spot Information

With your Gel Spots file loaded, you can find out basic information about your spots using the commands on the Spots menu.

Note: Spot information is stored in the Gel Spots file, but can be viewed in the Gel Image as well—provided the Gel Spots file is loaded in the same window.

Spot Quantity

Spot Quantity is the total intensity of a defined spot. It is calculated during spot detection and Gaussian fitting. The formula for calculating spot quantity in PDQUEST is:

Spot height * $\pi * \sigma_x * \sigma_y$, where:

Spot height is the peak of the Gaussian representation of the spot. It is measured in O.D.s or counts/IU², depending on the type of imaging device used.

Note: IU = Image Unit. For all Gel Spots images, 1 IU = 100 micrometers, even if the original images were scanned at different resolutions.

PDQUEST User Guide

σ_x is the standard deviation of the Gaussian distribution of the spot in the direction of the x axis, and σ_y is the standard deviation in the direction of the y axis. These are measured in IUs.

Note that this formula results in spot quantity units of O.D. * IU² for images measured in O.D.s, and spot quantity units of counts for images measured in counts.

Note: This formula results in more accurate quantitation of imaged spots than simply summing the intensities of the pixels in the spot. Measuring the area under the Gaussian compensates for peak intensity values that approach image saturation, as well as background intensity that submerges the intensity of a spot at its edges.

To display the quantities of the spots in a single image, select Spots > Spot Quantity (ALT+F6). To display the quantities of all the spots in all the subwindows of a selected window, select Spots > All Quantity.

Low Quantity Spots

You can define a set of spots that do not meet a threshold quantity value.

From the Low Quantity Spots submenu on the Spots menu, select Define Low Quan. Set (you must have a Gel Spots file loaded for this command to be available).

In the pop-up box, enter a quantity threshold value. All spots below this threshold will be flagged as low quantity.

When you click on OK, a pop-up box will list the total number of low quantity spots and their total quantity, and how they compare to the total number of spots in the image and their total quantity.

Detecting and Editing Spots

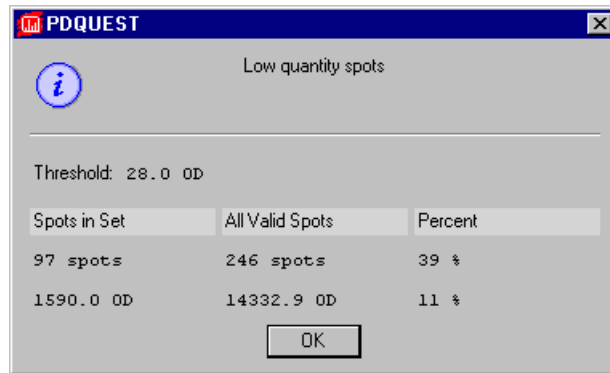


Fig.5-10. Pop-up box displaying low quantity spot information.

When you click on the OK button, the low quantity spots will appear highlighted in the image. You can also highlight the low quantity spots in your set by selecting Spots > Low Quantity Spots > Mark Low Quan. Spots.

To cancel all low quantity spots in the set, select Spots > Low Quantity Spots > Cancel Low Quan. Spots (see section 5.6 above for information about cancelling spots).

To restore the cancelled low quantity spots, select Spots > Low Quantity Spots > Restore Low Quan. Spots (see section 5.6 above for information about restoring canceled spots).

Spot Quality

Spot Quality on the Spots menu is a number ranging from 0 to 100 that is calculated based on five attributes of each spot:

- Gaussian fit. PDQUEST determines how well the spot fits the Gaussian model and assigns a value based on this fit.
- X streaking. The software determines how much of the spot is affected by gel streaking in the x direction.
- Y streaking. The software determines how much of the spot is affected by gel streaking in the y direction.

PDQUEST User Guide

- **Overlap.** The software determines the percentage of the spot that overlaps with another spot or spots.
- **Linear range of scanner.** PDQUEST determines whether the peak intensity value of the spot is within the linear range of the scanner.

Each spot attribute is evaluated and weighted to produce a single numerical value for Spot Quality. If a spot fits the Gaussian model perfectly, has no streaking in the x or y direction, does not overlap with any other spot, and has a peak intensity within the linear range of the scanner, PDQUEST will assign a value of 100 to the spot.

To display the quality of all the spots in all the subwindows of the selected window, select Spots > All Quality.

Low Quality Spots

You can define a set of spots that do not meet a threshold quality value.

From the Low Quality Spots submenu on the Spots menu, select Define Low Qual. Set (you must have a Gel Spots file loaded for this command to be available).

In the pop-up box, enter a quality threshold value. All spots below this threshold will be flagged as low quality.

When you click on the OK button, the low quality spots will appear highlighted in the image. You can also highlight the low quality spots in your image by selecting Spots > Low Quality Spots > Mark Low Qual. Spots.

To cancel all low quality spots in the set, select Spots > Low Quality Spots > Cancel Low Qual. Spots (see section 5.6 above for information about cancelling spots).

To restore the cancelled low quality spots, select Spots > Low Quality Spots > Restore Low Qual. Spots (see section 5.6 above for information about restoring canceled spots).

Detecting and Editing Spots

Spot Parameters

The Spot Parameters information box provides a variety of information about an individual spot, including its quantity, peak intensity, location in the gel, Mr pI data (if any), and the sigmas of its Gaussian curve.

Select Spots > Spot Parameters (F6) and click on a spot in a Gel Spots file or Gel Image to display the Spot Parameters data for that spot.

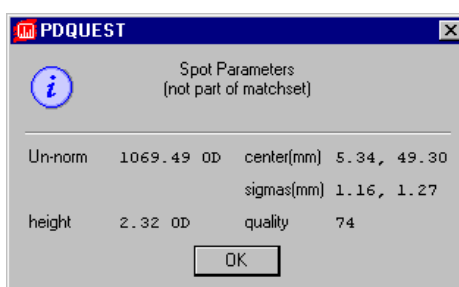


Fig.5-11. Spot Parameters information box.

Spot Details

Spots > Spot Details displays all available information about a spot. Select the command and then click on a spot. A pop-up box will display all the information associated with the spot.

5.6 Editing the Results of Spot Detection

Note: Usually, it is easier to correct problems with spot detection during matching (see the next chapter). This allows you to compare questionable spots on different gels. However, you can also correct problems before matching.

Spot detection errors include undetected spots, too many spots detected, misplaced spot centers in saturated regions, and spots that have not been properly fitted (crosshairs off-center, incorrect spot geometry, etc.).

PDQUEST User Guide

When checking the shape and position of spots, use Interchange Images (F11) to toggle between the Gel Image and the Gel Spots files.

Adding Spots

Spot detection will correctly identify most of the spots in your gel. However, some spots of low intensity or in clusters may be missed. Spots in the Gel Image that not marked with crosshairs will not be included in the Gel Spots file, and will not be available for analysis.

To manually mark a spot on your Gel Image, select Make Spot Crosshair from the Spots menu or toolbar, then click on the unmarked spot. The spot will be added to the Gel Spots file.

Make Spot Crosshair will invoke the Gaussian modeling algorithm for quantifying each new spot.

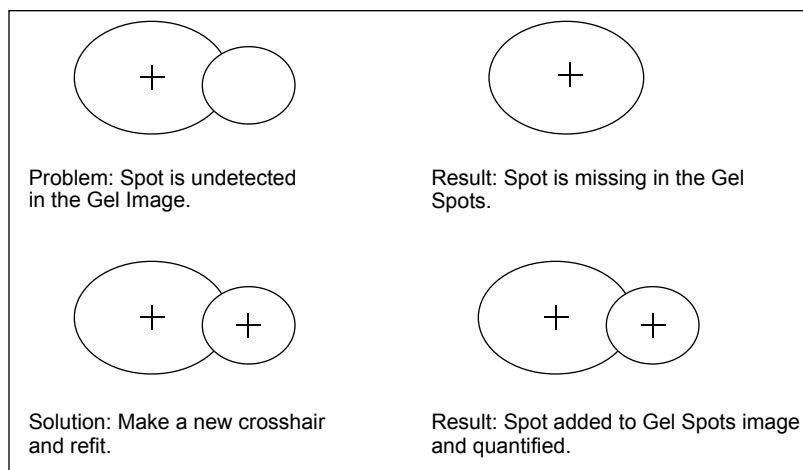


Fig.5-12. Adding a spot.

Use Plot Cross-section (View > View Density > Plot Cross-section) to aid in accurately positioning new spots (see section 3.5).

Detecting and Editing Spots

Erasing Spots

Occasionally, spot detection will identify too many spots, particularly in a cluster.

To remove an incorrect spot center, select Erase Spot Crosshair from the Spots menu or toolbar, then click on the wrongly marked spot on the Gel Image. The spot will be removed from the Gel Spots file.

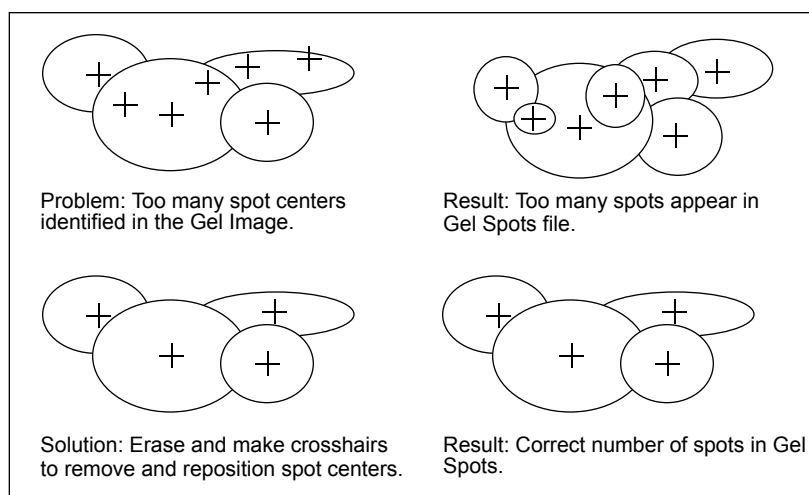


Fig.5-13. Erasing incorrectly identified spots in a cluster.

Note: You should use Make Spot Crosshair in conjunction with Erase Spot Crosshair to remove and reposition spot centers as needed.

Cancelling Spots

Cancel will preserve a marked spot in the Gel Spots file but not include it in matching or any other analysis.

Select Cancel from the Spots > More Edit Tools submenu or Edit Spot toolbar and click on the spot that you want to cancel.

PDQUEST User Guide

If you change your mind, select Restore from the Edit Spots submenu and click on the point of the cancelled spot. The spot will reappear.

Cancelled Spots (CTRL+F7) on the Spots > More Overlays submenu and Edit Spot toolbar highlights those spots that have been pushed into the background using the Cancel command.

Combining Spots

Combine will combine the quantitation of a cluster of spots in the Gel Image into one spot. The cluster will appear as individual spots in the Gel Image but as one spot in the Gel Spots file.

Select Combine from the Spots > More Edit Tools submenu or Edit Spots toolbar and draw a box around the cluster of spots that you want to combine. When you finish drawing the box, click on the point in the box where you want the combined spot to be placed.

If you change your mind, select Un-combine from the Edit Spots submenu and draw a box around the combined spot. The spot will be uncombined into the original multiple spots.

Combined Spots (CTRL+F8) on the Spots > More Overlays submenu or Edit Spots toolbar highlights those spots that have been combined from spot clusters using the Combine command.

Faint Spots

Faint Spots on the Spots > More Overlays submenu and Edit Spot toolbar highlights those spots whose signal falls within the bottom 5 percent of the intensity range of the image.

Saturated Spots

Because spot detection may not correctly fit spots in saturated regions, saturated spots are not quantified. However, you can still use spot crosshairs as markers for saturated spots on the image.

To highlight specific spots that have been saturated, select Spots > More Overlays > Saturated Spots (ALT+F7). Locate saturated regions in the Gel

Detecting and Editing Spots

Image using the Highlight Saturated Pixels option in the Transform window (Image > Transform).

After you have identified saturated spots, redisplay the spot crosshairs.

Look for crosshairs in the saturated spots that appear to be off-center. Correct the problem by erasing the off-center crosshairs (Erase Spot Crosshair) and adding new crosshairs (Make Spot Crosshair) in the correct locations.

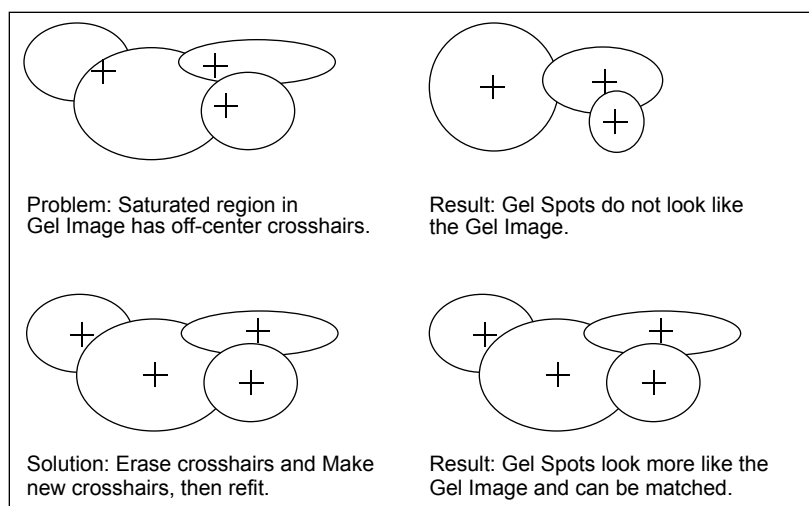


Fig.5-14. Correcting misplaced spot centers in saturated regions.

Remember that saturated spots are not quantified and centers in Gel Spot files are merely place holders for these spots.

PDQUEST User Guide

6. Creating and Editing Matchsets

After the protein spots in your gels have been detected, you are ready to compare spots across gels. To do so, you must create a matchset.

Matchsets are sets of gels that have been grouped together for the purpose of qualitative and quantitative comparison. Spots on one gel can be compared to spots on every other gel in the matchset.

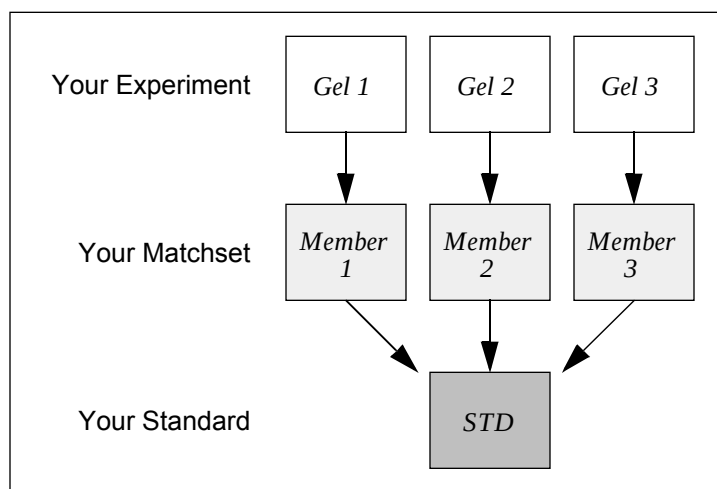


Fig.6-1. How experiments are organized into matchsets.

Note: While matchsets can contain gels with different quantitation units, quantitative spot comparisons can only be made between gels that have the same units (i.e., spot quantities measured in O.D.s cannot be compared to spot quantities measured in counts). Qualitative comparisons are still valid.

Matchsets can be created from any number of gels. Matching is coordinated through an image called the “**Standard**,” which is a synthetic image created

PDQUEST User Guide

from the members of your matchset. Ultimately, your Standard will contain a copy of every protein spot from every member gel in the matchset.

6.1 Creating a Matchset

Matchsets are created from Gel Spots files and can also contain Gel Image files.

Note: You should include Gel Images in your matchset unless disk space is a severe problem. Including Gel Images greatly facilitates spot matching and editing.

Matchsets can be created automatically after you perform spot detection. You will be prompted to create the matchset after spot detection is complete.

Otherwise, select Match > Create/Edit Matchset > Create Matchset to open the Create Matchset dialog.

Creating and Editing Matchsets

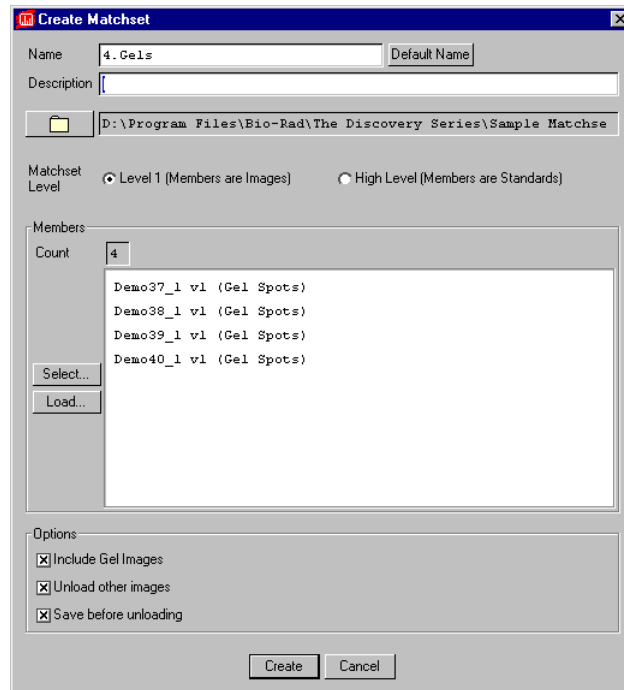


Fig.6-2. Dialog box for creating a matchset.

At the top of the dialog box, you are prompted to enter a name for your new matchset. To assign your matchset a default name generated by the program, click on the Default Name button.

You can also enter a description of the matchset in the text box next to the Description prompt.

Select the directory in which the matchset should be stored. The currently selected directory is displayed next to the Directory button. To save the matchset in a different directory, click on the Directory button; the Select Directory form will appear. Click on the Matchsets button in the Select Directory dialog to select the working Matchsets directory (specified in Edit > File Paths).

PDQUEST User Guide

Next, specify whether you are creating a “lower level” matchset or a “higher level” matchset by clicking on the appropriate button next to the Matchset Level prompt. (Lower level matchsets are created from original gel images; higher level matchsets are created from the Standards of lower level matchsets. For more information on higher level matchsets, see section 9.2)

Selecting Matchset Files

If you are creating a lower level matchset, select the Gel Spots files that should be included. If the files to be included are already loaded, click on the button labelled Select.

A list of currently loaded Gel Spots will be displayed. Select the files that you want include in your matchset and click on Done. The files you chose will be listed next to the Members prompt.

To include Gel Spots that are not currently open, click on the Load button. The Select New Member dialog box will be displayed. Select the files that you want to include in the matchset and click on Add Selected.

The names of the files you selected will appear in the Members list in the Create Matchset dialog box.

Other Options

Gel Images are included in the matchset by default. To not include Gel Images in the matchset, click in the Include Gel Images checkbox to deselect it.

Note: You should include Gel Images in your matchset unless disk space is a severe problem. Including Gel Images greatly facilitates spot matching and editing.

Any other loaded images will be unloaded before this matchset is created, unless you click in the Unload Other Images checkbox to deselect it.

Any loaded images will be saved before they are unloaded unless you click in the Save Before Unloading checkbox to deselect it.

Creating and Editing Matchsets

Creating the Matchset

Once you have finished entering information on the Create Matchset form, click on the Create button.

You will be prompted to select an image to serve as the template for the Standard. Generally, this will be an image with the most spots, or an image that is most representative of the entire matchset. (See the following section for more information.) If you are unsure which file you want to select for your Standard, click on Cancel—you can select the Standard template later. Otherwise, select the appropriate file name.

A new window will open displaying all the members of the newly created matchset.

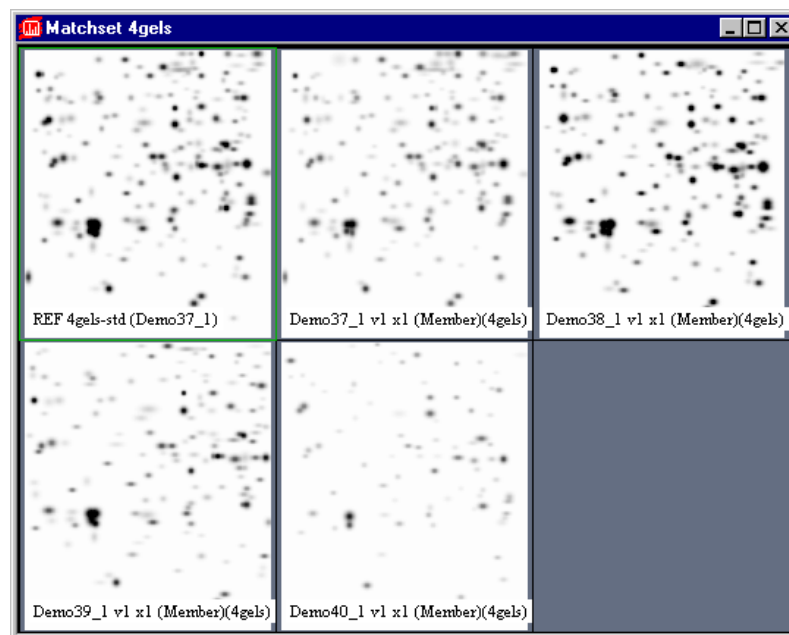


Fig.6-3. Example of a matchset. The Standard is in the upper left subwindow.

6.2 Matchset Standard

For each matchset, you must choose one gel that will serve as the template for the matchset Standard. A copy of that image will be created and designated the Standard.

The goal is to match all gels in the matchset to this Standard. The Standard represents all the data in the matchset and serves as the basis for comparing different experiments in higher level matchsets (see Chapter 8).

The Standard does not contain any image data, although many calculations, such as those that yield relative molecular weights and isoelectric points for protein spots, are performed in the Standard. Also, spot annotation occurs in the Standard.

6.2.a Criteria for Selecting the Standard Gel

Before choosing a gel to serve as the template for your Standard, you may want to review how many spots were detected in each Gel Spots file you are adding to the matchset. With the Gel Spots file open and its window selected, select File Info from the File menu. The info box will tell you how many spots were detected in the gel, as well as file size, total data in the spots, etc.

This information can help you decide which gel should serve as the Standard template.

The criteria for selecting the template for the Standard are listed below. Ideally, you want the Standard to satisfy all of the following criteria. If you intend to build a database, or analyze data from many different experiments, we recommend that you select a gel that is representative of your system and/or is a gel run from a representative control sample. Otherwise, choose the gel with the most spots that comes closest to satisfying the other criteria.

Number of Spots

In creating a Standard, you start with the spots in the template, then you manually add spots from the other gels in the matchset. Choosing a template with the most spots minimizes number of spots that you will need to add later.

Creating and Editing Matchsets

Representative Gel

The Standard is a representation of all the data in your matchset and is often used in reports. Therefore, you may not want to use an atypical gel as a template, even if it has the most spots. We suggest that you choose the most representative Gel Spots file (i.e., the one that looks like a typical gel run from samples in your biological system).

Control Gels

Choose a control gel as the Standard so that it can be included in every experiment you plan to analyze and used as the Standard in every matchset. This makes it easier to compare results between different experiments.

Spot Quality

Editing the Standard and matching other gels to the Standard are easier when the Standard consists of high quality spots.

6.2.b Creating a Standard

You will be prompted to select an image for the Standard when you first create the matchset. To create a Standard later, select Make Standard from the Match > Create / Edit Matchset submenu. A pop-up box will appear. Select the member image to serve as the Standard template from the list.

The Standard image will appear in the upper left subwindow of the matchset and will have the name “REF <matchset name>-std.”

At this point, we suggest that you save the matchset by clicking on Save.

6.3 Editing a Matchset

The commands for editing the structure of a matchset are located on the Match > Create / Edit Matchset submenu.

PDQUEST User Guide

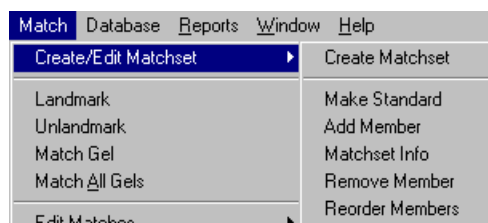


Fig.6-4. Create/Edit Matchset submenu.

You can add members to and remove members from an existing matchset. Selecting Add Member from the Match > Create / Edit Matchset submenu will open a version of the Create Matchset dialog box (see above). You can add a new Gel Spots file to the matchset using this dialog as if you were creating the matchset for the first time.

To remove a matchset member, select Remove Member and select the member to be removed from the pop-up list. You will receive a warning that this command cannot be undone.

Matchset Info opens a dialog box in which you can enter a description and comments about your matchset.

Reorder Matchset opens a dialog box in which you can select a member name and move it up or down in the matchset hierarchy. This hierarchy affects the order in which the matchset members are displayed on the screen (top to bottom, left to right).

Editing Spots in the Matchset Members

If you have not already done so, review the spots in the Gel Images of your matchset members and, if necessary, add, erase, cancel, or combine spots using the tools described in section 5.6, Editing the Results of Spot Detection.

Note: Spots can only be added to and erased from Gel Images, not Gel Spots files. Therefore, these commands will not work on the Standard, which does not have a Gel Image file.

Creating and Editing Matchsets

As always, when editing spots, it is useful to magnify selected portions of the images using the zoom tools and Imitate Zoom command described in Chapter 3.

Note: Spots edited in the Standard or any matchset member are not automatically updated in the other members or Standard. You must make the same changes to the same spots/clusters in all the images of the matchset.

6.4 Landmarking

Landmarks are reference spots used by PDQUEST to align and position matchset members for matching. They are used to compensate for slight differences and distortions in the member gels.

When choosing landmark spots in your member images, you should look for spots that are:

- Present in all members of the matchset.
- Well-resolved.
- Near a corner of the image.
- Isolated—i.e., very few spots in its immediate vicinity.
- Located in or near regions of high spot concentration.

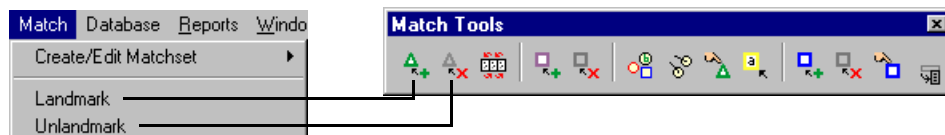


Fig.6-5. Landmarking commands.

We suggest that you landmark approximately five spots per section of the gel (unless you are matching colony plates—see section 6.8). As a rule of thumb, you should landmark approximately 10% of the spots on your gel.

PDQUEST User Guide

Entering Landmarks

To begin landmarking spots, select Landmark from the Match menu or toolbar and click on a spot on the Standard.

Note: Always initiate landmarking in the Standard.

The cursor will jump to the same area in the next member.

If the landmarking routine can identify the same spot in that member, it will appear marked.

If the routine is unsure about placing the landmark in that member, it will appear unmarked, and the area in which the routine expected to find the spot will be circled. Click on the appropriate spot manually.

If a spot is missing in one of the gels, skip that gel and continue landmarking the other gels.

Continue for all the members in the matchset.

To highlight all the landmarks in a matchset, select Mark Landmarks (F8) from the Match menu.

Unlandmarking

If you accidentally landmark the wrong spot, select Unlandmark from the Match menu or toolbar and click on the incorrect landmark.

Note: Unlandmarking a member image will only remove the landmark from that member. Removing a landmark from the Standard unlandmarks that landmark from all the members of the matchset.

6.5 Auto-Matching

Now you are ready to match the members of your matchset to the Standard.

First, display each member of your matchset, including the Standard. Make sure that you have created sufficient (more than 3) landmarks in the gels you want to match to the Standard.

Creating and Editing Matchsets

Select Match All Gels from the Match menu to match all the members in one step.

To match a single member, select Match Gel from the Match menu and click on the member to match.

Spots that have been matched in all the images will be marked on the Standard with a letter.

6.6 Checking and Editing Matches

The accuracy of the PDQUEST matching algorithm depends upon several factors, not the least of which is the number of landmarks you have entered. When the first round of matching is complete, we suggest that you review the results and make any corrections necessary to the matches called by PDQUEST.

The general procedure involves making as few corrections as possible, adding new landmarks where necessary, and rematching the entire gel to the Standard.

The commands for checking and editing matches are located on the Match menu, Match > Edit Matches submenu, and Match Tools toolbar.

PDQUEST User Guide

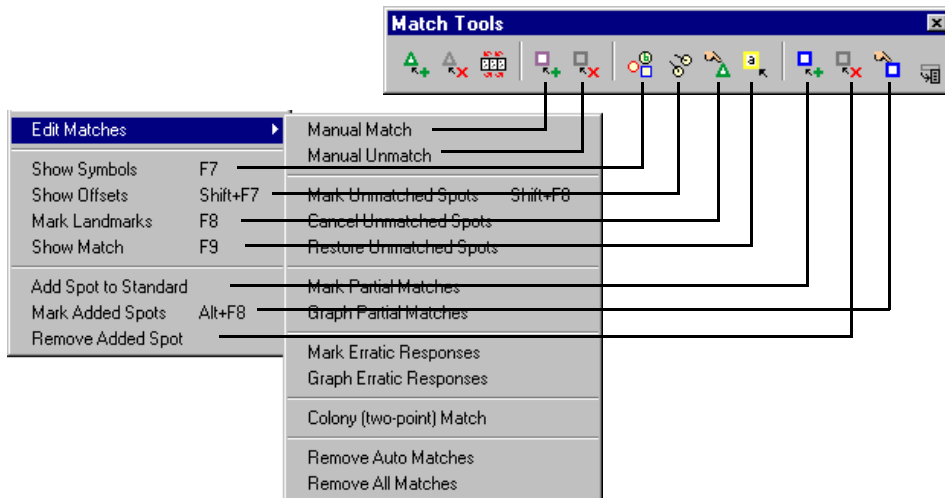


Fig.6-6. Edit Matches commands.

Note: To check and edit matches, configure two vertical subwindows (using the Configure Subwindows command from the main toolbar). Display the Standard image in the left window and the member image you want to edit in the right window. When you are done editing one member image, replace it with the next member.

6.6.a Checking Matches

PDQUEST has a number of tools for displaying and identifying your matched and unmatched spots.

Show Symbols

To label all the spots that were not matched in the automatic matching step, select Show Symbols (F7) from the Match menu or toolbar.

- Spots matched in all the images will be labeled with letters in both the Standard and member images.
- Spots that are in the Standard but not in any member image (except the image used to create the Standard) will be circled in red in the Standard.

Creating and Editing Matchsets

- Spots that are in a member image but not in the Standard will be circled in red in the member image.

Show Match

To look at individual matches, select Show Match (F9) and click on a spot in the Standard. The matched spots in the member images will be highlighted. If a spot is not matched in the member image, a circle will be drawn in its estimated location.

Show Offsets

Show Offsets displays the vector offsets from spot center locations in the member gels to the equivalent (matched) spot center locations in the Standard. Select Show Offsets (SHIFT+F7) from the Match menu or toolbar to display the offset lines.

Look for the following conditions that indicate possible matching problems:

- Offset lines that cross each other (this is a key problem).
- Grossly nonparallel offset lines in local regions.
- Offset lines of unequal lengths.
- Areas of the gel with few offsets; this indicates unmatched regions.

PDQUEST User Guide

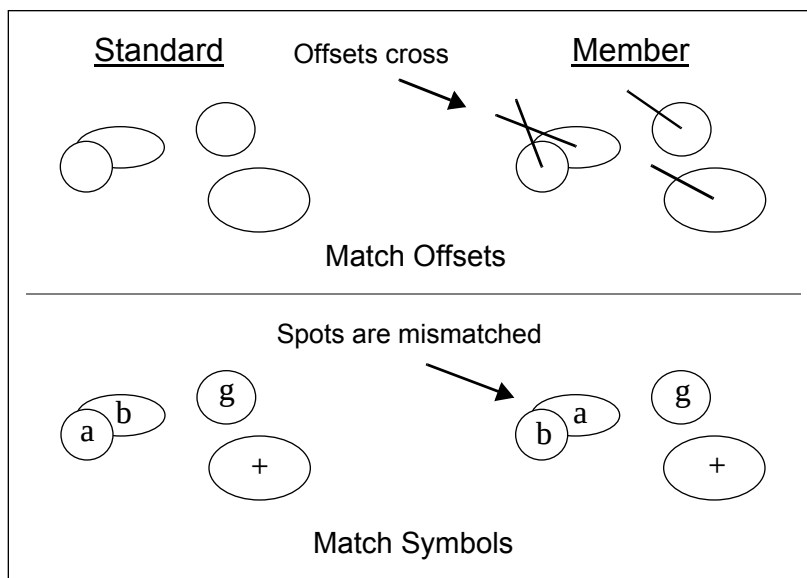


Fig.6-7. Example of incorrect matches.

Partial Matches

Partial matches are spots in the Standard that are not matched in one or more member images. This may indicate spots that either are not present in the gel or gels in question, in which case there is nothing to match, or spots that were missed during matching for one reason or another (e.g., the spot may have been cancelled in the gel).

Note: These commands only work on the Standard. In addition, Graph Partial Matches will only work in single-window mode.

Mark Partial Matches (Match > Edit Matches) highlights spots that are partially matched in the Standard.

Graph Partial Matches displays histograms connected to each partially matched spot in the Standard. These histograms show whether a spot is matched or unmatched in each member image. (For more information on graphs, see section 10.1)

Creating and Editing Matchsets

Each bar in a histogram represents the quantitation of the spot in each member of the matchset. The bars are ordered from left to the right based on the order in which the members were loaded. If a spot is not matched in a member, there will be no bar for that member in the histogram.

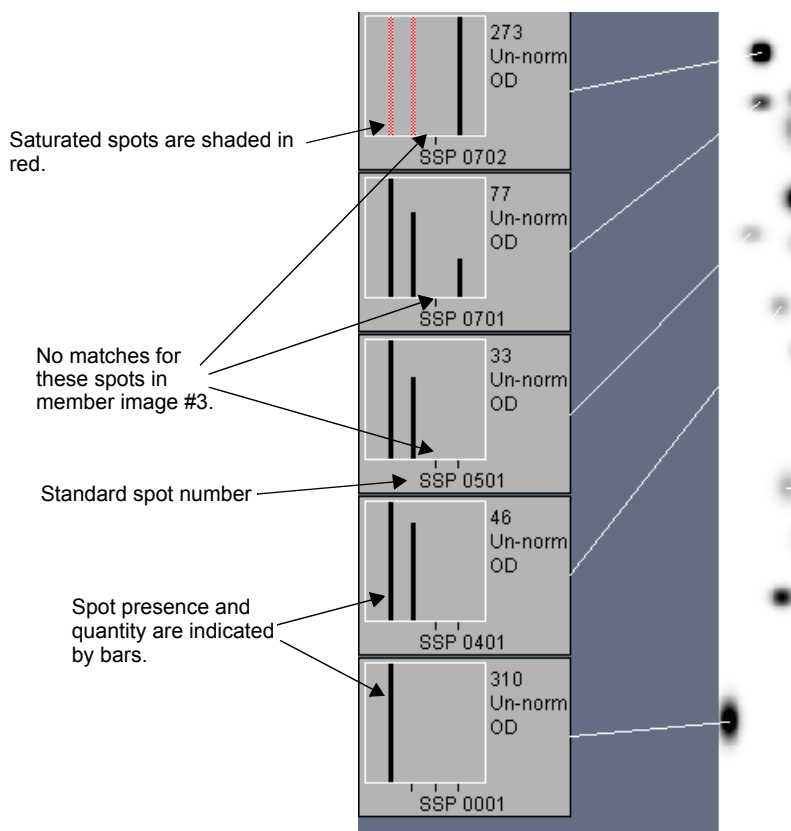


Fig.6-8. Examples of histograms of partially matched spots.

Histograms can only be displayed for one region of an image at a time. To shift regions, select Reports > Graph > Graph Advance (ALT+.) or Graph Retreat (ALT+,.). Graphs will be displayed for spots in the next region of the image.

PDQUEST User Guide

The histogram also displays the intensity and Standard Spot (SSP) number for each spot. (For more information on SSP numbers, see section 8.1) Decide which spots need to be checked and make a list of SSP numbers and gel numbers in which to reexamine the matching.

Be sure to check all the spots in one gel before moving on to the next.

To correct partial matches on a member gel, first configure two vertical windows, then display the Standard in the left and the member image in the right.

Use the Spots > Find Spot (SHIFT+F9) function on the Standard to locate a spot whose SSP number is included on the list of spots needing to be re-examined. Then use the Zoom Box function followed by Imitate Zoom to get an enlarged view of the area around the spot in both images.

Erratic Spots

Erratic spots are matched spots in the Standard whose quantity in the member images varies by a significant factor. To display these, select Mark Erratic Responses from the Match > Edit Matches submenu, and enter a fold change factor in the pop-up box. (For example, a change factor of 3 will highlight spots whose quantity varies by a factor of three or more among member images.)

Graph Erratic Responses on the Match > Edit Matches submenu displays histograms of those matched spots whose quantity in the member images varies by a significant factor. Enter a fold change factor in the pop-up box and click on Apply.

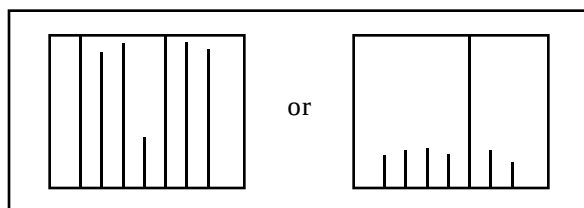


Fig.6-9. Example histograms of spots showing erratic responses.

Creating and Editing Matchsets

6.6.b Editing Matches

Additional Landmarking

Landmarking spots that have been mismatched correctly matches the spots in question and also improves matching in the region when you rematch the gel to the Standard. Select Landmark from the Match menu, click on the spot in the Standard, then click on the spot in the member to landmark it.

Note: Landmarking will improve the modeling of the Standard, facilitating additional matching.

Enter additional landmarks in regions that have few matched spots.

After you have added several new landmarks, rematch the entire gel with the Match All Gels command.

Manual Matching

If a spot is present in all the member images but is not matched in the Standard, select Manual Match from the Match > Edit Matches submenu or Match toolbar, click on the spot with the red circle in the Standard, then click on the same spot in all the member images. The circled spot in the Standard will changed to matched.

To unmatch a manually matched spot, select Manual Unmatch from the Match > Edit Matches submenu or Match toolbar and click on the manually matched spot.

Note: Manual matching will not change the modeling of the Standard, and therefore will not facilitate additional matching.

Adding Spots to Member Images

It is possible that a spot is in a Gel Image but was not detected and therefore was not included in the matchset member. If this is the case, use Make Spot Crosshair to add it to the member (see section 5.6). After that, you will want to match it to the Standard using Manual Match.

PDQUEST User Guide

If you cannot match a spot because it is not in the Standard, then you may need to add the spot to the Standard (see below).

Adding Spots to the Standard

In this step you look for spots in each member that are not included in the Standard. These may be spots induced by the experimental protocol for a particular gel.

It is important to locate these spots and add them to the Standard. Only then will you be able to match them to equivalent spots in other members.

Once again, configure your matchset into two vertical windows. Assign the Standard to the left window and the desired member to the right window.

Select Show Symbols (F7) and look for spots with red circles in the member image. These spots are unmatched. You can also use the function Mark Unmatched Spots (SHIFT+F8) to find these spots.

There are usually four situations that result in unmatched spots; only the last requires that you add spots to the Standard.

1. A spot near the unmatched spot may be mismatched to the Standard.

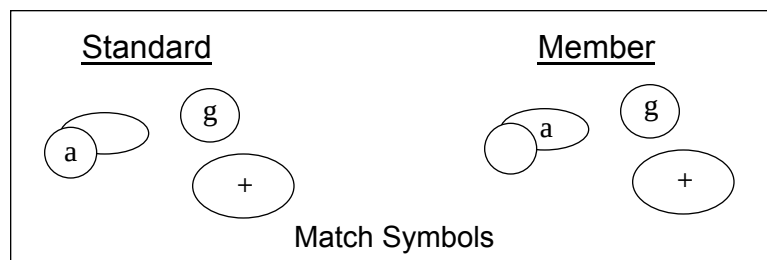


Fig.6-10. Mismatched spots.

2. The spot may be cancelled in the Standard. If so, restore the spot and match.

Creating and Editing Matchsets

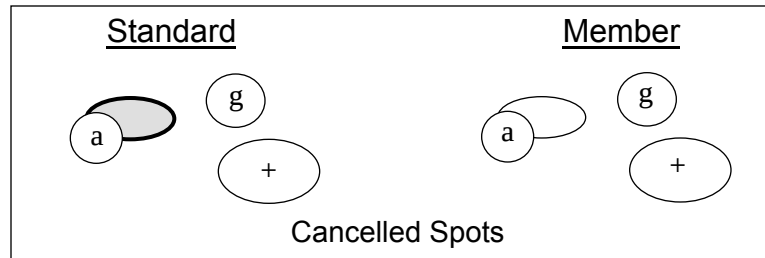


Fig.6-11. Restoring a spot that was cancelled in the Standard.

3. If the spot is in the Standard, but not matched in the member, match the spot.

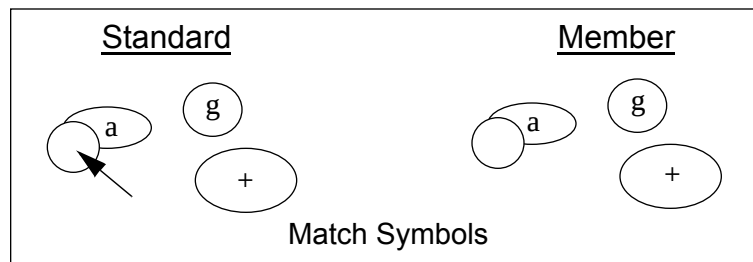


Fig.6-12. Matching unmatched spots.

4. The spot is not in the gel from which the Standard was made. This spot may be added to the Standard as described below.

PDQUEST User Guide

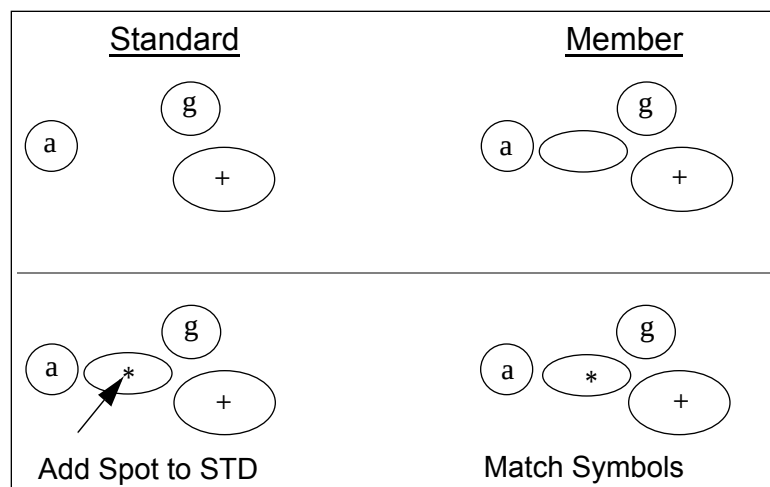


Fig.6-13. Adding a spot to the Standard.

Select Add Spot to Standard from the Match menu. Click on the center of the unmatched spot in the member. The cursor jumps into the Standard window to the position the spot is expected to be, based on the surrounding matches.

When you click in the Standard, the spot will appear and a match is automatically made to the spot in the member.

Repeat this procedure for the entire member image, then for all the members of the matchset.

To remove a spot that you have added, select Remove Added Spot from the Match menu and click on the added spot in the Standard.

Mark Added Spots (ALT+F8) highlights the spots that you've added to the Standard. They are marked by blue squares.

After adding a spot from one member to the Standard you may find that an unmatched spot in another member can be matched to the added spot. If so, make the appropriate matches.

Unless your system induces and/or represses many proteins, you may need to add no more than a dozen or so spots from all members to the Standard.

6.7 Removing All Matches from an Image

To remove all automatic matches from a member image, select Remove Auto Matches from the Match > Edit Matches submenu and click on the member image.

Manual matches and landmarks will not be removed. Because this command is irreversible, you will be prompted to complete the operation.

To remove all matches (automatic and manual), as well as landmarks, from a member image, select Remove All Matches from the Match > Edit Matches submenu and click on the member image. Because this command is irreversible, you will be prompted to complete the operation.

6.8 Colony Matching

If you have a matchset containing colony plate images, you can match colonies quickly and easily using the Colony (two-point) Match command on the Match > Edit Matches submenu.

Unlike the regular matching function, you only need to define two landmarks on your member images to use this command.

Note: Also unlike the regular matching function, this command will automatically correct for any rotational errors in the imaged plates (common in colony plates). You do not need to perform additional image rotation.

After you have created a matchset and Standard and defined two landmarks in the member images (see above), select Colony (two-point) Match from the submenu and click on a member image to automatically match the colonies.

PDQUEST User Guide

7. Mrpl Values and Spot Quantitation

PDQUEST has various tools for calculating and displaying the molecular weights/ isoelectric points and quantities of the spots in your matchsets.

Known molecular weights and isoelectric points can be entered for proteins using the Database > MrpI Standards commands. Using the values that you provide, PDQUEST calculates the molecular weight (Mr) and isoelectric point (pI) values for all the other proteins in your Standard.

Since most experiments are carried out in duplicate, data analysis can be performed on each gel individually or on sets of duplicate gels grouped together. Functions in the Database > Replicate Groups submenu allow you to group sets of duplicate gels and determine the average quantitation of their protein spots.

The Normalize function will normalize spot quantitation among gels in a matchset.

Graphs (histograms) of spot quantitation and MrpI data can be displayed using the functions in the Reports > Graph submenu. Examining the histograms can give you a general sense of the trends in protein expression.

7.1 Mrpl Standards

You can enter the molecular weight and isoelectric point values for your known protein spots. With this data, PDQUEST can calculate the Mr and pI values for all the spots in your Standard.

The functions needed to enter and display Mr and pI data are found in the Database > MrpI Standards submenu and MrpI Tools toolbar.

PDQUEST User Guide

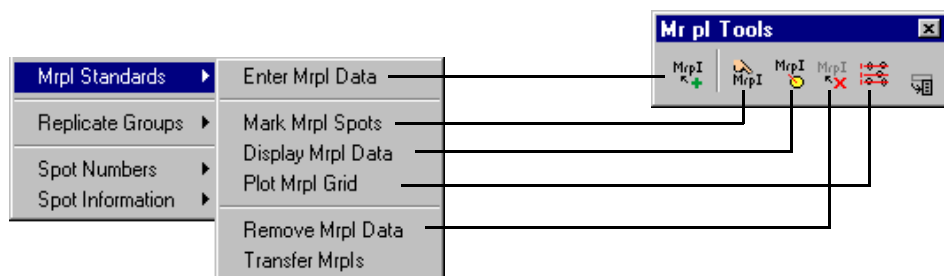


Fig.7-1. Mrpl commands.

Note: PDQUEST needs to know the orientation of the gel in order to calculate pI values. To enter the correct orientation, go to the Preferences dialog box (Edit > Preferences) and click on either Acid-Base or Base-Acid next to the pH Gradient prompt.

Entering Mrpl Values for Known Spots

To determine Mr and pI values for all the spots in your Standard, you first need to enter Mr and pI data for a few known spots.

Select Enter Mrpl Data from the Database > Mrpl Standards submenu or toolbar, then click on a spot in the Standard whose Mr and /or pI values are known. A dialog box will ask you to enter the Mr and pI data for that spot.

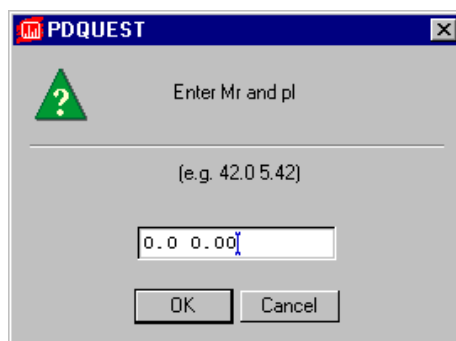


Fig.7-2. Pop-up box in which Mr and pI data are entered.

MrpI Values/Spot Quantitation

Type the Mr values first, then the pI data, and be sure to separate them with a space. Click OK when you are done.

If you make a mistake entering the data, you can select Enter MrpI Data again and click on the spot. The values you entered will appear in the dialog box and you can make the necessary corrections.

If you want to remove the Mr and pI data for a spot, select Remove MrpI Data from the submenu and toolbar and click on the spot.

Displaying Mr and pI Data

To see the proteins for which you have entered Mr and pI data, you can use either Mark MrpI Spots or Display MrpI Data. The first will highlight the spots for which you entered data. The second will display Mr and pI values next to the spots.

For a graphical display of the overall trend of Mr and pI values across the entire gel, first make sure the Standard subwindow is active, then select Plot MrpI Grid from the Database > MrpI Standards submenu or toolbar.

A grid of Mr and pI axes will appear over the Standard. You can use the grid to verify that the values that you entered make sense; you can also use it to estimate the Mr and pI values of proteins of interest.

To see the exact Mr and pI values of your favorite spot, select Spot Parameters from the Spots menu or toolbar and click on the spot.

In the pop-up box, the Mr and pI values will be displayed. These values are calculated based upon the known Mr and pI data that you entered.

Note: If the pI of the spot of interest falls outside of the highest and lowest known pI values, a pI value for that spot will not be calculated. Empirical results have shown that extrapolated pI values are unreliable, due to gel distortions and other factors.

7.2 Replicate Groups

Because we recommend that you perform your experiments more than once, PDQUEST offers you the option of looking at each gel individually or looking at the average of the data of duplicate gels.

A replicate group is the name given to a set of gels that are duplicates of each other. For example, three separate gels of a control cell lysate theoretically should produce the same quantitation for all spots on a 2-D gel. In practice, however, slight variations in the quantitation will probably be seen. Instead of choosing the one that you think is best, you can take the average and use those values for spot quantitation.

Creating a Replicate Group

Note: All your duplicate gels must be in the same matchset before you can create a replicate group.

To make a replicate group, open the Replicate Groups submenu (Database > Replicate Groups).

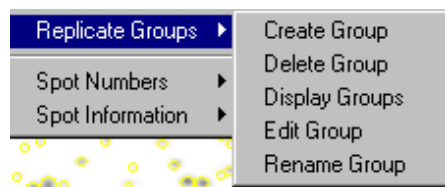


Fig.7-3. Replicate Groups submenu.

Select Create Group. In the dialog box, enter a name for the replicate group next to the Name prompt.

Mrpl Values/Spot Quantitation

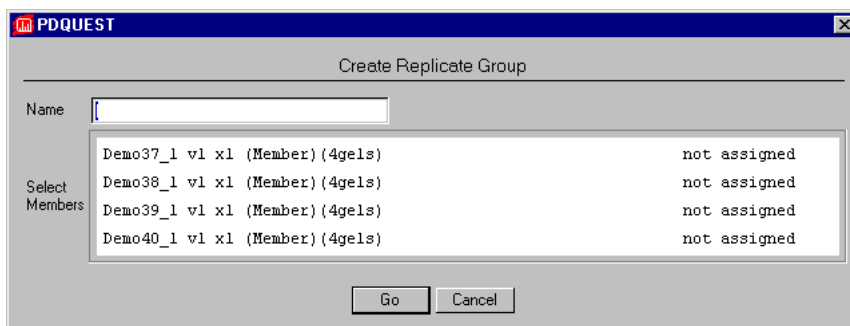


Fig.7-4. Enter a name for the replicate group in the dialog box.

Next you will be asked to select the gels that you want to include in your replicate group. Click on each gel that should be included.

Click on Go when you are finished. In the upper left message bar, you will see the total number of gels in this group.

Replicate groups can be altered using the commands Rename Group, Delete Group, and Edit Group.

- Rename Group allows you to enter a new name for the replicate group.
- Delete Group deletes the selected replicate group.
- Edit Group lets you add gels to or remove gels from a replicate group.

Quantitating a Replicate Group

After you have created a replicate group, you can calculate the average quantity of the members of the group by selecting Quantitation Mode under the Edit menu (CTRL+F6) and clicking on the Replicate Group Quantitation button in the pop-up form.

PDQUEST User Guide

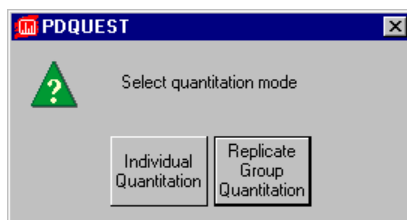


Fig.7-5. Quantitation Mode pop-up.

7.2.a Displaying Replicate Group Quantitation

Replicate group quantitation can be displayed using the Quantity Table function on the Spots menu (SHIFT+F6).

Click on the spot whose quantitation you want to see. Depending upon the quantitation mode you selected, one of two kinds of pop-up boxes will appear.

If you chose individual quantitation, the box will display the quantitation values found for that spot in each of the gels.


A screenshot of a Windows-style dialog box titled "PDQUEST". The dialog displays a table titled "Quantitation Table SSP 1401 DSN {}". The table has two columns: "Member Name" and "Quantity". There are four rows of data. Below the table is an "OK" button.

Quantitation Table SSP 1401 DSN {}	
Member Name	Quantity
Demo37_1	338.9
Demo38_1	314.1
Demo39_1	311.3
Demo40_1	130.1

Fig.7-6. Quantity table displaying individual spot quantitation.

If you selected Replicate Group Quantitation, the Average value of that spot's quantitation will be displayed, as will the coefficient of variance (CV) and the number of gels in the replicate group (Count).

MrpI Values/Spot Quantitation

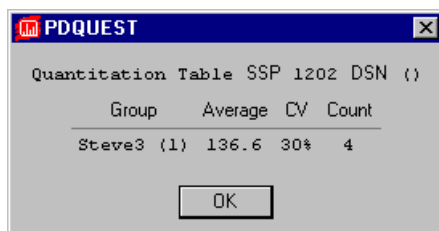

The image shows a software window titled 'PDQUEST' with a table of quantitation data. The table has four columns: Group, Average, CV, and Count. The data row shows 'Steve3 (1)' with an average of 136.6, a CV of 30%, and a count of 4. An 'OK' button is located below the table.

Group	Average	CV	Count
Steve3 (1)	136.6	30%	4

Fig.7-7. Quantity table displaying replicate group quantitation.

7.3 Normalization

When comparing 2-D gels in a matchset, you can look for on / off regulation or up / down regulation of proteins between samples. On / off regulatory changes are relatively easy to detect (a spot is either present or it isn't). However, up / down regulatory changes can be subtler and more difficult to detect.

Also, there is often some variation in spot intensity between gels that is not due to differential protein expression. This type of variation makes same-spot comparisons between gels difficult. The process of compensating for this nonexpression-related variation is called normalization.

Such variation in spot intensity can be caused by a number of factors:

- Inconsistent cell numbers (cell densities) in the sample prep,
- Inconsistent sample prep efficiency due to variations in reagents, protocols, etc.,
- Handling errors resulting in some sample loss during prep,
- Pipetting errors during sample prep,
- Sample loss during gel loading,
- Inconsistent staining times between gels,
- Inconsistent labeling efficiency with radioisotopes,
- Inconsistent excitation / emission efficiencies between different fluorescent tags,

PDQUEST User Guide


- Inconsistent detection energy sources between gels during image acquisition, and
- Inconsistent exposure times between gels during image acquisition.

All of the protein spots on a given 2-D gel are affected equally by these factors. However, the representation of the same spot on different gels can vary dramatically due to these factors. To accurately compare the quantity of any spot across two or more gels, it is essential to first establish that the value reported for a spot has been corrected to account for any variation between gels.

7.3.a Normalizing a Matchset

Note: Normalization affects histograms of spots and the quantitative data printed in reports.

With any image in the matchset active, select Normalize from the Edit menu. The Normalization dialog box will open.

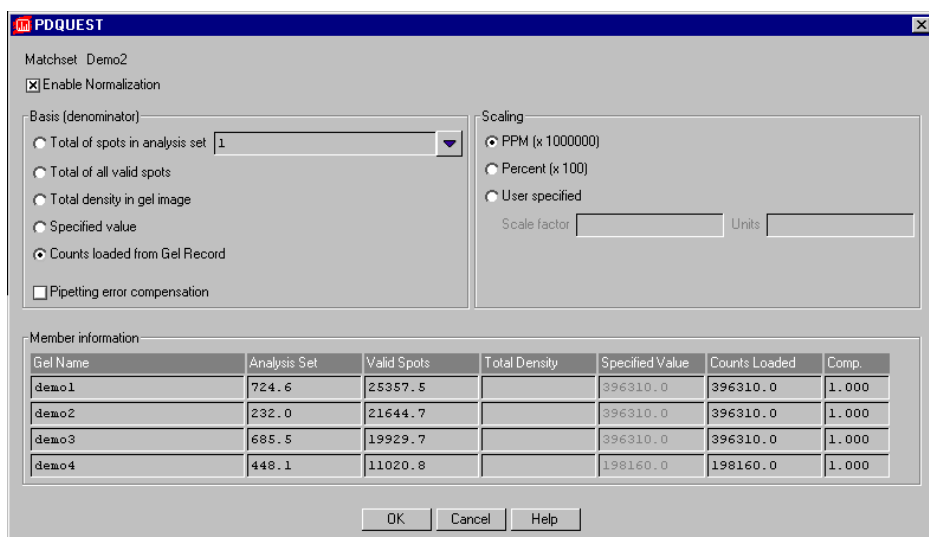

PDQUEST

Matchset Demo2

☒ Enable Normalization

Basis (denominator)

☐ Total of spots in analysis set 1

☐ Total of all valid spots

☐ Total density in gel image

☐ Specified value

☒ Counts loaded from Gel Record

☐ Pipetting error compensation

Scaling

☒ PPM (x 1000000)

☐ Percent (x 100)

☐ User specified

Scale factor Units

Member information

Gel Name	Analysis Set	Valid Spots	Total Density	Specified Value	Counts Loaded	Comp.
demo1	724.6	25357.5		396310.0	396310.0	1.000
demo2	232.0	21644.7		396310.0	396310.0	1.000
demo3	685.5	19929.7		396310.0	396310.0	1.000
demo4	448.1	11020.8		198160.0	198160.0	1.000

OK Cancel Help

Fig.7-8. Normalization dialog box for member images.

MrpI Values/Spot Quantitation

At the top of the form, click on the Enable Normalization checkbox to activate normalization for the matchset.

Each normalization option is described in detail in the following sections.

Note: Which normalization method you select will depend in part on how much you know about your gels and how they were prepared. For example, if you can readily identify housekeeping proteins on your gels, the Total of Spots in Analysis Set method is preferable. If you have known variations in your samples, you can select Specified Value and enter your known variation factor for each gel.

When you are satisfied with the settings in the dialog box, click on the OK button.

Normalization Formula

Normalization in PDQUEST is performed for each gel in a matchset using the following formula:

Normalized spot quantity = Raw spot quantity * Scaling factor * Pipetting error compensation factor (if any) / Normalization factor

1. Raw spot quantity is the unnormalized quantity of each spot,
2. Scaling factor is a factor you select to give a meaningful value (see below), and
3. Normalization factor is calculated for each gel based on the normalization method you select (Total of Spots in Analysis Set, Total of Valid Spots, Total Density, etc.)

Note: The pipetting error compensation factor is only available if you select Counts Loaded from Gel Record. See below.

7.3.b Scaling Factor

The scaling factor is a constant by which the normalized quantity is multiplied to give a more meaningful value. A common scaling factor is 10^6 , giving units of parts per million (PPM).

PDQUEST User Guide

Under Scaling, specify the scaling factor that you would like to use by clicking on PPM (x 1000000), Percent (x 100), or User Specified. If you select User Specified, enter your own scaling factor and units in the appropriate fields.

7.3.c Normalization Methods

Total of Spots in Analysis Set

This normalization method involves identifying a group of protein spots that are present in all matchset members and are conserved (i.e., their quantities should not vary between samples). These are proteins that you know are expressed in a very consistent manner at all times within a cell, with little or no regulation. These are commonly referred to as "housekeeping" proteins.

First, you create an Arbitrary Analysis Set containing several of these housekeeping proteins from different regions of the matchset standard.

Note: To create an Arbitrary Analysis Set, make sure your matchset standard is selected, then go to the Database menu and select Analysis Sets > Create Analysis Set. In the dialog box, click on the Type button and select Arbitrary from the list. In the Arbitrary Analysis Set dialog, enter a name for the set and select the No Spots option. Then click on Go. This creates an empty analysis set. To add spots to this set, choose Select Analysis Set from the Database > Analysis Sets submenu and select the name of the set you just created from the list. With the set selected, select Add Spot To Set from the Database > Analysis Sets submenu and click on an appropriate (i.e., housekeeping protein) spot. The spot will be highlighted, indicating that it has been added to the set. If you want to remove a spot from the analysis set, select the set, then select Remove Spot From Set from the Database > Analysis Sets submenu and click on a highlighted spot.

See section 9.1 for more information on creating analysis sets.

In the Normalization dialog box, select the Total of Spots in Analysis Set option, and select your set of housekeeping proteins from the pulldown list.

The total quantity of the spots in this analysis set will be displayed for each matchset gel in the table at the bottom of the dialog box.

Mrpl Values/Spot Quantitation

The raw quantity of each spot in a gel will be divided by the total quantity of the spots in the analysis set in that gel, according to the normalization formula.

Note: The advantage of selecting multiple housekeeping proteins is that occasional minor variations in the size and / or intensity of a particular spot will be averaged out when the group of spots is considered is a whole.

Total of All Valid Spots

In this normalization method, the raw quantity of each spot in a member gel is divided by the total quantity of all the spots in that gel that have been included in the Standard (see the normalization formula).

Select the Total of All Valid Spots option button to activate this normalization method. This method assumes that few protein spots change within the experiment, and that the changes average out across the whole gel. It can be useful if very little information is known about the possible source of sample variation.

Total Density in Gel Image

In this method, the raw quantity of each spot in a member gel is divided by the total intensity value of all the pixels in the image (see the normalization formula).

Select the Total Density in Gel Image option button to activate this normalization method. This model assumes that the total density of an image (i.e., background density plus spot density) will be relatively consistent from gel to gel. It can be useful if very little information is known about the possible source of sample variation.

Specified Value

In this method, the raw quantity of each spot in a member gel is divided by a user-definable normalization factor for that gel.

PDQUEST User Guide

Select the Specified Value option and enter a normalization value for each gel in the table at the bottom of the dialog box. This value will be used as the denominator in the normalization formula.

This model is useful if you want to correct for known variations in samples across gels (e.g., the amount of cells that went into an experiment, the amount of extract loaded, the length of time the images were detected, etc.).

Counts Loaded from Gel Record

In this method, the raw quantity of each spot in a member gel is divided by the total counts loaded into that gel, as specified in the Gel Record when you performed calibration.

Note: This method is only available if you have performed calibration on your matchset gels. See Chapter 4.

The counts loaded in each gel are listed in the table at the bottom of the Normalization dialog; this value will be used as the denominator in the normalization formula. Select the Counts Loaded from Gel Record option button to activate this normalization method.

This method assumes that few protein spots change within the experiment, and that the changes average out across the whole gel. This model can be useful if you want to correct for known variations in the amount of sample loaded across gels, as measured prior to loading with a scintillation counter.

If you select this method, you can compensate for pipetting error in all matchset members by selecting the Pipetting Error Compensation checkbox.

When this checkbox is selected, the software will attempt to correct for minor systematic errors inherent in scintillation counting, pipetting, and exposing films. The normalization process is based on the assumption that these minor errors affect all proteins equally. It also assumes that few protein spots change within the experiment or that when changes do occur, they do not alter the random distribution from one sample to the next. (See section J.4. for additional information about this function.)

When you select this checkbox, a pop-up box will display the calculated compensation value for each matchset member. These values will also appear in the in the table at the bottom of the Normalize dialog.

8. Spot Numbering and Annotations

PDQUEST has features for numbering and annotating your spots. These are located on the Database menu.

Every matched spot is assigned a Standard Spot (SSP) number. SSP numbers are increasingly important when you begin comparing spots between different matchsets using higher level matchsets (see the next chapter).

Besides the standard SSP numbering system, an alternative database spot number (DSN) system is also available. The advantage of the DSN system is that it indicates the general location of each spot, based on the MrpI grid.

Notes, comments, and other pertinent information can be kept in a file associated with each protein spot. The functions needed to enter and edit this data are found in the Annotations menu.

8.1 Standard Spot Numbers

Standard spot (SSP) numbers are unique numbers assigned to each spot in the matchset Standard. In general, spots in the lower left region of the image have low SSP numbers, while spots in the upper right region have high SSP numbers.

The commands for displaying the SSP numbers are located on the Database > Spot Numbers submenu.

PDQUEST User Guide

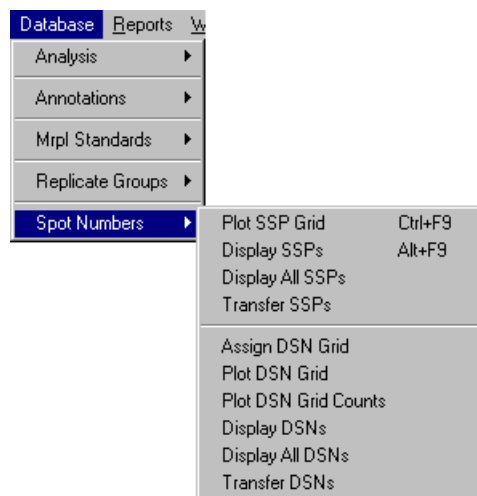


Fig.8-1. Spot Numbers submenu.

8.1.a Assigning SSP Numbers

SSP numbers are automatically assigned when you create a matchset Standard. To create SSP numbers, PDQUEST first divides the Standard into rows and columns, placed so that there are approximately the same number of spots in each section and that no section has over 100 spots. The first two digits of a spot's SSP number correspond to its section's X and Y coordinates. Spots within a section are numbered sequentially.

8.1.b Displaying SSP Numbers

Select Plot SSP Grid from the Database > Spot Numbers submenu to display the SSP grid on the Standard. Numbered vertical and horizontal grid lines will appear on the Standard.

The grid is helpful for locating spots by SSP number. However, it is an approximation and should only be used as a guide. If enough spots are added to a single grid block to use up the 99 available numbers, numbers will be assigned from other blocks. SSP numbers are always four digits.

Numbering and Annotations

Display SSPs displays SSP numbers next to every spot in the active window. You may need to zoom in on the image to be able to read them.

Display All SSPs displays SSP numbers next to every spot in every open window.

Transfer SSPs is used to copy the SSP numbers from one higher level matchset to another.

8.2 Database Spot Numbers

Database spot numbers (DSNs) are an alternative to the SSP numbers. The DSN system provides a more flexible way to number protein spots in a Standard. A coordinate grid can be established for each axis: Mr and pI. The grid coordinates are identified with letters (A–R, excluding I and O because they look like numbers). The letters are intended to represent specific Mr and pI ranges so that the numbering between different gels can be comparable. These ranges are as follows:

Mr values:	0-10=A	87-117=J
	11-14=B	118-158=K
	15-19=C	159-215=L
	20-25=D	216-293=M
	26-34=E	294-398=N
	35-46=F	399-541=P
	47-63=G	542-736=Q
	64-86=H	737-1000=R

pI values:	0-2=A	7.0-7.5=J
	2-3=B	7.5-8.0=K
	3-4=C	8.0-8.5=L
	4-5=D	8.5-9.0=M
	5.0-5.5=E	9-10=N
	5.5-6.0=F	10-11=P
	6.0-6.5=G	11-12=Q
	6.5-7.0=H	12-14=R

PDQUEST User Guide

DSNs have the format ABnn. “A” is the Mr grid letter, “B” is the pI grid letter, and “nn” is the number of the spot within that section of the grid.

8.2.a Assigning DSNs

If you have already entered Mr and pI values for known spots, you can create a DSN grid using this data by selecting Assign DSN Grid from the Database > Spot Numbers submenu. This automatically assigns a DSN number to each spot in the Standard.

To display the grid (this does not occur automatically when the grid is created), click in the Standard window and select Plot DSN Grid. A grid of vertical and horizontal axes will appear, with the letter assignments of each axis displayed in the window margins.

To see the spot assignments, select Display DSNs. To read the numbers, you will probably need to zoom in on a smaller region of image.

Transfer DSNs is used to copy the DSN numbers from one higher level matchset to another (see Chapter 9).

8.3 Finding Spots

You can find spots by SSP or DSN number.

Select Find Spot from the Spots menu (SHIFT+F9), click in the member or Standard you want to search, and enter the SSP or DSN number in the pop-up box.

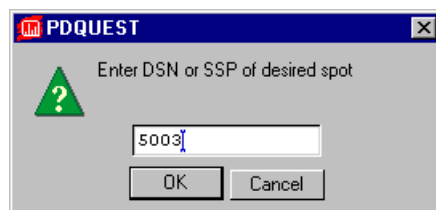


Fig.8-2. Find Spot pop-up box.

Numbering and Annotations

Click on OK to complete the search. The spot will appear highlighted and numbered on the image.

8.4 Annotations

You can add notes, comments, and other pertinent information to a file associated with each protein spot. As the size and complexity of your spot database increases, annotations can be an essential tool for keeping track of experimental data.

The functions for creating, displaying, and printing annotations are located on the Annotations submenu of the Database menu and on the Annotations toolbar.

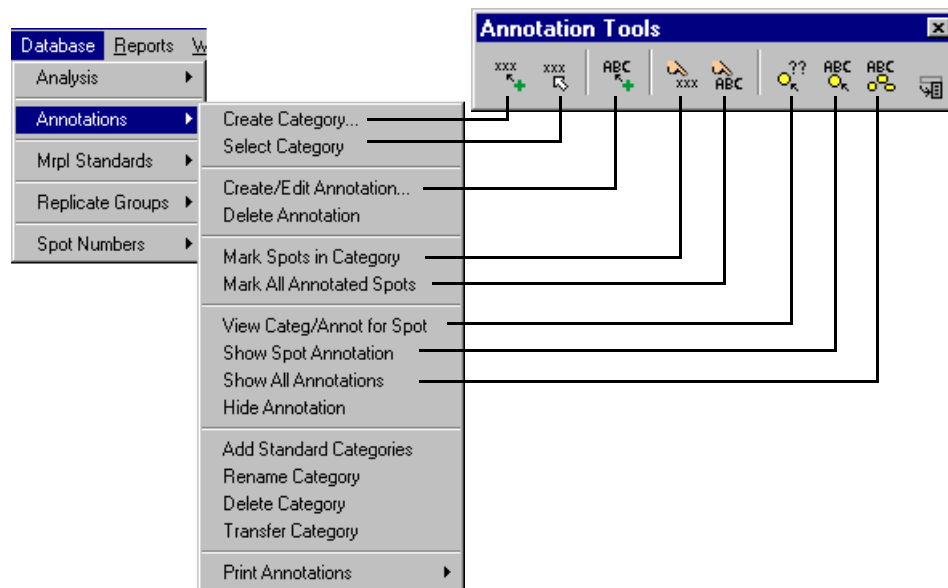


Fig.8-3. Annotation tools.

PDQUEST User Guide

8.4.a Creating Annotations

Before you can begin creating annotations for individual protein spots, you first need to create the category of information to be entered.

Some examples of annotation categories are: Protein Name, Amino Acid Composition, Binding Properties, Translational Regulation, etc.

Creating a Category

To create a category, select Create Category from the Database > Annotations submenu or toolbar. A dialog box will ask you to enter a name and a description for the category.

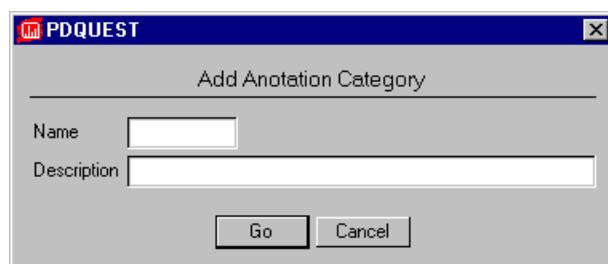


Fig.8-4. Enter a name for the new annotation category.

A description of the category is optional. Click on the Go button when you are done.

Once a category has been created, it is automatically active. You can enter spot annotations for the category without having to select it first.

Selecting an Existing Category

To enter information in an existing category, first open the category by choosing Select Category from the Database > Annotations submenu or toolbar and selecting the category name.

If annotations have been entered for spots in this category, the spots will be highlighted when the category is selected.

Numbering and Annotations

Creating, Editing, and Deleting Annotations

Once a category is active, you can enter annotations for individual spots. Select Create/Edit Annotation and click on a spot to open the annotations form.

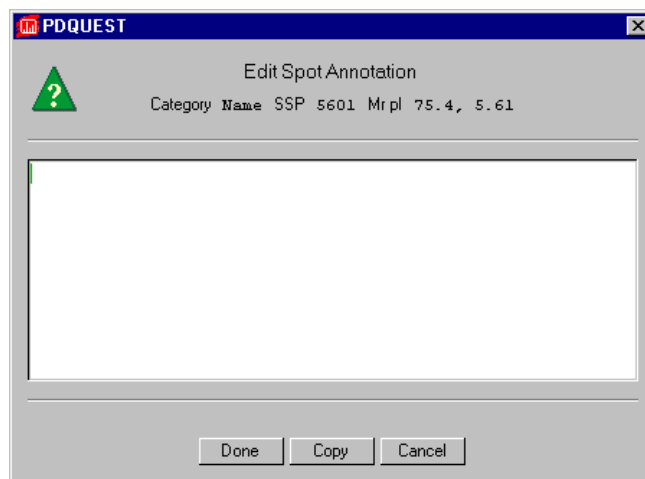


Fig.8-5. Edit Annotation dialog box.

When you have finished entering the information for that spot, click on the Done button at the bottom of the box. The box will disappear and the spot will become highlighted.

Repeat these steps for all the spots you want to annotate in the selected category.

To change an annotation, simply click on the annotated spot with the Create/Edit Annotation tool and the existing annotation will be displayed in the edit box.

To delete an annotation, select the category to delete from, then select Delete Annotation from the Database > Annotations submenu. You will be prompted to confirm the deletion.

PDQUEST User Guide

Deleting/Renaming a Category

To delete a category, select Delete Category from the Database > Annotations submenu and select the category to delete from the list.

To rename a category, select Rename Category and select the category name from the list. Enter a new category name in the pop-up box.

Creating Standard Categories

Create Standard Categories is a function that adds 38 predefined categories to your annotation database. These categories are then available for you to enter and edit annotation information about spots in your matchset.

Select Create Standard Categories from the Database > Annotations submenu. The categories will be added to the matchset. To view the categories, use the Select Category command.

Saving Annotations and Categories

Annotations and categories for a matchset are saved when you save the matchset. No separate operation is necessary.

8.4.b Displaying Annotations

Several functions are available for displaying annotation data. These are located on the Database > Annotations submenu and the Annotations toolbar.

Select Mark Spots in Category to highlight all the spots that have information entered in the selected category.

Select Mark All Annotated Spots to highlight all the spots in the Standard that have data entered in any category, regardless of whether the category is selected.

Show Spot Annotation will display the annotation information in the selected category for an individual spot that you click on. Select Show Spot Annotation, then click on a spot.

Numbering and Annotations

Show All Annotations will display the annotation information in a selected category for every spot in the Standard.

View Categ / Annot for Spot displays the annotation categories that have information for a particular spot, as well as the specific annotations in those categories. Select this command from the submenu or toolbar, then click on a spot of interest. A list of annotation categories for that spot will be displayed. Select a category from the list to display the annotation information in a pop-up box next to the spot.

To conceal an annotation, select Hide Annotation and click on the spot displaying the annotation.

8.4.c Printing Annotations

You can print your spot annotations using the functions in the Print Annotations submenu of the Database > Annotations submenu.

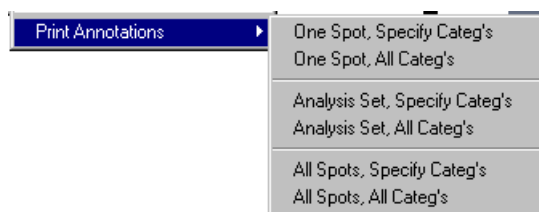


Fig.8-6. Print Annotations submenu.

Using the commands in this submenu, you can print annotations for individual spots, the protein spots in an Analysis Set, or all the spots in the Standard. You can also print annotations for particular categories or all categories.

Individual Spots

To print the annotations for a particular spot, select One Spot, Specify Categ's or One Spot, All Categ's, then click on the spot in the Standard.

If you select One Spot, Specify Categ's, a pop-up list will prompt you to select annotation categories to print. Once you have made your selections, click on

PDQUEST User Guide

the Done button below the list. The print dialog box will appear, in which you can select your printer settings.

If you select One Spot, All Categ's, the print dialog box will appear immediately.

Analysis Sets

To print annotations for spots in a specific Analysis Set, select Analysis Set, Specify Categ's or Analysis Set, All Categ's from the Annotation Report submenu.

A list of Analysis Sets will be displayed. Click on the name of the Analysis Set whose annotations you want to print.

If you selected Analysis Set, Specify Categ's, a pop-up list of categories will appear. Make your selection(s), then click on Done. The print dialog box will appear.

If you selected Analysis Set, All Categ's, the print dialog box will appear as soon as you have chosen an Analysis Set.

All Spots

All Spots, Specify Categ's prints selected annotation categories for all the spots in the Standard.

All Spots, All Categ's prints the annotations in all available categories for all annotated spots in the Standard.

9. Databasing Spots

As you accumulate more information about your proteins, it becomes increasingly important to store the data in an organized and accessible format. PDQUEST provides several functions to build and query such protein spot databases.

You can create Analysis Sets containing statistically and biologically related spots. The spots in an Analysis Set meet quantitative or qualitative criteria that you specify. There are six different types of Analysis Sets that you can create: Quantitative, Qualitative, Arbitrary, Boolean, and Statistic.

An important part of databasing is the ability to compare data between experiments. PDQUEST allows you to create a database from different experiments through the use of higher level matchsets.

9.1 Analysis Sets

The Analysis Set functions in PDQUEST are designed to help you create and study sets of proteins that are statistically and scientifically meaningful.

It is virtually impossible to analyze the gel or experiment as a whole, so PDQUEST provides a way to break down all the information on the gel into smaller, more manageable pieces. This is done by creating Analysis Sets (i.e., groups of spots that satisfy criteria that you specify).

There are six different kinds of Analysis Sets that you can create:

Quantitative Analysis Sets: Include spots whose quantitation has increased or decreased “x” fold (where x is a user-determined value) or whose quantitation has changed beyond or within the fold change bounds that you specify.

Qualitative Analysis Sets: Include spots whose quantitation is “x” times above the minimum level of detection in one gel but below the minimum level of detection in

PDQUEST User Guide

another gel—i.e., spots that have been turned “on” or “off.”

Arbitrary Analysis Sets: Are composed of spots that you choose. They begin as empty sets and you select the proteins that are to be members of them.

Boolean Analysis Sets: Contain a group of spots formed by the combination of two other sets. For example, Analysis Set C could be made up of the spots that are present in both Analysis Set A and Analysis Set B.

Statistical Analysis Sets: Include spots that are found to be significant according to the statistical test that you indicate.

The tools for creating and editing these sets are located on the Database > Analysis submenu and Analysis toolbar.

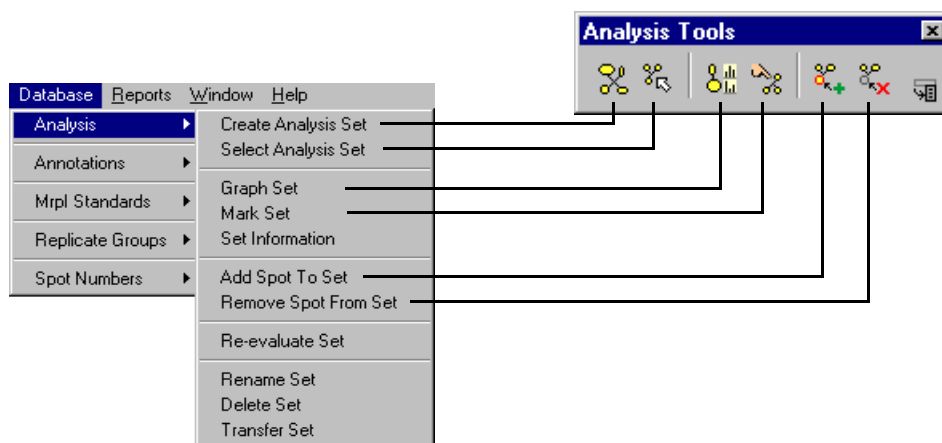


Fig.9-1. Analysis Set tools.

9.1.a Creating an Analysis Set

With a matchset open, select Create Analysis Set from the Database > Analysis submenu.

Databasing Spots

The Create Analysis Set dialog box will open.

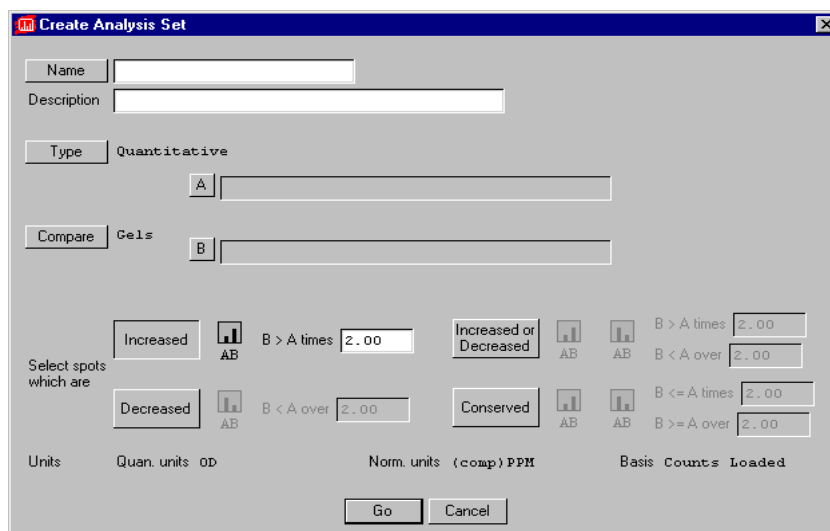


Fig.9-2. Example of a Create Analysis Set dialog.

If you want to update an Analysis Set that already exists, click on the Name button at the top of the Analysis Set dialog box. A list of existing sets will be displayed. Select the name of the set you want to update by clicking on it. The dialog box will be reconfigured and the information about that Analysis Set will be displayed.

Selecting the Type of Analysis Set

The first step in creating the Analysis Set is selecting the type you want to create. In the Analysis Set dialog box, click on the Type button.

A list of Analysis Set types will appear. Click on a type to select it. (See the following sections for descriptions of the different types of Analysis Sets.)

The list will disappear and the Analysis Set dialog box will reconfigure to display the test you chose. You will see the name of your selection displayed next to the Type button.

PDQUEST User Guide

Selecting the Operands to Compare

All Analysis Sets (except Arbitrary) involve comparisons: How does a protein or group of proteins (in a gel or replicate group) behave under one set of conditions versus another? To answer this question, you need to specify which sample should be compared with which other sample.

First, indicate whether you want to look at the quantitation of individual gels or the average quantitation of groups of gels. Click on the Compare button and select either Gels or Groups from the list displayed. Your choice will appear to the right of the Compare button.

Next, click on the A button and choose a sample from the list displayed. (If you selected Gels to compare, a list of all the gels will be displayed. If you selected Groups, a list of replicate groups will be displayed.)

After clicking on your choice, you will see it displayed next to the A button.

Repeat this procedure with the B button to enter the second sample.

Displaying the Spots that Are Included in a Set

After you have entered the settings for an Analysis Set, click on the Go button to create the set. The dialog box will disappear. On the Standard, the spots that are included in the Analysis Set will be highlighted.

Once an Analysis Set is active, use the function Mark Set (in the Analysis Sets submenu) to highlight all the spots contained in that Analysis Set.

Another way to display the proteins in an Analysis Set is with the Graph Set function (in the Database > Analysis submenu). This function will display histograms of the spot quantitation along the left and right margins of the screen. (For more information on histograms, see section 10.1)

To clear the Analysis Set display, click on Hide Overlays on the main toolbar.

Opening/Deleting/Editing an Existing Analysis Set

To open an existing set, choose Select Analysis Set from the Database > Analysis submenu and select the Analysis Set from the list.

Databasing Spots

You can also open an existing Analysis Set from within the Analysis Set dialog box. Click on the Name button at the top of the dialog and select the name of the set you want to load from the list. The settings and information for that set will be displayed in the dialog box.

To delete a saved Analysis Set, select Delete Set from the Database > Analysis submenu and select the set to delete from the list of names.

To rename an Analysis Set, select Rename Set from the Database > Analysis submenu and select the set to rename from the list. Enter a new name for the set in the pop-up box.

To recalculate the members of an Analysis Set (e.g., if you have added or matched more spots), select Re-evaluate Set from the Database > Analysis submenu. The set will be recalculated and the members will be displayed on the Standard.

9.1.b Quantitative Analysis Sets

Quantitative Analysis Sets allow you to perform fold change analyses. To create a Quantitative Analysis Set, select Database > Analysis > Create Analysis Set. Quantitative Analysis Set is the default configuration of the Analysis Set dialog box.

PDQUEST User Guide

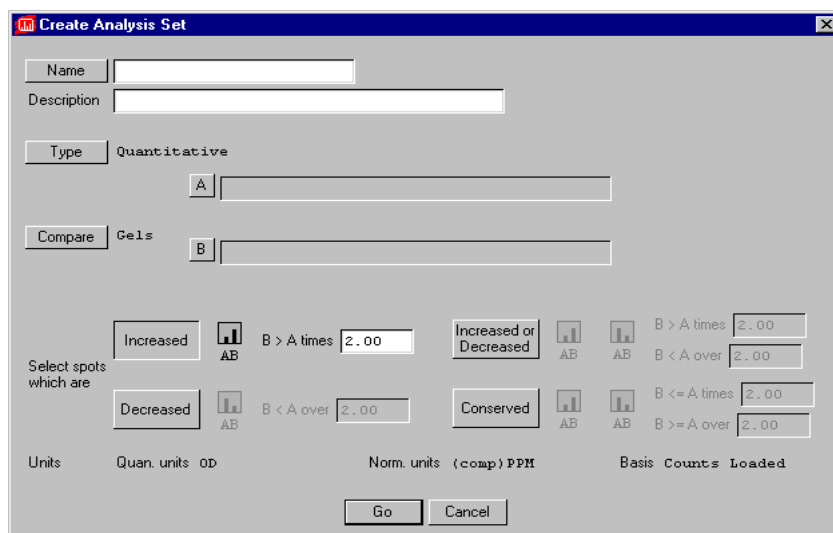


Fig.9-3. Quantitative Analysis Set dialog.

After entering your selections for comparison (e.g., control gels versus experimental gels), you need to decide what kind of quantitative changes you want to examine. You can select Increased, Decreased, Increased or Decreased, or Conserved.

Next to each choice is a diagrammatic representation of the type of change the test will find and an explanation of the type of quantitative criteria a spot must meet to be included in the Analysis Set.

Do you want to look for a protein whose quantitation went up under the experimental conditions as compared to the control? If so, click on the Increased button next to the “Select spots which are” prompt.

“A” represents all matched spots in sample A and “B” represents the corresponding matched spots in sample B.

Comparisons are made using spot quantity only. Positional comparisons, etc., are not possible.

Databasing Spots

The comparison is performed as follows: A threshold value is established by multiplying and/or dividing A by the factor you specify. B is then compared to the threshold.

For example:

1. If you ask to create an Analysis Set by clicking on the Increased button and enter a factor of 2.0, the set will include spots whose quantity in B is at least twice that of the corresponding spot in A.
2. If you click on the Increased button and enter a factor of 0.5, the Analysis Set will include spots whose quantity in B is greater than half the quantity of the corresponding spot in A.
3. Clicking on the Decreased button and entering a factor of 2.0 will find spots whose quantity in B is less than half the quantity of the corresponding spot in A.
4. Creating an Analysis Set by clicking on Decreased and entering a factor of 3.0 finds all the spots whose quantity in B is less than one-third the quantity of the same spot in A.
5. If you are interested in finding the spots whose quantitation has significantly changed either way (either increased or decreased), select the Increased or Decreased option. The proteins found in this test have increased in quantitation above the upper bound or decreased below the lower bound. If you enter an upper bound of 2.0 and a lower bound of 3.0 you will find spots whose quantity in B is more than twice the quantity in A or less than one-third the quantity in A.
6. To find spots whose quantitation in B has not changed significantly from A, click on the Conserved button. Enter factors to define the upper and the lower bounds. If, for example, you enter an upper bound factor of 3.0 and a lower bound factor of 2.0, the Analysis Set will include spots whose quantity in B is between three times and one-half times the quantity of those spots in A.

After entering the fold-change factor(s), the corresponding diagram will automatically update to show you the quantitative relationship these criteria define.

PDQUEST User Guide

Click on the Go button at the bottom of the Analysis Set dialog box to create the set you have specified.

Be sure to enter a name for the Analysis Set. If you forget, a pop-up box will remind you to do so.

9.1.c Qualitative Analysis Sets

Qualitative Analysis Sets contain spots that were not detected in one gel but were detected in the other. This test would detect a protein that was expressed under experimental conditions but was not expressed in the control sample, or visa versa.

Note: We strongly recommended that all the spots in such an Analysis Set be double-checked for accuracy. It is possible that a protein might actually be present in both gels but was missed during matching.

Creating a Qualitative Analysis Set

To create a Qualitative Analysis Set, select Database > Analysis > Create Analysis Set, click on the Type button in the dialog box, and select Qualitative from the list displayed. The dialog box will change to a Qualitative Analysis Set configuration.

Databasing Spots

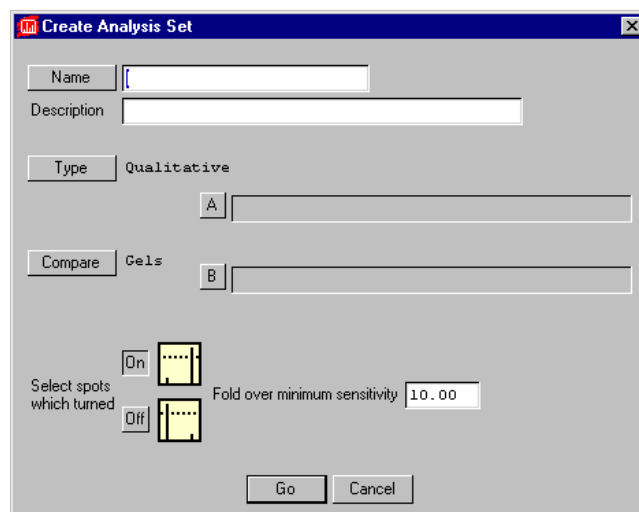
The image shows a software dialog box titled "Create Analysis Set". It has a blue title bar with a red icon on the left and a close button on the right. The main area is light gray. At the top, there are two text input fields: "Name" and "Description". Below these is a "Type" button, which is currently set to "Qualitative". Underneath, there are two more input fields labeled "A" and "B", preceded by a "Compare" button. The "Compare" button is currently set to "Gels". At the bottom left, there are two radio buttons labeled "On" and "Off", each with a small gel electrophoresis diagram next to it. To the right of these is a text input field labeled "Fold over minimum sensitivity" with the value "10.00" entered. At the very bottom are two buttons: "Go" and "Cancel".

Fig.9-4. Qualitative Analysis Set dialog.

As described above, choose your control and experimental gels.

To find proteins that were not present in the control sample (Gel A) but were expressed under experimental conditions (Gel B), click the On button next to the “Select spots which turned” prompt.

For spots present only in the controls (Gel A) and not in the experimental samples (Gel B), click on the Off button.

Finally, you need to tell PDQUEST what “turned on” or “turned off” means. How many times greater than the minimum detectable spot should the quantitation be for a spot to be included in this Analysis Set? Should it be 5-fold greater than the minimum? If so, enter 5.00 next to the Fold over minimum sensitivity prompt.

Click on the Go button to find all the spots that fit the criteria that you specified. They will be highlighted on the Standard.

PDQUEST User Guide

9.1.d Arbitrary Analysis Sets

Arbitrary Analysis Sets are composed of any group of spots that you choose. To create an Arbitrary Analysis Set, select Database > Analysis > Create Analysis Set, click on the Type button in the dialog box, and select Arbitrary from the list displayed.

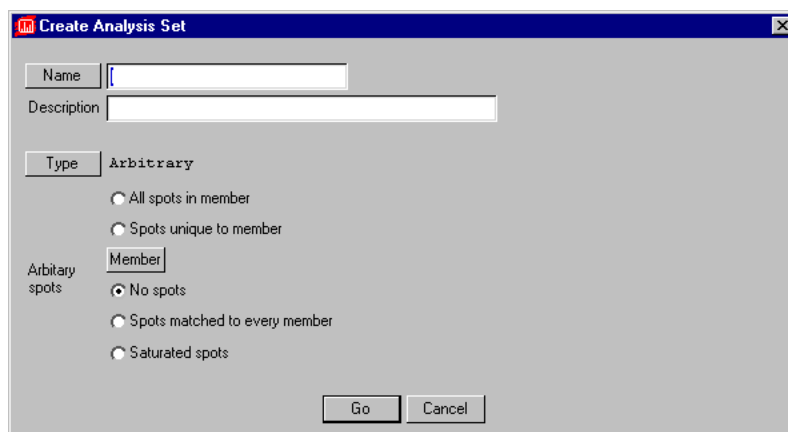


Fig.9-5. Arbitrary Analysis Set dialog.

The first option (“All spots in member”) will create an Analysis Set containing all the spots in a selected member that are matched to the Standard. Click on the Member button to display a list of the gels from which to choose. Click on your choice and it will be displayed next to the Member button. The Analysis Set will include all the spots in the member you selected that are matched to the Standard.

“Spots unique to member” will create an Analysis Set of those spots found only in the member that you specify and not in any of the other members of the matchset. Click on the Member button to display a list of gels from which to choose. Click on your selection and it will be displayed on the form next to the Member button.

“Spots matched to every member” will create an Analysis Set of spots that are matched in all the members.

Databasing Spots

“Saturated spots” will create an Analysis Set of all spots that are saturated in any member.

The “No spots” option creates an empty Analysis Set to which you can add any spots you choose. See the following section for detailed instructions.

Adding Spots to/Removing Spots from an Arbitrary Analysis Set

To add spots to an empty Analysis Set, you first select the Arbitrary Analysis Set you want to add to, then you click on the spots you want to add to it.

Note: Spots may be added on a spot-by-spot basis only to Arbitrary Analysis Sets.

Choose Select Analysis Set from the Database > Analysis submenu to display a list of Analysis Sets. Select the Arbitrary Analysis Set you want to add spots to.

With the set selected, select Add Spot To Set from the Database > Analysis submenu and click on a spot. The spot will be highlighted, indicating that it has been added to the Analysis Set. The name of the Analysis Set to which it has been added is displayed in the upper left corner of the window.

If you want to remove a spot from the Analysis Set, select the set, then select Remove Spot From Set from the Database > Analysis submenu and click on an added spot. The spot will be unhighlighted, indicating that it has been removed.

9.1.e Boolean Analysis Sets

Boolean Analysis Sets are formed by combining two previously defined Analysis Sets using Boolean operators. You can create Boolean Sets that include:

1. The intersection of two sets of spots (proteins present in both Analysis Set A and Analysis Set B).
2. The union of two Analysis Sets (spots in set A or in set B).
3. Spots that are unique to set A plus spots that are unique to set B.

PDQUEST User Guide

4. Spots that are found only in set A.
5. Spots that are found only in set B.

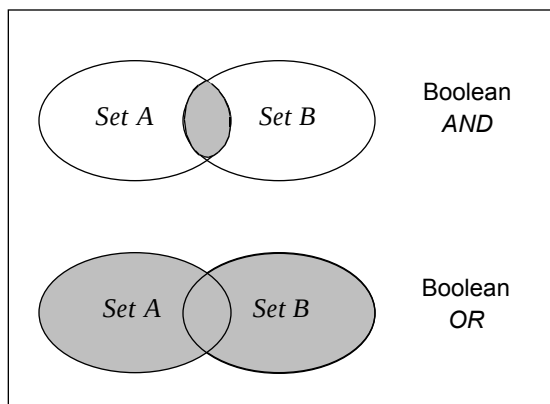


Fig.9-6. Venn diagrams depicting the AND and OR Boolean operations.

To create a Boolean Analysis Set, select Database > Analysis > Create Analysis Set, click on the Type button in the dialog box, and select Boolean from the list displayed.

Databasing Spots

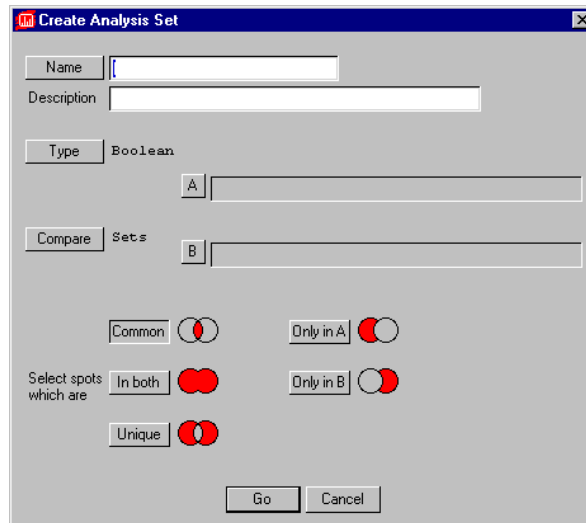


Fig.9-7. Boolean Analysis Set dialog.

To make a Boolean Analysis Set, select the two Analysis Sets you want to compare to one another from the lists displayed when you click the A and B buttons.

Enter the name and description to be given to the Analysis Set.

Select the Boolean operator to be used to create the set by clicking on one of the buttons next to the “Select spots which are” prompt.

When you click on the Go button, the dialog box will disappear and the spots that meet the criteria that you specified will be highlighted in the Standard image.

9.1.f Statistical Analysis Sets

Statistical Analysis Sets are composed of spots from replicate groups found to be significant according to the statistical test that you specify.

To create a Statistical Analysis Set, select Database > Analysis > Create Analysis Set, click on the Type button in the dialog box, and select Statistic

PDQUEST User Guide

from the list displayed. The dialog box will be reconfigured to let you select the replicate groups to be compared and the statistical test to be performed.

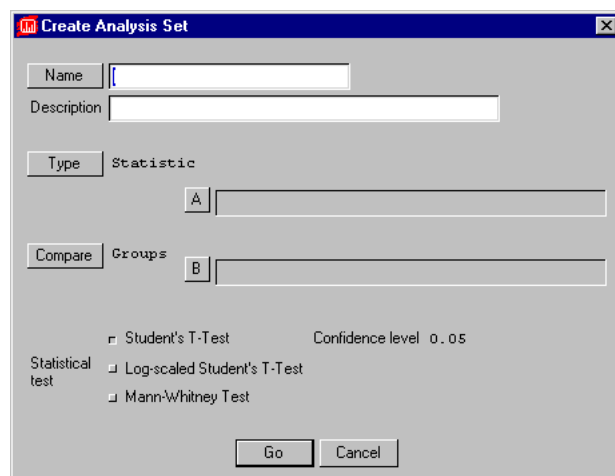


Fig.9-8. Statistic Analysis Set dialog.

Statistical Analysis Sets only work with groups, because the tests require multiple gels. The Student's T Test and Log Student's T Test require at least two gels per group. The Mann-Whitney Test requires at least four gels per group.

Choose the groups by clicking on the A and B buttons as described in previous sections.

Indicate the statistical test you want to perform by clicking on the button next to the name of the test.

Enter a name for the Analysis Set at the top of the dialog box and click the Go button.

Significant protein spots will be highlighted on the Standard.

9.2 Higher Level Matchsets

PDQUEST allows you to create matchsets from matchsets, so you can build a database by comparing and summarizing the results of a series of experiments.

“Higher level” matchsets are matchsets created from the Standard images of other (“lower level”) matchsets. A lower level matchset might be a matchset created from your original gels (a “level-1” matchset), or it could be another higher level matchset (that is nevertheless “lower” than the one you are creating). Many experimental protocols can be designed to take advantage of higher level matching.

PDQUEST allows you to edit spots in the member images of lower level matchsets and have those changes take immediate effect in the higher level matchset. Having lower level members accessible from the higher level Standard helps you to visualize the original gels that contain the data analyzed in the higher level Standard.

For example, when working with a higher level matchset, you may notice that a spot is present in one member and not the other. Your next step would be to go back to the level-1 matchset from which that member came, to see if any of the gels in the level-1 matchset contain the spot of interest. If you found the spot, you would add it to the Standard of the level-1 matchset so that it could be subsequently added to the Standard of the higher level matchset. Being able to access and edit the members of the level-1 matchset from the higher level matchset eliminates the need to unload the higher level matchset, load the level-1 matchset to edit the spot, unload the level-1 matchset and then reload the higher level matchset.

One experiment may find all the proteins that change during cell growth in culture, while another may identify nuclear proteins by cell fractionation. These two experiments can be analyzed individually in level-1 matchsets, each with their own Standard images. These two Standards can then be combined into a higher level matchset with its own Standard. Data from the two experiments are now linked through a third Standard. This higher level matchset may shed light on growth-related proteins found in the nucleus.

PDQUEST User Guide

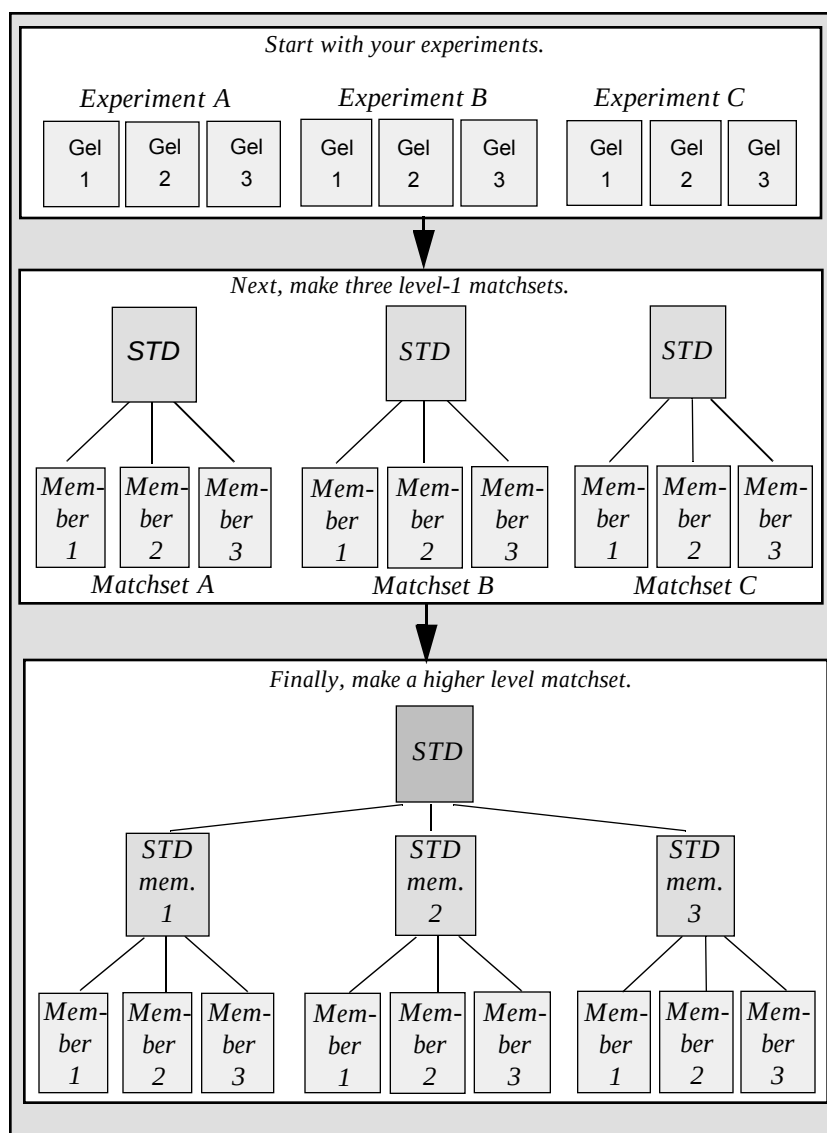


Fig.9-9. Creating a higher level matchset: In this example, the higher level matchset is made from 3 lower level matchsets, each containing 3 members.

9.2.a Creating a Higher Level Matchset

As mentioned previously, higher level matchsets consist of lower level matchsets. To make a higher level matchset, select Create Matchset from the Create/Edit Matchset submenu of the Match menu. The Create Matchset dialog box will open.

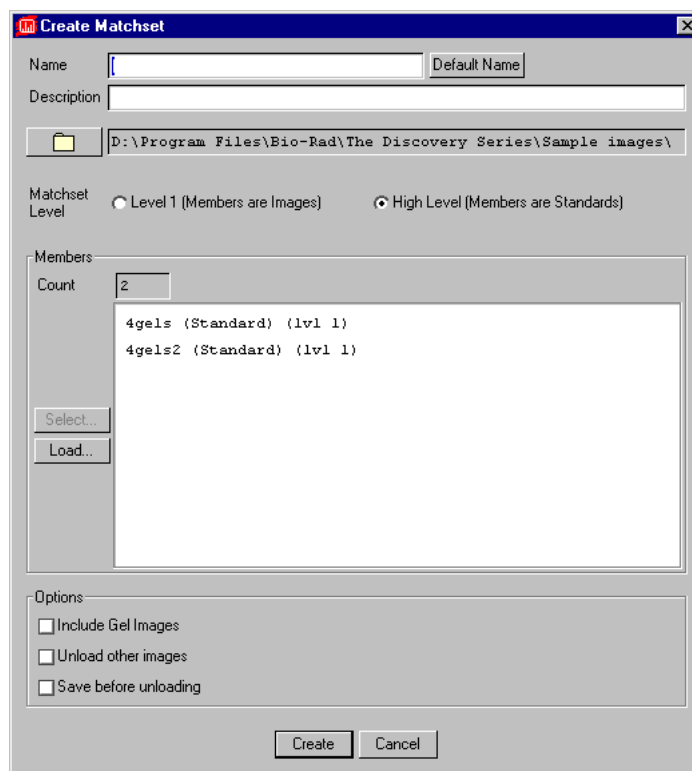


Fig.9-10. Dialog box for creating a higher level matchset.

Higher level matchsets are created in the same way as lower level matchsets (see Chapter 6), with the following exceptions:

- Next, to the Matchset Type prompt, select High Level instead of Low Level.

PDQUEST User Guide

- You cannot include Gel Images in a higher level matchset, because higher level matchsets are made up of Standards, which do not have Gel Images.
- The members that you select for a higher level matchset are lower level matchsets rather than Gel Spots files.

To select the matchsets for inclusion in your higher level matchset, click on the Load button. This will open the Select New Member dialog.

Select the matchsets you want to include by clicking on them. Once you have chosen all the matchsets that you want to include, click on the Add Selected button.

The names of those matchsets will appear in the Members list in the Create Matchset dialog box. When you are satisfied with your choices, click on the Create button.

You will be asked to select a Standard for the higher level matchset (see below). This is similar to the Standard in a lower level matchset.

The screen will be configured into enough windows to display each member of the matchset. Each member will be assigned to a window.

9.2.b Editing Higher Level Matchsets

You can add members to, remove members from, and reorder members in higher level matchsets just like lower level matchsets. These commands are described in section 6.3

9.2.c Creating a Standard for a Higher Level Matchset

As with lower level matchsets, each higher level matchset has a Standard. You can select the member to use as the template for the Standard when you create the matchset, or you can select it later using the Make Standard function.

Select Make Standard from the Create/Edit Matchset submenu of the Match menu, then click on the matchset member that will be the template for the Standard.

Databasing Spots

To assign the same SSP numbers to the higher level matchset Standard as exist in the Standard template, click on Yes in the pop-up box. For more information about retaining and transferring SSP numbers, see below.

9.2.d Matching and Editing Higher Level Matchsets

Once your higher level matchset has been created, you need to edit and match its members. The same protocols used with lower level matchsets can be applied to higher level matchsets.

Landmarks need to be entered before member matchsets can be matched to each other. Members should be carefully reviewed for spots that need to be added to the higher level matchset.

For more information about matching and editing, see Chapter 6.

9.2.e Transferring DSNs or SSP Numbers

Each matchset has unique DSN or SSP numbers assigned to proteins that may not correspond to those in other matchsets. So how can you correlate data between two or more matchsets?

Begin by creating a higher level matchset consisting of the two lower level matchsets and creating its Standard.

Apply the numbering scheme from one of the members to all other members of the matchset (SSPs or DSNs).

Depending on the numbering system, select Transfer SSPs or Transfer DSNs from the Database > Spot Numbers submenu.

When you use either of these functions, you will first be asked to select the image from which the numbers should be copied. Select the image, then indicate the matchset member to which the numbers should be applied.

Confirm the transfer of numbers by clicking on the Yes button.

If you transfer numbers between all the members, the spot numbers for the matched spots in the lower level matchsets will correspond to each other and to the spots in the higher level matchset as well. This makes it easier to correlate data from the reports of the member matchsets.

PDQUEST User Guide

9.2.f Transferring Mr and pI Data

The function Transfer MrpIs (Database > MrpI Standards > Transfer MrpIs) allows you to apply Mr and pI values from one Standard to the other members of a higher level matchset.

When you select this function, you will first be asked to select the matchset member the MrpI data should be transferred from.

Next, choose the matchset members to which the data should be applied.

Confirm the transfer of the MrpI data by clicking on the Yes button.

9.2.g Transferring Annotations

If annotations are entered in the member matchsets of a higher level matchset, you can transfer them all to the higher level matchset Standard. This makes them all available from a single file. You can do this with the Transfer Category function.

A selected annotation category can be transferred between members in a higher level matchset, but only for spots that are matched between both members. If no spots are matched between the two, only the category name is transferred (no data is transferred).

Before you can transfer annotations, you need to select the annotation category. From the Database > Annotations submenu, choose Select Category, then click on the Standard with the annotations you want to transfer. A list of the available annotation categories for that Standard will open.

Select the category you want to transfer. This category will become active and ready to be transferred.

From the Annotations submenu, select Transfer Category. From the pop-up list, select the Standard from which the annotations should be transferred. Next, select the member to which the annotations should be transferred.

Complete the annotation transfer by clicking on the Transfer OK button.

A message in the status bar will indicate the number of spots with annotations in that category and how many were transferred.

Databasing Spots

To see to which protein spots the annotations were transferred, select the category in the window displaying the “recipient” member and use any of the annotation display functions described in section 8.4

PDQUEST User Guide

10. Graphs and Reports

The graphs and reports described in this section are all available from the Reports menu.

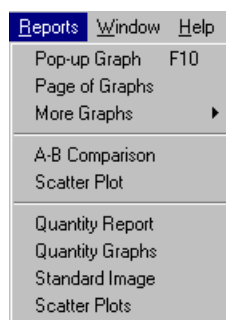


Fig.10-1. Reports menu.

10.1 Histogram Graphs

Histogram graphs can be used to quickly compare spot quantitation in each member of a matchset. A brief look at a histogram can give you a sense of the general trends in spot quantitation.

Histograms are also useful as tools for detecting problems such as spots that are not matched and spots whose quantitation is suspiciously off the average. Reviewing the histograms of spots in an Analysis Set is a quick and easy way to confirm that the members of the Analysis Set meet the criteria that was specified when the set was created.

10.1.a Data Displayed in the Histogram

Histograms are graphs consisting of bars whose lengths represent the data being plotted. Each bar represents the quantitation of the spot in each member of the matchset.

PDQUEST User Guide

The bars represent the gels from the left to the right in the order in which they were loaded. To check the order, use Assign (ALT+A). If you would like to change the order, use the function Reorder Members (Match > Create/Edit Matchset > Reorder Members).

As long as the spot is matched to the Standard, you can display a histogram for it.

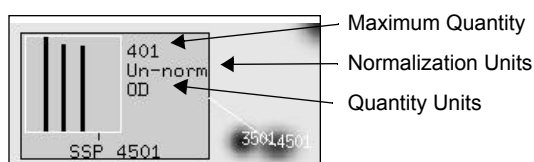


Fig.10-2. Sample histogram with individual quantitation.

The Standard Spot (SSP) number is displayed beneath the histogram.

The number at the upper right of the histogram is the quantitation of the maximum bar in the graph. The other bars are drawn proportional to the highest bar.

Beneath the maximum quantitation, the normalization units are displayed. Normalization is a process by which quantitative data from different gels are adjusted so that you can compare different samples to one another (for more information on normalization, see section 7.3).

If a spot is saturated, the bar representing it on the histogram will be shaded red instead of black.

If replicate groups are defined for the matchset (see section 7.2), you can look at the quantitation of individual spots or the average data of duplicate gels. You can specify which way you want to review spot quantitation by selecting Quantitation Mode from the Edit menu.

If you specify Replicate Group Quantitation, the bars on the histograms will display not only the maximum spot quantitation for each replicate group, but also the average spot quantitation and the standard deviation. Each bar represents one replicate group. The order of the bars corresponds to the order in which the replicate groups were made.

Reports

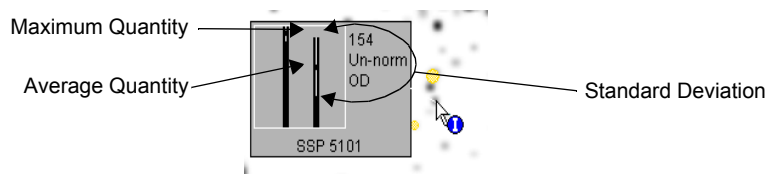


Fig.10-3. Sample histogram with Replicate Group Quantitation.

10.1.b Types of Graphs

Once an Analysis Set has been selected, use the Graph Set function to display the histograms of the spots in the set (Database > Analysis Sets > Graph Set).

The other tools for displaying graphs are located on the Reports menu, More Graphs submenu, and Graph Tools toolbar. The majority of these functions will only work on the Standard image.

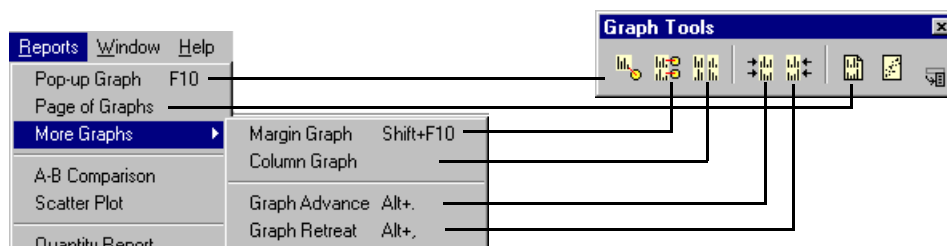


Fig.10-4. Graph commands.

Pop-up Graph

Pop-up Graph (F10) shows a histogram next to a selected spot. Select Pop-up Graph from the menu or toolbar, then click on a spot.

With Pop-up Graph, only one histogram can be displayed in a window at a time. If you select another spot, the first histogram will disappear and a new one will be displayed.

PDQUEST User Guide

Margin Graph

Margin Graph (SHIFT+F10) on the Reports > More Graphs submenu and Graph toolbar will display histograms for several spots along the left and the right margins of the screen. Select Margin Graph, then click on several spots. Their histograms will be displayed in the left and right margins.

Depending on the size of your window, you will be able to display up to 10 margin graphs at once. If you keep clicking on spots after you've displayed 10 graphs, you will begin replacing graphs with your more recent selections.

Column Graphs

Column Graph on the Reports > More Graphs submenu and Graph toolbar will fill the right and left margins with histograms of spots in the Standard. Select Column Graph; the spots whose graphs are displayed are all located in a narrow column of the image, centered around the location of the cursor.

To change the column of spots, use the Graph Advance/Graph Retreat tools (see below).

Graph Advance/Graph Retreat

If you are displaying Column Graphs or graphing spots using the Graph Partial Matches or Graph Analysis Set commands, you will need to shift the display to graph spots in different regions of the image. The commands for doing this are located on the Reports > More Graphs submenu or the Graph toolbar.

Graph Advance (ALT+.) will shift the column of spots being graphed to the right. Graph Retreat (ALT+,) will shift the column of spots being graphed to the left.

Page of Graphs

Page of Graphs displays the histograms of all the spots in your Standard as they would appear in a printed report. Select Page of Graphs from the Reports menu or Graph toolbar to open the window.

Reports

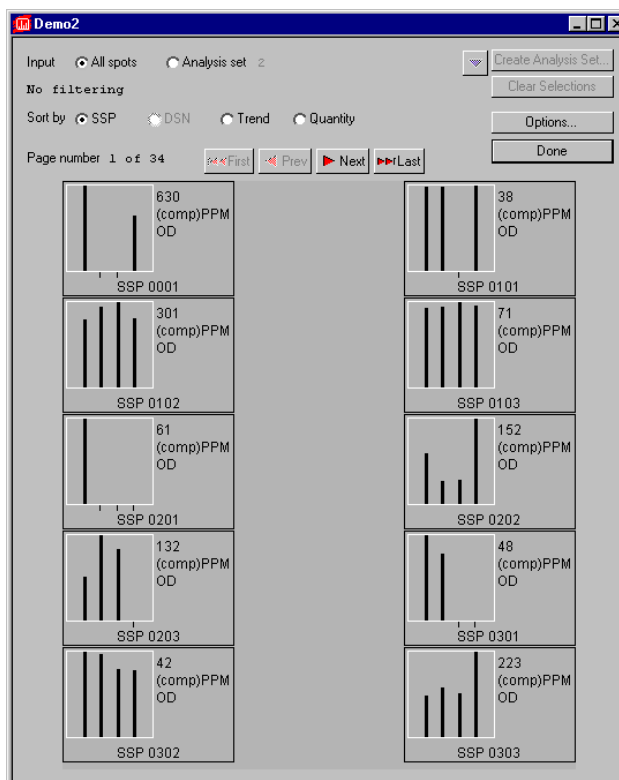


Fig.10-5. Page of Graphs window.

Note: Your displayed matchset will be reconfigured alongside the Page of Graphs report to allow you to compare the spots in the image with the spots in the report. When you close the Page of Graphs, your matchset will return to normal configuration. You can turn off this automatic reconfiguration under the report options (see below).

To peruse the Page of Graphs, click on the arrow buttons above the graphs.

You can display all spots in the matchset (All Spots) or specify an Analysis Set of spots to display (Analysis Set). Select the Analysis Set option button, then select the Analysis Set to display from the pop-up list.

PDQUEST User Guide

You can rearrange the graphs by clicking on the “Sort by” buttons. You can sort by SSP number or DSN number (if any are defined). If you select Trend, the graphs will be arranged based on a linear regression model of spot presence and quantity in each member. If you select Quantity, the graphs will be arranged starting with spots of greatest quantity to spots of least quantity.

You can create an Analysis Set containing some or all of the graphed spots. Click on the graph of each spot that you want to include in the Analysis Set. As you click, the border around each graph will become highlighted, indicating that it is selected.

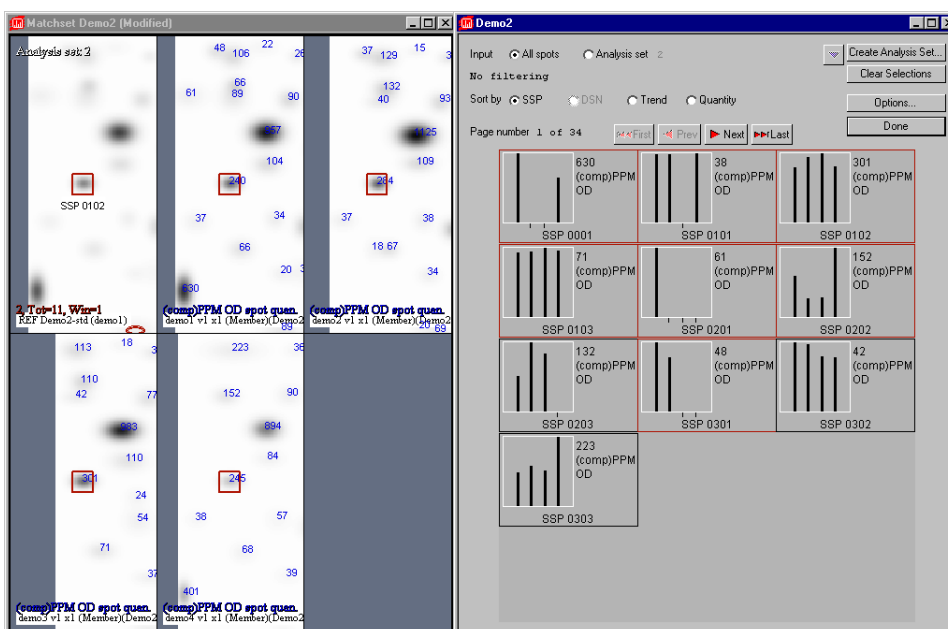


Fig.10-6. Creating an Analysis Set using Page of Graphs.

If you make a mistake, click on a selected graph to deselect it or click on the Clear Selections button to deselect all the graphs.

After selecting all the graphs of interest, click on the Create Analysis Set button and enter a name for the set and any descriptive information in the pop-up box. Click on Create to create the set.

Reports

See section 9.1 for more information on Analysis Sets.

Click on the Options button to change the display options of the report

To close the Page of Graphs window, click on Done.

10.2 A-B Comparison

The A-B Comparison function is used on matchset Standards to highlight and compare the spots that appear in two different matchset members.

Select A-B Comparison from the Reports menu and click on the Standard of the matchset in which you want to perform the comparison. A pop-box will prompt you to select the first member you want to compare (Member A). Click on the member name. A second pop-up box will prompt you to select the second member you want to compare (Member B).

When you select the second member, the box will disappear and the spots in the Standard will appear highlighted.

Spots that are marked in both members will be highlighted in green, spots that appear only in Member A will appear highlighted in red, and spots that appear only in Member B will appear highlighted in blue.

Spots that appear in neither member will not be highlighted.

PDQUEST User Guide

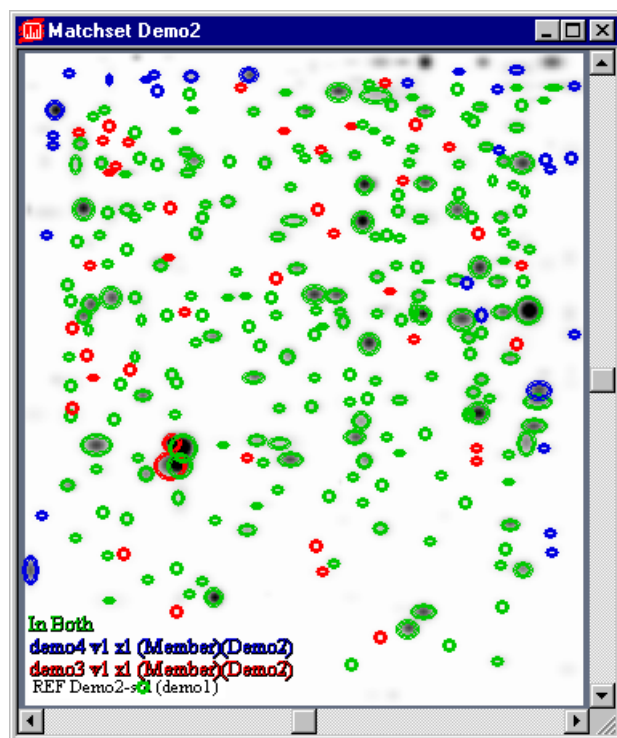


Fig.10-7. Example of an A-B Comparison display.

A color key will appear at the bottom of the Standard.

10.3 Scatter Plots

Scatter plots are graphs that show the relatedness of two gels. Each matched spot is plotted on a log scale as a point according to its quantitation in the first gel (X-axis) and its quantitation in the second gel (Y-axis). If the quantitation value of a spot is identical in both gels, the point representing that spot falls on the line with a slope of 1.00 (45° line). If the values are not equal, the point falls to the left or to the right of the line.

Reports

You can display a scatter plot for any two gels in a matchset, or print out the scatter plots for all the members of a matchset.

Displaying the Scatter Plot for Any Two Matchset Members

Select Scatter Plot from the Reports menu. The Scatter Plot window will open.

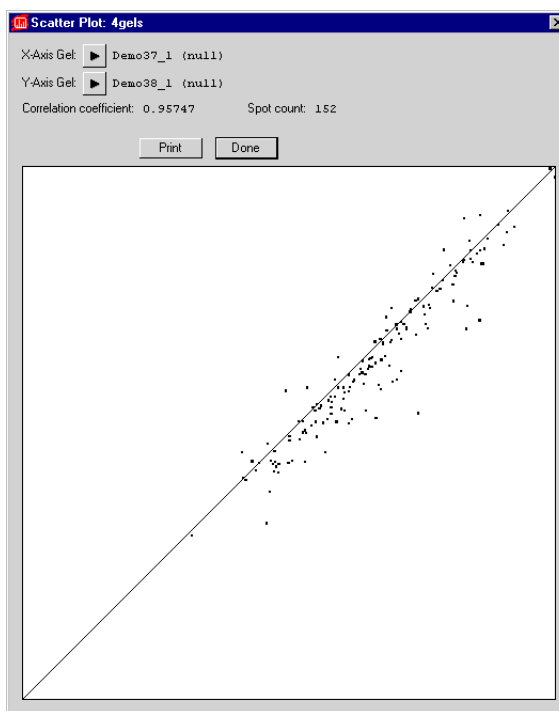


Fig.10-8. Scatter Plot window.

Select the members to compare by clicking on the arrows buttons next to X-Axis Gel and Y-Axis Gel.

The Spot Count is the number of matched spots in the two members, and the Correlation Coefficient is a measure of quantitative similarity between gels. A coefficient value of 1.00 indicates that the two gels are perfectly related (i.e.,

PDQUEST User Guide

equivalent), while a low coefficient value (e.g., 0.40) indicates that the two gels are not very similar.

To print this scatter plot, click on the Print button in the window. The standard print report dialog will open, in which you can select the settings for your print and proceed to print.

Printing Scatter Plots of All Members of a Matchset

To print the scatter plots for all members of a matchset, select Scatter Plots from the Reports menu.

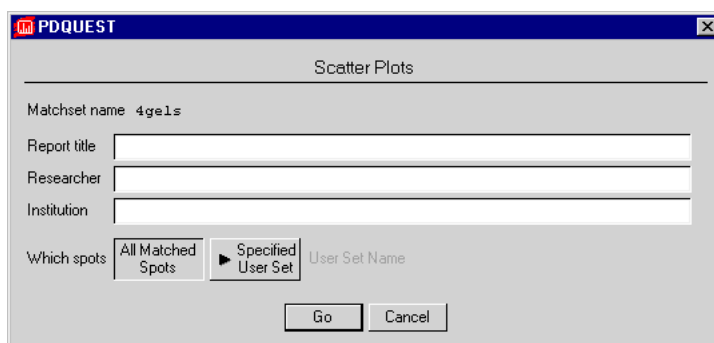
The image shows a screenshot of the 'Scatter Plots' dialog box within the PDQUEST software. The window has a title bar with the PDQUEST logo and the text 'PDQUEST'. The main title of the dialog is 'Scatter Plots'. Below the title, there is a label 'Matchset name' followed by the text '4gels'. There are three text input fields: 'Report title', 'Researcher', and 'Institution'. Below these fields, there is a section labeled 'Which spots' with two radio buttons: 'All Matched Spots' and 'Specified User Set'. The 'Specified User Set' radio button is selected. To the right of the 'Specified User Set' radio button is a text field labeled 'User Set Name'. At the bottom of the dialog are two buttons: 'Go' and 'Cancel'.

Fig.10-9. Scatter Plots dialog box.


You can type in a report title, the name of the researcher, and the name of the research institution in the appropriate fields.

You can specify the spots that will be included in your report: All Matched Spots or the spots in a Specified Analysis Set. If you click on the Specified Analysis Set button, a list of available Analysis Sets will be displayed. The name of the set you choose will appear next to the Analysis Set Name prompt.

Once you have selected your options, click on the Go button at the bottom of the form. This will open a smaller version of the standard print dialog, in which you can select your printer settings and proceed to print.

10.4 Quantity Reports

If you select Quantity Report from the Reports menu, the general report information dialog box will open.

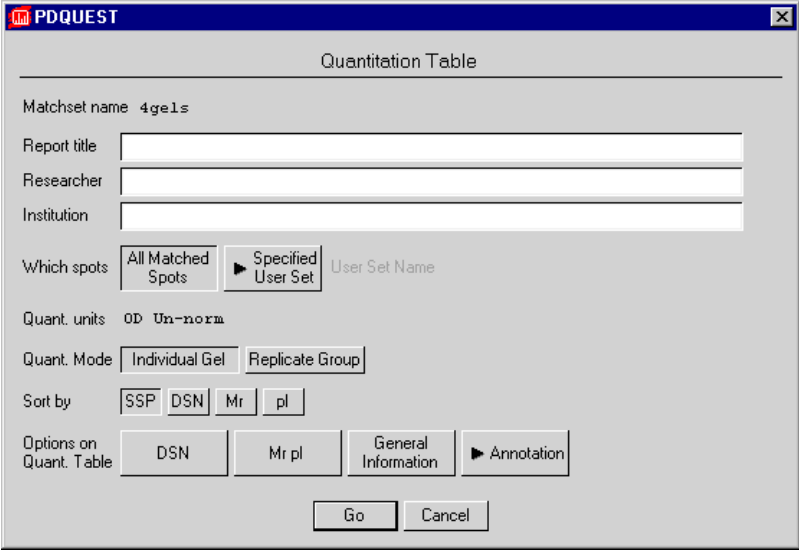

The image shows a software dialog box titled "PDQUEST" with a subtitle "Quantitation Table". It contains several input fields and buttons. The "Matchset name" field is pre-filled with "4gels". Below it are three empty text boxes for "Report title", "Researcher", and "Institution". The "Which spots" section has two buttons: "All Matched Spots" and "Specified User Set", with a "User Set Name" label to the right. The "Quant. units" field is set to "OD Un-norm". The "Quant. Mode" section has two buttons: "Individual Gel" and "Replicate Group". The "Sort by" section has four buttons: "SSP", "DSN", "Mr", and "pl". The "Options on Quant. Table" section has four buttons: "DSN", "Mr pl", "General Information", and "Annotation". At the bottom are "Go" and "Cancel" buttons.

Fig.10-10. Quantitation Table dialog box.

You can type in a report title, the name of the researcher, and the name of the research institution in the appropriate fields.

You can specify the spots that will be included in your report: All Matched Spots or the spots in a Specified Analysis Set. If you click on the Specified Analysis Set button, a list of available Analysis Sets will be displayed. The name of the set you choose will appear next to the Analysis Set Name prompt.

If you want to see spot quantitation for each individual member of the matchset, select Individual Gel next to the Quant. Mode prompt. If you want to see the average quantitation for a Replicate Group, click on the Replicate Group button.

PDQUEST User Guide

You can sort your spots by SSP number, DSN, molecular weight, or isoelectric point by clicking on the appropriate “Sort by” button.

Other Options

If DSN numbers are assigned to the protein spots and you would like to see them on the report, click on the DSN button.

Mr and pI data will be reported if you click on the Mr pI button.

Clicking on the General Information button will include a title page, a list of symbols and abbreviations, and information about the matchset, member gels, and their calibration data in the report.

You can print spot annotations by clicking on the Annotation button and selecting from the list of annotation categories. The report will include annotation information for the spots that have annotation data in the selected category.

Once you have selected your options, click on the Go button at the bottom of the form. This will open a smaller version of the standard print dialog, in which you can select your printer settings and proceed to print.

10.5 Quantity Graphs

If you select Quantity Graphs from the Reports menu, a dialog box will open in which you can specify the format of your quantitation histograms.

Reports

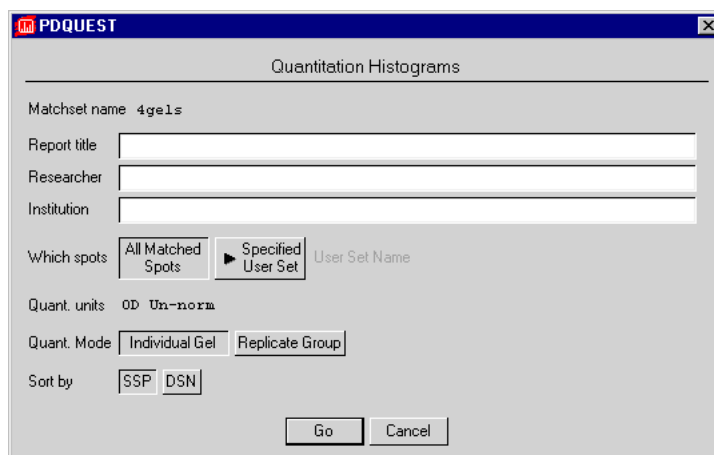
The image shows a software dialog box titled "PDQUEST" with a subtitle "Quantitation Histograms". It contains several input fields and buttons. The "Matchset name" field is pre-filled with "4gels". Below it are three empty text boxes for "Report title", "Researcher", and "Institution". The "Which spots" section has two buttons: "All Matched Spots" and "Specified User Set", with the latter being selected. To the right of these buttons is a text field labeled "User Set Name". The "Quant. units" field is set to "OD Un-norm". The "Quant. Mode" section has two buttons: "Individual Gel" and "Replicate Group", with the latter being selected. The "Sort by" section has two buttons: "SSP" and "DSN", with "SSP" being selected. At the bottom right are "Go" and "Cancel" buttons.

Fig.10-11. Quantitation Histograms dialog box.

You can type in a report title, the name of the researcher, and the name of the research institution in the appropriate fields.

You can specify the spots that will be included in your report: All Matched Spots or the spots in a Specified Analysis Set. If you click on the Specified Analysis Set button, a list of available Analysis Sets will be displayed. The name of the set you choose will appear next to the Analysis Set Name prompt.

If you want to see spot quantitation for each individual member of the matchset, select Individual Gel next to the Quant. Mode prompt. If you want to see the average quantitation for a Replicate Group, click on the Replicate Group button.

You can sort your spots by SSP number or DSN by clicking on the appropriate "Sort by" button.

Once you have selected your options, click on the Go button at the bottom of the form. This will open a smaller version of the standard print dialog, in which you can select your printer settings and proceed to print.

10.6 Standard Image Reports

Selecting Standard Image from the Reports menu allows you to print a map of up to 300 spots from the matchset Standard.

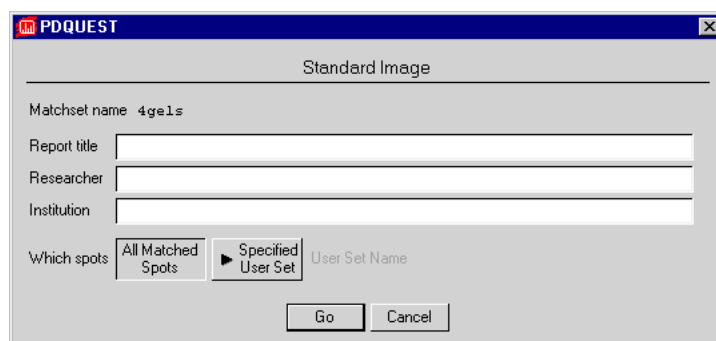
The image shows a software dialog box titled "PDQUEST" with a sub-header "Standard Image". Inside the dialog, there is a label "Matchset name" followed by the text "4ge1s". Below this are three text input fields labeled "Report title", "Researcher", and "Institution". At the bottom left, there is a label "Which spots" followed by two buttons: "All Matched Spots" and "Specified User Set". The "Specified User Set" button is currently selected, indicated by a small black triangle. To the right of these buttons is a text field labeled "User Set Name". At the bottom center of the dialog are two buttons: "Go" and "Cancel".

Fig.10-12. Standard Image dialog box.

You can type in a report title, the name of the researcher, and the name of the research institution in the appropriate fields.

You can specify the spots that will be included in your report: All Matched Spots or the spots in a Specified Analysis Set. If you click on the Specified Analysis Set button, a list of available Analysis Sets will be displayed. The name of the set you choose will appear next to the Analysis Set Name prompt.

Once you have selected your options, click on the Go button at the bottom of the form. This will open a smaller version of the standard print dialog, in which you can select your printer settings and proceed to print.

11. Printing and Exporting

This chapter describes the commands for printing and exporting images and data. These commands are located on the File menu.

11.1 Printing

The commands for printing images, actual-size images, image reports, and video printouts are located on the File > Print submenu.

1. Print Image prints a copy of the image and any overlays that are currently displayed.
2. Print Actual Size prints an actual-size copy of the image.
3. Image Report prints the image and information about its scan history, number of pixels, data range, etc.
4. Video Print prints images and reports to a video printer.*

*Video printing requires installation of the video board and cable that come with the Gel Doc 1000/2000 or Chemi Doc gel documentation system. The video board and cable can also be ordered separately.

11.1.a Printing Images

Print Image prints a copy of the active image window and any image overlays that are displayed.

File > Print > Print Image opens the Print Image dialog box.

PDQUEST User Guide

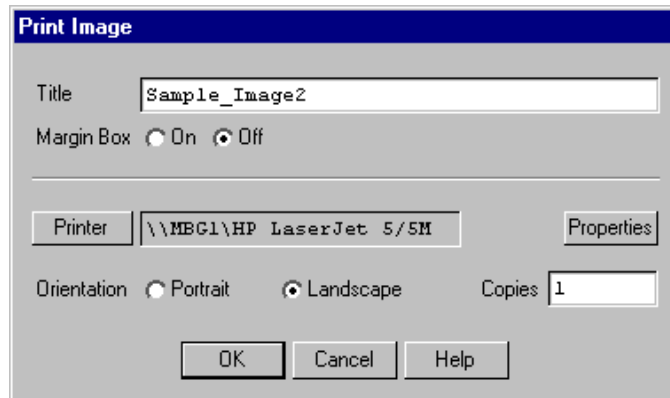


Fig.11-1. Print Image dialog box.

At the top of the dialog box, enter the title that you want to appear above the image when you print.

Choose whether or not a margin border should appear around the edge of the image by clicking either the On or Off button next to the Margin Box prompt.

Click on Printer to display a list of your available printers. After you have selected a printer, you can click on Properties to select the standard Windows print settings for that printer.

Note: If you are using Windows 95 with certain printer drivers, you may experience difficulty printing some overlays and overlay symbols on images. While there is no single solution to this problem, you should try adjusting some settings in the Printer Properties dialog box, under the Graphics and Advanced tabs—specifically, the Dithering, DPI, Rasterizing, and RET settings.

Select the paper orientation by clicking on Portrait or Landscape next to the Layout prompt. (This can also be set in the Properties dialog box.)

Indicate how many copies should be printed in the text box next to the Copies prompt.

When you are satisfied with all the print parameters, click on the OK button to send the image to the printer.

Printing and Exporting

11.1.b Print Actual Size

You can print your images at their actual size using the Print Actual Size command.

Note: If you are using the Gel Doc, Chemi Doc, Fluor-S, or Fluor-S MAX, you must specify the correct image area size when capturing your images to ensure accurate 1:1 printing. You can specify the image area size in the acquisition window for the instrument. See the chapter on each imaging device for more information.

Select Print Actual Size from the File > Print menu. The standard Print Image dialog box will open (see previous section).

11.1.c Image Report

Image Report allows you to print out a single-page report of an image and its associated information. The format of the report is designed to provide a concise yet thorough summary of the most relevant features of an image for documentation purposes.

The report includes the following information:

- The image.
- The report date.
- The image title and a brief description.
- The directory location.
- The date the image was scanned.
- The type of imaging device.
- The imaging area and number of pixels.
- Pixel size.
- The intensity range, image color, and memory size.
- Image background information.
- Relevant lane and band information.

PDQUEST User Guide

To print a image report of a particular gel, select Image Report from the File > Print submenu. This will open a smaller version of the standard print dialog box (see Print Image, above).

Note: TIFF images do not contain all the tagged information that would normally be included in an image file (for example, imaging device, scan date, image color, etc.). For this reason, the Image Report may list this information as “Unknown” for imported TIFF files.

11.1.d Video Printing

The Video Print command allows you to print images and reports on a video printer. To create a video printout of the active window, select Video Print from the File > Print submenu.

Note: Video printing requires installation of the video board and cable that come with the Gel Doc 1000/2000 and Chemi Doc gel documentation systems. The video board and cable can also be ordered separately.

Settings for the Mitsubishi P90W/P91W Video Printer

There are three settings for the Mitsubishi P90W / P91W video printer. Set Contrast to 0, Brightness to 0, and Gamma to 5.

The dip switches should stay in the orientation in which they are shipped: Pin 1 is up (on), and Pins 2–10 are down (off).

11.2 Exporting

Matchset data can be exported for analysis by spreadsheet programs, scans can be exported as TIFF images for analysis by other image analysis software, and annotations can be exported as text files for inclusion in documents and spreadsheets.

11.2.a Exporting Matchset Data

Data from a matchset can be exported in ASCII format to your hard disk, file server, or other storage medium for analysis by spreadsheet programs.

Printing and Exporting

With the matchset loaded, select Export > Matchset Data from the File menu. This will open the Export Matchset Data dialog box.

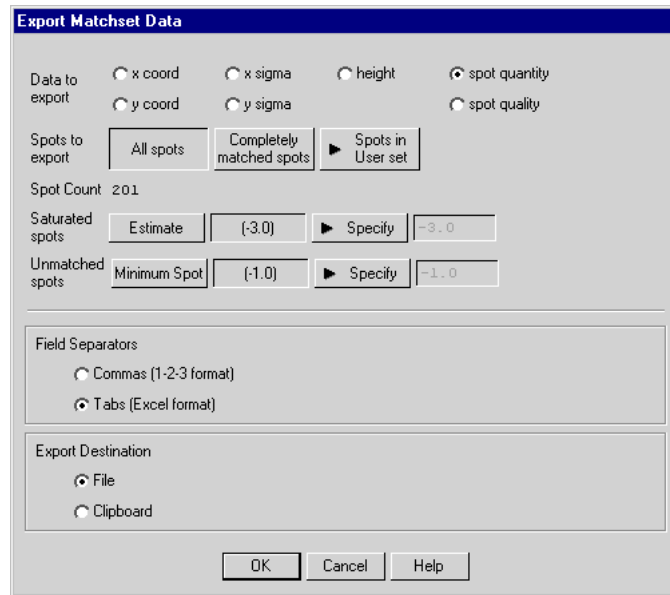
The image shows a software dialog box titled "Export Matchset Data". It contains several sections for configuring export options. The "Data to export" section has radio buttons for "x coord", "x sigma", "height", "spot quantity" (selected), "y coord", "y sigma", and "spot quality". The "Spots to export" section has buttons for "All spots", "Completely matched spots", and "Spots in User set". The "Spot Count" is displayed as "201". The "Saturated spots" section has an "Estimate" button, a text field with "(-3.0)", a "Specify" button, and another text field with "-3.0". The "Unmatched spots" section has a "Minimum Spot" button, a text field with "(-1.0)", a "Specify" button, and another text field with "-1.0". The "Field Separators" section has radio buttons for "Commas (1-2-3 format)" and "Tabs (Excel format)" (selected). The "Export Destination" section has radio buttons for "File" (selected) and "Clipboard". At the bottom are "OK", "Cancel", and "Help" buttons.

Fig.11-2. Export Matchset Data dialog box.

Next to the “Data to export” prompt, select the characteristic of the spot that you want to export.

Specify the spots whose data should be exported. You can choose either all the spots in the matchset, those spots that are matched in every member of the matchset or the spots in a particular Analysis Set.

The Spot Count will automatically adjust to reflect the total number of spots selected.

If you are exporting spot quantity data, select the way in which saturated spots and unmatched spots are reported.

PDQUEST User Guide

The quantitation of saturated spots can be estimated, reported as (-3.0), or reported in some other way that you specify. If you opt for the latter, click on the Specify button and enter the value that should be reported.

If a spot is unmatched in one or more members of a matchset, it can be reported as the minimum spot quantitation value (the minimum quantity needed for a spot to be detected for the gel). Alternatively, it can be reported as (-1.0) or as a value that you enter in the text box next to the Specify button.

Next, select the separators you want to use to separate columns of data in your spreadsheets: commas or tabs.

Your data can be exported as a file or to the clipboard.

If you choose the File option button, clicking on OK will open a version of the Save As dialog box, in which you can enter a name and export location for the file.

If you choose the Clipboard option button, clicking on OK will save your data to the clipboard.

11.2.b Exporting an Image

To export a 2-D scan or calstrip scan as a TIFF image, select Export > Export to TIFF Image from the File menu. The Export Scan dialog will open.

Printing and Exporting

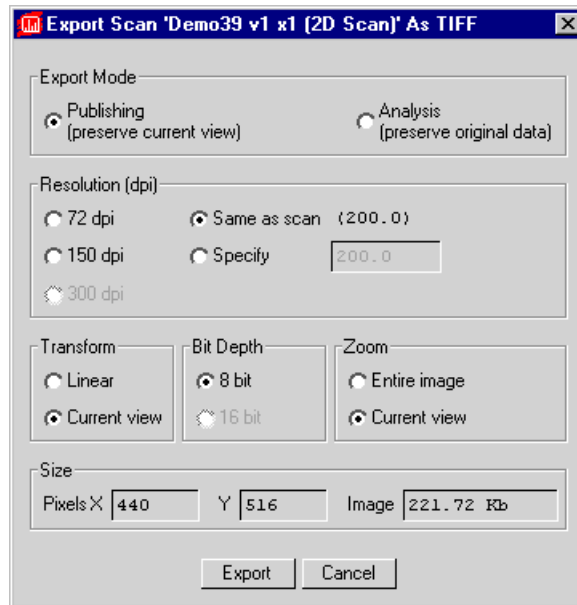


Fig.11-3. Exporting TIFF Images.

Note: Only 2-D scans and calstrip scans can be exported as TIFF images.

Analysis Export Mode

Analysis mode exports the scan data, unmodified by PDQUEST's viewing controls.

If you select Analysis mode, the other controls in this dialog box will become inactive.

Publishing Export Mode

Publishing mode exports a TIFF image that looks like the image as it is currently displayed in PDQUEST (with whatever Transform adjustments you may have made).

PDQUEST User Guide

In Publishing mode, you can specify a resolution for the TIFF image by selecting 72 dpi (typical computer screen resolution), 150 or 300 dpi (standard printing resolutions), the same resolution as the scan file, or any resolution you want to specify (up to the resolution of the scan).

If you have log transformed the image, you can specify a linear transform, or preserve the current view.

If your image is 16 bits, you can compress it to 8 bits for export to TIFF.

Finally, if you are only displaying part of the image due to magnification or repositioning, you can preserve the current view or export the entire image.

Exporting the Image

The size of the pixels in the image and the file size of the image are listed at the bottom of the dialog. When you are ready to export, click on the Export button.

A Save As dialog box will open. The default file name will have a .tif extension, and the file type will indicate that this is a TIFF image. You can change the file name or select a different directory to save in.

PDQUEST will generate a grayscale TIFF image compatible with Melanie software.

11.2.c Exporting Annotations

You can export any matchset annotations to a text file for inclusion in a document or spreadsheet.

Select Annotations from the File > Export submenu. A box will open prompting you to select the annotation category or categories to export.

Printing and Exporting

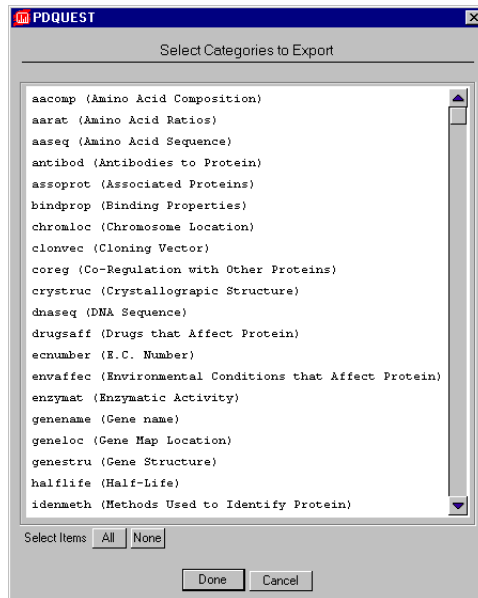


Fig.11-4. Export Annotations pop-up box.

Select from the list of available categories or select all available categories by clicking on the All button. Deselect any categories you select by clicking on the None button.

When you have made your selection, click on Done. A Save dialog box will open allowing you to name your text file and select the directory to save to.

Text files can be imported into standard Windows and Macintosh documents and spreadsheets.

PDQUEST User Guide

Appendix A

Gel Doc 1000/2000



Fig.A-1. Gel Doc 1000/2000.

Before you can begin acquiring images using the Gel Doc[®] 1000/2000 gel documentation system, the system must be properly installed and connected with the host computer. See the Gel Doc hardware manual for installation, startup, and operating instructions.

PDQUEST User Guide

To use the Gel Doc, you will need to have the Bio-Rad-supplied acquisition board installed in your PC or Macintosh. The drivers for this board will be installed when you install PDQUEST.

Make sure that your Gel Doc camera is turned on. If the camera is not turned on, the Gel Doc acquisition window will open but the screen will be black, and you will be unable to Video Print.

Simulation Mode

Any of the imaging device acquisition windows can be opened in “simulation mode.” In this mode, an acquisition window will open and the controls will appear active, but instead of capturing real images, the window will create “dummy” images of manufactured data.

You do not need to be connected to an imaging device to open a simulated acquisition window. This is useful for demonstration purposes or practice scans.

To enter simulation mode, hold down the CTRL key and select the name of the device from the File menu. The title of the acquisition window will indicate that it is simulated.

A.1 Gel Doc Acquisition Window

To acquire images using the Gel Doc, go to the File menu and select Gel Doc.... The acquisition window for the instrument will open, displaying a control panel and a video display window.

Appendix A. Gel Doc

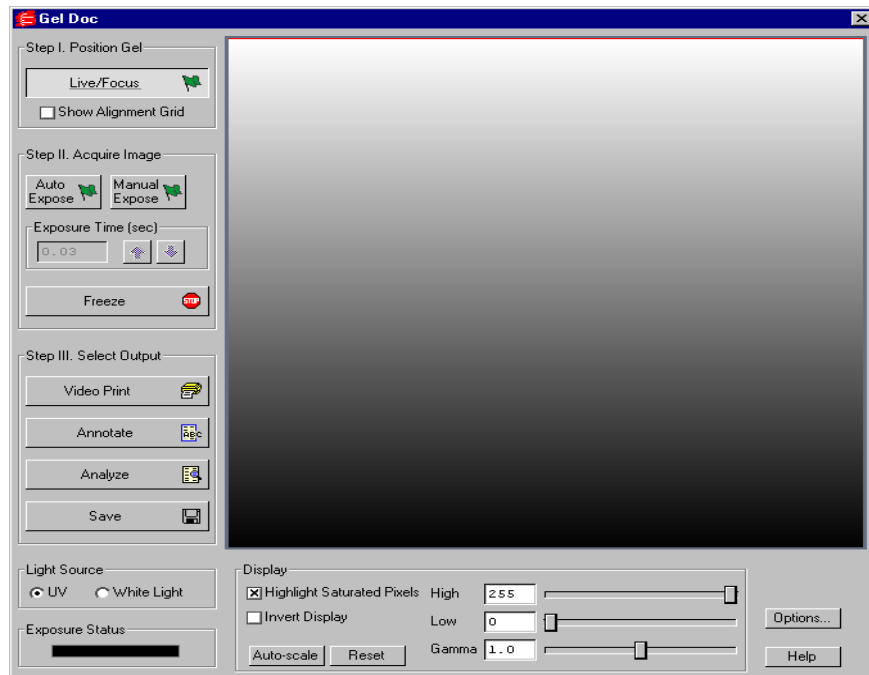


Fig.A-2. Gel Doc acquisition window

The Gel Doc video display window will open in “live” mode, giving you a live video display of your sample. If no image is visible, make sure the camera is on, check the cable connections, make sure the f-stop on the camera is not closed, and make sure that the protective cap is off the camera lens. Also check to see that the transilluminator is on and working.

The control panel has been arranged from top to bottom to guide you through the acquisition procedure. There are three basic steps to scanning an image using the Gel Doc:

1. Position and focus the gel or other object to be imaged.
2. Acquire the image.
3. Select the output.

PDQUEST User Guide

A.2 Step I. Position Gel

The Gel Doc window will open in “live” mode, giving you a live video display of your sample. (Frames will be captured and displayed at about 10 frames per second, depending on the speed of your computer.) In this mode, the Live/Focus button will appear selected.

You can use live mode to zoom and focus the camera on your sample area and position your sample within the area. You can also select the Show Alignment Grid checkbox to facilitate positioning.

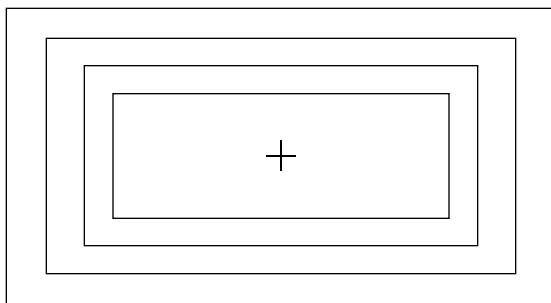


Fig.A-3. Gel Doc alignment grid.

The image should stay in focus while zooming. If focus is not maintained, refer to your Gel Doc hardware manual.

Note: After positioning your sample, you should check the Imaging Area dimensions under Options (see below) to make sure that they conform to the size of the area you are focusing on. To determine the size of the area you are focusing on, you can place a ruler in the Gel Doc box so that it is visible by the camera.

A.3 Step II. Acquire Image

For many white light applications, you can skip this step and save and print images directly from Live/Focus mode.

Appendix A. Gel Doc

For UV light or chemiluminescent applications, or for faint samples, you will need to select an exposure time for the image.

“Exposure” refers to the integration of the image on the camera CCD over a set period of time. The effect is analogous to exposing photographic film to light over a period of time. The Gel Doc has two exposure modes: Auto and Manual.

Auto Expose

Clicking on the Auto Expose button will cancel Live/Focus mode and begin an automatic exposure. (The Auto Expose button will appear selected throughout the exposure.) If you know the approximate exposure time you want (± 3 seconds), you can skip this step and go directly to Manual Expose.

In Auto Expose mode, the image is continuously integrated on the camera CCD until it has reached a certain percentage of saturated pixels. This percentage is set in the Options dialog box. (Default = 0.75 percent. See Options below.)

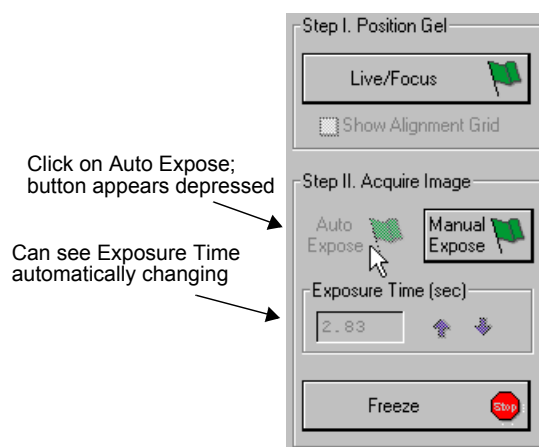


Fig.A-4. Acquiring an image, step 1: Select Auto Expose.

Once an image has reached the specified percent of saturated pixels, it is captured and displayed in the video display window, Auto Expose is

PDQUEST User Guide

automatically deactivated, the exposure time appears active in the Exposure Time field, and Manual Expose is activated.

Note: If you are having difficulty auto-exposing your sample, you can use Manual Expose to adjust your exposure time directly. Most applications only require an exposure time of a few seconds, which can be quickly adjusted using Manual Expose.

Manual Expose

Manual Expose is automatically activated after Auto Expose has deactivated. Or, if you know the approximate exposure time you want, you can simply click on the Manual Expose button.

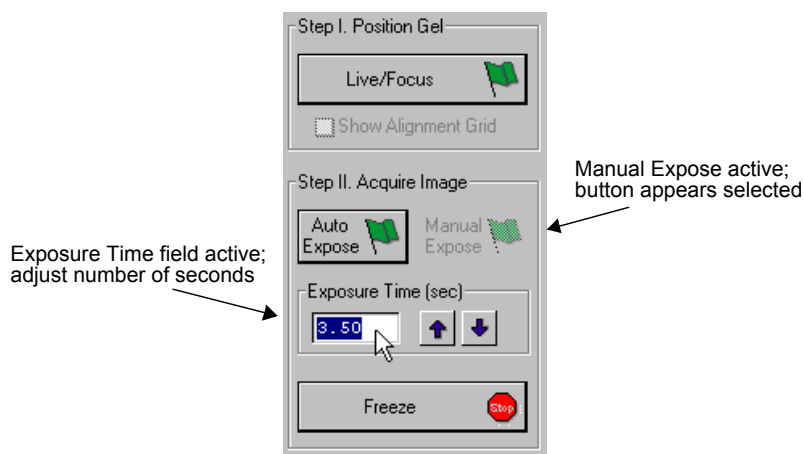


Fig.A-5. Acquiring an image, step 2: Adjust exposure time.

When Manual Expose is active (the button will appear selected), the last captured image will be displayed in the Gel Doc image window. The camera continues to integrate the image on the CCD, updating the display whenever the specified number of seconds is reached.

You can adjust the exposure time by changing the number of seconds in the Exposure Time field. To do so, enter a number directly or use the arrow buttons next to the field.

Appendix A. Gel Doc

Freeze

Once you are satisfied with the quality of the displayed image, click on the Freeze button to stop the exposure process. The last full exposure will be displayed in the image window.

Note: Freeze is automatically activated if you adjust any of the subsequent controls (e.g., Video Print, Image Mode, Display controls, etc.).

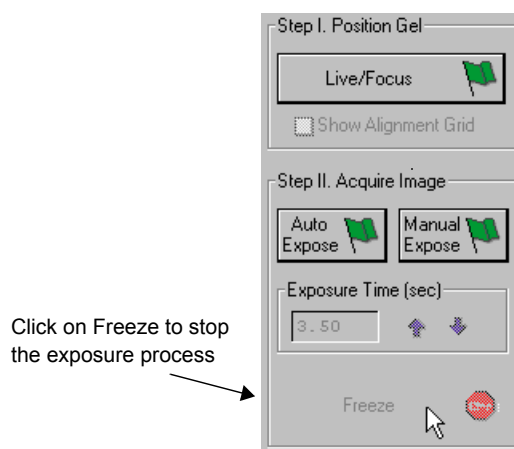


Fig.A-6. Acquiring an image, Step 3: Freeze the exposure.

A.4 Step III. Select Output

The Gel Doc window has several output options.

Video Print

Clicking on Video Print will automatically send the currently displayed frame (either live or integrated) to a video printer. You can add information about your image to the bottom of the printout by selecting the appropriate checkboxes in the Options dialog box. (See Options, below.)

PDQUEST User Guide

Annotate

Clicking on Annotate will open a separate image window displaying the captured image. The default name for the image will include the date, time, and user (if known).

The Annotation Tools toolbar will also pop up to allow you to annotate your image.

The image will not be saved until you select Save or Save As from the File menu.

Analyze

Clicking on Analyze will open a separate image window displaying the captured image. The default name for the image will include the date, time, and user (if known).

You can then analyze the image using the menu and toolbar functions.

The image will not be saved until you select Save or Save As from the File menu.

Save

Clicking on Save will open a separate image window displaying the captured image. A Save As dialog box will automatically open displaying the default file name for the image, which will include the date, time, and user (if known). You can then change the file name and storage directory.

You can also save your image as a TIFF image for export to other applications.

A.5 Image Mode

The Image Mode option buttons allow you to set the type and scale of your data.

Appendix A. Gel Doc

UV

Select this mode for fluorescent and chemiluminescent samples. With this mode selected, the data will be measured in linear intensity units.

White Light

Select this mode for reflective and transmissive samples. With this mode selected, the data will be measured in uncalibrated optical density (uOD) units.

A.6 Exposure Status

The Exposure Status bar shows the progress of your exposure. If your exposure time is greater than 1 second, the status bar display will give you a graphical representation of the remaining time before exposure is complete.

If the exposure time is less than 1 second, the status bar will not refresh itself for each exposure; it will remain at 100 percent.

A.7 Display

The Display controls are useful for quickly adjusting the appearance of your image for output to a video printer. Adjusting these controls will automatically freeze the video display and allow you to alter the image within the Gel Doc window.

These controls are similar to those in the Transform dialog box.

Note: The Display controls will only change the appearance of the image. They will not change the underlying data.

Highlight Saturated Pixels

When this box is checked, any saturated pixels in the image will appear highlighted in red in the scan window and in the pop-up image window. To

PDQUEST User Guide

view / hide saturated pixels in the pop-up image window, use the Image > Transform command.

Invert Display

This checkbox will switch light spots on a dark background to dark spots on a light background, and visa versa. This will only affect how the image is displayed on the screen, not the actual image data.

Auto-scale

Clicking on Auto-scale will adjust your displayed image automatically. The lightest part of the image will be set to the minimum intensity (e.g., white), and the darkest will be set to the maximum intensity (e.g., black). You can then “fine tune” the display using the High, Low, and Gamma sliders described below.

High/Low Sliders

If Auto-scale doesn't give you the appearance you want, you can use the High and Low sliders to redraw the image yourself. Dragging the High slider handle to the left will make weak signals appear darker. Dragging the Low slider handle to the right will reduce background noise.

You can also type specific High and Low values in the text boxes next to the sliders. Clicking anywhere on the slider bars will move the sliders incrementally.

Gamma Slider

Some images may be more effectively visualized if their data are mapped to the computer screen in a nonlinear fashion. Adjusting the Gamma slider handle changes the light and dark contrast nonlinearly.

Reset

Reset will return the image to its original, unmodified appearance.

A.8 Options

Clicking on the Options button opens the Options dialog box. Here you can specify certain settings for your Gel Doc system.

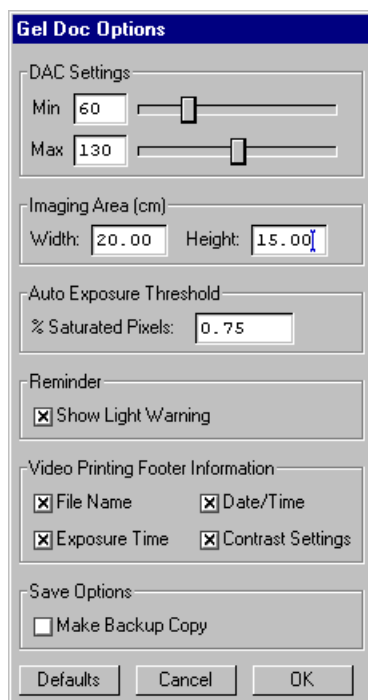


Fig.A-7. Available options in the Gel Doc acquisition window.

Clicking on OK implements any changes you make in this box. Clicking on Defaults restores the settings to the factory defaults.

DAC Settings

Note: The default DAC settings are highly recommended and should be changed with caution.

PDQUEST User Guide

These sliders may be used to adjust the minimum and maximum voltage settings of your video capture board. The minimum slider defines the pixel value that will appear as white in the image, while the maximum slider defines the pixel value that will appear as black. The slider scale is 0–255, with the defaults set to 60 minimum and 130 maximum.

Imaging Area

These fields are used to specify the size of your imaging area in centimeters, which in turns determines the size of the pixels in your image (i.e., resolution). When you adjust one imaging area dimension, the other dimension will change to maintain the aspect ratio of the camera lens.

Note: Your imaging area settings must be correct if you want to do 1:1 printing. These are also important if you are comparing the size of objects (e.g., using the Volume Tools) between images.

Auto Exposure Threshold

When you click on Auto Expose, the exposure time is determined by the percentage of saturated pixels you want in your image. This field allows you to specify that percentage.

Typically, you will want less than 1 percent of the pixels in your image saturated. Consequently, the default value for this field is 0.75 percent.

Reminder

When this checkbox is selected, the software will warn you to turn off your transilluminator light when you exit the Gel Doc acquisition window or when your system is “idle” for more than 5 minutes.

Note: If you are performing experiments that are longer than 5 minutes (e.g., chemiluminescence), this should be deselected.

Video Printing Footer Information

The checkboxes in this group allow you to specify the information that will appear at the bottom of your video printer printouts.

Appendix A. Gel Doc

Save Options

To automatically create a backup copy of any scan you create, select the Make Backup Copy checkbox.

With this checkbox selected, when you save a scan, a backup copy will be placed in the same directory as the scanned image. Windows backup files will have an “.sbk” extension. Macintosh backup files will have the word “backup” after the file name.

This backup copy will be read-only, which means that you cannot make changes to it. You can open it like a normal file, but you must save it under a different file name before editing the image or performing analysis.

PDQUEST User Guide

Appendix B

Chemi Doc



Fig.B-1. Chemi Doc.

Before you can begin acquiring images using the Chemi Doc gel documentation system, the system must be properly installed and connected with the host computer. See the Chemi Doc hardware manual for installation, startup, and operating instructions.

PDQUEST User Guide

To use the Chemi Doc, you will need to have the Bio-Rad-supplied acquisition board installed in your PC or Macintosh. The drivers for this board will be installed when you install PDQUEST.

Make sure that your Chemi Doc camera is turned on. If the camera is not turned on, the Chemi Doc acquisition window will open but the screen will be black, and you will be unable to Video Print.

Simulation Mode

Any of the imaging device acquisition windows can be opened in “simulation mode.” In this mode, an acquisition window will open and the controls will appear active, but instead of capturing real images, the window will create “dummy” images of manufactured data.

You do not need to be connected to an imaging device to open a simulated acquisition window. This is useful for demonstration purposes or practice scans.

To enter simulation mode, hold down the CTRL key and select the name of the device from the File menu. The title of the acquisition window will indicate that it is simulated.

B.1 Chemi Doc Acquisition Window

To acquire images using the Chemi Doc, go to the File menu and select Chemi Doc.... The acquisition window for the instrument will open, displaying a control panel and a video display window.

Appendix B. Chemi Doc

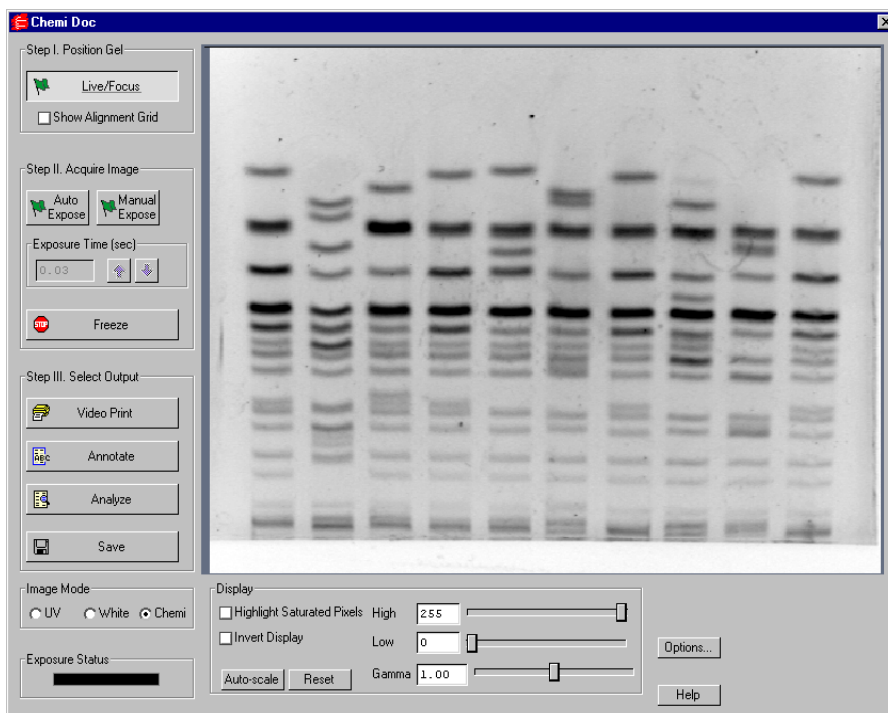


Fig.B-2. Chemi Doc acquisition window

The Chemi Doc video display window will open in “live” mode, giving you a live video display of your sample. If no image is visible, make sure the camera is on, check the cable connections, make sure the f-stop on the camera is not closed, and make sure that the protective cap is off the camera lens. Also check to see that the transilluminator is on and working.

The control panel has been arranged from top to bottom to guide you through the acquisition procedure. There are three basic steps to scanning an image using the Chemi Doc:

1. Position and focus the gel or other object to be imaged.
2. Acquire the image.

PDQUEST User Guide

3. Select the output.

B.2 Step I. Position Gel

The Chemi Doc window will open in “live” mode, giving you a live video display of your sample. (Frames will be captured and displayed at about 10 frames per second, depending on the speed of your computer.) In this mode, the Live/Focus button will appear selected.

You can use live mode to zoom and focus the camera on your sample area and position your sample within the area. You can also select the Show Alignment Grid checkbox to facilitate positioning.

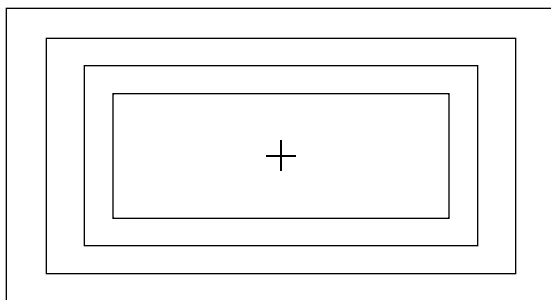


Fig.B-3. Chemi Doc alignment grid.

The image should stay in focus while zooming. If focus is not maintained, refer to your Chemi Doc hardware manual.

Note: After positioning your sample, you should check the Imaging Area dimensions under Options (see below) to make sure that they conform to the size of the area you are focusing on. To determine the size of the area you are focusing on, you can place a ruler in the Chemi Doc box so that it is visible by the camera.

B.3 Step II. Acquire Image

For many white light applications, you can skip this step and save and print images directly from Live/Focus mode.

For UV light, chemiluminescent applications, or faint samples, you will need to select an exposure time for the image.

“Exposure” refers to the integration of the image on the camera CCD over a set period of time. The effect is analogous to exposing photographic film to light over a period of time. The Chemi Doc has two exposure modes: Auto and Manual.

Auto Expose

Clicking on the Auto Expose button will cancel Live/Focus mode and begin an automatic exposure. (The Auto Expose button will appear selected throughout the exposure.) If you know the approximate exposure time you want (± 3 seconds), you can skip this step and go directly to Manual Expose.

In Auto Expose mode, the image is continuously integrated on the camera CCD until it has reached a certain percentage of saturated pixels. This percentage is set in the Options dialog box. (Default = 0.75 percent. See Options below.)

PDQUEST User Guide

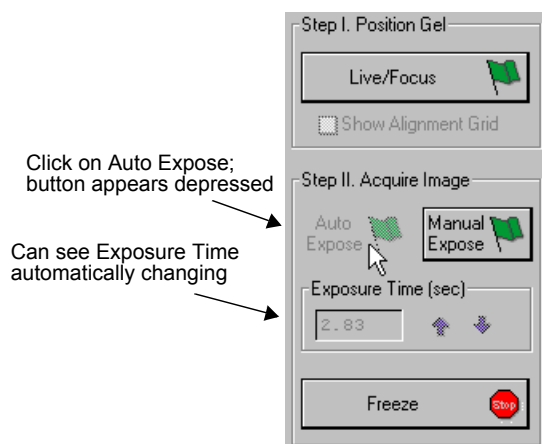


Fig.B-4. Acquiring an image, step 1: Select Auto Expose.

Once an image has reached the specified percent of saturated pixels, it is captured and displayed in the video display window, Auto Expose is automatically deactivated, and the exposure time appears active in the Exposure Time field. At this point, if you are in UV or White image mode, Manual Expose will be automatically activated. If you are in Chemi mode, the image will be frozen after Auto Expose is complete.

Note: If you are having difficulty auto-exposing your sample, you can use Manual Expose to adjust your exposure time directly. Most non-chemi applications only require an exposure time of a few seconds, which can be quickly adjusted using Manual Expose.

Manual Expose

If you know the approximate exposure time you want, you can click on the Manual Expose button. In UV or White image mode, Manual Expose is automatically activated after Auto Expose is complete.

Appendix B. Chemi Doc

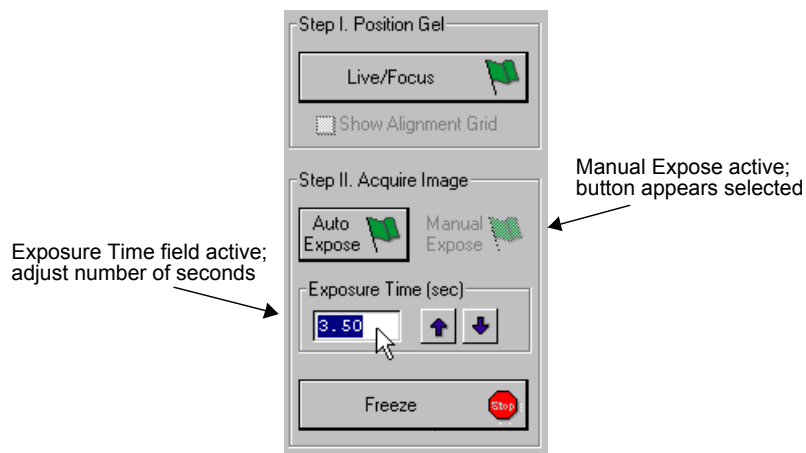


Fig.B-5. Acquiring an image, step 2: Adjust exposure time.

With Manual Expose activated, you can adjust the exposure time directly by changing the number of seconds in the Exposure Time field. Type in a number or use the arrow buttons next to the field.

In UV or White image mode, when Manual Expose is active, the last captured image will be displayed in the Chemi Doc image window. The camera will continue to integrate the image on the CCD, updating the display whenever the specified number of seconds is reached.

In Chemi mode, Manual Expose will expose an image over the specified exposure time and then stop.

Freeze

In UV and White image mode, once you are satisfied with the quality of the displayed image, click on the Freeze button to stop the exposure process. The last full exposure will be displayed in the image window.

In Chemi image mode, Freeze will be automatically activated after Auto or Manual Expose is complete.

PDQUEST User Guide

Note: Freeze is automatically activated if you adjust any of the subsequent controls (e.g., Video Print, Image Mode, Display controls, etc.).

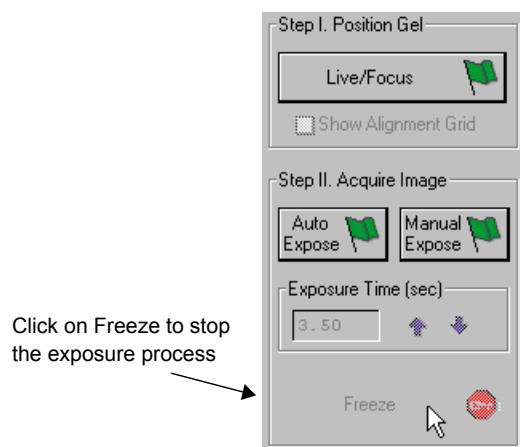


Fig.B-6. Acquiring an image, Step 3: Freeze the exposure.

B.4 Step III. Select Output

The Chemi Doc window has several output options.

Video Print

Clicking on Video Print will automatically send the currently displayed frame (either live or integrated) to a video printer. You can add information about your image to the bottom of the printout by selecting the appropriate checkboxes in the Options dialog box. (See Options, below.)

Annotate

Clicking on Annotate will open a separate image window displaying the captured image. The default name for the image will include the date, time, and user (if known).

Appendix B. Chemi Doc

The Annotation Tools toolbar will also pop up to allow you to annotate your image.

The image will not be saved until you select Save or Save As from the File menu.

Analyze

Clicking on Analyze will open a separate image window displaying the captured image. The default name for the image will include the date, time, and user (if known).

You can then save and analyze the image using the standard menu and toolbar commands.

Save

Clicking on Save will open a separate image window displaying the captured image. A Save As dialog box will automatically open displaying the default file name for the image, which will include the date, time, and user (if known). You can then change the file name and storage directory.

You can also save your image as a TIFF image for export to other applications.

B.5 Image Mode

The Image Mode option buttons change the type and scale of the data, as well as the behavior of the Chemi Doc when acquiring an image.

UV

Select this mode for fluorescent and chemiluminescent samples. With this mode selected, the data will be measured in linear intensity units.

PDQUEST User Guide

White Light

Select this mode for reflective and transmissive samples. With this mode selected, the data will be measured in uncalibrated optical density (uOD) units.

Chemi

This mode is designed for chemiluminescent samples. With this mode selected, the data is measured in linear intensity units; however, the data is inverted, so that samples will appear dark on a light background.

Note: This setting will invert the data in your image, so that a pixel with an intensity of 255 will be changed to 0, and visa versa.

Also, Chemi mode changes the behavior of the Auto and Manual Expose functions, as described above.

B.6 Exposure Status

The Exposure Status bar shows the progress of your exposure. If your exposure time is greater than 1 second, the status bar display will give you a graphical representation of the remaining time before exposure is complete.

If the exposure time is less than 1 second, the status bar will not refresh itself for each exposure; it will remain at 100 percent.

B.7 Display

The Display controls are useful for quickly adjusting the appearance of your image for output to a video printer. Adjusting these controls will automatically freeze the video display and allow you to alter the image within the Chemi Doc window.

These controls are similar to those in the Transform dialog box.

Appendix B. Chemi Doc

Note: The Display controls will only change the appearance of the image. They will not change the underlying data.

Highlight Saturated Pixels

When this box is checked, any saturated pixels in the image will appear highlighted in red in the scan window and in the pop-up image window. To view / hide saturated pixels in the pop-up image window, use the Image > Transform command.

Invert Display

This checkbox will switch light spots on a dark background to dark spots on a light background, and visa versa.

Note: This will only affect how the image appears on your screen; it will not change the image data.

Auto-scale

Clicking on Auto-scale will adjust your displayed image automatically. The lightest part of the image will be set to the minimum intensity (e.g., white), and the darkest will be set to the maximum intensity (e.g., black). You can then “fine tune” the display using the High, Low, and Gamma sliders described below.

High/Low Sliders

If Auto-scale doesn't give you the appearance you want, you can use the High and Low sliders to redraw the image yourself. Dragging the High slider handle to the left will make weak signals appear darker. Dragging the Low slider handle to the right will reduce background noise.

You can also type specific High and Low values in the text boxes next to the sliders. Clicking anywhere on the slider bars will move the sliders incrementally.

PDQUEST User Guide

Gamma Slider

Some images may be more effectively visualized if their data are mapped to the computer screen in a nonlinear fashion. Adjusting the Gamma slider handle changes the light and dark contrast nonlinearly.

Reset

Reset will return the image to its original, unmodified appearance.

B.8 Options

Clicking on the Options button opens the Options dialog box. Here you can specify certain settings for your Chemi Doc system.

Appendix B. Chemi Doc

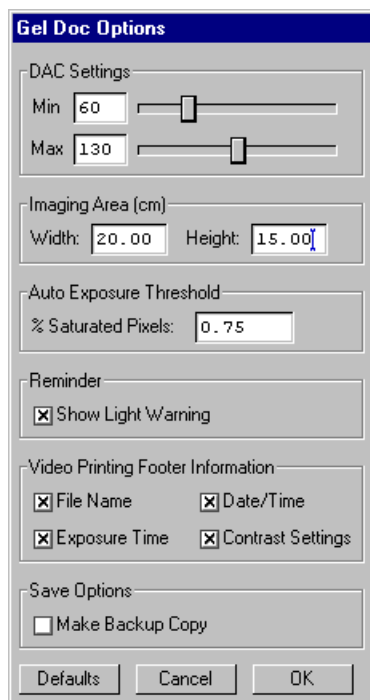


Fig.B-7. Available options in the Chemi Doc acquisition window.

Clicking on OK implements any changes you make in this box. Clicking on Defaults restores the settings to the factory defaults.

DAC Settings

Note: The default DAC settings are highly recommended and should be changed with caution.

These sliders may be used to adjust the minimum and maximum voltage settings of your video capture board. The minimum slider defines the pixel value that will appear as white in the image, while the maximum slider defines the pixel value that will appear as black. The slider scale is 0–255, with the defaults set to 60 minimum and 130 maximum.

PDQUEST User Guide

Imaging Area

These fields are used to specify the size of your imaging area in centimeters, which in turns determines the size of the pixels in your image (i.e., resolution). When you adjust one imaging area dimension, the other dimension will change to maintain the aspect ratio of the camera lens.

Note: Your imaging area settings must be correct if you want to do 1:1 printing. These are also important if you are comparing the size of objects (e.g., using the Volume Tools) between images.

Auto Exposure Threshold

When you click on Auto Expose, the exposure time is determined by the percentage of saturated pixels you want in your image. This field allows you to specify that percentage.

Typically, you will want less than 1 percent of the pixels in your image saturated. Consequently, the default value for this field is 0.75 percent.

Reminder

When this checkbox is selected, the software will warn you to turn off your transilluminator light when you exit the Chemi Doc acquisition window or when your system is “idle” for more than 5 minutes.

Note: If you are performing experiments that are longer than 5 minutes (e.g., chemiluminescence), this should be deselected.

Video Printing Footer Information

The checkboxes in this group allow you to specify the information that will appear at the bottom of your video printer printouts.

Save Options

To automatically create a backup copy of any scan you create, select the Make Backup Copy checkbox.

Appendix B. Chemi Doc

With this checkbox selected, when you save a scan, a backup copy will be placed in the same directory as the scanned image. Windows backup files will have an “.sbk” extension. Macintosh backup files will have the word “backup” after the file name.

This backup copy will be read-only, which means that you cannot make changes to it. You can open it like a normal file, but you must save it under a different file name before editing the image or performing analysis.

PDQUEST User Guide

Appendix C

GS-700 Imaging Densitometer



Fig.C-1. GS-700 Imaging Densitometer

Before you can begin scanning images with the GS-700 Imaging Densitometer®, your instrument must be properly installed and connected with the host computer. See the hardware manual for installation, startup, and operating instructions.

PC Only: A Note About SCSI Cards

The GS-700 is connected to your computer by a Small Computer System Interface (SCSI) cable. To use the GS-700, you must have a SCSI card installed in your PC. If you have an older PC, you may also need to load the SCSI and WinASPI drivers that came with your card.

Simulation Mode

Any of the imaging device acquisition windows can be opened in “simulation mode.” In this mode, an acquisition window will open and the controls will

PDQUEST User Guide

appear active, but instead of capturing real images, the window will create “dummy” images of manufactured data.

You do not need to be connected to an imaging device to open a simulated acquisition window. This is useful for demonstration purposes or practice scans.

To enter simulation mode, hold down the CTRL key and select the name of the device from the File menu. The title of the acquisition window will indicate that it is simulated.

Note: There is no simulated calibration for the GS-700 and GS-710 Imaging Densitometers.

C.1 GS-700 Acquisition Window

To acquire images using the GS-700, go to the File menu and select GS-700.... The acquisition window for the densitometer will open, displaying a control panel and a scanning window.

Appendix C. GS-700

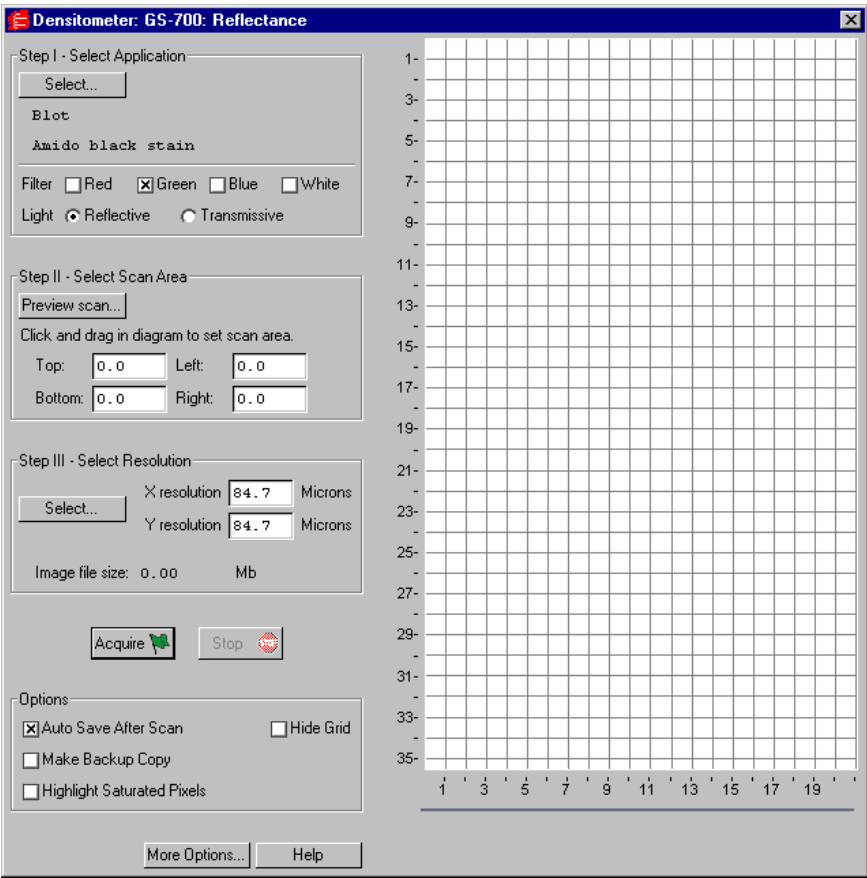


Fig.C-2. GS-700 acquisition window

The scanning window is marked by grid lines that divide the area into square centimeters. These are numbered 1–35 top to bottom and 1–21 left to right if the light source is uncalibrated reflective, and 1–25 top to bottom and 1–20 left to right if the light source is uncalibrated transmissive (see below for details). To hide the gridlines, click on the Hide Grid checkbox under Options.

PDQUEST User Guide

The control panel has been arranged from top to bottom to guide you through the acquisition procedure. There are four basic steps to scanning an image using the GS-700:

1. Select the application.
2. Select the scan area.
3. Select the resolution of the scan.
4. Acquire the image.

C.2 Step I. Select Application

To set the parameters for your particular scan, you can:

1. Select from a list of possible applications, or
2. Choose your own filter and light source settings.

Selecting from the List of Applications

When you select from the list of applications, the software automatically sets the appropriate filter(s) and other parameters for that particular application.

To select from the list of applications, click on the Select button under Select Application.

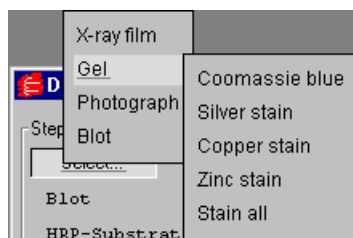


Fig.C-3. Example of the application tree in the GS-700 dialog box.

Appendix C. GS-700

The applications are listed in a tree that expands from left to right. First you select the category of your application, then you select your particular application.

To exit the tree without selecting, press the ESC key.

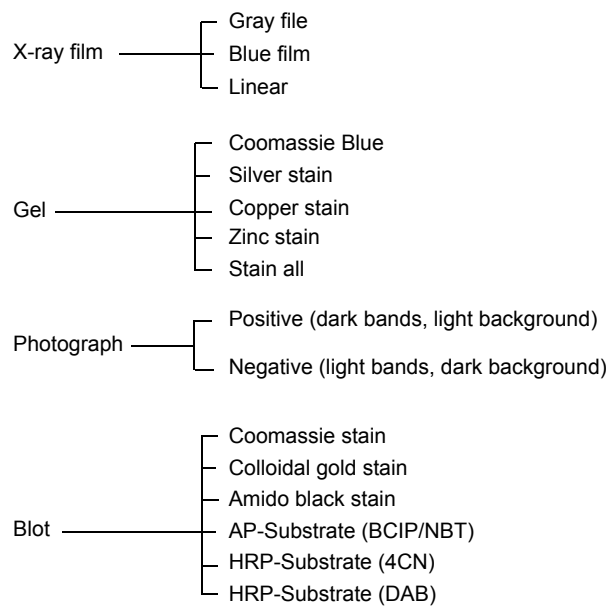


Fig.C-4. Applications available in the GS-700 acquisition window.

Choosing Your Own Settings

If you know the filter and light source settings you want, or want to experiment with different settings, you can choose them yourself.

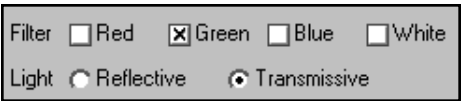


Fig.C-5. GS-700 Custom Application controls.

PDQUEST User Guide

Next to Filter, click on either the Red, Green, Blue, or White checkbox, or a combination of two of the first three (Red-Green, Green-Blue, Red-Blue).

Filter Color	Wavelength	Application Examples
Red	595–750 nm	Coomassie G-250, BCIP/NBT, Fast Green FCF, Methylene Blue
Green	520–570 nm	Coomassie R-250, Basic Fuchsin
Blue	400–530 nm	Crocein Scarlet
Gray Scale	400–750 nm	Silver Stains, Copper Stains, Film

Fig.C-6. Examples of filter colors and applications.

Next to Light, select the appropriate light source.

- Reflective mode is for scanning opaque mediums such as dried gels on filter paper, arrays, TLC plates, and photographs. The scanning dimensions in this mode are 21 cm x 35 cm (uncalibrated).
- Transmissive mode is for scanning transparent mediums such as films, gels, and slides. The scanning dimensions in this mode are 20 cm x 25 cm (uncalibrated).

C.3 Step II. Select Scan Area

Preview Scan

Before selecting the particular area to scan, you can preview the entire scanning area to determine the exact position of your sample.

Click on Preview Scan. A quick, low-resolution scan of the entire scanning area will appear in the scanning window. Your sample should be visible within a portion of this scan.

Appendix C. GS-700

Selecting an Area

Using the preview scan as a guide, select your scan area by dragging your mouse within the scanning window. The border of the scan area you are selecting will be marked by a frame.

Note: The scan area you select must be at least 1 cm wide.

When you release the mouse button, the border changes to a dashed blue line, indicating a selected area.

- To *reposition* the scanning area box you have created, position your cursor inside the box and drag. The entire box will move.
- To *resize* the box, position your cursor on a box side and drag. The side you have selected will move.
- To *redo* the box entirely, position your cursor outside the box and drag. The old box will disappear and a new box will be created.

You can also select the scanning area by entering coordinates in the appropriate fields (Top, Bottom, Left, Right). As you enter a coordinate, the position of the scanning area box will change accordingly.

When selecting, be sure to include the entire area of interest, and be generous with borders. You can always crop the image later.

C.4 Step III. Select Resolution

To select from a list of possible scanner resolutions, click on the Select button under Select Resolution. This will open the Select Scan Resolution dialog box.

PDQUEST User Guide

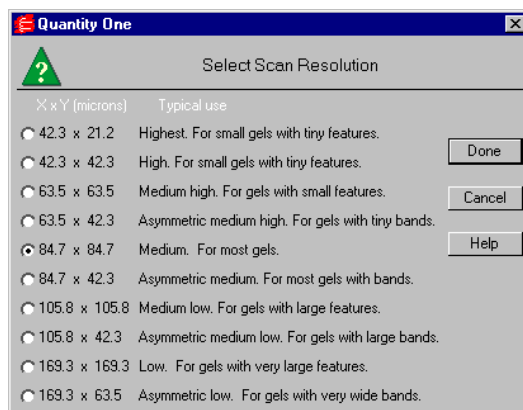


Fig.C-7. Select Scan Resolution dialog box.

Available resolutions are listed from highest to lowest in terms of the dimensions of the resulting pixels (in microns). Smaller pixels equal higher resolution. Each resolution is listed with its typical use.

In general, the size of your pixels should be one-tenth the height of your smallest object.

Some of the resolutions are asymmetrical, meaning that the resolution is higher in the vertical dimension (i.e., the pixels in the resulting image are larger in the horizontal dimension than in the vertical dimension). This is useful for gels with bands, where you are more interested in resolving in the y dimension to determine the size of bands and the spacing between them.

Specifying Your Own Resolution

If you select Oversample under More Options (see section C.7 for details), you can specify your own resolution within the range of 43–169 microns (micrometers). With Oversample selected, enter values directly in the fields next to X resolution and Y resolution in the main acquisition window.

Appendix C. GS-700

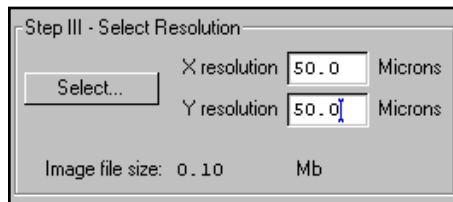


Fig.C-8. Entering a custom resolution (with Oversample selected).

Image File Size

Image File Size (under Select Resolution) shows the size of the scan file you are about to create. If you do not have enough computer memory for the specified file size, an error message will appear when you attempt to scan. If this happens, select a lower resolution or decrease the size of the area to be scanned. (Macintosh users can also increase the application memory partition. See your Macintosh computer documentation for guidance.)

C.5 Acquire the Image

If you want to calibrate your scans, read the following section on calibration before scanning.

To begin to scan, click on the Acquire button. The scanned image will begin to appear in the scanning window, line by line.

To interrupt a scan, click on the Stop button. A message will ask you to confirm the interrupt, and then you will be asked if you want to keep the partial scan. This feature is useful if you overestimated the size of the area you selected.

After the scan is complete, a window will open displaying the scan image, at which point you can analyze and save it.

Note: The image will open with a default file name that includes the date, time, and (if applicable) user name. However, unless you have selected Auto

PDQUEST User Guide

Save After Scan, the file will not be saved until you select Save or Save As from the File menu.

C.6 Calibration

If you have installed a calibration overlay, you can automatically calibrate your transmissive and reflective scans. (Calibration overlays for the GS-700 can be ordered from Bio-Rad.)

To set the automatic calibration settings, click on the More Options button in the GS-700 acquisition window. This will open the Densitometer Options dialog box.

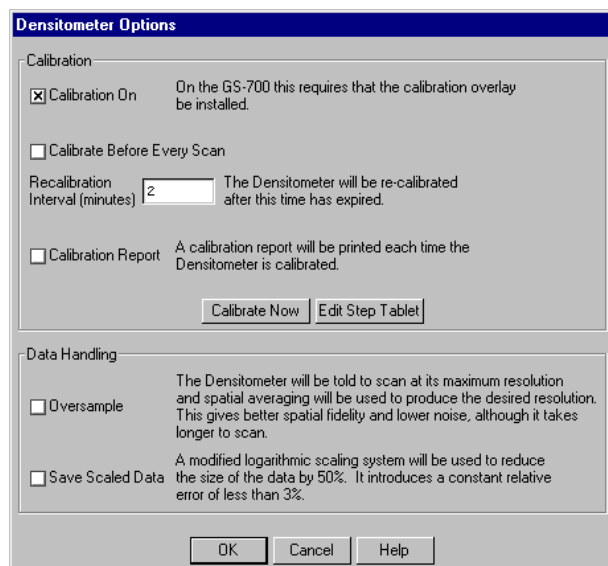


Fig.C-9. More Options in the GS-700 acquisition window.

To enable automatic calibration, click on the Calibration On checkbox.

Note: Before you can calibrate, you must enter the correct step tablet values for your transmissive calibration strip into the Step Tablet Values form, as

Appendix C. GS-700

described below. You do not need to change the values in the reflective step tablet form; you can use the default values.

With Calibration On selected, the other calibration settings become active.

Calibration Strip Window

When calibration is turned on, a calibration strip window will appear below the main scanning window and the length of the main scanning window will be reduced to 29 cm reflective and 23 cm transmissive.

Every time the densitometer calibrates, an image of the calibration strip will appear in the calibration strip window.

C.6.a Editing the Step Tablet


Before you can calibrate, you must enter the step tablet values for your transmissive calibration strip into the Step Tablet Values form.

In the package with your calibration strip, you will find a printout of the diffuse density values for each segment of your transmissive calibration strip. Those exact density values must be entered into the computer for the software to associate a correct density value with each step on the step tablet.

Note: Scanning in transmissive mode with incorrect step tablet values entered into the form can cause significant errors in the reported densities of your scans.

To enter the density values for your step tablet, click on the Edit Step Tablet button in the Densitometer Options box. The Step Tablet Values dialog box for the GS-700 will open.

PDQUEST User Guide

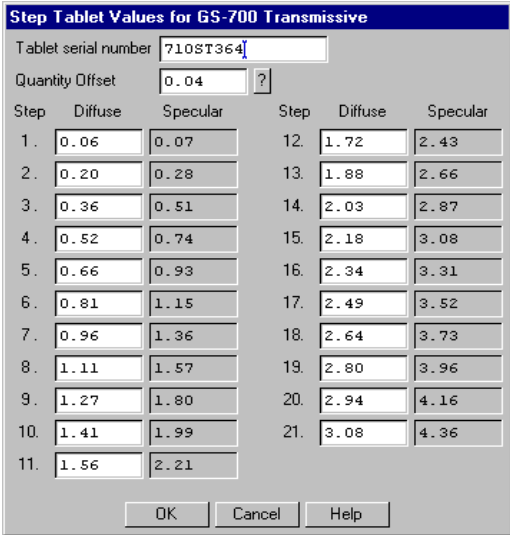

The dialog box is titled "Step Tablet Values for GS-700 Transmissive". It contains a "Tablet serial number" field with the value "710ST364" and a "Quantity Offset" field with the value "0.04" and a question mark icon. Below these fields is a table with two columns of data, each with three sub-columns: "Step", "Diffuse", and "Specular". The table contains 21 rows of data. At the bottom of the dialog box are three buttons: "OK", "Cancel", and "Help".

Step	Diffuse	Specular	Step	Diffuse	Specular
1.	0.06	0.07	12.	1.72	2.43
2.	0.20	0.28	13.	1.88	2.66
3.	0.36	0.51	14.	2.03	2.87
4.	0.52	0.74	15.	2.18	3.08
5.	0.66	0.93	16.	2.34	3.31
6.	0.81	1.15	17.	2.49	3.52
7.	0.96	1.36	18.	2.64	3.73
8.	1.11	1.57	19.	2.80	3.96
9.	1.27	1.80	20.	2.94	4.16
10.	1.41	1.99	21.	3.08	4.36
11.	1.56	2.21			

Fig.C-10. Step Tablet dialog box.

The dialog box will indicate whether the displayed form is for the transmissive or reflective step tablet, depending on whether you are in transmissive or reflective mode.

Note: In transmissive scanning, the scanner is calibrated against a transmissive step tablet. In reflective scanning, the scanner is calibrated against a reflective step tablet, unless the reflective calibration strip is uncalibrated. If the reflective calibration strip does not have step tablet values associated with it, you can use the default values in the form.

In the Tablet Serial Number field, you can type in the serial number for the step tablet values you are entering. (The reflective step tablet will probably not have a serial number.)

Your calibration strip will have a clear plastic cover. If the quantity offset value for the clear cover is included in your calibration strip package, you can enter this number in the Quantity Offset field. Otherwise, use a value of 0.045.

Appendix C. GS-700

Finally, enter the values in the appropriate fields under the Diffuse column. After the step tablet is scanned, the software will associate each density value with its corresponding segment on the step tablet.

The step tablet density values do not need to be entered each time you calibrate. Once they have been entered and saved, they will be automatically recalled when the calibration strip is scanned. You only need to enter new values if you use a new step tablet.

When you finished filling out the Step Tablet Values form, click on OK.

Diffuse Versus Specular O.D.

In the step tablet form, you enter O.D. as diffuse density, and then the software automatically calculates the specular density.

Specular density is a measure of the light that passes directly through a medium. Diffuse density includes light that is scattered as it passes through the medium. Step tablet values are given in diffuse density, but are measured by the scanner in specular density, and therefore must be converted according to the specular / diffuse O.D. ratio. This conversion does not affect quantitation.

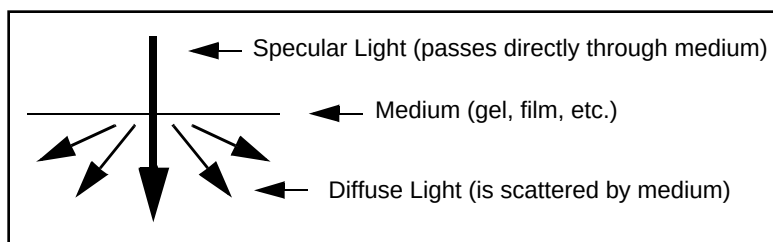


Fig.C-11. Specular and diffuse density

Diffuse density values are converted to specular optical density units according to the following formula:

$$\text{Specular OD} = 1.4 \text{ Diffuse OD}$$

PDQUEST User Guide

C.6.b Calibration Settings

After you have entered the step tablet values, you can immediately calibrate by clicking on the Calibrate Now button (in the Densitometer Options dialog box).

You can also specify how often you want the densitometer to automatically recalibrate. Either click on the Calibrate Before Every Scan checkbox or enter a Recalibration Interval (in minutes) in the appropriate field.

Note: With calibration turned on, the scanner will automatically recalibrate each time you change your filter or your reflective/transmissive setting. (If you select a different application with the same filter and light settings, it will not auto recalibrate.)

Calibration Report

To print out a calibration report each time the densitometer calibrates, click on the Calibration Report checkbox.

C.7 Other Options

Oversample

This feature allows you to scan at the maximum resolution of the GS-700 (42.3 x 42.3 microns) and then use spatial averaging to create an image with lower resolution. This can result in better images at lower resolution—however, it takes longer to scan.

To turn on oversampling, click on the More Options button in the acquisition window and click on the Oversample checkbox.

With oversampling on, you can specify your own resolution within the range of 43–169 microns by entering values directly in the fields next to X resolution and Y resolution in the main acquisition window.

Appendix C. GS-700

Auto Save After Scan

To automatically save any scan you create, click on the Auto Save After Scan checkbox.

With this checkbox selected, when you click on Acquire, a Save As dialog box will open asking you to specify a file name and location for the image you are about to create. The scan will begin when you click on the Save button.

Make Backup Copy

If you have checked Auto Save After Scan, you can also automatically create a backup copy of any scan you create.

Click on the Make Backup Copy checkbox. With this checkbox selected, when a scan is created and saved, a backup copy will be placed in the same directory as the scanned image. Windows backup files will have an “.sbk” extension. Macintosh backup files will have the word “backup” after the file name.

This backup copy will be read-only, which means that you cannot make changes to it. You can open it like a normal file, but you must save it under a different file name before editing the image or performing analysis.

Highlight Saturated Pixels

When this box is checked, any saturated pixels in the image will appear highlighted in red in the scan window and in the pop-up image window. To view / hide saturated pixels in the pop-up image window, use the Image > Transform command.

Hide Grid

To hide the gridlines in the scanning area window, click on the Hide Grid checkbox.

PDQUEST User Guide

Appendix D

GS-710 Imaging Densitometer



Fig.D-1. GS-710 Imaging Densitometer

Before you can begin scanning images with the GS-710 Imaging Densitometer[®], your instrument must be properly installed and connected with the host computer. See the hardware manual for installation, startup, and operating instructions.

PC Only: A Note About SCSI Cards

The GS-710 is connected to your computer by a Small Computer System Interface (SCSI) cable. To use the GS-710, you must have a SCSI card installed in your PC. If you have an older PC, you may also need to load the SCSI and WinASPI drivers that came with your card.

Simulation Mode

Any of the imaging device acquisition windows can be opened in “simulation mode.” In this mode, an acquisition window will open and the controls will

PDQUEST User Guide

appear active, but instead of capturing real images, the window will create “dummy” images of manufactured data.

You do not need to be connected to an imaging device to open a simulated acquisition window. This is useful for demonstration purposes or practice scans.

To enter simulation mode, hold down the CTRL key and select the name of the device from the File menu. The title of the acquisition window will indicate that it is simulated.

Note: There is no simulated calibration for the GS-700 and GS-710 Imaging Densitometers.

D.1 GS-710 Acquisition Window

To acquire images using the GS-710, go to the File menu and select GS-710.... The acquisition window for the densitometer will open, displaying a control panel and a scanning window.

Appendix D. GS-710

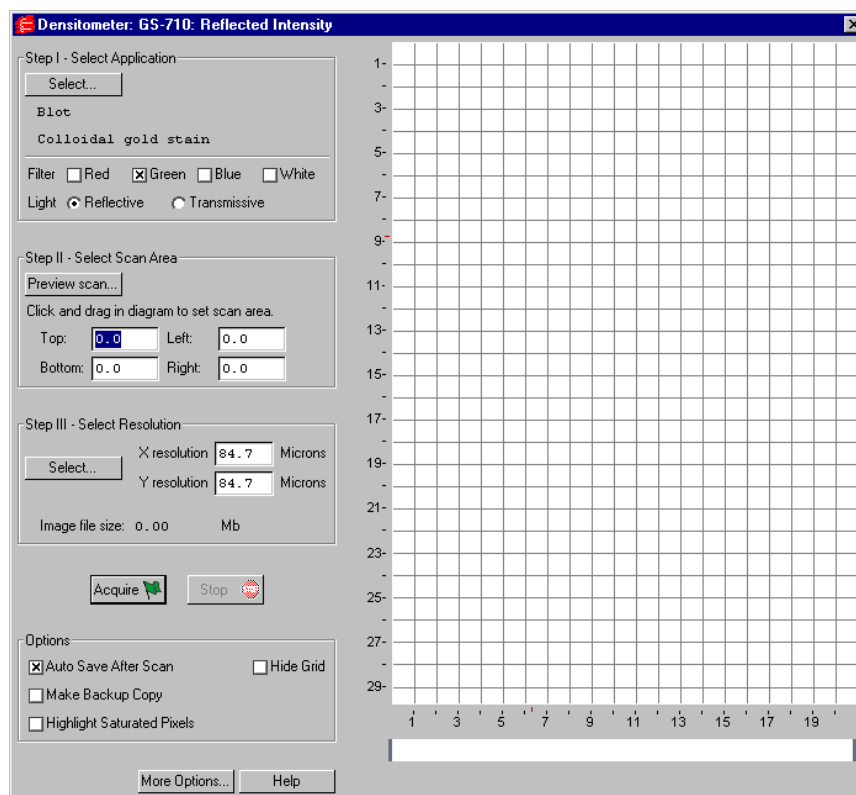


Fig.D-2. GS-710 acquisition window

The scanning window is marked by grid lines that divide the area into square centimeters. These are numbered 1–29 top to bottom and 1–21 left to right if the light source is reflective, and 1–23 top to bottom and 1–20 left to right if the light source is transmissive (see below).

To hide the gridlines, click on the Hide Grid checkbox under Options.

Below the main scanning window is the calibration strip window. Every time the densitometer calibrates, an image of the calibration strip will appear in this window.

PDQUEST User Guide

The control panel has been arranged from top to bottom to guide you through the acquisition procedure. There are five basic steps to scanning an image using the GS-710:

1. Select the application.
2. Select the scan area.
3. Select the resolution of the scan.
4. Calibrate the instrument. (This is automatic, after you enter the step tablet values before you first scan after installation.)
5. Acquire the image.

D.2 Step 1. Select Application

To set the parameters for your particular scan, you can:

1. Select from a list of possible applications, or
2. Choose your own filter and light source settings.

Selecting from the List of Applications

When you select from the list of applications, the software automatically sets the appropriate filter(s) and other parameters for that particular application.

To select from the list of applications, click on the Select button under Select Application.

Appendix D. GS-710

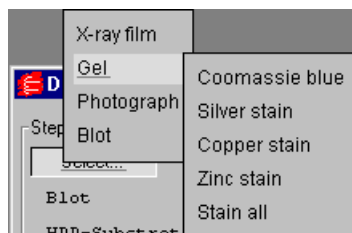


Fig.D-3. Example of the application tree in the GS-710 dialog box.

The applications are listed in a tree that expands from left to right. First you select the category of your application, then you select your particular application.

To exit the tree without selecting, press the ESC key.

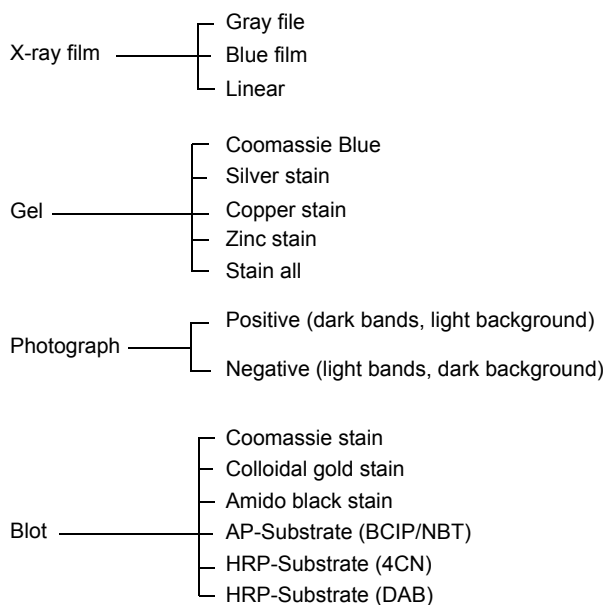


Fig.D-4. Applications available in the GS-710 acquisition window.

PDQUEST User Guide

Choosing Your Own Settings

If you know the filter and light source settings you want, or want to experiment with different settings, you can choose them yourself.

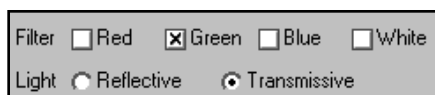


Fig.D-5. GS-710 Custom Application controls.

Next to Filter, click on either the Red, Green, Blue, or White checkbox, or a combination of two of the first three (Red-Green, Green-Blue, Red-Blue).

Filter Color	Wavelength	Application Examples
Red	595–750 nm	Coomassie G-250, BCIP/NBT, Fast Green FCF, Methylene Blue
Green	520–570 nm	Coomassie R-250, Basic Fuchsin
Blue	400–530 nm	Crocein Scarlet
Gray Scale	400–750 nm	Silver Stains, Copper Stains, Film

Fig.D-6. Examples of filter colors and applications.

Next to Light, select the appropriate light source.

- Reflective mode is for scanning opaque mediums such as dried gels on filter paper, arrays, or photographs. The scanning dimensions in this mode are 21 cm x 29 cm.
- Transmissive mode is for scanning transparent mediums such as films, gels, or slides. The scanning dimensions in this mode are 20 cm x 23 cm.

D.3 Step II. Select Scan Area

Preview Scan

Before selecting the particular area to scan, you can preview the entire scanning area to determine the exact position of your sample.

Click on Preview Scan. A quick, low-resolution scan of the entire scanning area will appear in the scanning window. Your sample should be visible within a portion of this scan.

Selecting an Area

Using the preview scan as a guide, select your scan area by dragging your mouse within the scanning window. The border of the scan area you are selecting will be marked by a frame.

Note: The scan area you select must be at least 1 cm wide.

When you release the mouse button, the border changes to a dashed blue line, indicating a selected area.

- To *reposition* the scanning area box you have created, position your cursor inside the box and drag. The entire box will move.
- To *resize* the box, position your cursor on a box side and drag. The side you have selected will move.
- To *redo* the box entirely, position your cursor outside the box and drag. The old box will disappear and a new box will be created.

You can also select the scanning area by entering coordinates in the appropriate fields (Top, Bottom, Left, Right). As you enter a coordinate, the position of the scanning area box will change accordingly.

When selecting, be sure to include the entire area of interest, and be generous with borders. You can always crop the image later.

PDQUEST User Guide

D.4 Step III. Select Resolution

To select from a list of possible scanner resolutions, click on the Select button under Select Resolution. This will open the Select Scan Resolution dialog box.

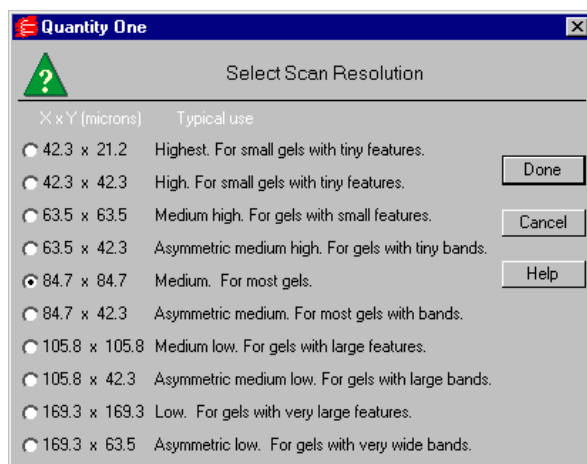


Fig.D-7. Select Scan Resolution dialog box.

Available resolutions are listed from highest to lowest in terms of the dimensions of the resulting pixels (in microns). Smaller pixels equal higher resolution. Each resolution is listed with its typical use.

In general, the size of your pixels should be one-tenth the height of your smallest object.

Some of the resolutions are asymmetrical, meaning that the resolution is higher in the vertical dimension (i.e., the pixels in the resulting image are larger in the horizontal dimension than in the vertical dimension). This is useful for gels with bands, where you are more interested in resolving in the y dimension to determine the size of bands and the spacing between them.

Appendix D. GS-710

Specifying Your Own Resolution

If you select Oversample under More Options, you can specify your own resolution within the range of 43–169 microns (micrometers). With Oversample selected, enter values directly in the fields next to X resolution and Y resolution in the main acquisition window.

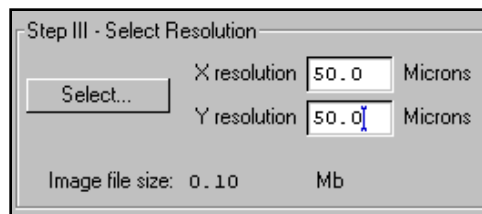


Fig.D-8. Entering a custom resolution (with Oversample selected).

Image File Size

Image File Size (under Select Resolution) shows the size of the scan file you are about to create. If you do not have enough computer memory for the specified file size, an error message will appear when you attempt to scan. If this happens, select a lower resolution or decrease the size of the area to be scanned. (Macintosh users can also increase the application memory partition. See your Macintosh computer documentation for guidance.)

D.5 Calibration

Images scanned using the GS-710 are automatically calibrated. The scanner plate has a built-in calibration overlay for both transmissive and reflective scanning.

- The transmissive calibration strip calibrates the densitometer from 0 to 3.0 O.D.
- The reflective calibration strip calibrates the densitometer to 2.0 O.D.

PDQUEST User Guide

The first time you use the GS-710, you must select some calibration settings to ensure that your calibration is accurate.

To set the automatic calibration settings, click on the More Options button in the GS-710 acquisition window. This will open the Densitometer Options dialog box.

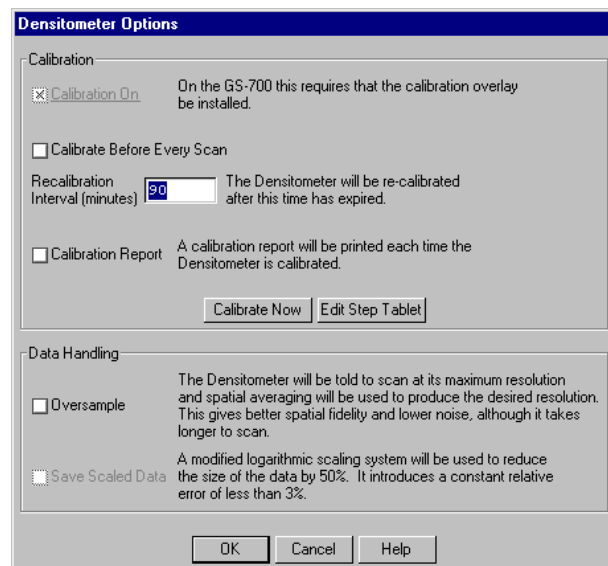


Fig.D-9. Densitometer Options dialog box.

D.5.a Editing the Step Tablet

Before you scan for the first time, you must enter the step tablet values for your transmissive calibration strip into the Step Tablet Values form.

Note: In reflective scanning, the scanner is calibrated against the reflective step tablet. However, you do not need to enter the values for this step tablet; you can use the default values in the Step Tablet Values form.

Attached to the outside of your GS-710, you will find a copy of the manufacturer's printout of the diffuse density values for each segment of your transmissive step tablet, as well as a serial number that identifies the

Appendix D. GS-710

step tablet to which those density values belong. Those exact density values must be entered into the computer for the software to associate a correct density value with each step on the step tablet.

Note: Scanning in transmissive mode with incorrect step tablet values entered into the computer can cause significant errors in the reported densities of your scans.

To enter the density values for your step tablet, click on the Edit Step Tablet button in the Densitometer Options box. The Step Tablet Values dialog box for the GS-710 will open.

Step	Diffuse	Specular	Step	Diffuse	Specular
1.	0.05	0.07	12.	1.72	2.43
2.	0.20	0.28	13.	1.88	2.66
3.	0.36	0.51	14.	2.03	2.87
4.	0.52	0.74	15.	2.18	3.08
5.	0.66	0.93	16.	2.34	3.31
6.	0.81	1.15	17.	2.49	3.52
7.	0.96	1.36	18.	2.64	3.73
8.	1.11	1.57	19.	2.80	3.96
9.	1.27	1.80	20.	2.94	4.16
10.	1.41	1.99	21.	3.08	4.36
11.	1.56	2.21			

Fig.D-10. Step Tablet dialog box.

The dialog box will indicate whether the values are for the transmissive or reflective step tablet, depending on whether you are in transmissive or reflective mode.

If you are in reflective mode, you do not need to change the default step tablet values.

PDQUEST User Guide

If you are in transmissive mode, you do need to enter the correct step tablet values that were shipped with your GS-710. You can type the serial number for these values into the Tablet Serial Number field.

The Quantity Offset field does not apply in the GS-710. This value should be set to zero.

Next, enter the values for the transmissive step tablet in the appropriate fields under the Diffuse column. After the step tablet is scanned, the software will associate each density value with its corresponding segment on the step tablet.

The step tablet density values do not need to be entered each time you calibrate. Once they have been entered and saved, they will be automatically recalled when the calibration strip is scanned.

When you are finished filling out the Step Tablet Values form, click on OK.

Diffuse Versus Specular O.D.

In the step tablet form, you enter O.D. as diffuse density, and then the software automatically calculates the specular density.

Specular density is a measure of the light that passes directly through a medium. Diffuse density includes light that is scattered as it passes through the medium. Step tablet values are given in diffuse density, but are measured by the scanner in specular density, and therefore must be converted according to the specular / diffuse O.D. ratio. This conversion does not affect quantitation.

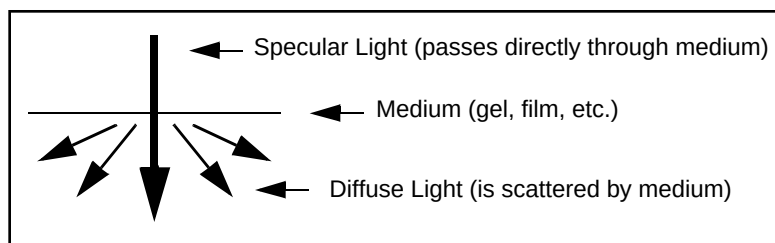


Fig.D-11. Specular and diffuse density

Appendix D. GS-710

Diffuse density values are converted to specular optical density units according to the following formula:

$$\text{Specular OD} = 1.4 \text{ Diffuse OD}$$

D.5.b Calibration Settings

After you have entered the step tablet values, you can immediately calibrate by clicking on the Calibrate Now button (in the Densitometer Options dialog box).

You can also specify how often you want the GS-710 to automatically recalibrate. Either click on the Calibrate Before Every Scan checkbox or enter a Recalibration Interval (in minutes) in the appropriate field.

Note: The scanner will automatically recalibrate each time you change your filter or your reflective/transmissive setting. (If you select a different application with the same filter and light settings, it will not auto recalibrate.)

Calibration Report

To print out a calibration report each time the densitometer calibrates, click on the Calibration Report checkbox.

D.6 Acquire the Image

To begin to scan, click on the Acquire button. The scanned image will begin to appear in the scanning window, line by line.

To interrupt a scan, click on the Stop button. A message will ask you to confirm the interrupt, and then you will be asked if you want to keep the partial scan. This feature is useful if you overestimated the size of the area you selected.

After the scan is complete, a window will open displaying the scan image, at which point you can analyze and save it.

PDQUEST User Guide

Note: The image will open with a default file name that includes the date, time, and (if applicable) user name. However, unless you have selected Auto Save After Scan, the file will not be saved until you select Save or Save As from the File menu.

D.7 Other Options

Oversample

This feature allows you to scan at the maximum resolution of the GS-710 (42.3 x 42.3 microns) and then use spatial averaging to create an image with lower resolution. This can result in better images at lower resolution—however, it takes longer to scan.

To turn on oversampling, click on the More Options button in the acquisition window and click on the Oversample checkbox.

With oversampling on, you can specify your own resolution within the range of 43–169 microns by entering values directly in the fields next to X resolution and Y resolution in the main acquisition window.

Auto Save After Scan

To automatically save any scan you create, click on the Auto Save After Scan checkbox.

With this checkbox selected, when you click on Acquire, a Save As dialog box will open asking you to specify a file name and location for the image you are about to create. The scan will begin when you click on the Save button.

Make Backup Copy

If you have checked Auto Save After Scan, you can also automatically create a backup copy of any scan you create.

Click on the Make Backup Copy checkbox. With this checkbox selected, when a scan is created and saved, a backup copy will be placed in the same

Appendix D. GS-710

directory as the scanned image. Windows backup files will have an “.sbk” extension. Macintosh backup files will have the word “backup” after the file name.

This backup copy will be read-only, which means that you cannot make changes to it. You can open it like a normal file, but you must save it under a different file name before editing the image or performing analysis.

Highlight Saturated Pixels

When this box is checked, any saturated pixels in the image will appear highlighted in red in the scan window and in the pop-up image window. To view/hide saturated pixels in the pop-up image window, use the Image > Transform command.

Hide Grid

To hide the gridlines in the scanning area window, click on the Hide Grid checkbox.

PDQUEST User Guide

Appendix E

Fluor-S Multimager



Fig.E-1. Fluor-S Multimager.

Before you can begin acquiring images using the Fluor-S[®] MultiImager, the imaging system must be properly installed and connected with the host computer. See the Fluor-S hardware manual for installation, startup, and operating instructions.

Note: Make sure that the temperature light on the Fluor-S is green before attempting to capture an image. If you are using a PC, the Fluor-S should be turned on and the initialization sequence completed *before* the host computer is turned on. See the hardware manual for more details.

PDQUEST User Guide

PC Only: A Note About SCSI Cards

The Fluor-S is connected to your computer by a Small Computer System Interface (SCSI) cable. To use the Fluor-S, you must have a SCSI card installed in your PC. If you have an older PC, you may also need to load the SCSI and WinASPI drivers that came with your card.

Simulation Mode

Any of the imaging device acquisition windows can be opened in “simulation mode.” In this mode, an acquisition window will open and the controls will appear active, but instead of capturing real images, the window will create “dummy” images of manufactured data.

You do not need to be connected to an imaging device to open a simulated acquisition window. This is useful for demonstration purposes or practice scans.

To enter simulation mode, hold down the CTRL key and select the name of the device from the File menu. The title of the acquisition window will indicate that it is simulated.

E.1 Fluor-S Acquisition Window

To acquire images using the Fluor-S, go to the File menu and select Fluor-S.... The acquisition window for the imager will open, displaying a control panel and an image display window.

Appendix E. Fluor-S



Fig.E-2. Fluor-S acquisition window

When the Fluor-S window first opens, no image will be displayed.

The control panel has been arranged from top to bottom to guide you through the acquisition procedure. There are four basic steps to acquiring an image using the Fluor-S:

1. Select the application.
2. Position and focus the gel or other object to be imaged.
3. Set the exposure time.
4. Acquire the image.

PDQUEST User Guide

E.2 Step 1. Select Application

To set the appropriate filter and other parameters for the type of object you are imaging, click on the Select button under Select Application. The available applications and their associated settings are listed in a tree that expands from left to right.

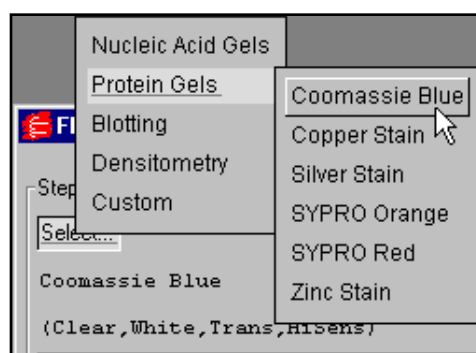


Fig.E-3. The application tree in the Fluor-S acquisition window.

First select your general application, then select the particular stain or medium you are using. If you select the chemiluminescent application under “Blotting,” you also have the option of selecting High Resolution or High Sensitivity (see below).

When you select an application, the software automatically sets the appropriate standard filter (520LP, 530DF60, 610LP, clear, or none), light type (UV, white, or none), and light source (trans, epi, or neither) in the Fluor-S for that particular application.

For applications involving trans illumination, you must also specify a scan dimension (see below).

Your selection will be displayed below the Select button. To exit the tree without selecting, press the ESC key.

Appendix E. Fluor-S

Custom Applications

If your application is not listed, if you want to use a user-installed filter, or if you want to access High Sensitivity mode (see below), you can create and save your own custom application.

From the application tree, select Custom, then Create. This will open a dialog box in which you can name your application and select your settings.

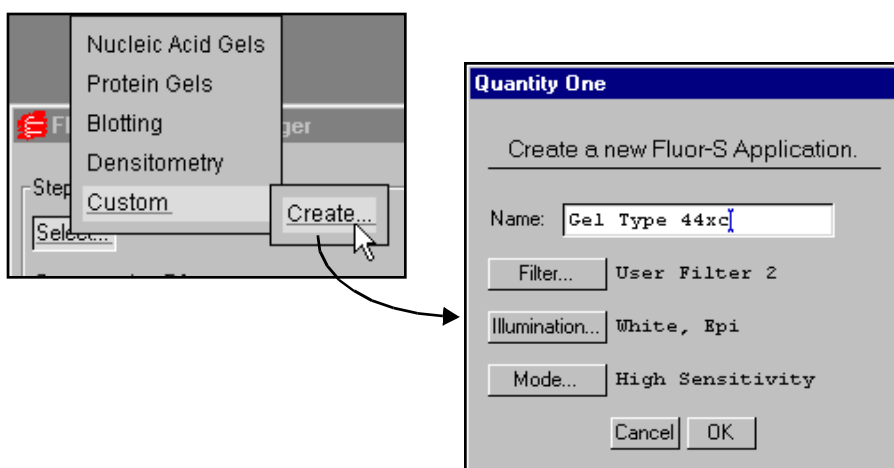


Fig.E-4. Creating a new custom application.

To select the filter (including user-defined), type of illumination, and camera mode, click on the appropriate buttons and select from the pop-up lists.

Note: Under Illumination, there is a listing for a spare UV light source. This selects the spare UV bulb in the Fluor-S. Select this light source if your main UV bulb fails.

Enter a name for your application in the Name field. Click on OK to implement your changes.

After you have created an application, you can select it from the application tree by selecting Custom and the name you created. You can delete the application by selecting Custom, Delete, and the name of the application.

PDQUEST User Guide

Scan Dimension

If an application uses trans illumination, the Scan Dimension buttons become active. The scan dimension is the distance traveled by the transilluminating light source as it scans horizontally across the platen.



Fig.E-5. Scan dimension settings for trans illumination.

The full scanning range is 300 mm. Select a smaller range if your sample is small and you do not want to wait while the light source travels over the maximum scan width.

High Resolution Versus High Sensitivity

High Resolution and High Sensitivity are mutually exclusive options. High Resolution is the normal operating mode for the Fluor-S. High Sensitivity provides optimal sensitivity for low-light applications. It is the default selection for the chemiluminescence application, and may be selected for a custom application.

In High Resolution mode, images are captured at the maximum resolution of the camera.

In High Sensitivity mode, the pixels in the camera are “binned” (i.e., four pixels are combined into one) to increase the amount of signal per pixel. However, combining the pixels results in a reduction in the resolution of the image.

E.3 Step II. Position/Focus

Note: When you click on the Position or Focus button, the light inside the Fluor-S box automatically turns on. To turn the light off while positioning or focusing, hold down the SHIFT key when clicking on the button.

Appendix E. Fluor-S

Position

After you have selected your application, you are ready to center your gel or other object within the camera frame. To do so, click on the Position button. The Fluor-S will begin capturing a “live” image and updating it every second.

With the Position button selected, look at the image in the acquisition window while you position your object in the center of the platen. If you have a zoom lens on the camera, you can adjust the magnification while you position.

You can select the Show Alignment Grid checkbox to facilitate positioning.

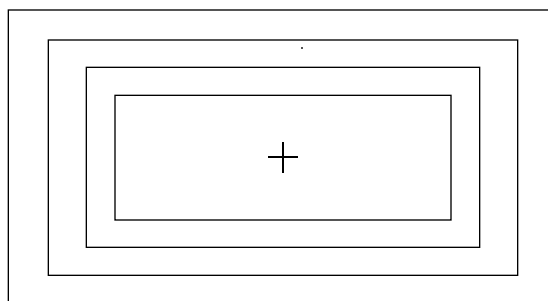


Fig.E-6. Fluor-S alignment grid.

When you are finished positioning, click on the Stop button.

Focus

Note: Before focusing, you should adjust the f-stop on the camera to the lowest setting (i.e., the maximum aperture). This reduces the depth of field, allowing you to more accurately focus the camera. Then, after focusing, increase the f-stop to the desired setting. (See the following table on recommended f-stops.)

After you have positioned your sample, click on the Focus button and look at the image in the acquisition window while aligning the two focusing arrows on the camera lens. While focusing, the camera will limit its focus to a small

PDQUEST User Guide

portion of the sample (this will not affect any zoom lens adjustments you may have made.)

When you are finished focusing, click on the Stop button.

E.4 Step III. Set Exposure Time

When you are ready to capture an image, you will need to select an exposure time. “Exposure” refers to the integration of image captures on the CCD over a set period of time. The effect is analogous to exposing photographic film to light.

The exposure time you select should be based on your application and your “best guess” as to what exposure will give you the best image.

Note: The minimum exposure time in trans illuminated mode is 1 second. The minimum exposure time in epi illuminated mode is 0.1 second.

You can enter an exposure time (in seconds) directly in the field, or use the Arrow buttons to adjust the exposure time in 10 percent increments.



Fig.E-7. Selecting an exposure time.

The following table provides recommended exposure times for various applications.

Appendix E. Fluor-S

Recommended Exposure Times and Lenses

Sample	Recommended Exposure ¹	Lens & Filter ²	Accessories Used
Fluorescent Stain Gel	3–40 sec.	Zoom/IR	None
Fluorescence End-Label Gel	30 sec.–3 min.	Zoom/IR	None
Fluorescent Blot	0.5–5 sec.	Zoom/IR	None
Chemifluorescent Blot	0.5–5 sec.	Zoom/IR	None
Colorimetric Gel	1–10 sec.	Zoom/IR	White Diffusion Plate
Colorimetric Blot	0.5–20 sec.	Zoom/IR	None
X-ray film	1–10 sec.	Zoom/IR	White Diffusion Plate
Weak Chemiluminescence ³	5–20 min.	50 mm	Chemi Tray (if sample is small)
Strong Chemiluminescence ³	30 sec.–2 min.	50 mm	Chemi Tray (if sample is small)
¹ Increase exposure time two fold for each step increase in f-stop. ² For sharper focusing, close the f-stop down 1–2 stops from full open while focusing. ³ For chemi applications, the 50mm lens is recommended. Always remove the 660 filter.			

For most applications, you can select an exposure time, capture an image, study it, then adjust the exposure time accordingly. Repeat this procedure as many times as necessary to obtain a good image.

For chemiluminescent samples, which degrade over time and emit low levels of light, you can select a high exposure time initially and use Live Acquire mode to save intermediate exposures (see following section).

Preview

For shorter exposures, you can use Preview to test different exposure times. Click on the Preview button create a preview exposure and display it in the acquisition window.

PDQUEST User Guide

A preview scan takes only half as long to create as a real scan, because the preview scan does not capture a “dark” image (see the following section on Options). The progress of the exposure will be displayed in the Exposure Status bar at the bottom of the dialog box.

You cannot save preview scans.

If you want to stop a preview scan that is in progress, click on the Stop button.

E.5 Acquire the Image

The Fluor-S gives you the option of simply acquiring and displaying a fully exposed image, or preserving intermediate exposures.

Acquire

To acquire and display a fully exposed image, click on the Acquire button. An exposure will be taken based on the time selected in Step III. This is appropriate for most short exposures.

The progress of the exposure will be displayed in the Exposure Status bar at the bottom of the acquisition window.

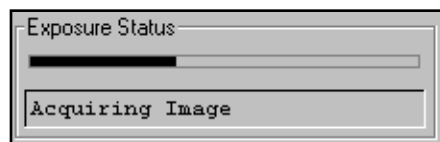


Fig.E-8. Exposure Status bar when acquiring an image.

Depending on which dark subtraction type you have selected, a dark count may be acquired immediately following image acquisition. See Dark Subtraction Type under Options, below.

If you want to stop a scan that is in progress, click on the Stop button. The acquisition will be terminated.

Appendix E. Fluor-S

After an image has been acquired, a separate window will pop up containing the new image. The window will have a default file name that includes the date, time, and user (if known). To save the image, select Save or Save As from the File menu.

You can then analyze the image using the analysis functions.

Live Acquire

Live Acquire allows you to view and preserve intermediate exposures leading up to a full exposure. This is useful for longer exposure times, such as those required by chemiluminescent samples, where there is the potential for image saturation.

When you click on the Live Acquire button, the exposure time you selected is divided by the number of exposure counts set in the Fluor-S Options dialog box (see Options, below). For example, if you enter an exposure time of 10 minutes and an exposure count of 20, then 20 intermediate exposures will be produced at 30-second intervals. Each intermediate exposure will be displayed in the scan window. The final, full exposure will be displayed in a separate image window.

Note: The first intermediate exposure will take longer than 30 seconds to display if dark subtraction is performed.

You can automatically save your intermediate exposures as separate files for later review using the Auto Save After Scan option. See Options, below.

If you see an intermediate exposure that you like in the scan window, click on the Stop button. Live Acquire mode will end and the last intermediate exposure to be completed will open in a separate image window. You can then save it for analysis.

E.6 Options

Click on the Options button to open the Fluor-S Options dialog box.

PDQUEST User Guide

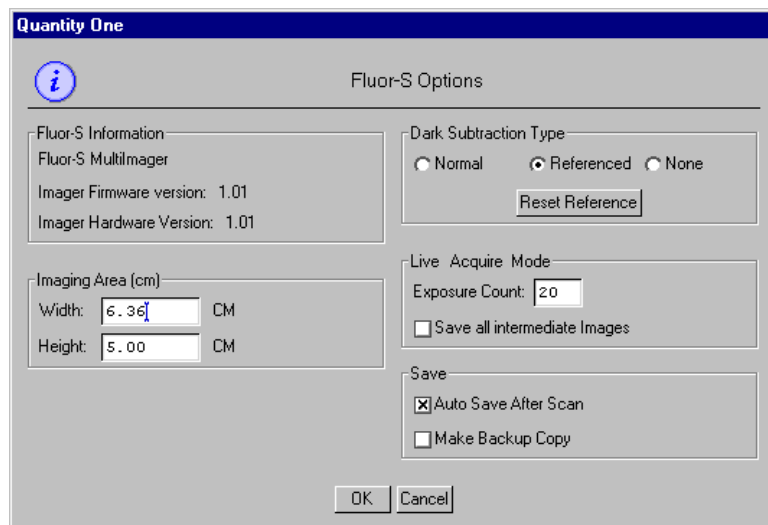


Fig.E-9. Options dialog box.

E.6.a Dark Subtraction Type

All CCD cameras accumulate electrons that produce a signal that is indistinguishable from light. This “dark count” adds to the noise in your images. In most cases, you will want to subtract this dark count from your images.

Normal

The Normal option button selects the default dark subtraction type. In this mode, after you acquire an image, a “dark” image of the same exposure length will be taken, and this will be subtracted from your image.

The progress of the dark exposure will be displayed in the Exposure Status bar following the regular image exposure.

Note: In Normal mode, a dark image is only acquired the first time you perform a scan with particular application and exposure settings. If you perform subsequent scans with the same settings, no dark exposure will be taken.

Appendix E. Fluor-S

Referenced

If you do not want to perform a dark exposure with each acquisition, you can take a “reference” dark exposure that will be saved and subtracted from all subsequent acquisitions. Click on the Referenced button to activate this feature.

The first time you acquire an image after selecting this option, the Fluor-S will take a 60-second dark exposure that will be saved and used to subtract the dark count from all subsequent acquisitions.

For image exposures of greater or less than 60 seconds, the reference dark will be scaled accordingly and then subtracted. You can change the default reference dark exposure time using the Reset Reference button (see below).

If you deselect the Referenced button and then reselect it, the old reference dark exposure will still be available.

Note: Separate reference dark exposures will be taken for High Resolution mode and High Sensitivity mode. Once you have created a reference dark in each of these modes, each reference dark will be used according to the mode you are in.

Reset Reference

If you would like a reference dark with an exposure time that more closely matches that of your typical scans, click on the Reset Reference button.



Fig.E-10. Reset Reference Dark pop-up box.

PDQUEST User Guide

A pop-up box will prompt you to enter a new reference dark exposure time in seconds. Click on OK to implement your change. The new reference dark will be created when you acquire your next image.

Note: Because of the high sensitivity of the CCD, fluctuations in background radiation and/or temperature in the room can affect the level of dark count. If you feel that radiation/temperature conditions have changed in the room since your last reference dark was created, use the Reset Reference button to delete your old reference and create a new one under current conditions.

None

If you do not want to perform dark subtraction, select None. No dark exposure will be acquired or subtracted.

E.6.b Live Exposure Mode

Exposure Count

If you are using the Live Acquire function (see previous section), you need to specify how many intermediate exposures you want to view/save during acquisition. Enter this number in the Exposure Count field.

The total exposure time will be divided by the number you enter in the Exposure Count field. If you enter an exposure time of 10 minutes and a count of 10, you will create 10 intermediate exposures at 1 minute intervals.

Note: Do not enter a count that will result in an intermediate exposure time that is less than the minimum exposure time for the mode you are in. The minimum exposure time in trans illuminated mode is 1 second, and the minimum exposure time in epi illuminated mode is 0.1 second. (Example: For a trans illuminated application, an exposure time of 20 seconds and an exposure count of 21 would result in an error.)

Save All Intermediate Images

If Auto Save After Scan is selected (see following section), the Save All Intermediate Images checkbox will become active. If you select this checkbox,

Appendix E. Fluor-S

all your intermediate exposures will be saved as separate files. These files will have the same root name appended by a number indicating the exposure sequence. The final, full exposure will have the root name only, with no exposure number.

E.6.c Save

Auto Save After Scan

To automatically save any image you create, click on the Auto Save After Scan checkbox.

With this checkbox selected, when you click on Acquire or Live Acquire, a Save As dialog box will open asking you to specify a file name and location for the image you are about to create. The scan will begin when you click on the Save button.

Note that in Live Acquire mode you can save your intermediate exposures by selecting Auto Save After Scan and then Save All Intermediate Images.

Make Backup Copy

You can automatically create a backup copy of any scan you create. To do so, first select Auto Save After Scan (see above), then select the Make Backup Copy checkbox.

With this checkbox selected, when you save a scan, a backup copy will be placed in the same directory as the scanned image. Windows backup files will have an “.sbk” extension. Macintosh backup files will have the word “backup” after the file name.

This backup copy will be read-only, which means that you cannot make changes to it. You can open it like a normal file, but you must save it under a different file name before editing the image or performing analysis.

PDQUEST User Guide

E.6.d Imaging Area Size

The imaging area is the area of the sample that is captured by the camera and displayed in the scan window. To specify the size of this area, enter a dimension in the appropriate field under Imaging Area.

When you change one imaging area dimension, the other will change to maintain the aspect ratio of the camera lens.

The imaging area will change depending on your zoom factor. For example, if you have zoomed in on a area that is 4.5 x 3.5 cm, then you would enter 4.5 for the width (3.5 for the height would be calculated automatically).

Note: Your imaging area settings must be correct if you want to do 1:1 printing. They must also be correct if you want to compare the quantities of objects (e.g., using the Volume Tools) in different images.

The imaging area dimensions also determine the size of the pixels in your image (i.e., resolution). A smaller imaging area will result in a higher resolution.

E.7 Other Features

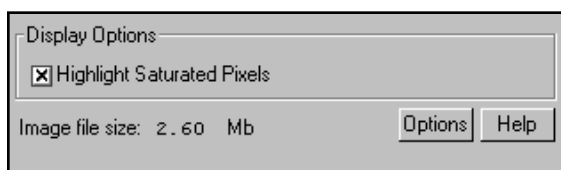


Fig.E-11. Other Fluor-S acquisition window features.

Highlight Saturated Pixels

When this box is checked, any saturated pixels in the image will appear highlighted in red in the scan window and in the pop-up image window. To view /hide saturated pixels in the pop-up image window, use the Image > Transform command.

Appendix E. Fluor-S

File Size of Images

Image File Size shows the size of the image file you are about to create. This size is determined by whether the image was created in High Resolution or High Sensitivity mode.

If you do not have enough computer memory for the specified file size, an error message will appear when you attempt to acquire an image. (Macintosh users can increase the application memory partition. See your Macintosh computer documentation for guidance.)

PDQUEST User Guide

Appendix F

Fluor-S MAX Multimager



Fig.F-1. Fluor-S MAX Multimager.

Before you can begin acquiring images using the Fluor-S[®] MAX MultiImager, the imaging system must be properly installed and connected with the host computer. See the Fluor-S MAX hardware manual for installation, startup, and operating instructions.

Note: Make sure that the temperature light on the Fluor-S MAX is green before attempting to capture an image. If you are using a PC, the Fluor-S MAX should be turned on and the initialization sequence completed *before* the host computer is turned on. See the hardware manual for more details.

PDQUEST User Guide

PC Only: A Note About SCSI Cards

The Fluor-S MAX is connected to your computer by a Small Computer System Interface (SCSI) cable. To use the Fluor-S MAX, you must have a SCSI card installed in your PC. If you have an older PC, you may also need to load the SCSI and WinASPI drivers that came with your card.

Simulation Mode

Any of the imaging device acquisition windows can be opened in “simulation mode.” In this mode, an acquisition window will open and the controls will appear active, but instead of capturing real images, the window will create “dummy” images of manufactured data.

You do not need to be connected to an imaging device to open a simulated acquisition window. This is useful for demonstration purposes or practice scans.

To enter simulation mode, hold down the CTRL key and select the name of the device from the File menu. The title of the acquisition window will indicate that it is simulated.

F.1 Fluor-S MAX Acquisition Window

To acquire images using the Fluor-S MAX, go to the File menu and select Fluor-S MAX.... The acquisition window for the imager will open, displaying a control panel and an image display window.

Appendix F. Fluor-S MAX

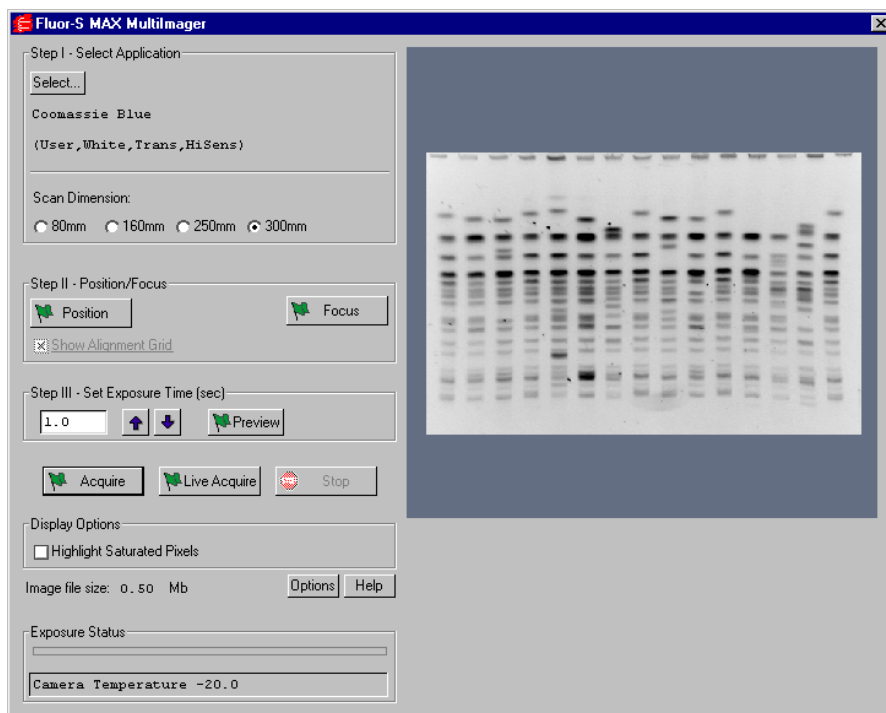


Fig.F-2. Fluor-S MAX acquisition window.

When the Fluor-S MAX window first opens, no image will be displayed.

The control panel has been arranged from top to bottom to guide you through the acquisition procedure. There are four basic steps to acquiring an image using the Fluor-S MAX:

1. Select the application.
2. Position and focus the gel or other object to be imaged.
3. Set the exposure time.
4. Acquire the image.

F.2 Step I. Select Application

To set the appropriate filter and other parameters for the type of object you are imaging, click on the Select button under Select Application. The available applications and their associated settings are listed in a tree that expands from left to right.

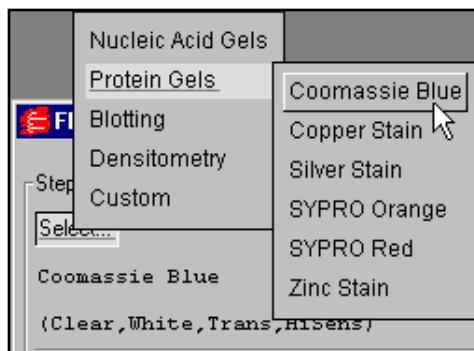


Fig.F-3. The application tree in the Fluor-S MAX acquisition window.

First select your general application, then select the particular stain or medium you are using. If you select the chemiluminescent application under “Blotting,” you also have the option of selecting High Sensitivity or Ultra Sensitivity (see below).

When you select an application, the software automatically sets the appropriate standard filter (520LP, 530DF60, 610LP, clear, or none), light type (UV, white, or none), and light source (trans, epi, or neither) in the Fluor-S MAX for that particular application.

For applications involving trans illumination, you must also specify a scan dimension (see below).

Your selection will be displayed below the Select button. To exit the tree without selecting, press the ESC key.

Appendix F. Fluor-S MAX

Custom Applications

If your application is not listed, if you want to use a user-installed filter, or if you want to access Ultra Sensitivity mode (see below), you can create and save your own custom application.

From the application tree, select Custom, then Create. This will open a dialog box in which you can name your application and select your settings.

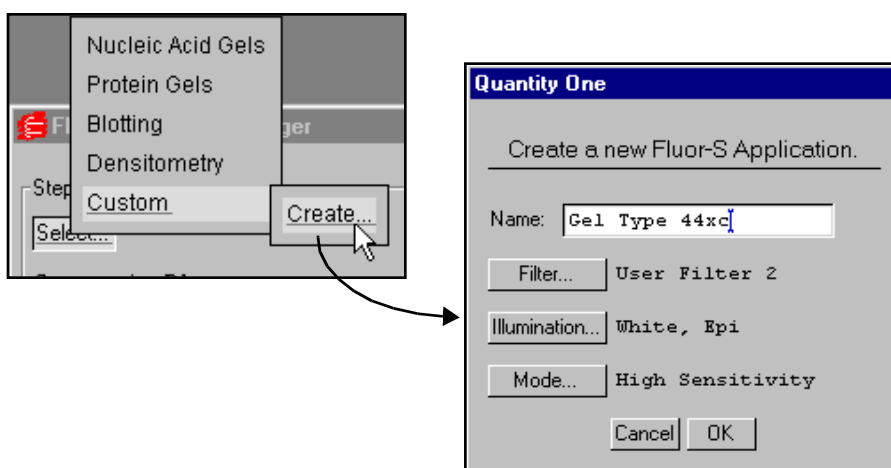


Fig.F-4. Creating a new custom application.

To select the filter (including user-defined), type of illumination, and camera mode, click on the appropriate buttons.

Note: Under Illumination, there is a listing for a spare UV light source. This selects the spare UV bulb in the Fluor-S MAX. Select this light source if your main UV bulb fails.

Enter a name for your application in the Name field. Click on OK to implement your changes.

After you have created an application, you can select it from the application tree by selecting Custom and the name you created. You can delete the application by selecting Custom, Delete, and the name of the application.

PDQUEST User Guide

Scan Dimension

If an application uses trans illumination, the Scan Dimension buttons become active. The scan dimension is the distance traveled by the transilluminating light source as it scans horizontally across the platen.



Fig.F-5. Scan dimension settings for trans illumination.

The full scanning range is 300 mm. Select a smaller range if your sample is small and you do not want to wait while the light source travels over the maximum scan width.

High Sensitivity and Ultra Sensitivity

High Sensitivity and Ultra Sensitivity are different camera modes. In Ultra Sensitivity mode, the Fluor-S MAX camera is cooled to -35.0 degrees C. In High Sensitivity mode, the camera is cooled to -20.0 degrees C. The current camera temperature is displayed at the bottom of the acquisition window.

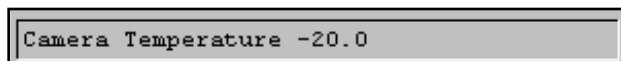


Fig.F-6. Camera Temperature display.

High Sensitivity is the normal operating mode of the Fluor-S MAX. Ultra Sensitivity provides optimal sensitivity for low-light applications. It is the default selection for the chemiluminescence application, and may be selected for a custom application.

Note: When you change from High to Ultra Sensitivity or visa versa, there will be a delay of several minutes while the Fluor-S MAX camera cools down or warms up. If you attempt to acquire an image during this period, you will be notified of the changing temperature. If you do not want to wait, you can cancel the mode change.

F.3 Step II. Position/Focus

Note: When you click on the Position or Focus button, the light inside the Fluor-S MAX box automatically turns off. This is because if you focus with the camera at maximum aperture (see note below), leaving the light on would make it difficult to view the image. To turn the light on while positioning or focusing, hold down the SHIFT key when clicking on the button.

Position

The next step in acquiring an image is centering your gel or other object within the camera frame. To do this, click on the Position button. The Fluor-S MAX will begin capturing a “live” image and updating it every second.

With the Position button selected, look at the image in the acquisition window while you position your object in the center of the platen. If you have a zoom lens on the camera, you can adjust the magnification while you position.

You can select the Show Alignment Grid checkbox to facilitate positioning.

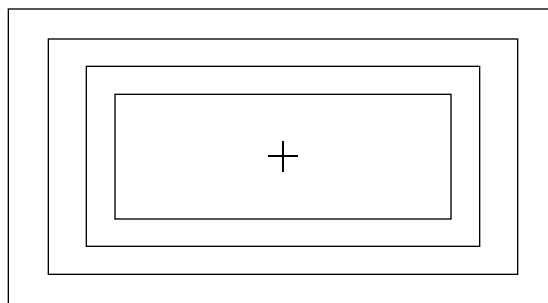


Fig.F-7. Fluor-S MAX alignment grid.

When you are finished positioning, click on the Stop button.

PDQUEST User Guide

Focus

Note: Before focusing, you should adjust the f-stop on the camera to the lowest setting (i.e., the maximum aperture). This reduces the depth of field, allowing you to more accurately focus the camera. Then, after focusing, increase the f-stop to the desired setting.

After you have positioned your sample, click on the Focus button and look at the image in the acquisition window while aligning the two focusing arrows on the camera lens. While focusing, the camera will limit its focus to a small portion of the sample (this will not affect any zoom lens adjustments you may have made.)

When you are finished focusing, click on the Stop button.

F.4 Step III. Set Exposure Time

When you are ready to capture an image, you will need to select an exposure time. “Exposure” refers to the integration of image captures on the CCD over a set period of time. The effect is analogous to exposing photographic film to light.

The exposure time you select should be based on your application and your “best guess” as to what exposure will give you the best image.

Note: The minimum exposure time in trans illuminated mode is 1 second. The minimum exposure time in epi illuminated mode is 0.1 second.

You can enter an exposure time (in seconds) directly in the field, or use the Arrow buttons to adjust the exposure time in 10 percent increments.



Fig.F-8. Selecting an exposure time.

Appendix F. Fluor-S MAX

The following table provides recommended exposure times for various applications.

Recommended Exposure Times and Lenses

Sample	Recommended Exposure	Lens & Filter ¹	Accessories Used
Fluorescent Stain Gel	1–20 sec.	Zoom/IR	None
Fluorescence End-Label Gel	10 sec.–2 min.	Zoom/IR	None
Fluorescent Blot	0.1–3 sec.	Zoom/IR	None
Chemifluorescent Blot	0.1–3 sec.	Zoom/IR	None
Colorimetric Gel	1–5 sec.	Zoom/IR	White Diffusion Plate
Colorimetric Blot	0.2–10 sec.	Zoom/IR	None
X-ray film	0.1–5 sec.	Zoom/IR	White Diffusion Plate
Weak Chemiluminescence ²	2–10 min.	50 mm	Chemi Tray (if sample is small)
Strong Chemiluminescence ²	5 sec.–1 min.	50 mm	Chemi Tray (if sample is small)
¹ For sharper focusing, close the f-stop down 1–2 stops from full open while focusing. ² For chemi applications, the 50mm lens is recommended. Always remove the 660 filter.			

For most applications, you can select an exposure time, capture an image, study it, then adjust the exposure time accordingly. Repeat this procedure as many times as necessary to obtain a good image.

For chemiluminescent samples, which degrade over time and emit low levels of light, you can select a high exposure time initially and use Live Acquire mode to save intermediate exposures (see following section).

PDQUEST User Guide

Preview

For shorter exposures, you can use Preview to test different exposure times. Click on the Preview button create a preview exposure and display it in the acquisition window.

A preview scan takes only half as long to create as a real scan, because the preview scan does not capture a “dark” image (see below). The progress of the exposure will be displayed in the Exposure Status bar at the bottom of the dialog box.

You cannot save preview scans.

If you want to stop a preview scan that is in progress, click on the Stop button.

F.5 Acquire the Image

The Fluor-S MAX gives you the option of simply acquiring and displaying a fully exposed image, or preserving intermediate exposures.

Acquire

To acquire and display a fully exposed image, click on the Acquire button. An exposure will be taken based on the time selected in Step III. This is appropriate for most short exposures.

The progress of each exposure will be displayed in the Exposure Status bar at the bottom of the acquisition window.

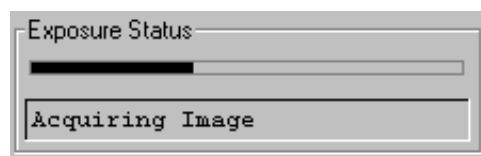


Fig.F-9. Exposure Status bar when acquiring an image.

Appendix F. Fluor-S MAX

Depending on which dark subtraction type you have selected, a dark count may be acquired immediately following image acquisition. See Dark Subtraction Type under Options, below.

If you want to stop a scan that is in progress, click on the Stop button. The acquisition will be terminated.

After an image has been acquired, a separate window will pop up containing the new image. The window will have a default file name that includes the date, time, and user (if known). To save the image, select Save or Save As from the File menu.

You can then analyze the image using the standard analysis functions.

Live Acquire

Live Acquire allows you to view and preserve intermediate exposures leading up to a full exposure. This is useful for longer exposure times, such as those required by chemiluminescent samples, where there is the potential for image saturation.

When you click on the Live Acquire button, the exposure time you selected is divided by the number of exposure counts set in the Options dialog box (see Options, below). For example, if you enter an Exposure Time of 10 minutes and an Exposure Count of 20, then 20 intermediate exposures will be produced at 30-second intervals. Each intermediate exposure will be displayed in the scan window. The final, full exposure will be displayed in a separate image window.

Note: The first intermediate exposure will take longer than 30 seconds to display if dark subtraction is performed.

You can automatically save your intermediate exposures as separate files for later review using the Auto Save After Scan option. See Options, below.

If you see an intermediate exposure that you like in the scan window, click on the Stop button. Live Acquire will end and the last intermediate exposure to be completed will open in a separate image window. You can then save it for analysis.

PDQUEST User Guide

F.6 Options

Click on the Options button to open the Options dialog box.

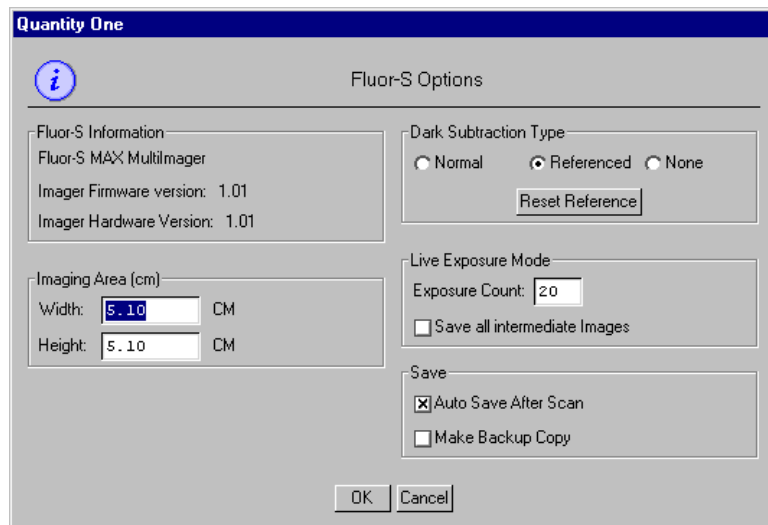


Fig.F-10. Options dialog box.

F.6.a Dark Subtraction Type

All CCD cameras accumulate electrons that produce a “signal” that is indistinguishable from light. This “dark count” adds to the noise in your images. In most cases, you will want to subtract this dark count from your images.

Normal

The Normal option button selects the default dark subtraction type. In this mode, after you acquire an image, a “dark” image of the same exposure length will be taken, and this will be subtracted from your image.

Appendix F. Fluor-S MAX

The progress of the dark exposure will be displayed in the Exposure Status bar following the regular image exposure.

Note: In Normal mode, a dark image is only acquired the first time you perform a scan with particular application and exposure settings. If you perform subsequent scans with the same settings, no dark exposure will be taken.

Referenced

If you do not want to perform a dark exposure with each acquisition, you can take a “reference” dark exposure that will be saved and subtracted from all subsequent acquisitions. Click on the Referenced button to activate this feature.

The first time you acquire an image after selecting this option, the Fluor-S MAX will take a 300-second dark exposure that will be saved and used to subtract the dark count from all subsequent acquisitions.

For image exposures of greater or less than 300 seconds, the reference dark will be scaled accordingly and then subtracted. You can change the default reference dark exposure time using the Reset Reference button (see below).

If you deselect the Referenced button and then reselect it, the old reference dark exposure will still be available.

Note: Separate reference dark exposures will be taken for High Sensitivity mode and Ultra Sensitivity mode. Once you have created a reference dark in each of these modes, each reference dark will be used according to the mode you are in.

Reset Reference Button

If you would like a reference dark with an exposure time that more closely matches that of your typical scans, click on the Reset Reference button.

PDQUEST User Guide



Fig.F-11. Reset Reference Dark pop-up box.

A pop-up box will prompt you to enter a new reference dark exposure time in seconds. Click on OK to implement your change. The new reference dark will be created when you acquire your next image.

Note: Because of the high sensitivity of the CCD, fluctuations in background radiation and / or temperature in the room can affect the level of dark count. If you feel that radiation / temperature conditions have changed in the room since your last reference dark was created, use the Reset Reference button to delete your old reference and create a new one under current conditions.

None

If you do not want to perform dark subtraction, select None. No dark exposure will be acquired or subtracted.

F.6.b Live Exposure Mode

Exposure Count

If you are using the Live Acquire function (see previous section), you need to specify how many intermediate exposures you want to view / save during acquisition. Enter this number in the Exposure Count field.

Appendix F. Fluor-S MAX

The total exposure time will be divided by the number you enter in the Exposure Count field. If you enter an exposure time of 10 minutes and a count of 10, you will create 10 intermediate exposures at 1 minute intervals.

Note: Do not enter a count that will result in an intermediate exposure time that is less than the minimum exposure time for the mode you are in. The minimum exposure time in trans illuminated mode is 1 second, and the minimum exposure time in epi illuminated mode is 0.1 second. (Example: For a trans illuminated application, an exposure time of 20 seconds and an exposure count of 21 would result in an error.)

Save All Intermediate Images

If Auto Save After Scan is selected (see following section), the Save All Intermediate Images checkbox will become active. If you select this checkbox, all your intermediate exposures will be saved as separate files. These files will have the same root name appended by a number indicating the exposure sequence. The final, full exposure will have the root name only, with no exposure number.

F.6.c Save

Auto Save After Scan

To automatically save any image you create, click on the Auto Save After Scan checkbox.

With this checkbox selected, when you click on Acquire or Live Acquire, a Save As dialog box will open asking you to specify a file name and location for the image you are about to create. The scan will begin when you click on the Save button.

Note that in Live Acquire mode you can save your intermediate exposures by selecting Auto Save After Scan and then Save All Intermediate Images.

PDQUEST User Guide

Make Backup Copy

You can automatically create a backup copy of any scan you create. To do so, first select Auto Save After Scan (see above), then select the Make Backup Copy checkbox.

With this checkbox selected, when you save a scan, a backup copy will be placed in the same directory as the scanned image. Windows backup files will have an “.sbk” extension. Macintosh backup files will have the word “backup” after the file name.

This backup copy will be read-only, which means that you cannot make changes to it. You can open it like a normal file, but you must save it under a different file name before editing the image or performing analysis.

F.6.d Imaging Area Size

The imaging area is the area of the sample that is captured by the camera and displayed in the scan window. To specify the size of your imaging area, enter a dimension in the appropriate field. When you change one imaging area dimension, the other will change to maintain the aspect ratio of the camera lens.

The imaging area will change depending on your zoom factor. For example, if you have zoomed in on a area that is 4.5 x 3.5 cm, then you would enter 4.5 for the width (3.5 for the height would be calculated automatically).

Note: Your imaging area settings must be correct if you want to do actual-size printing. They must also be correct if you want to compare the quantities of objects (e.g., using the Volume Tools) in different images.

The imaging area dimensions also determine the size of the pixels in your image (i.e., resolution). A smaller imaging area will result in a higher resolution.

F.7 Other Features

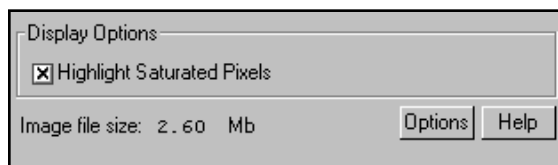


Fig.F-12. Other Fluor-S MAX acquisition window features.

Highlight Saturated Pixels

When this box is checked, any saturated pixels in the image will appear highlighted in red in the scan window and in the pop-up image window. To view / hide saturated pixels in the pop-up image window, use the Image > Transform command.

File Size of Images

Image File Size (below Options) shows the size of the image file you are about to create. This size is determined by whether the image was created in High Sensitivity or Ultra Sensitivity mode.

If you do not have enough computer memory for the specified file size, an error message will appear when you attempt to acquire an image. (Macintosh users can increase the application memory partition. See your Macintosh computer documentation for guidance.)

PDQUEST User Guide

Appendix G

Personal Molecular Imager FX

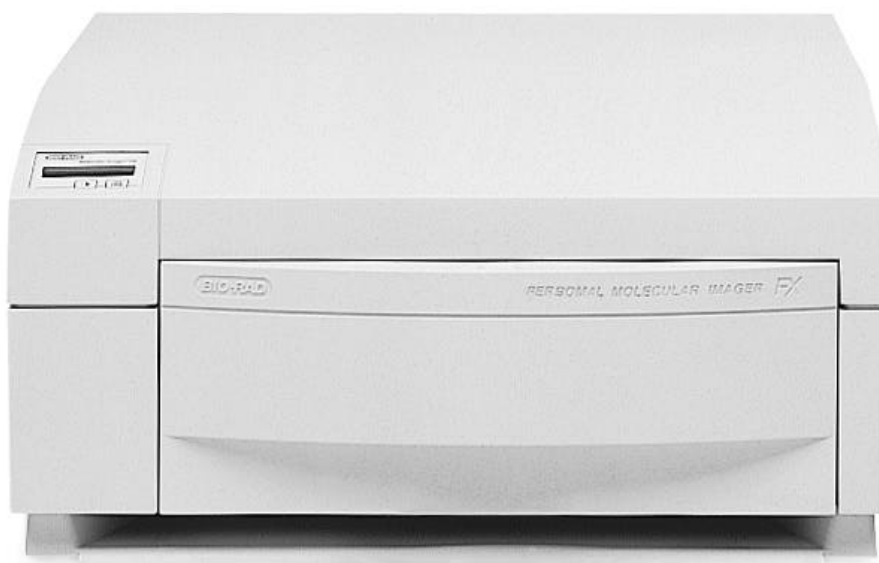


Fig.G-1. Personal Molecular Imager FX

Before you can begin acquiring images using the Personal Molecular Imager[®] FX, the instrument must be properly installed and connected with the host computer. See the Personal FX hardware manual for installation, startup, and operating instructions.

Note: The Personal FX should be turned on and the initialization sequence completed *before* the host computer is turned on (except in the case of certain Power Macintosh configurations). See the hardware manual for more details.

PDQUEST User Guide

PC Only: A Note About SCSI Cards

The Personal FX is connected to your computer by a Small Computer System Interface (SCSI) cable. To use the Personal FX, you must have a SCSI card installed in your PC. If you have an older PC, you may also need to load the SCSI and WinASPI drivers that came with your card.

Simulation Mode

Any of the imaging device acquisition windows can be opened in “simulation mode.” In this mode, an acquisition window will open and the controls will appear active, but instead of capturing real images, the window will create “dummy” images of manufactured data.

You do not need to be connected to an imaging device to open a simulated acquisition window. This is useful for demonstration purposes or practice scans.

To enter simulation mode, hold down the CTRL key and select the name of the device from the File menu. The title of the acquisition window will indicate that it is simulated.

G.1 Personal FX Acquisition Window

To acquire images using the Personal FX, go to the File menu and select Personal FX.... The acquisition window for the imager will open, displaying a control panel and the scanning area window.

Appendix G. Personal FX

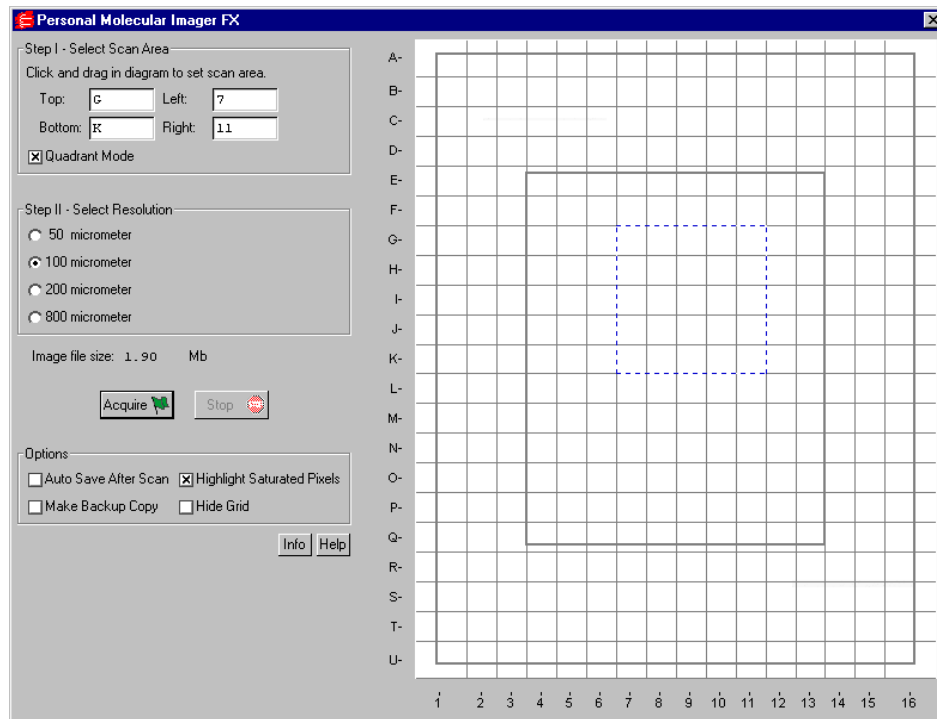


Fig.G-2. Personal FX acquisition window

The default scanning window is marked by grid lines that divide the area into quadrants. There is also an outer box and inner box marked by thicker lines. This conforms to the sample pad for the standard Bio-Rad Exposure Cassette that is supplied with the Personal FX. The quadrants are numbered 1–16 left to right and lettered A–U top to bottom.

If you prefer a scanning window measured in centimeters, deselect the Quadrant Mode checkbox in the control panel by clicking on it. To hide the gridlines, click on the Hide Grid checkbox under Options.

The control panel has been arranged from top to bottom to guide you through the acquisition procedure. There are three basic steps to scanning an image using the Personal FX:

PDQUEST User Guide

1. Select the scan area
2. Select the resolution
3. Acquire the image

G.2 Step I. Select Scan Area

To select a scan area, drag your mouse within the scanning window. (In the scanning window, your cursor appearance will change to a cross.) The border of the scan area you are selecting is marked by a frame.

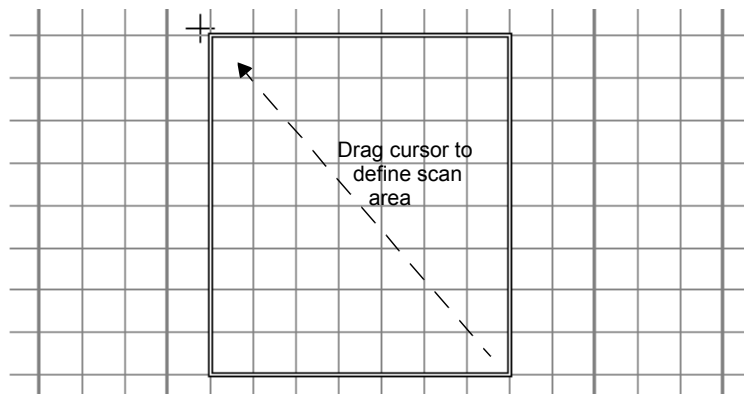


Fig.G-3. Selecting a scan area.

If you are in quadrant mode, note that the frame “locks” onto the next quadrant as you drag. When you release the mouse button, the border changes to a dashed blue line, indicating a selected area.

- To *reposition* the scanning box you have selected, position your cursor inside the box and drag. The entire box will move.
- To *resize* the box, position your cursor on a box side and drag. The side you have selected will move.
- To *redo* the box entirely, position your cursor outside the box and drag. The old box will disappear and a new box will be created.

Appendix G. Personal FX

You can also select the scanning area by entering coordinates in the appropriate fields (Top, Bottom, Left, Right). After you enter a coordinate, the position of the scanning area box will change accordingly.

When selecting, be sure to include the entire area of interest, and be generous with borders. You can always crop the image later.

G.3 Step II. Select Resolution

The Personal FX acquisition window allows you to scan at 50, 100, 200, or 800 micrometers. These resolutions are listed as option buttons in the control panel.

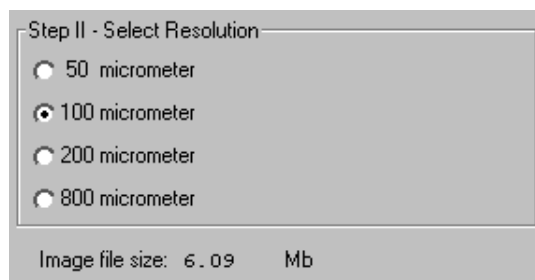


Fig.G-4. Resolution option buttons.

The resolution you select should be based on the size of the objects (e.g., bands, spots) you are interested in. For example:

- 50 micrometer resolution should be reserved for images requiring the highest level of detail, e.g., high density in situ samples, 1,536-well microplates, high density arrays, samples with very closely spaced bands. Files scanned at 50 micrometers can be very large.
- 100 micrometer resolution should be used for typical gels and arrays.
- 200 micrometer resolution is useful for gels with large bands and dot blots.
- 800 micrometer resolution should be reserved for very large objects, such as CAT assays.

PDQUEST User Guide

File Size of Images

Image File Size (below Select Resolution) shows the size of the scan file you are about to create. If you do not have enough computer memory for the specified file size, an error message will appear when you attempt to scan. If this happens, select a lower resolution or decrease the size of the area to be scanned. (Macintosh users can also increase the application memory partition. See your Macintosh computer documentation for guidance.)

G.4 Acquire the Image

Once you have selected your scan area and resolution, you are ready to acquire an image.

Click on the Acquire button. There may be a short delay while the image laser warms up; then the scanned image will begin to appear in the scanning window, line by line.

To interrupt a scan, click on the Stop button. A message will ask you to confirm the interrupt, and then you will be asked if you want to keep the partial scan. This feature is useful if you overestimated the size of the area you selected.

Note: If the image you are scanning has more than 10 saturated pixels, you will receive a warning message.

Saving the Image

After the scan is complete, a message will appear asking you if you want to keep the scan. If you select Yes, a separate window will pop up containing the new image.

You can then save and analyze the image using the standard menu and toolbar functions.

G.5 Options

Auto Save After Scan

To automatically save any scan you create, click on the Auto Save After Scan checkbox.

With this checkbox selected, when you click on Acquire, a Save As dialog box will open asking you to specify a file name and location for the image you are about to create. The scan will begin when you click on the Save button.

Make Backup Copy

You can automatically create a backup copy of any scan you create. To do so, first select Auto Save After Scan (see above), then select the Make Backup Copy checkbox.

With this checkbox selected, when you save a scan, a backup copy will be placed in the same directory as the scanned image. Windows backup files will have an “.sbk” extension. Macintosh backup files will have the word “backup” after the file name.

This backup copy will be read-only, which means that you cannot make changes to it. You can open it like a normal file, but you must save it under a different file name before editing the image or performing analysis.

Highlight Saturated Pixels

When this box is checked, any saturated pixels in the image will appear highlighted in red in the scan window and in the pop-up image window. To view / hide saturated pixels in the pop-up image window, use the Image > Transform command.

Hide Grid

To hide the gridlines in the scanning area window, click on the Hide Grid checkbox.

PDQUEST User Guide

Appendix H

Molecular Imager FX

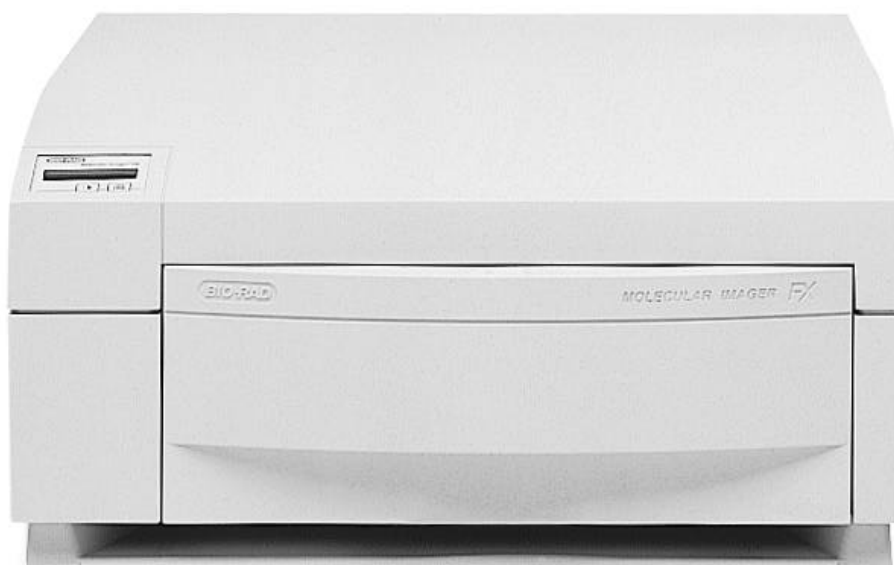


Fig.H-1. Molecular Imager FX

Before you can begin acquiring images using the Molecular Imager[®] FX, the instrument must be properly installed and connected with the host computer. See the FX hardware manual for installation, startup, and operating instructions.

Note: The FX should be turned on and the initialization sequence completed *before* the host computer is turned on (except in the case of certain Power Macintosh configurations). See the hardware manual for more details.

PDQUEST User Guide

PC Only: A Note About SCSI Cards

The FX is connected to your computer by a Small Computer System Interface (SCSI) cable. To use the FX, you must have a SCSI card installed in your PC. If you have an older PC, you may also need to load the SCSI and WinASPI drivers that came with your card.

Simulation Mode

Any of the imaging device acquisition windows can be opened in “simulation mode.” In this mode, an acquisition window will open and the controls will appear active, but instead of capturing real images, the window will create “dummy” images of manufactured data.

You do not need to be connected to an imaging device to open a simulated acquisition window. This is useful for demonstration purposes or practice scans.

To enter simulation mode, hold down the CTRL key and select the name of the device from the File menu. The title of the acquisition window will indicate that it is simulated.

H.1 FX Acquisition Window

To acquire images using the FX, go to the File menu and select FX.... The acquisition window for the imager will open, displaying the control panel for the imager and the scanning area window.

Appendix H. Molecular Imager FX

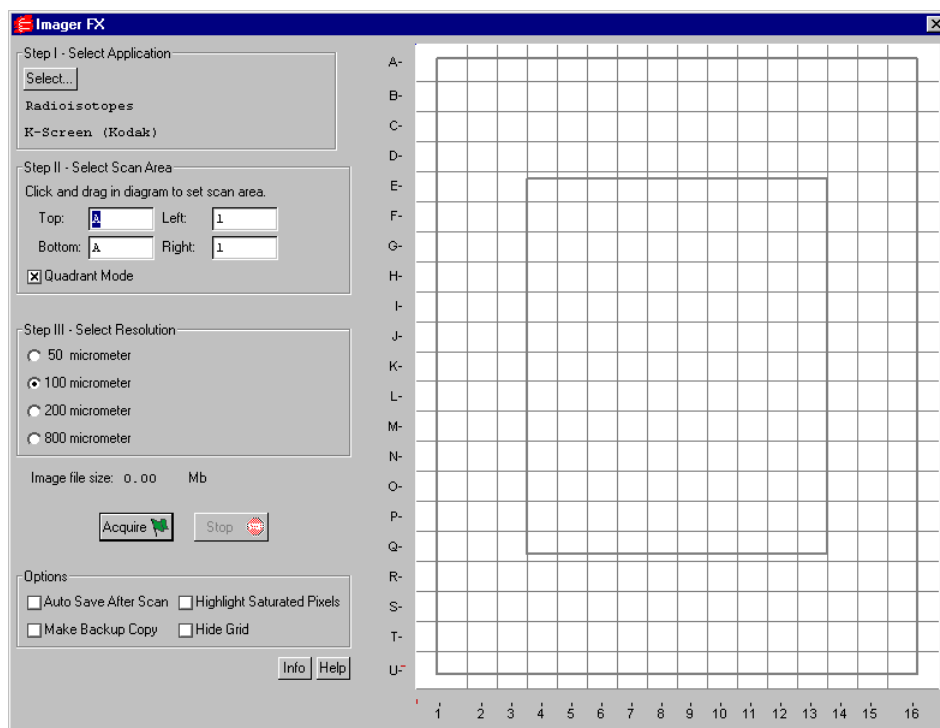


Fig.H-2. FX acquisition window

The default scanning window is marked by grid lines that divide the area into quadrants. There is also an outer box and inner box marked by thicker lines. This conforms to the sample pad for the standard Bio-Rad Exposure Cassette that is supplied with the FX. The quadrants are numbered 1–16 left to right and lettered A–U top to bottom.

If you prefer a scanning window measured in centimeters, deselect the Quadrant Mode checkbox in the control panel by clicking on it. To hide the gridlines, click on the Hide Grid checkbox under Options.

The control panel has been arranged from top to bottom to guide you through the acquisition procedure. There are four basic steps to scanning an image using the FX:

PDQUEST User Guide

1. Select the application.
2. Select the scan area.
3. Select the resolution.
4. Acquire the image.

H.2 Step I. Select Application

To select the appropriate filter and other parameters for the type of object you are imaging, click on the Select button under Select Application.

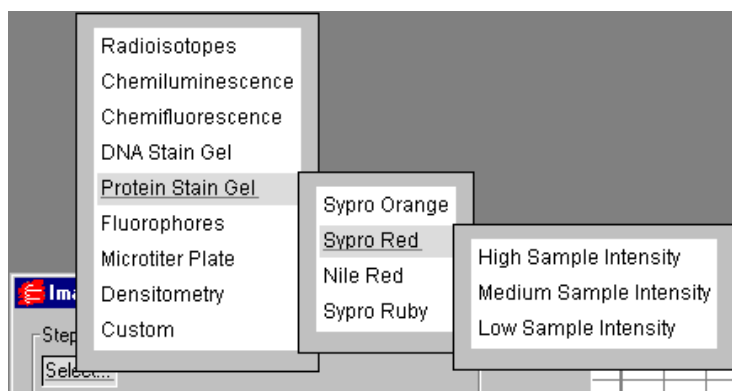


Fig.8-3. Example of an application tree: Ethidium Bromide gel.

Standard Applications

The standard applications and associated settings are listed in a tree that expands from left to right. When you select a standard application, the software automatically sets the appropriate filter(s) in the FX for that particular application.

Appendix H. Molecular Imager FX

Standard FX Applications

Category	Application
Radioisotopes	CS- or BI-Screen (Bio-Rad) K-Screen (Kodak) Fuji-Screen
Chemiluminescence	Chemi-Screen (Bio-Rad)
Chemifluorescence	ECL-Plus Attophos
DNA Stain Gel	Ethidium Bromide Sybr Green I & II Sybr Gold
Protein Stain Gel	Sypro Orange Sypro Red Nile Red Sypro Ruby
Fluorophores	Alexa 488 Alexa 532 Alexa 546 FITC FAM CY3 HEX R6G Texas Red
Microtiter Plate	DNA (Sybr Green I) Protein (Nano Orange) ssDNA (Oligreen) DNA (Picogreen) B-Gal (FDG) GUS (FDG)
Densitometry	Coomassie Blue Gel/Blot Copper Stain Gel/Blot Silver Stain Gel/Blot X-Ray Film (Grey Type)

PDQUEST User Guide

First select your general application, next select the particular stain or medium you are using, and finally (if appropriate) select the intensity of your samples.

Note: Some applications require an external laser. If you choose one of these but do not have an external laser attached, you will receive a warning.

To exit the tree without selecting, press the ESC key.

Your selection will be displayed below the Select button.

Sample Intensity

Most FX applications require that you specify a sample intensity. This is simply a rough estimate of how much sample is visible in your gel or other object.

If you are unsure of the level of intensity of your sample, you can always select a level, capture an image, then adjust the level and capture another image.

For example, if you select Low Sample Intensity and the resulting image has too many saturated pixels, you will receive a warning message. Simply change the setting to Medium Sample Intensity and rescan. If you select High Sample Intensity and the resulting image is too faint, select Medium or Low and rescan.

Custom Applications

If your application is not listed, if you want to use user-installed filters, or if you want to use an external laser, you can create and save your own custom application.

From the application tree, select Custom, then Create. This will open a dialog box in which you can name your application and select your settings.

Appendix H. Molecular Imager FX

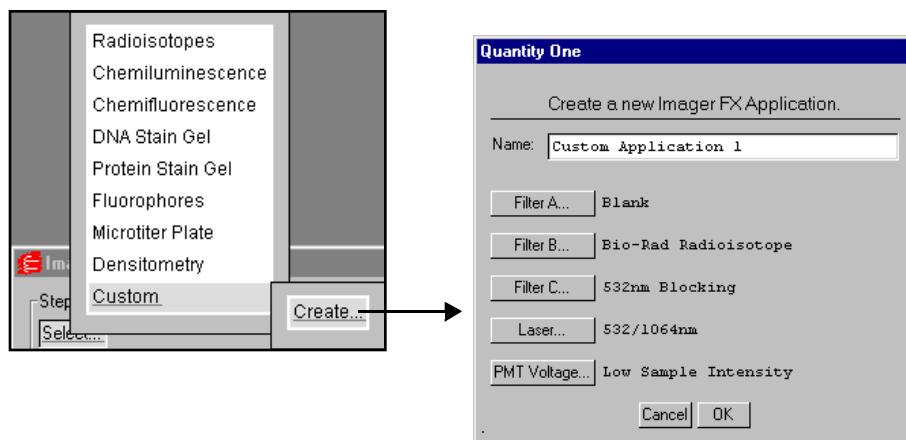


Fig.8-4. Creating a custom application.

To select a filter (including user-defined) or filter combination, click on the buttons for Filters A, B, and C, and make your choice from each pop-up list.

Note: The user-defined filters (User1, User2, etc.) cannot be renamed in the pop-up list, so be sure to remember which filter you insert into each position in the FX.

To use an external laser, click on the Laser button and select it from the pop-up list. Otherwise, use the default internal laser (532/1064nm).

Click on the PMT Voltage button to select a standard voltage for your custom application or create a custom PMT voltage.

To select a custom voltage, click on the Custom option. In the dialog box, adjust the slider to select a PMT voltage as a percentage of the maximum.

PDQUEST User Guide

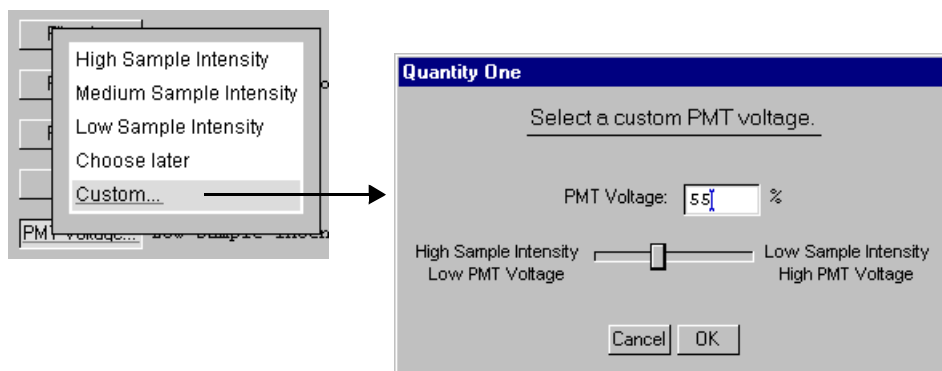


Fig.8-5. Selecting a custom PMT voltage.

Note: For voltages above 80% of maximum, you will receive a warning message that the high voltage could damage the PMT.

If you select Choose Later from the list of PMT voltages, the choices of sample intensity will be displayed when you select your custom application.

Finally, enter a name for your application in the Name field and click on OK to implement your changes.

After you have created an application, you can select it from the application tree by selecting Custom and the name you created. You can delete the application by selecting Custom, Delete, and the name of the application.

Appendix H. Molecular Imager FX

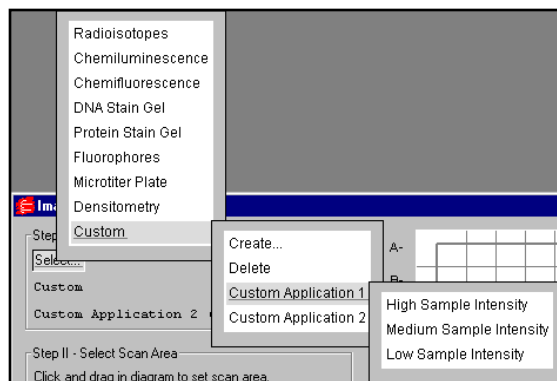


Fig.8-6. Selecting a custom application.

H.3 Step II. Select Scan Area

To select a scan area, drag your mouse within the scanning window. (In the scanning window, your cursor appearance will change to a cross.) The border of the scan area you are selecting is marked by a frame.

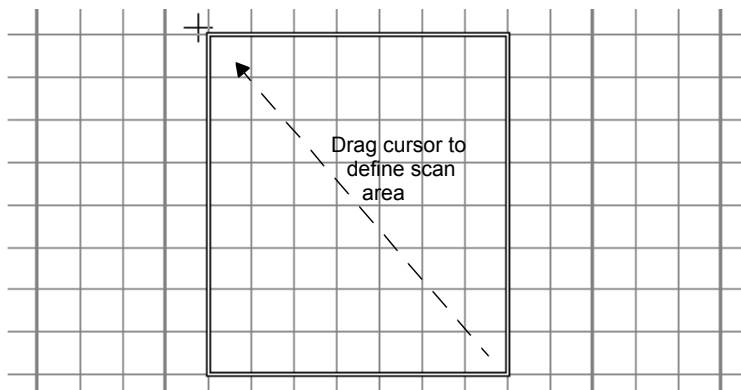


Fig.H-7. Selecting a scan area.

PDQUEST User Guide

If you are in quadrant mode, note that the frame “locks” onto the next quadrant as you drag. When you release the mouse button, the border changes to a dashed blue line, indicating a selected area.

- To *reposition* the scanning box you have selected, position your cursor inside the box and drag. The entire box will move.
- To *resize* the box, position your cursor on a box side and drag. The side you have selected will move.
- To *redo* the box entirely, position your cursor outside the box and drag. The old box will disappear and a new box will be created.

You can also select the scanning area by entering coordinates in the appropriate fields (Top, Bottom, Left, Right). After you enter a coordinate, the position of the scanning area box will change accordingly.

When selecting, be sure to include the entire area of interest, and be generous with borders. You can always crop the image later.

H.4 Step III. Select Resolution

The FX acquisition window allows you to scan at 50, 100, 200, or 800 micrometers. These resolutions are listed as option buttons in the control panel.

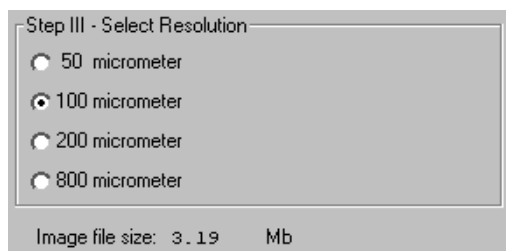


Fig.H-8. Resolution option buttons.

The resolution you select should be based on the size of the objects (e.g., bands, spots) you are interested in. For example:

Appendix H. Molecular Imager FX

- 50 micrometer resolution should be reserved for images requiring the highest level of detail, e.g., high density in situ samples, 1,536-well microplates, high density arrays, samples with very closely spaced bands. Files scanned at 50 micrometers can be very large.
- 100 micrometer resolution is useful for typical gels and arrays.
- 200 micrometer resolution is useful for gels with large bands and dot blots.
- 800 micrometer resolution should be reserved for very large objects, such as CAT assays.

File Size of Images

Image File Size (below Select Resolution) shows the size of the scan file you are about to create. If you do not have enough computer memory for the specified file size, an error message will appear when you attempt to scan. If this happens, select a lower resolution or decrease the size of the area to be scanned. (Macintosh users can also increase the application memory partition. See your Macintosh computer documentation for guidance.)

H.5 Acquire the Image

Once you have selected your application, scan area, and resolution, you are ready to acquire an image.

Click on the Acquire button. There may be a short delay while the image laser warms up; then the scanned image will begin to appear in the scanning window, line by line.

To interrupt a scan, click on the Stop button. A message will ask you to confirm the interrupt, and then you will be asked if you want to keep the partial scan. This feature is useful if you overestimated the size of the area you selected.

Note: If the image you are scanning has more than 10 saturated pixels, you will receive a warning message. If this happens, you can go back and select a higher sample intensity in the application tree.

PDQUEST User Guide

Saving the Image

After the scan is complete, a message will appear asking you if you want to keep the scan. If you select Yes, a separate window will pop up containing the new image.

You can then save and analyze the image using the standard menu and toolbar functions.

H.6 Options

Auto Save After Scan

To automatically save any scan you create, click on the Auto Save After Scan checkbox.

With this checkbox selected, when you click on Acquire, a Save As dialog box will open asking you to specify a file name and location for the image you are about to create. The scan will begin when you click on the Save button.

Make Backup Copy

You can automatically create a backup copy of any scan you create. To do so, first select Auto Save After Scan (see above), then select the Make Backup Copy checkbox.

With this checkbox selected, when you save a scan, a backup copy will be placed in the same directory as the scanned image. Windows backup files will have an “.sbk” extension. Macintosh backup files will have the word “backup” after the file name.

This backup copy will be read-only, which means that you cannot make changes to it. You can open it like a normal file, but you must save it under a different file name before editing the image or performing analysis.

Appendix H. Molecular Imager FX

Highlight Saturated Pixels

When this box is checked, any saturated pixels in the image will appear highlighted in red in the scan window and in the pop-up image window. To view/hide saturated pixels in the pop-up image window, use the Image > Transform command.

Hide Grid

To hide the gridlines in the scanning area window, click on the Hide Grid checkbox.

PDQUEST User Guide

Appendix I

Cross-Platform File Exchange

It is possible to move image data between Discovery Series software applications on different platforms. There are different protocols depending on the platform (PC or Macintosh) you are transferring from and to.

I.1 Macintosh to PC

To transfer a file from a Discovery Series application running on a Macintosh to a Discovery Series application running on a PC, you need to tag the file name with the suffix appropriate to the image file type (e.g., 2-D scan).

For gel scans, use the suffix: .gsc

For gel images, use the suffix: .gim

For gel spots files, use the suffix: .gsp

For calstrip scans, use the suffix: .csc

For matchsets, use the suffix: .ms

For higher-level matchsets, use the suffix: .hms

I.2 PC to Macintosh

PC versions of Discovery Series files can be read directly by Macintosh applications, with no required modifications.

PDQUEST User Guide

Appendix J

Techniques

The techniques described in this appendix apply to detecting spots and normalizing for quantity among matchset members.

J.1. Calibration

When using a scanner or imager, detection of spots can vary on the same 2-D scan depending on whether or not it is calibrated. The following is a discussion of calibration and spot detection's dependency on it.

J.1.a. Calibration

Calibration consists of the following steps: determine the table of values relating film density to counts, calculate a curve which maps every possible density to a count value, translate the 2-D scan using the calibration curve.

Conditions

Calibration is an attempt to compensate for the nonlinear response of the film to exposing radiation. The calstrip must undergo the same processing as its associated film since the film's response curve and background noise are dependent on the film processing. If processed correctly, X-ray film will not saturate until over 2 O.D., although the top of the linear range of its response is at approximately 2 O.D., which is why we recommend that as a ceiling.

Density patches on the calstrip must have uniform counts over their entire area.

Determine the table of values relating density to counts

Density is a value from 0–255 (8-bit binary number). When scanning, the user selects the O.D. value to which 255 corresponds. Typically this is 2 O.D. Each patch on the calibration strip is averaged within the specified box to obtain a

PDQUEST User Guide

density value for that patch. The counts are obtained from a database which the user must enter for each batch of calstrips.

Calculate a curve which maps density to counts

The lowest point of the curve relates background density to 0.0 counts. Background is obtained from the calstrip scan by finding the mode of density values which are below the mean of all density values.

A linear interpolation is performed between adjacent sample points to obtain values there. If the density value of any segment is within 1% of 255, it is considered saturated. The highest two nonsaturated samples are used to extrapolate the curve to saturation. This establishes the effective counts at saturation.

If the slope of a curve segment falls below 10% of the highest value ever seen, saturation is set at the sample point at the lower end of the offending segment. This is only necessary when the film response decreases below 2 O.D. This normally indicates a failure of film processing. The user should diagnose the problem in the film exposure/development process and reprocess the gel if possible. The linearity of the film response has been compromised.

Translate the 2-D scan

Every pixel in the 2-D scan file is mapped to a count value using the calibration table of interpolated and extrapolated values.

A scale factor is calculated which attempts to map 0–2 O.D. to the first 255 steps in the range. This is to accommodate the spot detection slope detector discussed below. This is obtained in the following way:

The average slope of the curve is obtained by a least squares regression of the nonsaturated calibration samples.

The scale factor is calculated from the following:

$$\text{hgt_to_den} = (255 * \text{avg_slope} * \text{max_OD}) / (2.0 * \text{pix_at_max_OD} - \text{bkgd_dens} * \text{max_OD}).$$

This scale factor is applied to the mapped values above.

Appendix J. Techniques

J.1.b. Spot detection slope detector

A slope threshold technique is used to recognize spot edges. Slope values are obtained as Y/X where Y is density (uncalibrated) or counts (calibrated) and X is distance between scanner pixels in 0.1 mm units.

When using uncalibrated data, the data are scaled so that 0–255 represents 0–2 O.D. to give uniform results when detecting different images.

When using calibrated data, the equivalent scaling is difficult to obtain due to the nonlinear response of the film. A linear regression through the density-counts sample points is used to approximate this relationship.

J.2. Spot Detection

Spot detection is a series of steps which results in a list of spots which have been located and quantified. Due to the imperfections in gel and imaging technology, it is impossible to construct a single series of steps which optimally performs this task on all images. Therefore, the facility has been added to PDQUEST to allow a user to specify what processing steps should be performed by autodetection.

An image is obtained from an imaging device. This image is referred to as a 2-D scan. If available, calibration information is applied to the 2-D scan or in the case of multiple exposures of the same gel, the 2-D scans. The 2-D scan or scans are then merged to produce the Gel Image. This image includes all information available from all the 2-D scans on a single gel, hence the name Gel Image. Image processing to remove streaks and background is performed on the Gel Image. Spot detection locates spots on the image and spot fitting quantifies them. The resulting quantified spots are stored in a Gel Spots file.

PDQUEST User Guide

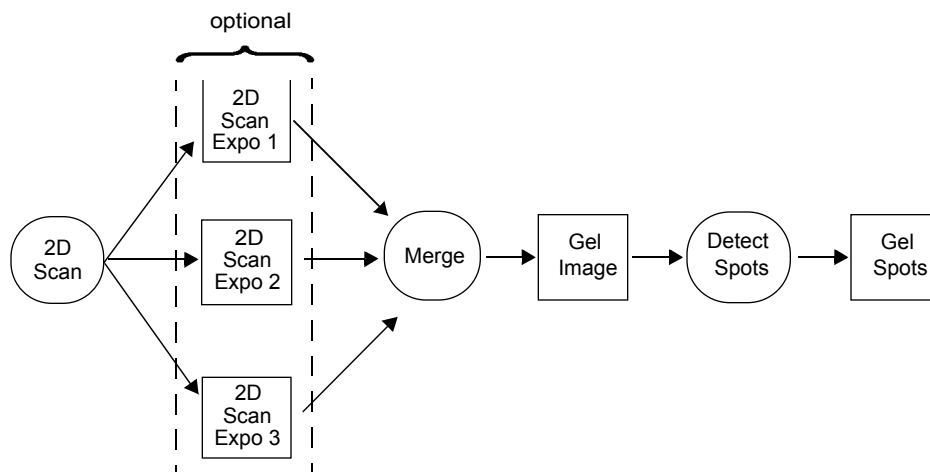


Fig. J-1. Auto-Detection procedure.

J.2.a. Auto-Detect Glossary

2-D scan —The image file produced by scanning.

Background —Film fog and unfocused protein on the image.

Calstrip —Either a calibration strip or the scan of a calibration strip

Gel Image —The image file produced by calibrating and merging scans.

Gel Spots —The list of spots from a single gel.

Radius —The radius of a disk used to determine background for the purpose of background subtraction or background normalization. Its size is a number of pixels in the image. Typical values are 50–150 for gels and over 800 for calstrips.

Spot Detection —The process that locates spots on the Gel Image.

Spot Fitting —The process that finds the optimal gaussian parameters in X and Y to represent the spots found by spot detection.

Streak —A protein that has focused only in one dimension.

J.3. Rolling Disk Streak Removal

The method we use for streak removal is based on a simple geometric model that is efficient enough to be performed in a relatively short time.

Consider the following diagram. The density trace represents O.D. values taken from a single vertical column of pixels in the image. (or several columns averaged together). A disk of radius ' r ' is used to obtain 'touch points' on the density trace. These 'touch points' are connected by straight line segments to obtain the background trace. The touch points are defined as those points on the trace where the topmost point of the circle can touch without any other part of the circle intersecting the density trace.

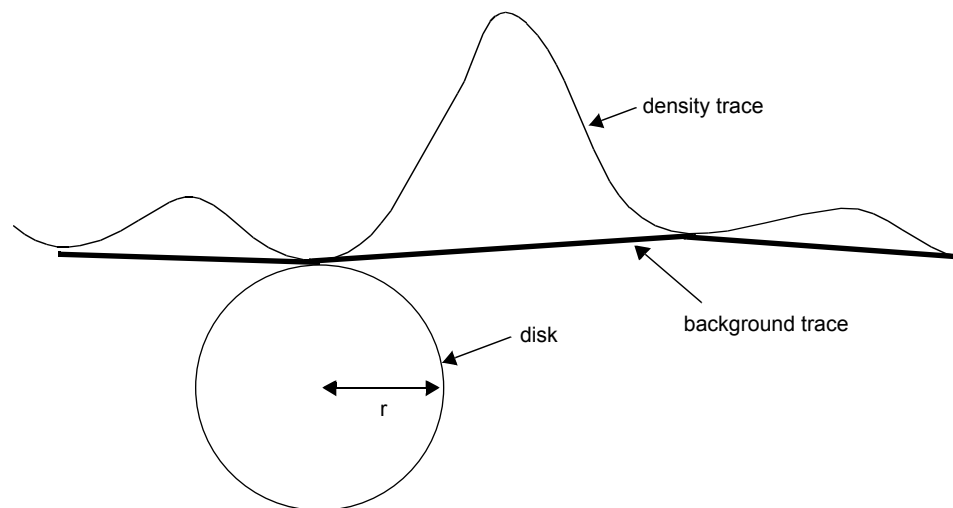


Fig. J-2. Rolling Disk background determination

Using a vertically-oriented disk results in vertical streak removal since the disk is able to contact the underside of the streak over its length.

PDQUEST User Guide

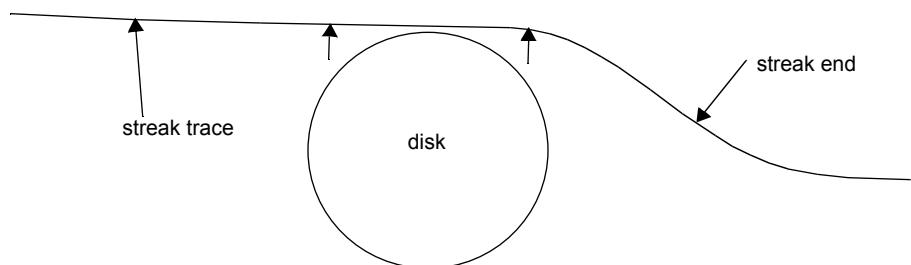


Fig. J-3. Disk contacting underside of streak trace

The background trace is nearly identical to the streak trace. There are some artifacts where the streak ends which usually do not introduce new spots into the Gel Spots file. The user must be aware that this may occur. If the streaks fade gradually the artifacts will not be present. If the streaks are saturated this often introduces pronounced artifacts around the saturated areas. The user must remove any spots detected in this region.

There is some disagreement about the validity of streak removal. Most researchers we have consulted share the opinion that the protein causing the streak and protein in spots embedded in the streak can co-exist in the gel and density measured at any point is the sum of the various contributing factors.

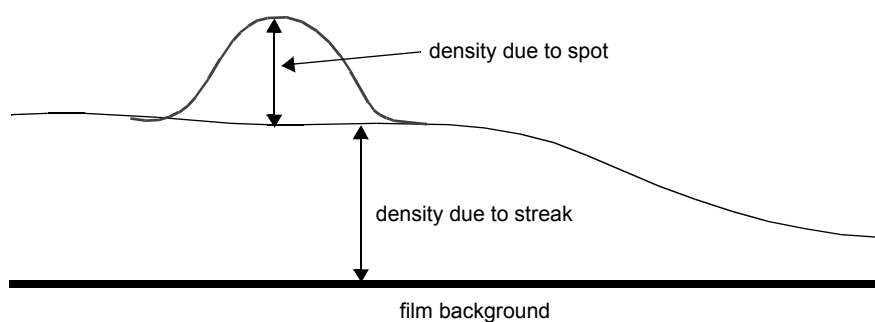


Fig. J-4. Summed streak and spot densities.

If you accept this model, then it is proper to remove the streak before quantitating the spot. This is the default method of spot detection.

Appendix J. Techniques

The contradictory view is that proteins in a gel are mutually exclusive: no two proteins can co-exist at the same location in the gel. The implication is that the density of the spot does not add to the density of the streak.

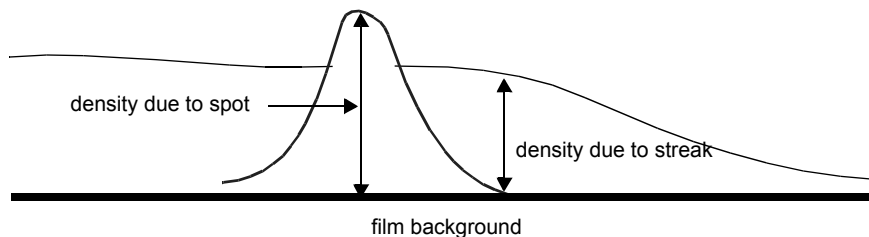


Fig. J-5. Contradictory view: spot and streak proteins do not co-exist.

If you accept this view, then change the auto-detect procedures to use Subtract Background (Ball) rather than Subtract Background (Vert. Disk).

Note that these considerations apply only to spots buried in streaks and does not affect the quantitation of most spots which are either isolated or overlap other spots.

J.4. Pipetting Error Compensation

Pipetting Error Compensation (in the Quantity Normalization dialog) can be used to correct systemic errors such as pipetting or counting errors which affect the quantitation of all spots in the gel equally. It results in a small correction factor which is applied to all spots when calculating quantity in PPM. When properly applied, it improves the accuracy of quantitative comparisons between gels.

However, it relies on two key assumptions: the log of spot quantity data is normally distributed, and changes to the spots are random.

In samples we have studied which are derived from whole cell lysates, the assumption of normal distribution has been true. The user must verify this for his own samples.

PDQUEST User Guide

Changes to the spots will not be random in experiments where overall protein synthesis is affected, e.g. the introduction of a toxic substance. In this case the Compensation will obscure the experimental result and apply a 'correction' factor which removes it.

The Compensation factors are obtained by calculating the average ratios of the log of the quantitation for all spots between all gels in a matchset. From the second assumption, the overall average ratio should be 1.0, indicating no net change. The correction factor for each gel is the reciprocal of its average change.

Index

A

A-B Comparison	10-7
Acid-Base	2-25
Add Spot To Set	9-11
Adding a spot to the standard	6-20
Advanced View	
All Windows	3-5
Advanced viewing tools	3-4
All Entire Images	3-6
All Same	3-6
Analysis sets	9-1
Arbitrary	9-10
Boolean	9-11
Creating	9-2
Displaying	9-4
Graphing	9-4
Highlighting spots in	9-4
Printing annotations for	8-10
Qualitative	9-8
Quantitative	9-5
Removing a spot from	9-11
Statistical	9-13
Annotations	8-5
Category, creating	8-6
Category, deleting or changing	8-8
Category, selecting	8-6
Category, standard	8-8
Creating	8-7
Deleting	8-7
Displaying	8-8
Editing	8-7

PDQUEST User Guide

Exporting	11-8
Hiding	8-9
Listing categories and annotations for a spot	8-9
Printing	8-9
Transferring	9-20
Arbitrary Analysis Sets	9-10
Arrow keys	3-3
Assigning images to windows	3-6
Authorization	
Unable to Obtain	1-15
Auto All Same	2-26

B

Base-Acid	2-25
Bio-Rad Technical Service	1-22
Boolean Analysis Sets	9-11

C

Calibrate Gels. See Calstrips.	
Calibration	
Overview	1-6
Calstrips	
Boxing segments	4-8
Calibration curve	4-9
Configuring in subwindows	4-2
Cropping out of images	4-1
Deactivating a segment	4-11
Entering data	4-2
Gel information	4-5
Range of values	4-4
Records	4-2
Removing segments	4-10

Index

Sampling segments	4-7
Categories of annotations	8-6
Centering an image	3-4
Change Version / Exposure	2-16
Chemi Doc	
Acquiring the image	B-5
Adjusting the display	B-11
Annotating images	B-8
Auto expose	B-5
Auto exposure settings	B-14
Backup images	B-14
Chemi mode	B-10
DAC settings	B-13
Display settings	B-10
Exposure Status bar	B-10
Exposure time	B-7
Freeze	B-7
Image mode	B-9
Imaging area	B-14
Invert display	B-11
Live image display	B-4
Manual expose	B-6
Options	B-12
Positioning the sample	B-4
Saturated pixels, highlighting	B-11
Saving images	B-9
Simulation mode	B-2
UV mode	B-9
Video card	B-2
Video display window	B-3
Video printing	B-8
Video printing, footer info	B-14
White light mode	B-10
Click to Focus Tile	2-25
Closing all files	2-15
Closing an image	2-15
Colony matching	6-21
Colors, setting	3-14
Column Graph	10-4

PDQUEST User Guide

Combined spots, viewing	2-25
Computer Requirements	1-10
Computer requirements	
Macintosh	1-11
PC	1-10
Configure Subwindows	3-12
Contrast, see Transform.	
Create	
Create Matchset	6-2
Cropping	3-23
Advanced crop commands	3-24
Drawing the crop box	3-24
Cursor, jump on alert	2-25
Custom Rotation	3-27

D

Database Node ID	2-24
Database spot numbers. See DSNs	
Delete	
Group	7-5
Density Tools	3-8
Density at Cursor	3-9
Density in Box	3-9
Plot Cross-section	3-9
Plot Density Distribution	3-9
Plot Vertical Trace	3-10
Display	
MrpI Data	7-3
Downloading from Internet	
Macintosh	1-13
Windows	1-12
DSNs	8-3
Assigning	8-4
Displaying	8-4
Transferring	9-19

E

Edit	
Replicate Group	7-5
Enter MrpI Data	7-2
Erase Spot Crosshair	5-23
Erratic spots	6-16
Exit	2-22
Exporting	
Annotations	11-8
Images	11-6
Matchset data	11-4
Exposures	
Merging	4-11

F

File	
Closing	2-15
Unloading	2-15
File information	2-17
File locations in Windows	1-12
File paths	2-26
File transfer between Windows and Macintosh	2-13
Filter Wizard	3-30
Filtering image noise	3-29
Find Spot	6-16
Finding spots	8-4
Flipping images	3-26
Fluor-S MAX MultiImager	
Acquiring the image	F-10
Control panel	F-3
Custom settings	F-5
Exposure time	F-8
Exposure times, recommended	F-9

PDQUEST User Guide

File size of images	F-17
Focusing	F-8
Options	F-12
Positioning	F-7
Preview scan	F-10
Saturated pixels, highlighting	F-17
Scan dimension	F-6
Selecting an application	F-4
Sensitivity, high and ultra	F-6
Simulation mode	F-2
Fluor-S MultiImager	
Acquiring the image	E-10
Control panel	E-3
Custom settings	E-5
Exposure time	E-8
Exposure times, recommended	E-9
File size of images	E-17
Focusing	E-7
High resolution / high sensitivity	E-6
Options	E-11
Positioning	E-6
Preview scan	E-9
Saturated pixels, highlighting	E-16
Scan dimension	E-6
Selecting an application	E-4
Simulation mode	E-2
FX, see Molecular Imager FX	

G

Gaussian modeling	1-7
Gel Doc	
Acquiring the image	A-4
Annotating images	A-8
Auto expose	A-5
Auto exposure settings	A-12

Index

Backup images	A-13
DAC settings	A-11
Display settings	A-9
Exposure Status bar	A-9
Freeze	A-7
Image mode	A-8
Imaging area	A-12
Invert display	A-10
Live image display	A-4
Manual expose	A-6
Options	A-11
Positioning the sample	A-4
Saturated pixels, highlighting	A-9
Saving images	A-8
Simulation mode	A-2
UV mode	A-9
Video card	A-2
Video display window	A-3
Video printing, footer info	A-12
White light mode	A-9
Gel Image file	5-14
Gel Record	4-5
Gel scans	
Calibrating	4-1
Gel Spots	
Background	2-25
Gel Spots file	5-14
Grab	3-3
Graph	
Advance	10-4
Column Graph	10-4
Erratic Responses	6-16
Margin Graph	10-4
Page of Graphs	10-4
Partial Matches	6-14
Pop-up Graph	10-3
Retreat	10-4
Graphs	10-1
Display order of histogram bars	10-2

PDQUEST User Guide

GS-700 Imaging Densitometer	
Applications, selecting	C-4
Calibration	C-10
Calibration settings	C-14
Filters and light source, selecting	C-5
Options	C-14
Oversampling	C-14
Preview scan	C-6
Resolution, selecting	C-7
Saturated pixels, highlighting	C-15
Scan area, selecting	C-7
Scanning an image	C-9
Scanning window	C-3
SCSI card	C-1
Simulation mode	C-1
Step tablet, editing	C-11
GS-710 Imaging Densitometer	
Applications, selecting	D-4
Calibration	D-9
Calibration settings	D-13
Filters and light source, selecting	D-6
Options	D-14
Oversampling	D-14
Preview scan	D-7
Resolution, selecting	D-8
Saturated pixels, highlighting	D-15
Scan area, selecting	D-7
Scanning an image	D-13
Scanning window	D-3
SCSI card	D-1
Simulation mode	D-1
Step tablet, editing	D-10

H

Hardware Protection Key

Index

Macintosh	1-17
PC	1-16
Hide Overlays	3-4
Higher Level Matchsets	
Creating	9-17
Diagram	9-16
Standard	9-18
Transferring annotations	9-20
Transferring MrpI data	9-20
Transferring numbers	9-19
Histograms	10-1

I

Image report, printing	11-3
Image size, changing	2-18
Images	
Assigning to windows	3-6
Exporting	11-6
Interchanging	3-8
Imitate Zoom	3-11
Importing TIFF files	2-13
Installation	1-12
Macintosh	1-13
PC	1-12
Interchanging images	3-8
Invert	
Image data	3-35
Image display	3-22

K

Keyboard shortcuts	2-10
--------------------------	------

PDQUEST User Guide

L

Landmarking	6-9
Auto-match after	2-24
License	1-14
Load	
Crop Settings	3-25
Loading files	2-11
Location in Gel	3-5
Low quality spots	5-20
Low quantity spots	5-18

M

Macintosh, memory assigned to Quantity One	1-11
Magnifying images	3-1
Make Spot Crosshair	5-22
Make Standard	6-7, 9-18
Manual Match	6-17
Manual Unmatch	6-17
Margin Graph	10-4
Mark	
Added Spots	6-20
Erratic Responses	6-16
MrpI Spots	7-3
Partial Matches	6-14
Unmatched Spots	6-18
Match	
Offsets	6-13
Overview	1-9
Symbols	6-12
Matches	
Editing	6-11
Removing all	6-21
Matching	6-10

Index

Adding a spot to the Standard	6-20
All gels	6-11
Colonies	6-21
Editing the results of	6-11
Landmarking	6-9
Manual	6-17
Removing a spot from the Standard	6-20
Unmatching	6-17
Matchset	
Partial Matches	6-14
Show Match	6-13
Standard	6-6
Matchsets	
Adding a member	6-8
Automatic matching	6-10
Creating	6-2
Creating a standard	6-7
Editing spots in	6-8
Exporting data	11-4
Info	6-8
Removing a member	6-8
Reordering	6-8
Menus	2-1
Merging	
Cropping images to be merged	3-24
Merging multiple exposures	4-11
Molecular Imager FX	
Acquiring the image	H-11
Control panel	H-2
File size of images	H-11
Options	H-12
Saturated pixels, highlighting	H-13
Saving the image	H-12
Scanning window	H-3
Selecting an application	H-4
Selecting resolution	H-10
Selecting scan area	H-9
Simulation mode	H-2
Mouse-assignable tools	2-9

PDQUEST User Guide

Behavior with keyboard commands	2-10
MrpI	
Display	7-3
Entering Data	7-2
Plot Grid	7-3
Remove Data	7-3
Standards	7-1
MrpIs	
Transferring	9-20

N

Normalization	7-7
---------------------	-----

O

Optical Density	1-7
Optimizing images, see Transform.	
Overlays, hiding	3-4

P

Pages of Graphs	10-4
Panning	3-5
Password, entering	1-21
Personal Molecular Imager FX	
Acquiring an image	G-6
Control panel	G-2
File size of images	G-6
Options	G-7
Saturated pixels, highlighting	G-7
Saving the image	G-6

Index

Scanning window	G-3
Selecting resolution	G-5
Selecting scan area	G-4
Simulation mode	G-2
pH Gradient	2-25
Plot	
DSN Grid	8-4
MrpI Grid	7-3
SSP Grid	8-2
Pop-up Graph	10-3
Preferences	
Display	2-25
General	2-22
Toolbars	2-23
Printing	
Actual-size images	11-3
Annotations	8-9
Image reports	11-3
Images	11-1
Video printouts	11-4

Q

Qualitative Analysis Sets	9-8
Quantitation	
Mode	7-5
Replicate Groups	7-4
Single Gels	7-4
Quantitation of spots	5-17
Quantitative Analysis Sets	9-5
Quantity Graphs	10-12
Quantity Table	7-6

PDQUEST User Guide

R

Reduce File Size	2-19
Refresh	3-5
Registration	
By fax / e-mail	1-21
By Internet	1-20
Form	1-19
Software license	1-14
Trial Period	1-17
Registration Form	1-19
Remove MrpI Data	7-3
Remove Segment	4-10
Removing a spot from the standard	6-20
Removing all matches	6-21
Replicate Group	
Creating	7-4
Edit Group	7-5
Rename	7-5
Reports	
A-B Comparison	10-7
Normalization	7-7
Quantitation Table	10-11
Quantity	10-11
Quantity Graphs	10-12
Scatter Plots	10-8
Standard Image	10-14
Revert to Saved	2-16
Rotating images	3-26

S

Sample images	2-13
Saturated pixels, highlighting	3-22
See also individual imaging devices	

Index

Saturated spots	5-24
Save All command	2-16
Save As command	2-15
Saving a file	2-15
Scanning	
Overview	1-6
Scans	
Printing	11-1
Scatter Plots	10-8
Show Match	6-13
Simulation mode, see individual imaging devices	
Size of images, changing	2-18
Software license	1-14
Spot Crosshairs	5-15
Spot Detection	
Overview	1-6
Reviewing	5-21
Steps	1-8
Spot Detection Wizard	5-1
Background subtraction	5-6
Finding spot centers	5-7
Fix faint spots	5-9
Loading images	5-2
Noise removal	5-6
Parameters, loading	5-3
Processing gels	5-10
Saving parameters	5-10
Sensitivity and size scale	5-9
Streak removal	5-7
Spots	
Adding crosshairs	5-22
Cancelling	5-23
Combining	5-24
Crosshairs, displaying	5-15
Detecting	5-10
Detection parameters, selecting	5-1
Detection report	5-13
Detection Wizard	5-1
Editing	5-21

PDQUEST User Guide

Ellipses	5-16
Finding	8-4
Information	5-17
Low quality sets	5-20
Low quantity sets	5-18
Parameters	5-21
Quality	5-19
Quantity	5-17
Restoring cancelled spots	5-24
Saturated	5-24
Uncombining combined spots	5-24
SSP numbers	8-1
Displaying	8-2
Transferring	9-19
Standard Image	10-14
Standard image	6-6
Creating	6-7
Selecting	6-6
Standard spot numbers. See SSP numbers	
Starting Quantity One	1-14
Statistical Analysis Sets	9-13
Status boxes	2-2
Subwindows	
Commands	3-10
Configuring	3-12
Definition	2-5
Focusing	2-8
Magnifying all the same	3-6
Selecting	2-8, 2-25

T

Technical Service, contacting	1-22
Temporary HPK License	1-18
Three by three (3x3) sections, viewing	3-5
TIFF files, importing	2-13

Index

Tile Windows	3-11
Tool Help	2-2
Toolbars	
Main	2-2
Secondary	2-3
Tools, mouse-assignable	2-9
Transfer	
DSN's	9-19
SSP's	9-19
Transform	3-17
All subwindows	3-23
Auto-scale	3-19
Controls	3-19
Gamma slider	3-21
High/Low Sliders	3-20
Highlight Saturated Pixels	3-22
Histogram	3-19
Invert Display	3-22

U

Unable to Obtain Authorization message	1-15
--	------

V

Video printing	11-4
View Entire Image	3-4
Viewing images	
Centering on the cursor	3-4
Undoing view commands	3-5

PDQUEST User Guide

W

Windows	
Assigning/interchanging images	3-6
Imitate Zoom	3-11
Multiple subwindow viewing commands	3-5
Tiling	3-11
Windows and subwindows, overview	2-5
Windows file locations	1-12

Z

Zoom	
Imitate in all subwindows	2-26
Zoom Box	3-2
Zoom In/Zoom Out	3-3
Zoom tools	
Changing behavior	2-25