



Multispot HIV-1/HIV-2 Rapid Test

The Clear Choice to See HIV-1 and HIV-2 — Separately

BIO-RAD



Multispot HIV-1/HIV-2 Rapid Test

Trusted Interpretation


- HIV-1 Positive
- HIV-2 Positive
- HIV-1 Indeterminate
(Applicable to Diagnostic Testing Algorithm
only and not to Rapid Testing Criteria)
- HIV Positive
(Undifferentiated)

**The only HIV antibody
differentiation test
for use in a Diagnostic
Testing Algorithm***

*CLSI M53-A (Vol 31 No 13)

Test Result Appearance and Interpretation

Spot Identity

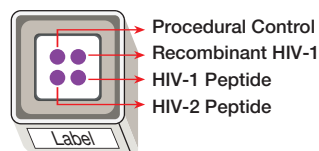

Negative



HIV-2 Reactive:




Diagnostic Testing Algorithm - Positive
Rapid Testing - Preliminary Positive

HIV-1 Reactive: Diagnostic Algorithm

HIV-1 Positive


Purple color development for both HIV-1 spots

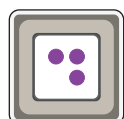

HIV-1 Reactive: Diagnostic Algorithm

HIV-1 Indeterminate*

Purple color development for one but not both HIV-1 spots

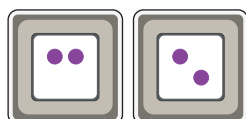


* HIV-1 Indeterminate is a category for diagnostic testing algorithm only. There are no indeterminate results when test is run as a rapid test screen.

HIV-1 Reactive: Rapid Testing

HIV-1 Preliminary Positive


Purple color development for one or both HIV-1 spots

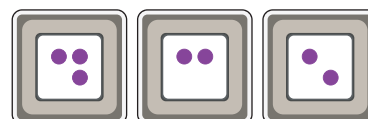

HIV Positive (Undifferentiated)

Diagnostic Testing Algorithm - Positive

Rapid Testing - Preliminary Positive

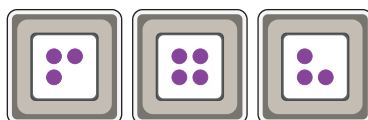

A dilutional procedure is used to differentiate samples that have purple color development in both HIV-1 and HIV-2 spots. After dilutional testing, if the dual reactivity does not disappear or both spots become nonreactive at the same dilution, the specimen should be reported as Undifferentiated. Refer to the Package Insert for detailed instructions.

Invalid

If no color develops in the Procedural Control Spot, the results are invalid



Multispot HIV-1/HIV-2 Rapid Test

INTENDED USE:

The Multispot HIV-1/HIV-2 Rapid Test is a single use qualitative immunoassay to detect and differentiate circulating antibodies to Human Immunodeficiency Virus Types 1 and 2 (HIV-1, HIV-2) in fresh or frozen human serum and plasma. This rapid HIV-1/HIV-2 test kit is intended as an aid in the diagnosis of infection with HIV-1 and/or HIV-2 in fresh or frozen human serum or plasma. This test is suitable for use in multi-test algorithms designed for statistical validation of rapid HIV test results or as part of an HIV-1/HIV-2 diagnostic testing algorithm that includes differentiation of HIV-1 and HIV-2 antibodies.


RESTRICTIONS:

- Sale of Multispot HIV-1/HIV-2 Rapid Test is restricted to clinical laboratories that have an adequate quality assurance program, including planned systematic activities to provide adequate confidence that requirements for quality will be met and where there is assurance that operators will receive and use the instructional materials.
- The Multispot HIV-1/HIV-2 Rapid Test is approved for use only by an agent of a clinical laboratory.
- Test subjects must receive the "Subject Information Notice" prior to specimen collection, and appropriate information when test results are provided, unless this test is used as part of a multi-test diagnostic algorithm.
- The Multispot HIV-1/HIV-2 Rapid Test is not approved for use to screen blood, plasma, cell, or tissue donors.


Multispot HIV-1/HIV-2 Rapid Test

Procedure

(Refer to the Package Insert for detailed procedure instructions.)

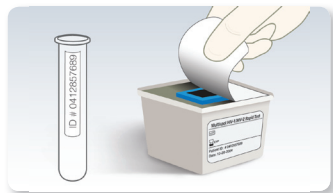

Remove foil. Label cartridge and specimen or control test tubes.

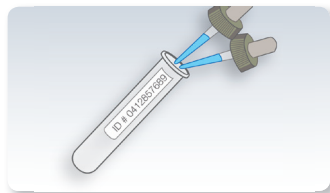

Add 2 dropperfuls (300 µL) of specimen diluent to each test tube. Add 1 drop (30 µL) specimen or control. Mix well.



Pour specimen into the prefilter. Wait 2 min. Remove and discard prefilter.



Fill cartridge with wash solution and let absorb. Add 3 drops conjugate. Wait 2 min.



Fill with wash solution and let absorb. Repeat. Add 3 drops development reagent. Wait 5 min.



Fill with stop solution. Allow to absorb and read results.

Ordering Information

Catalog No. **Description**

25228 Multispot HIV-1/HIV-2 Rapid Test* 50 tests

* Only Available in the US.



**Bio-Rad
Laboratories**

For further information, please contact the Bio-Rad office nearest you or visit our website at www.bio-rad.com/diagnostics

*Clinical
Diagnostics Group*

Website www.bio-rad.com/diagnostics **U.S.** 1-800-224-6723 **Australia** 61-2-9914-2800 **Austria** 43-1-877-8901 **Belgium** 32-9-385-5511 **Brazil** 5521-3237-9400
Canada 1-514-334-4372 **China** 86-21-61698500 **Czech Republic** 420-241-430-532 **Denmark** +45-4452-1000 **Finland** 358-9-804-22-00 **France** 33-1-47-95-60-00
Germany +49 (0)89-318-840 **Greece** 30-210-7774396 **Hong Kong** 852-2789-3300 **Hungary** +36-1-459-6100 **India** 1800-180-1224 **Israel** 972-3-9636050 **Italy** +39-02-216091 **Japan** 81-3-6361-7070 **Korea** 82-2-3473-4460 **Mexico** +52 (55)5488-7670 **The Netherlands** +31-318-540666 **New Zealand** 64-9-415-2280 **Norway** 47-23-38-41-30 **Poland** 48-22-3319999 **Portugal** 351-21-472-7700 **Russia** 7-495-721-14-04 **Singapore** 65-6415-3170 **South Africa** 27-11-442-85-08 **Spain** 34-91-530-5200
Sweden 46-8-555-127-00 **Switzerland** +41 (0)26-674-52-01/02/03 **Taiwan** 886-2-2578-7189 **Thailand** 662-651-8311 **United Kingdom** +44 (0)20-8328-2000