



VARIANT™ II TURBO Link Hemoglobin Testing System

**TO USE AS AN AID IN
THE DIAGNOSIS OF DIABETES**



BIO-RAD



VARIANT™ II TURBO Link Hemoglobin Testing System

Efficient
streamlined process
Powerful
integrated system
Automated
saves time and resources

With so many patients counting on your results each day, delivering the accuracy they need as quickly as possible is a top priority. We designed the VARIANT™ II TURBO Link Hemoglobin Testing System with your needs—and high standards—in mind.

The VARIANT™ II TURBO Link Hemoglobin Testing System integrates our proven HbA_{1c} assay performance with the Sysmex® XN-9000 Automated Hematology System to handle over 90% of whole blood EDTA samples while delivering maximum flexibility and accuracy. Confidently report HbA_{1c} results to aid physicians with diagnosing those who have diabetes, identifying those who may be at risk of developing the disease, or monitoring the progress of those living with diabetes.

Let the consolidated VARIANT™ II TURBO Link Hemoglobin Testing System help you see the full picture and make a difference in the lives of your patients.

Bio-Rad VARIANT II TURBO Link HbA1c
CDM 5.2 Software

PATIENT REPORT
V2TURBO_STL

Patient Data

Sample ID:
Patient ID:
Name:
Physician:
Sex:
DOB:

Analysis Data

Analysis Performed: 04/27/2017 10:06:12
Injection Number: 36
Run Number: 3
Rack ID: 18
Tube Number: 8
Report Generated: 04/27/2017 10:06:12
Operator ID:


Comments:

Peak Name	IFCC mmol/mol	NGSP %	Area %	Retention Time (min)
Unknown	—	—	0.4	0.106
A1a	—	—	1.1	0.163
A1b	—	—	0.9	0.225
F	—	—	1.3	0.287
LA1c	—	—	1.0	0.424
A1c	62*	7.8*	—	0.519
P3	—	—	2.5	0.773
P4	—	—	0.9	0.869
Ao	—	—	49.3	1.018
Variant Window	—	—	38.7	1.126

*Values outside of expected range

Total Area

HbA1c (IFCC) = 62* mmol/mol HbA1c (NGSP) = 7.8*%

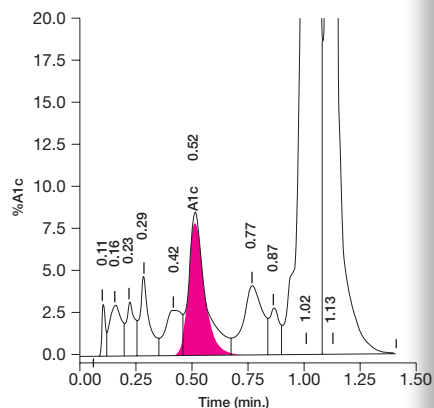

HbA1c reportable in the presence of hemoglobin S, C, D, and E traits*

Bio-Rad VARIANT II TURBO Link HbA1c
CDM 5.2 Software

PATIENT REPORT
V2TURBO_STL

Patient Data

Sample ID:
Patient ID:
Name:
Physician:
Sex:
DOB:

Analysis Data


Analysis Performed: 02/14/2017 10:30:39
Injection Number: 4513
Run Number: 101
Rack ID: 0011
Tube Number: 1
Report Generated: 02/14/2017 10:44:15
Operator ID:

Comments:

Peak Name	NGSP%	Area %	Retention Time (min)	Peak Area
Unknown	—	0.4	0.112	9414
A1a	—	8.4	0.176	218589
F	—	23.9	0.292	619799
A1c	5.9	—	0.510	98995
P3	—	3.6	0.777	93134
P4	—	1.0	0.875	26310
Unknown	—	0.5	0.967	13930
Ao	—	58.3	1.035	1507942

Total Area: 2,588,114

HbA1c (NGSP) = 5.9%

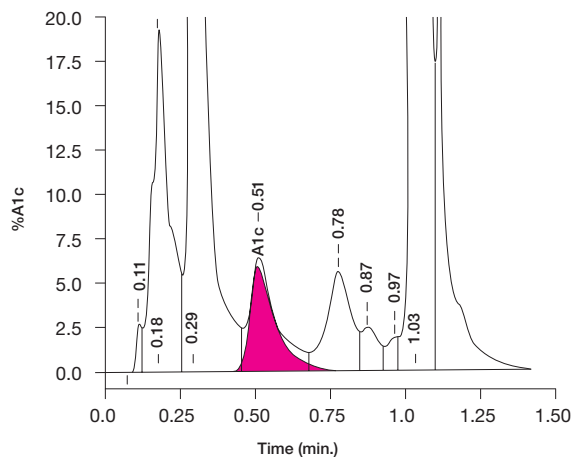

HbA1c reportable in the presence of HbF up to 25%*

* See VARIANT™ II TURBO HbA1c Kit - 2.0 Instructions for Use (Catalog # 12000447)

Every Patient Deserves **Quality Results**



Turn greater efficiency into better patient care

With the VARIANT™ II TURBO Link Hemoglobin Testing System, we're combining Bio-Rad's decades of leadership in HbA_{1c} with Sysmex's proven expertise in automation to improve your laboratory's workflow. The result is a fully-automated solution that helps you be the difference with timely, accurate results.

Streamlined Process

- On-demand testing 24/7
- Sysmex Lavender Top Management®
- Workflow optimization
- Reduced sample handling errors and delays

Save Time and Labor

- Specimen retrieval and archiving using a sorter
- Standardized automated sample data management using Sysmex WAM™
- Automated reflex and repeat testing using Sysmex WAM™

More than just an HbA1c Result

Bio-Rad A1c provides clinically significant information by

- Detecting common hemoglobin variants
- Detecting potentially interfering substances
- Identifying those who may be at risk of developing diabetes
- Monitoring patients living with diabetes
- Aiding in the diagnosis of diabetes

Powerful Fully-integrated Lavender Top Management® solution

- CBC testing
- Integrated slide maker/stainer
- HbA_{1c} testing
- Automated sample sorting and archiving
- Sysmex WAM decision logic rules database

Lavender Top Management is a registered trademark of Sysmex

*Every Patient Deserves Quality Results. You Deserve **operational efficiency.***




Experience superior HbA_{1c} quality, accuracy, and efficiency

- Gold standard ion-exchange High Performance Liquid Chromatography (HPLC) methodology
- No interference in the presence of:
 - Hemoglobin S, C, D, or E traits*
 - Hemoglobin F (HbF) concentrations up to 25%*
 - Carbamylated (CHb) and Labile A_{1c} (LA1c)*
 - Acetylated hemoglobin*
- Anchored to the IFCC reference method
- NGSP certified and traceable to the Diabetes Control and Complications Trial (DCCT)
- Sample results in 97 seconds
- Cartridge lifetime of 2,500 injections
- Hemoglobin Capillary Collection System (HCCS) allows remote finger stick sample collection

* See VARIANT™ II TURBO HbA_{1c} Kit - 2.0 Instructions for Use (Catalog # 12000447)

VARIANT™ II TURBO Link Hemoglobin Testing System

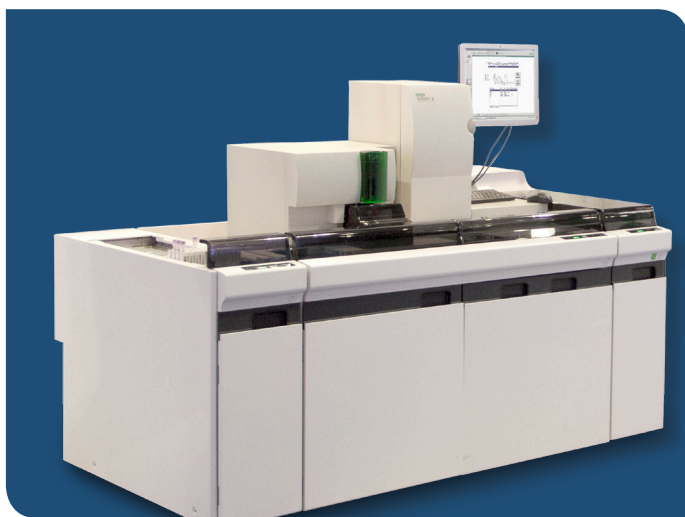

Unparalleled Support

Bio-Rad's world class service and support team is backed by our decades-long leadership in the field of HbA_{1c} testing. The VARIANT™ II TURBO Link Hemoglobin Testing System represents our continued commitment to helping you be the difference for the patients who count on you.

Ordering Information

Catalog #	Description
270-2700	VARIANT™ II TURBO Link Hemoglobin Testing System (110/220V)..... 1 unit
12000447	VARIANT™ II TURBO HbA _{1c} Kit – 2.0..... 2500 tests
270-2730	VARIANT™ II TURBO Hemoglobin A _{1c} Wash/Diluent Set..... 4 x 2500 mL
740	Lyphochek® Diabetes Control, Bilevel (3 each of 2 levels)..... 6 x 0.5 mL
12000070	Lyphochek® Hemoglobin A _{1c} Linearity Set (1 each of 6 levels)..... 6 x 0.5 mL
BC80	EQAS Hemoglobin Program (3 analytes) 12 x 0.5 mL
171	Liquichek™ Diabetes Control, Level 1 6 x 1 mL
172	Liquichek™ Diabetes Control, Level 2 6 x 1 mL
173	Liquichek™ Diabetes Control, Level 3 6 x 1 mL
196-2051	Hemoglobin Capillary Collection System (HCCS) 5 tests
196-2052	Hemoglobin Capillary Collection System (HCCS) 100 tests
196-2053	Hemoglobin Capillary Collection System (HCCS) 5000 tests
196-2100	Plastic Capillaries, 100 per dispenser 20 dispensers

BIO-RAD

**Bio-Rad
Laboratories**

For further information, please contact the Bio-Rad office nearest you or visit our website at www.bio-rad.com/diagnostics

*Clinical
Diagnostics Group*

Website www.bio-rad.com/diagnostics **U.S.** 1-800-224-6723 **Australia** 61-2-9914-2800 **Austria** 43-1-877-8901 **Belgium** +32 (3)710-53-00 **Brazil** +55 (31)3689-6600 **Canada** 1-514-334-4372 **China** 86-21-61698500 **Czech Republic** 420-241-430-532 **Denmark** +45-4452-1000 **Finland** 358-9-804-22-00 **France** 33-1-47-95-60-00 **Germany** +49 (0)89-318-840 **Greece** 30-210-7774396 **Hong Kong** 852-2789-3300 **Hungary** +36-1-459-6100 **India** 1800-180-1224 **Israel** 972-3-9636050 **Italy** +39-02-216091 **Japan** 81-3-6361-7070 **Korea** 82-2-3473-4460 **Mexico** +52 (55)5488-7670 **The Netherlands** +31-318-540666 **New Zealand** 64-9-415-2280 **Norway** +47-23-38-41-30 **Poland** 48-22-3319999 **Portugal** 351-21-472-7700 **Russia** +7-495-721-1404 **Singapore** 65-6415-3170 **South Africa** 27-11-442-85-08 **Spain** 34-91-590-5200 **Sweden** +46-8-555-127-00 **Switzerland** +41 (0)26-674-55-05/06 **Taiwan** 886-2-2578-7189 **Thailand** 662-651-8311 **United Kingdom** +44 (0)20-8328-2000