

VARIANT™ II Hemoglobin A1c Reorder Pack

You step up to provide the best care for your patients everyday. And now we're stepping up our commitment to helping you deliver reliable, accurate results you require. The enhanced VARIANT™ II HbA_{1c} Reorder Pack lets you be the difference for your laboratory's patients.

Take Advantage of Expanded Reportability Claims

With an enhanced cartridge and new prefilter, our VARIANT™ II HbA_{1c} assay delivers better peak separation resulting in an A1C reportability in the presence of:

- Hemoglobin E, D traits in addition to S or C traits
- Carbamylated Hemoglobin (CHb)
- Labile A1c (LA1c)


Manage your VARIANT™ II with ease

The enhanced VARIANT™ II HbA_{1c} Reorder Pack offers:

- Doubled cartridge lifetime to 1000 injections minimizing hands-on time
- Reduced calibration frequency while maintaining consistently precise HbA_{1c} results
- Interchangeability of different lots of Buffers A and B within an analytical cartridge resin lot simplifying workflow
- Extended kit shelf life providing streamlined inventory management

Accuracy


The following graph shows the correlation of HbA_{1c} values on NGSP samples with that obtained on the new VARIANT™ II HbA_{1c} assay.

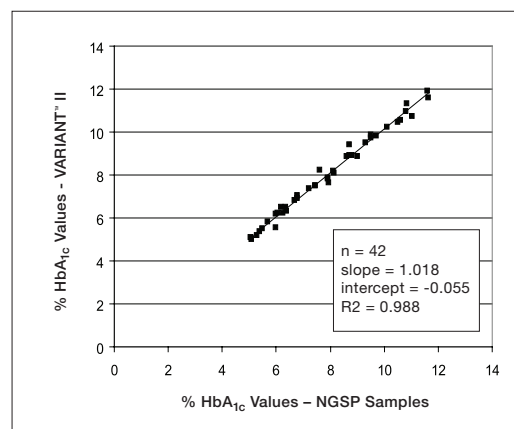**BIO-RAD**

VARIANT™ II Hemoglobin A1c Reorder Pack


Sample Chromatograms

The enhanced resolution of peaks facilitates reporting of A1C in presence of CHb and labile A1c and variant traits such as HbE and D thereby eliminating the need to review chromatograms before reporting A1C values.

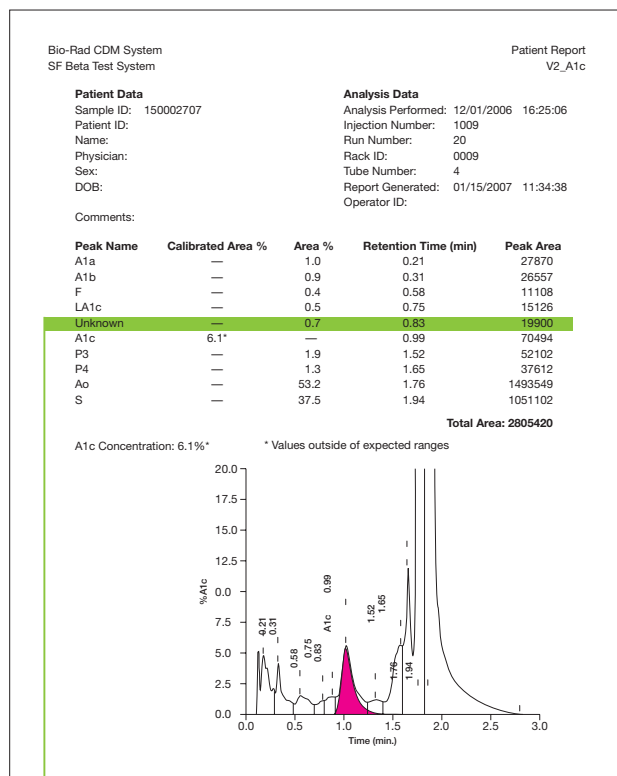

A1C reportable in presence of CHb

CHb identified as unknown peak in between A1C and LA1c

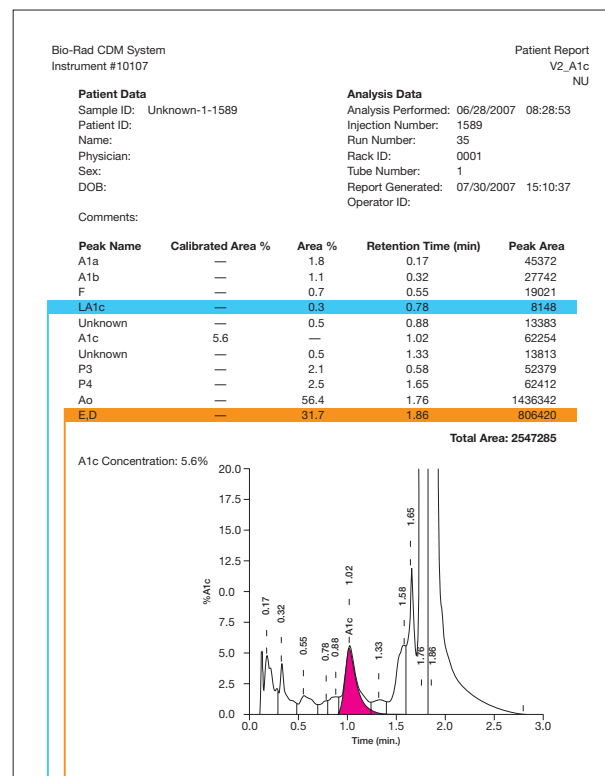

A1C reportable in presence of HbE trait

A1C reportable in presence of labile A1c

Ordering Information

Catalog No. Description

270-2000	VARIANT™ II Hemoglobin Testing System	1 unit
270-2101NU	VARIANT™ II Hemoglobin A1c Reorder Pack	1000 tests
740	Lyphochek® Diabetes Control, Bilevel (3 each of 2 levels).....	6 x 0.5 mL
120	Lyphochek® Hemoglobin A1c Linearity Set, (1 each of 4 levels)	4 x 0.5 mL
BC80	EQAS™ Hemoglobin Program (3 analytes)	12 x 0.5 mL

Bio-Rad A1c • Be the difference



Bio-Rad
Laboratories

For further information, please contact the Bio-Rad office nearest you or visit our website at www.bio-rad.com/diagnostics.

Clinical
Diagnostics Group

Website www.bio-rad.com/diagnostics U.S. 1-800-2BIO-RAD Australia 61-2-9914-2800 Austria 43-1-877-8901 Belgium 32-9-385-5511 Brazil 5521-3237-9400
Canada 1-514-334-4372 China 86-21-64260808 Czech Republic 420-241-430-532 Denmark +45-4452-1000 Finland 358-9-804-22-00 France 33-1-47-95-60-00
Germany +49-(0)89-318-840 Greece 30-210-7774396 Hong Kong 852-2789-3300 Hungary 36-1-455-8800 India 91-124-4029300 Israel 972-3-9636050 Italy +39-02-216091
Japan 81-3-6361-6902 Korea 82-2-3473-4460 Mexico +52(55)5488-7670 The Netherlands +31-318-540666 New Zealand 64-9-415-2280 Norway 47-23-38-41-30
Poland 48-22-3319999 Portugal 351-21-472-7700 Russia 7-495-721-14-04 Singapore 65-6415-3188 South Africa 27-11-442-85-08 Spain 34-91-590-5200
Sweden 46-8-555-127-00 Switzerland 41-61-717-95-55 Thailand 662-651-8311 United Kingdom +44-(0)20-8328-2000