



ANALYZER




PRINTER



KEYPAD




CARTRIDGE

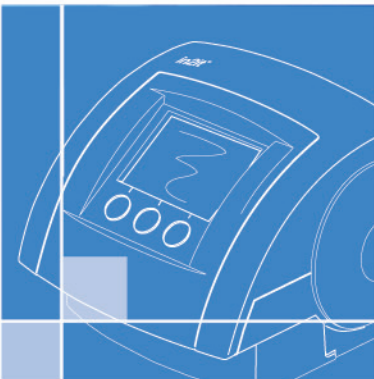


in2it™ A1C Test Procedure

1. Getting started


- Test cartridges must be warmed to room temperature
- Consult the user manual for instructions on collecting control or venous samples

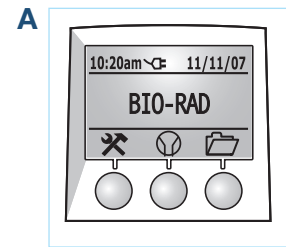


When the in2it™ is ready to run a test, the cartridge symbol  will flash. Press the middle button (under the cartridge symbol), to start a test.

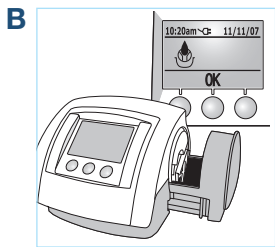


The in2it™ door will open. It is ready to accept a cartridge.

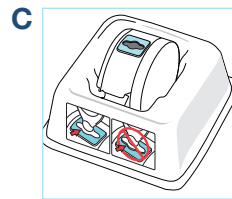

Place a cartridge in the workstand with the blue blood port facing up.

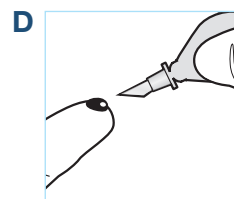


With a sterile lancet, prick the patient's finger and gently squeeze to obtain a generous drop of blood.

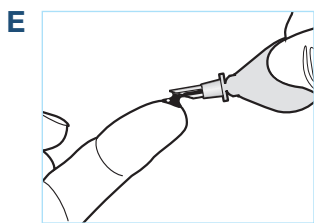


Holding the blood key by the round handle, touch the flat end to the blood drop. Blood will automatically flow into the key. **Take care not to get blood on the outside of the blood key.**

2. Insert blood key into cartridge


- Proper insertion of blood key into the test cartridge is essential. See Quick Reference Guide

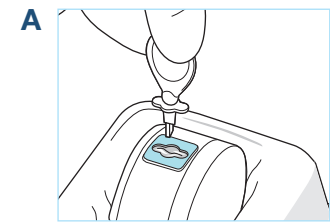


Align
Properly align the blood key prior to pressing the key through the blue blood port. The winged flanges on the blood key must face the same direction as the notch in the blood port.

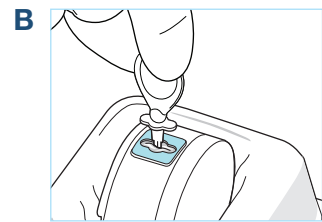

Pierce
Pierce the film covering the blue blood port. The blood key will drop slightly into the port.

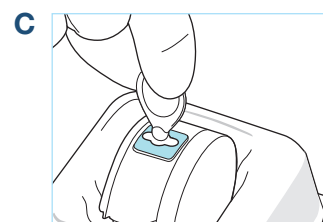

Push
Press **FIRMLY** down on the blood key until it locks in position. **You may feel two separate and distinct clicks.**


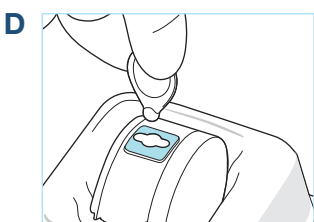


Twist
Twist the handle of the blood key until it snaps off.


3. Place cartridge into instrument




Immediately (within 10 seconds), place the cartridge in the door with the white plate facing inward and the blue blood port facing up.

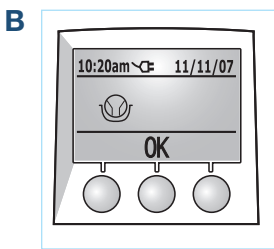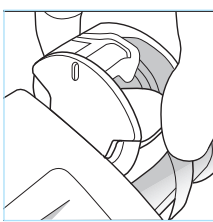

Immediately (within 10 seconds), press the OK button to close the door. Testing will start automatically.

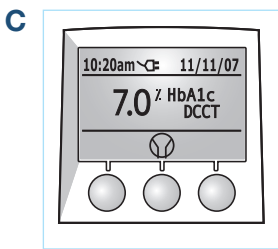

When the test is completed, the result will be displayed. Record the result.

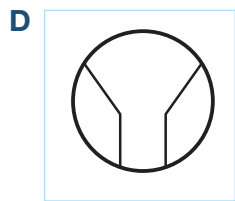

Press the cartridge symbol button to open the door and remove the cartridge. **After removing the cartridge, press the OK button to close the door prior to running the next test.**

4. Results displayed in 10 minutes

Contact Information

Customer/Technical Service
• 1-800-224-6723

Get in2it™!

BIO-RAD