

SSCP (Single-Strand Conformation Polymorphism)

The Sequi-Gen® II Electrophoresis Cells Offer the Optimal System for SSCP Analysis

Single-strand conformation polymorphism is a popular technique for detection of sequence variations in PCR*- amplified fragments of DNA. SSCP analysis requires uniform, low temperature, non-denaturing electrophoresis to maintain single-stranded DNA secondary structure.

The long, thin gel format of the Sequi-Gen II cell yields:

- Single base pair resolution
- Smile free nucleic acid separations

The patented upper buffer chamber of the Sequi-Gen II cell provides uniform temperature distribution throughout the entire gel surface.

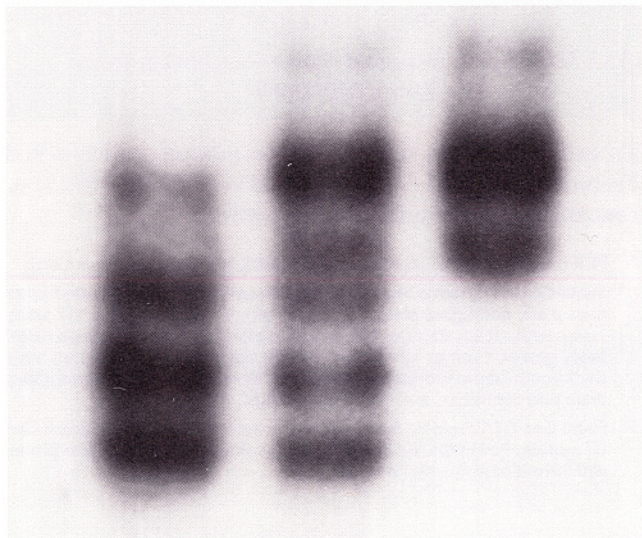


Photo contributed by Richard A. Gatti, M.D., Department of Pathology, UCLA School of Medicine, Los Angeles, California.

A technical review and detailed protocol for SSCP using the Sequi-Gen II cell is available (Bulletin 1857).

To order, or to receive more information, call Bio-Rad Technical Service at 1-800-4BIORAD (1-800-424-6723) or contact your local Bio-Rad Representative.

New Product Update

The Sequi-Gen II cell has been replaced with the new Sequi-Gen GT cell. All of the features that make the Sequi-Gen ideal for microsatellite and SSCP analysis have been kept in the new model. In addition, the new models offer a more convenient horizontal casting method and a size developed specifically for microsatellite and SSCP analysis (38 cm wide x 30 cm long). Request bulletin 2006 for more information.

Ordering Information

Catalog Number	Product Description
Sequi-Gen GT Sequencing Systems	
165-3860	Sequi-Gen GT System, 21 x 40 cm**
165-3861	Sequi-Gen GT System, 21 x 50 cm**
165-3862	Sequi-Gen GT System, 38 x 30 cm**
165-3863	Sequi-Gen GT System, 38 x 50 cm**
165-3802	Sequi-Gen GT/PowerPac 3000 System, 21 x 40 cm, 100/120 V†
165-3805	Sequi-Gen GT/PowerPac 3000 System, 21 x 40 cm, 220/240 V†
165-3803	Sequi-Gen GT/PowerPac 3000 System, 21 x 50 cm, 100/120 V†
165-3806	Sequi-Gen GT/PowerPac 3000 System, 21 x 50 cm, 220/240 V†
165-3810	Sequi-Gen GT/PowerPac 3000 System, 38 x 30 cm, 100/120 V†
165-3811	Sequi-Gen GT/PowerPac 3000 System, 38 x 30 cm, 220/240 V†
165-3804	Sequi-Gen GT/PowerPac 3000 System, 38 x 50 cm, 100/120 V†
165-3807	Sequi-Gen GT/PowerPac 3000 System, 38 x 50 cm, 220/240 V†

Microsatellite Mapping (Short Tandem Repeat Analysis)

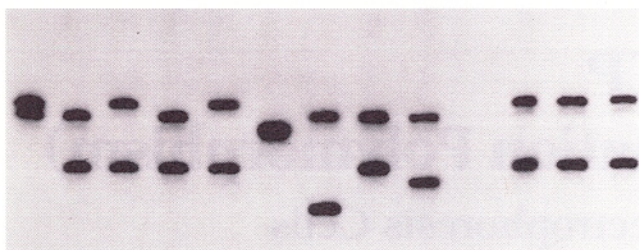


Photo contributed by Catriona MacGeoch, David P. Kelsell, and Nigel K. Spurr, Imperial Cancer Research Fund, U.K.

Microsatellite mapping is a very useful technique for genetic linkage analysis and identification of DNA polymorphisms. DNA is amplified by PCR* to detect variations in simple-sequence tandem repeats (microsatellites). Differences of only one tandem repeat (two base pairs) are possible, thus nucleic acid sequencing gels are required for product separation and analysis.

The Sequi-Gen II electrophoresis cell routinely provides:

- Single base pair resolution
- Smile free nucleic acid separation
- Multiple sample loading capabilities (1-96 samples)
- Three popular gel sizes (21 x 40, 21 x 50, and 38 x 50 cm)
- Wide range of single-thickness and wedge spacers, and sharks-tooth and well-forming combs for microsatellite mapping.

In addition, the patented upper buffer chamber of the Sequi-Gen II cell provides uniform temperature distribution throughout the entire gel surface.

Technical review and detailed protocol for microsatellite mapping using the Sequi-Gen II cell is available (Bulletin 1858).

New Product Update

The Sequi-Gen II has been replaced with the new Sequi-Gen GT. All of the features that make the Sequi-Gen ideal for microsatellite and SSCP analysis have been kept in the new model. In addition, the new models offer a more convenient horizontal casting method and a size developed specifically for microsatellite and SSCP analysis (38 cm wide x 30 cm long). Request bulletin 2006 for more information.

Ordering Information

Catalog Number	Product Description
Sequi-Gen GT Sequencing Systems	
165-3860	Sequi-Gen GT System, 21 x 40 cm**
165-3861	Sequi-Gen GT System, 21 x 50 cm**
165-3862	Sequi-Gen GT System, 38 x 30 cm**
165-3863	Sequi-Gen GT System, 38 x 50 cm**
165-3802	Sequi-Gen GT/PowerPac 3000 System, 21 x 40 cm, 100/120 V†
165-3805	Sequi-Gen GT/PowerPac 3000 System, 21 x 40 cm, 220/240 V†
165-3803	Sequi-Gen GT/PowerPac 3000 System, 21 x 50 cm, 100/120 V†
165-3806	Sequi-Gen GT/PowerPac 3000 System, 21 x 50 cm, 220/240 V†
165-3810	Sequi-Gen GT/PowerPac 3000 System, 38 x 30 cm, 100/120 V†
165-3811	Sequi-Gen GT/PowerPac 3000 System, 38 x 30 cm, 220/240 V†
165-3804	Sequi-Gen GT/PowerPac 3000 System, 38 x 50 cm, 100/120 V†
165-3807	Sequi-Gen GT/PowerPac 3000 System, 38 x 50 cm, 220/240 V†

To order, or to receive more information, call Bio-Rad Technical Service at 1-800-4BIORAD (1-800-424-6723) or contact your local Bio-Rad Representative.

* PCR is covered by U.S. patent number 4,683,202 issued to Cetus Corp.

** Sequi-Gen GT systems include GT IPC assembly (IPC and bonded inner glass plate, outer glass plate, and clamp set), GT universal base, GT safety covers with cables, stabilizer bar, precision caster assembly (precision caster base, gasket, tubing, luer tapers, tubing, and syringe), 0.40 mm vinyl sharktooth comb and spacers, gel temperature indicator, leveling bubble, drain port connector, and instruction manual.

† Sequi-Gen GT/PowerPac 3000 systems include the appropriate Sequi-Gen GT system, PowerPac 3000 Power Supply, PowerPac temperature probe, and PowerPac instruction manual.

BIO-RAD

Bio-Rad
Laboratories

Molecular
Bioscience Group

Bio-Rad Laboratories Main Office, 2000 Alfred Nobel Drive, Hercules, California 94547, Ph. (510) 741-1000, Fx. (510) 741-5800
Also in: Regents Park, Australia, Ph. 02-9914-2800, Fx. 02-9914-2889 Wien, Austria, Ph. (1) 877 89 01, Fx. (1) 876 56 29 Nazareth, Belgium, Ph. 09-385 55 11, Fx. 09-385 65 54
Mississauga, Canada, Ph. (905) 712-2771, Fx. (905) 712-2990 Beijing, China, Ph. (01) 2046622, Fx. (01) 2051876 Kybehavn, Denmark, Ph. 39 17 9947, Fx. 39 27 1698
Espoo, Finland, Ph. 90 804 2200, Fx. 90 804 1100 Ivry sur Seine Cedex, France, Ph. (1) 49 60 68 34, Fx. (1) 46 71 24 67 München, Germany, Ph. 089 31884-0, Fx. 089 31884-100
New Delhi, India, Ph. 91-11-461-0103, Fx. 91-11-461-0765 Milano, Italy, Ph. 02-21609.1, Fx. 02-21609.399 Tokyo, Japan, Ph. 03-5811-6270, Fx. 03-5811-6272 Veenendaal,
The Netherlands, Ph. 0318-540666, Fx. 0318-542218 Auckland, New Zealand, Ph. 09-443 3099, Fx. 09-443 3097 Kowloon, Hong Kong, Ph. 7893300, Fx. 7891257 Singapore,
Ph. (65) 272-9877, Fx. (65) 273-4835 Sundbyberg, Sweden, Ph. 46 (0) 8 627 50 00, Fx. 46 (0) 8 627 54 00 Madrid, Spain, Ph. (91) 661 70 85, Fx. (91) 661 96 98 Glattpburg, Switzerland,
Ph. 01/809 55 55, Fx. 01/809 55 00 Hemel Hempstead, United Kingdom, Ph. 0800 181134, Fx. 01442 259118