

Salieri
or
Mozart?

BIO-RAD

The choice between
very good



and genius.

At Bio-Rad,
we've composed a line
of IPG strips that feature
impressive range
and brilliant performance.



Discover the difference
for yourself.



pH 3

pH 4

pH 5

pH 6

pH 7


pH 8

pH 9

pH 10

The range of a true virtuoso

Bio-Rad micro-range ReadyStrip™ IPG strips deliver high-resolution results over the broadest pH range with fewer gels than those of the competition.

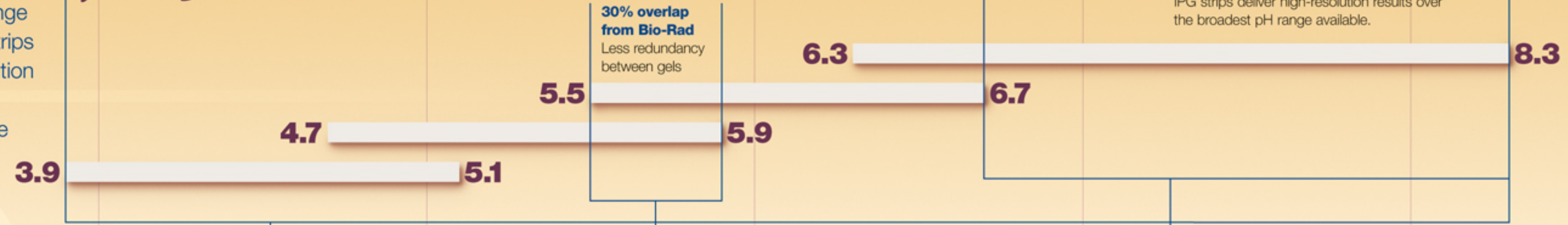

● The broadest total pH range

● Innovative gradient overlaps

● More data from 4 gels

Nobody else has this range!
The new Bio-Rad micro-range ReadyStrip IPG strips deliver high-resolution results over the broadest pH range available.

Get more data!

Streamline your analysis!


Save time and expense!

3.5

6.7

The competition's 50% overlap
Data repeated on each gel

The competition
Less data from more gels

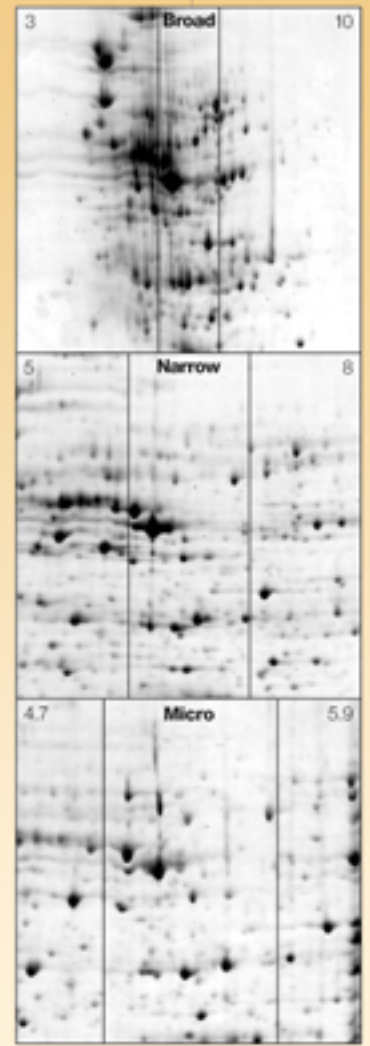

E. coli extract separated on 17 cm ReadyStrip IPG strips and PROTEAN® II SDS-PAGE gels. The increase in resolution of the center region of the pH 3–10 gel is shown with a narrow-range (pH 5–8) separation and micro-range (pH 4.7–5.9) separation.

	pH										Focusing Power			ReadyStrip IEF Buffer						
	3	4	5	6	7	8	9	10	7 cm	11 cm	17 cm	3-10	7-10	3.9-5.1	4.7-5.9	5.5-6.7	6.3-8.3			
Broad Range																				
3-10											1x	1.6x	2.4x	○						
3-10 nonlinear (NL)											1x	1.6x	2.4x	○						
Narrow Range																				
3-6											2.3x	3.7x	5.7x	○						
5-8											2.3x	3.7x	5.7x	○						
7-10											2.3x	3.7x	5.7x		○					
4-7											2.3x	3.7x	5.7x	○						
Micro Range																				
3.9-5.1											5.8x	9.2x	14.2x			○				
4.7-5.9											5.8x	9.2x	14.2x				○			
5.5-6.7											5.8x	9.2x	14.2x					○		
6.3-8.3											3.5x	5.5x	8.5x						○	

Bio-Rad Advantages

Bio-Rad's ReadyStrip IPG strips produce more spots than you've ever seen before on broad-, narrow-, and micro-range strips. You can take advantage of Bio-Rad's wide selection of pH ranges no matter the size of your second-dimension gels. Every pH range is offered in 3 strip lengths: 7 cm for mini gels, 11 cm for Criterion gels, and 17 cm for large gels. Each strip is also permanently labeled with the pH range at the anode end, making it easy to handle and track. All these features and more are built into the micro-range strips.

A Broader Total pH Range

Bio-Rad's micro-range ReadyStrip IPG strips span a total pH range of 3.9–8.3. The competitor's narrow-range strips span a total pH range of 3.5–6.7. This difference is very important because a significant portion of the proteome is found in the pH 6–8 range. However, a single-pH-unit strip near pH 7.0 is considered by 2-D experts to be problematic. Bio-Rad has solved this limitation by offering the pH 6.3–8.3 strip that spans 2 pH units and includes pH 7.0. No one else offers a zoom strip in this important pH range.

Streamlined Data Analysis

Each micro-range gradient overlaps the next by 30%. This is sufficient to allow easy gel matching. The competition's overlap is 50% from strip to strip. What does this mean to you? With the competitor's strips, up to 100% of each gel's data is repeated on the overlapping gels, meaning that you collect and process excessive redundant data. Data processing is time-consuming. With Bio-Rad's ReadyStrip IPG strips the amount of redundancy is limited to the amount needed for gel matching, streamlining data analysis.

More Data with Fewer Gels

An additional benefit of the innovative gradient overlap is that high resolution over the total pH range from 3.9 to 8.3 is achieved with just 4 gels. A customer using other strips runs 5 gels and only resolves a total pH range of 3.5–6.7. Bio-Rad strips generate more data with fewer gels, saving time and expense and increasing throughput.

Now that's genius!

Ordering Information

ReadyStrip IPG Strips, 12 per Package

pH Range	7 cm	11 cm	17 cm
3–10	163-2000	163-2014	163-2007
3–10 nonlinear (NL)	163-2002	163-2016	163-2009
4–7	163-2001	163-2015	163-2008
3–6	163-2003	163-2017	163-2010
5–8	163-2004	163-2018	163-2011
7–10	163-2005	163-2019	163-2012
3.9–5.1	163-2028	163-2024	163-2020
4.7–5.9	163-2029	163-2025	163-2021
5.5–6.7	163-2030	163-2026	163-2022
6.3–8.3	163-2031	163-2027	163-2023

Catalog # Description

Bio-Lyte® IEF Buffers*

163-2094	100x Bio-Lyte 3/10 Ampholyte, 1 ml
163-2093	100x ReadyStrip 7–10 Buffer, 1 ml
163-2098	100x ReadyStrip 3.9–5.1 Buffer, 1 ml
163-2097	100x ReadyStrip 4.7–5.9 Buffer, 1 ml
163-2096	100x ReadyStrip 5.5–6.7 Buffer, 1 ml
163-2095	100x ReadyStrip 6.3–8.3 Buffer, 1 ml
163-2099	ReadyStrip Instruction Manual, free upon request with ReadyStrip purchase

* Dilute ReadyStrip buffers to 1x final in each sample to equal 0.2% Bio-Lyte ampholyte.

BIO-RAD

**Bio-Rad
Laboratories, Inc.**

Life Science
Group

Web site www.bio-rad.com USA (800) 481-0RAD Australia 02 9914 2800 Austria (01) 877 89 01 Belgium 09-385 55 11 Brazil 55 21 507 6191
Canada (905) 712-2771 China (86-21) 63052255 Denmark 45 44 52-1000 Finland 358 (09 804 2200 France 01 47 95 69 65 Germany 089 318 84-177
Hong Kong 852-2789-3300 India (91-124) 6398112/113/114, 6450092/93 Israel 03 951 4127 Italy 39 02 216091 Japan 03-5811-6270
Korea 82-2-3473-4460 Latin America 305-894-5950 Mexico 52 5 534 2552 to 54 The Netherlands 0318-540666 New Zealand 64-9-4152280
Norway 47-23-38-41-30 Portugal 351-21-472-7700 Russia 7 095 721 1404 Singapore 65-2729677 South Africa 00 27 11 4428508 Spain 590 5200
Sweden 46 (08-55 51 27 00 Switzerland 061 717-9555 United Kingdom 0800-181134