


SaSelect™ Procedure and Interpretation

Procedure

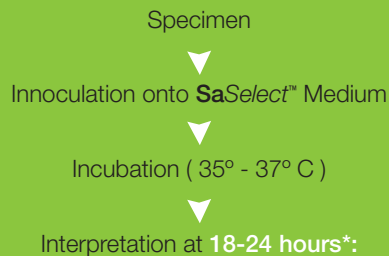

* If incubated > 24 hours, confirm any colonies


* Interpretation of the test results should be made taking into consideration the patient's history, the source of the specimen, colonial and microscopic morphology and if necessary, the results of any other tests performed.

SaSelect™ Direct Identification

No Additional testing required



Staphylococcus aureus
Pink to orange colonies

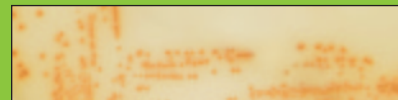

Staphylococcus aureus
Pink to orange colonies

▶ Presumptive & Mixed Colony ID on reverse side


For more information, visit us on the web at www.bio-rad.com/diagnostics
Call toll free at 1-800-2BIO-RAD (1-800-224-6723, option 2, option 3);
Outside the US, contact your local sales office.

BIO-RAD

SaSelect™ Interpretation

SaSelect™ Presumptive Identification

Additional testing required for species identification

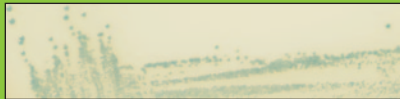


Staphylococcus cohnii
Blue to turquoise colonies




Staphylococcus saprophyticus
Blue to turquoise colonies




Staphylococcus intermedius
Purple or grey colonies

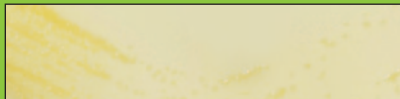

Staphylococcus lugdunensis
White to yellow colonies

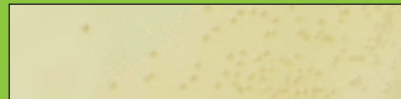

Staphylococcus haemolyticus
White to yellow colonies

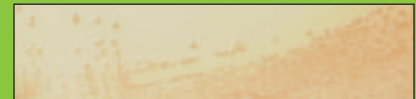

*Staphylococcus schleiferi**
Faint Pink colonies


SaSelect™ Mixed Colony Identification



Staphylococcus aureus
Pink to orange colonies
*Staphylococcus saprophyticus***
Blue to turquoise colonies



Staphylococcus aureus
Pink to orange colonies
*Staphylococcus schleiferi**,**
Faint Pink colonies

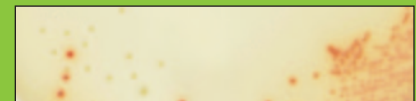

Staphylococcus aureus
Pink to orange colonies
*Staphylococcus epidermidis***
White to faint pink colonies

* This organism may appear after 18 – 24 hours of incubation

** Additional testing required for species identification

► Procedure and Direct ID on reverse side

BIO-RAD