



eFlex™ Software for the BioPlex™ 2200

Software for the first and only fully-automated,
random access, multiplex testing platform

BIO-RAD






Status Workspace

- Intuitive touchscreen interface
- Icon-driven menu
- Continuous access to supply status, system messages and control functions

- 1** Bi-directional LIS connectivity
- 2** One button maintenance tasks
- 3** On-board reagent kit specifics and levels
- 4** "Bread-crumbs" error tracking system to alert operator

The instrument status screen keeps the user informed on system performance, while offering simple icon-driven access to all quality monitoring methods utilized by the BioPlex™ 2200 system

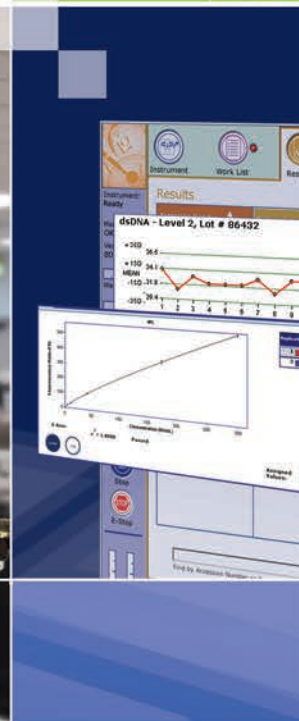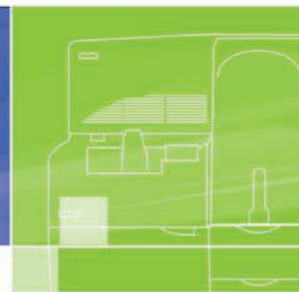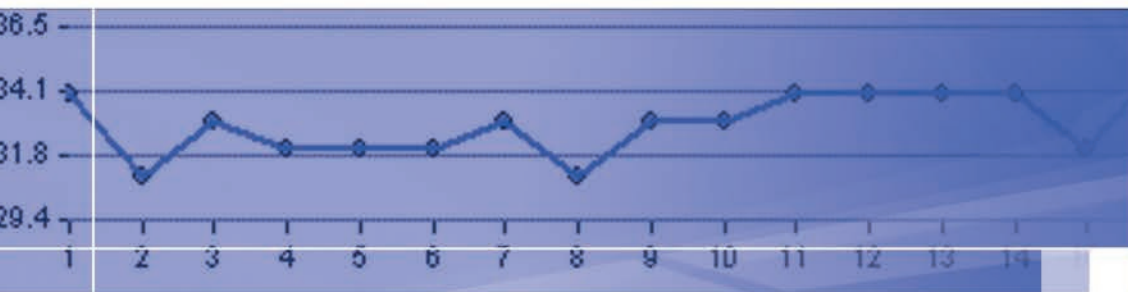

eFlex™ Software Features

Powerful Software

*The system features
a powerful package
of software tools that
provides optimum
flexibility to the
management, set-up,
and sample processing
of the BioPlex™ 2200*

- Simple intuitive design places complete system and sample status at your fingertips
- Add-on test orders can be obtained electronically or by processing new patient samples
- Create custom test groups to order commonly-used assays at the touch of a button
- Fine-tune quality control rules to accommodate your facility's unique needs, including customized Westgard QC rules.
- Take advantage of direct connectivity solutions for troubleshooting, reportable billing and Internet-based QC programs (i.e. Unity™)
- Bi-directional LIS interface streamlines data management
- Analyze results using pattern-recognition algorithms to aid in disease characterization

Operator: Maria Gonzalez
Date: 8/25/2004
Time: 3:16:16 PM

Instrument: Ready
Maintenance: OK
Vessels: 800
Wash buffer
DI water
Sheath fluid
Run
Pause
Stop
E-Stop
Liquid Solid Waste Waste

Work List Row 1 of 38 Go to...

Accession Number	Patient Info.	Assay	Status	Time Left
0975		(ENA6x)	Ordered	
0976		(SLE7x)	Ordered	
0977		(ANA11x)	Ordered	
0978		(ANADem) dsDNAx	Ordered	
0979		(ANADem) Sd70x	Ordered	
0980		(ENA6x)	Ordered	
0981		(ENA6x)	Ordered	
0982		(ENA6x)	Ordered	
0983		(ANA11x)	Ordered	
0984		(ANA)	Ordered	
0985		(SLE7x)	Ordered	
0986		(ENA6x)	Ordered	
0987		(ANA11x)	Ordered	
0988		(ANA11x)	Ordered	

Find by Accession Number or Name: Find View Change Delete Add

Messages:

Set Up System Print CSN LIS

Work List Workspace

- 1 Order all-inclusive panels, individual analytes, or custom test groups
- 2 Minute by minute updates on sample status
- 3 Use **Find** to locate samples quickly
- 4 Add orders manually (or via LIS)

The Work List Workspace provides constant updates on the status of all calibrator, QC and patient orders.

Operator: Maria Gonzalez
Date: 8/25/2004
Time: 3:16:16 PM

Instrument Work List Results

Results Row 1 of 32 Go to...

Accession Number	Patient Info.	Result Date/Time	Assay	Result	Status
0503		06/28/2005 4:46PM	(ANA)	Positive	
0506		06/28/2005 4:46PM	(ANA)	Positive	
0508 (STAT)		06/28/2005 4:46PM	(SLE7)	(multiple)	(multiple)
0509		06/28/2005 4:46PM	(ANA 11)	(multiple)	(multiple)
0511		06/28/2005 4:46PM	(ANA)	Negative	
0514		06/28/2005 4:46PM	(ANA)	Positive	
0533		06/28/2005 4:47PM	(ANA +)	(multiple)	(multiple)
0543		06/28/2005 4:47PM	(ANA 11)	(multiple)	(multiple)
0551		06/28/2005 4:47PM	(ANA +)	Positive	
0554		06/28/2005 4:47PM	(ENA6)	(multiple)	(multiple)
0556		06/28/2005 4:47PM	(SLE7)	1. IU/mL	
			dsDNA	1.4 , Pos	
			Chrom	>8.0 , Pos	
			Ribo P	>8.0 , Pos	
			SS-A	>8.0 , Pos	
			SS-B	>8.0 , Pos	
			Sm	0.4 , Neg	
			RNP	<0.2 , Neg	

Find by Accession Number or Patient Name Find Filter.. View Export Add-On Release Review

Messages:


Set Up System Print CSN LIS

Results Workspace

- 1 Use **Add-On** to place new orders or get additional test results without reprocessing samples
- 2 Easily review results, order dilutions and repeats, or release results
- 3 Expanded view shows complete results
- 4 Access results to pattern recognition algorithms for specific disease families
- 5 Icons identify special sample conditions

The Results Workspace allows real-time updates on patient results and data. Release information and errors are easily identified on any given result order.

eFlex™ Software Functionality Overview

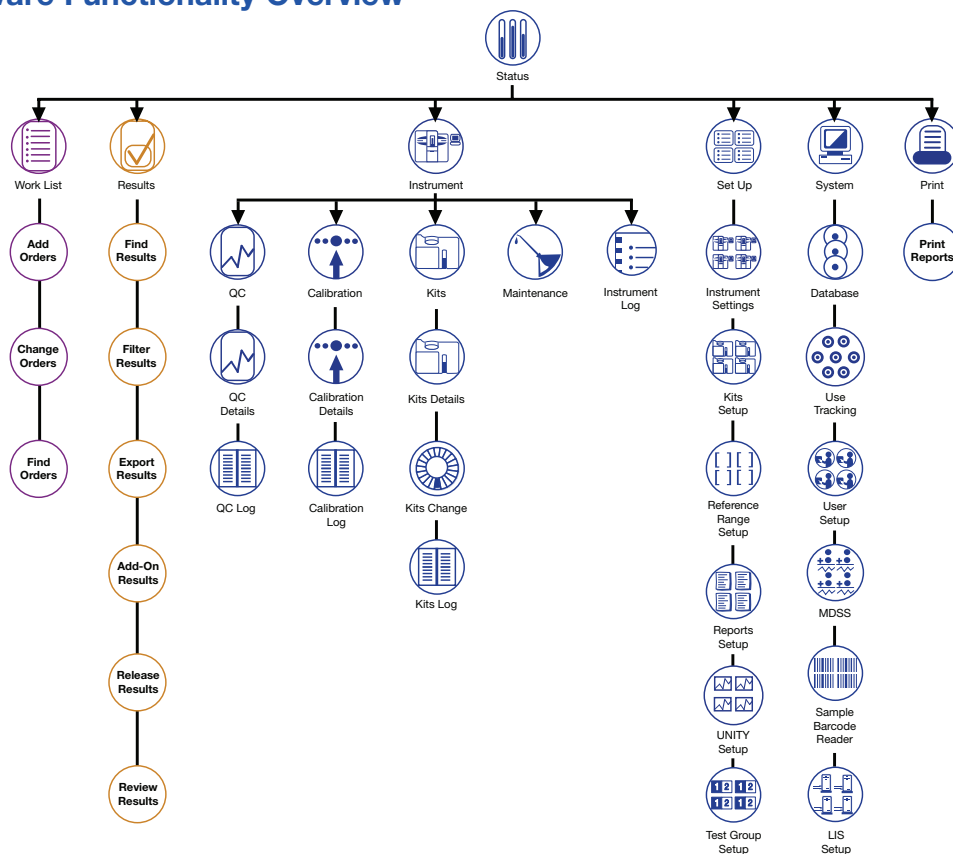

Host Computer

Dual 2.40 GHz CPU, Intel Pentium XEON or better, Memory: 1 GB or better, Hard Drive: 80 GB or better

Ordering Information

Catalog No.	Description	
660-0000	BioPlex™ 2200	1 instrument
660-0800	Sample Racks (30 per box)	1 box
660-0801	Sample Tray	1 tray
660-0817	Instrument Sheath Fluid (2 per box)	1 box
660-0818	Instrument Wash Buffer	1 box
660-2003	Reaction Vessels (1000 per bag)	1 bag
666-0001	Instrument Detector Calibration Pack	1 pack
666-0002	Instrument Detector Clean Pack	1 pack
666-0003	Instrument Probe Cleaning Solution (4 per box)	1 box

BioPlex 2200 System software is only available as part of the complete BioPlex 2200 System.



**Bio-Rad
Laboratories**

*For further information, please contact the Bio-Rad office nearest
you or visit our website at www.bio-rad.com/diagnostics*

*Clinical
Diagnostics Group*

Website www.bio-rad.com/diagnostics **U.S.** 1-800-2BIO-RAD **Australia** 61-2-9914-2800 **Austria** 43-1-877-8901 **Belgium** 32-9-385-5511 **Brazil** 5521-3237-9400
Canada 1-514-334-4372 **China** 86-21-64260808 **Czech Republic** 420-241-430-532 **Denmark** +45-4452-1000 **Finland** 358-9-804-22-00 **France** 33-1-47-95-60-00
Germany +49-(0)89-318-840 **Greece** 30-210-7774396 **Hong Kong** 852-2789-3300 **Hungary** +36-1-459-6100 **India** 91-124-4029300 **Israel** 972-3-9636050 **Italy** +39-02-216091
Japan 81-3-6361-7070 **Korea** 82-2-3473-4460 **Mexico** 52(55)5200-0520 **The Netherlands** +31-318-540666 **New Zealand** 64-9-415-2280 **Norway** 47-23-38-41-30
Poland 48-22-3319999 **Portugal** 351-21-472-7700 **Russia** 7-495-721-14-04 **Singapore** 65-6415-3188 **South Africa** 27-11-442-85-08 **Spain** 34-91-590-5200
Sweden 46-8-555-127-00 **Switzerland** 41-61-717-95-55 **Thailand** 662-651-8311 **United Kingdom** +44-(0)20-8328-2000