



Test VINEO™ Unstable Proteins

Quality diagnostics for enology testing

- Specific: Targets unstable proteins
- Precise: Allows exact dosing of bentonite
- Reliable: Preserves the wine's organoleptic characteristics
- Fast: Results in only 2 hours



Kit VINEO™ Unstable Proteins

For more information, visit us on the web at www.foodscience.bio-rad.com

BIO-RAD

The spoilage of proteins

The appearance of turbidity in white and rose wines after bottling does not affect the flavor of the wine in any way, but its presence can bother consumers.

This turbidity, which was described more than a century ago, is caused by the precipitation of unstable proteins (protein haze) (Laborde, 1904). The wine's stability and clarity over time is guaranteed if these proteins are not present.

Colloid formation is a complex, multifactorial phenomenon. However, experts agree that grape proteins (chitinases, thaumatin-like proteins) are involved (Manteau et al. 2010). Bentonite fining treatment eliminates these undesirable proteins through flocculation (Ribereau-Gayon and Penaud, 1961).

Although adding bentonite is the accepted solution to this problem, its use must be fully controlled. Excessive amounts of bentonite can reduce the wine's organoleptic quality.

The VINEO™ Unstable Proteins kit relies on an immunology-based technique to specifically target the proteins involved in protein haze. It detects the presence of unstable proteins in wine and determines the exact amount of bentonite needed to eliminate all of these proteins. Because of the kit's accuracy, the wine preserves its organoleptic assets.

Key Benefits

- Ready-to-use kits
- Easy protocol
- Easy handling
- High specificity
- Precision of immunoassays
- Easy reading and interpretation of results


Ordering Information

Catalog #	Description
VINEO™ Unstable Proteins	
354-8121	100 tests

The VINEO™ Unstable Proteins kit is manufactured by Bio-Rad under the know-how license of SOFRALAB.

Please contact Bio-Rad for all orders.

Protocol

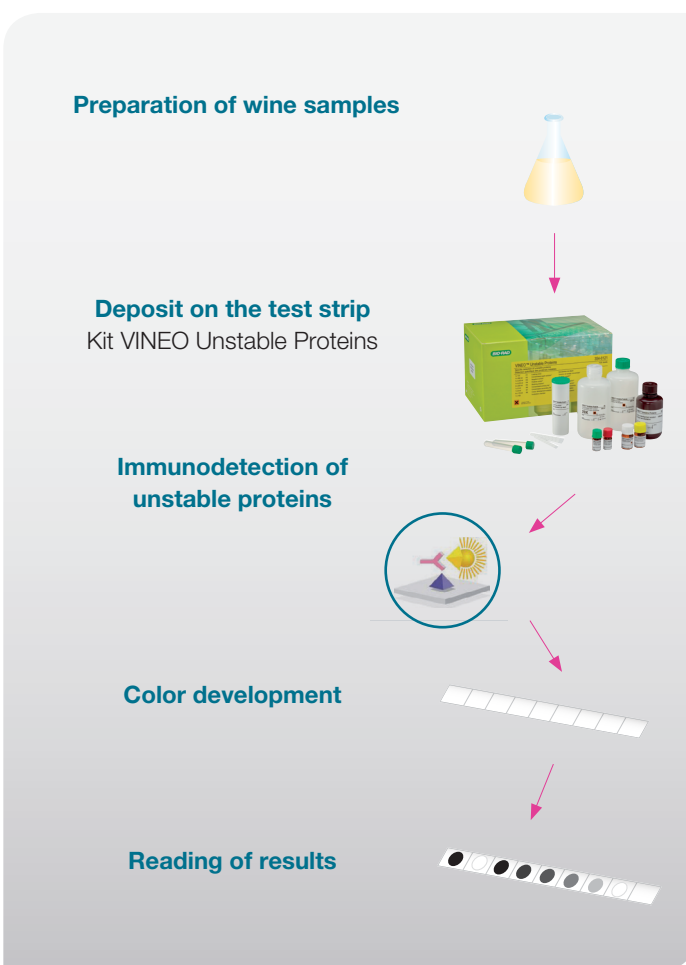

**Bio-Rad
Laboratories, Inc.**

Life Science
Group

Web site www.bio-rad.com **USA** 800 424 6723 **Australia** 61 2 9914 2800 **Austria** 01 877 89 01 **Belgium** 09 385 55 11 **Brazil** 55 31 3689 6600
Canada 905 364 3435 **China** 86 21 6169 8500 **Czech Republic** 420 241 430 532 **Denmark** 44 52 10 00 **Finland** 09 804 22 00
France 01 47 95 69 65 **Germany** 089 31 884 0 **Greece** 30 210 777 4396 **Hong Kong** 852 2789 3300 **Hungary** 36 1 459 6100 **India** 91 124 4029300
Israel 03 963 6050 **Italy** 39 02 216091 **Japan** 03 6361 7000 **Korea** 82 2 3473 4460 **Malaysia** 60 3 2117 5260 **Mexico** 52 555 488 7670
The Netherlands 0318 540666 **New Zealand** 64 9 415 2280 **Norway** 23 38 41 30 **Poland** 48 22 331 99 99 **Portugal** 351 21 472 7700
Russia 7 495 721 14 04 **Singapore** 65 6415 3170 **South Africa** 27 861 246 723 **Spain** 34 91 590 5200 **Sweden** 08 555 12700
Switzerland 061 717 95 55 **Taiwan** 886 2 2578 7189 **Thailand** 66 2 6518311 **United Kingdom** 020 8328 2000