



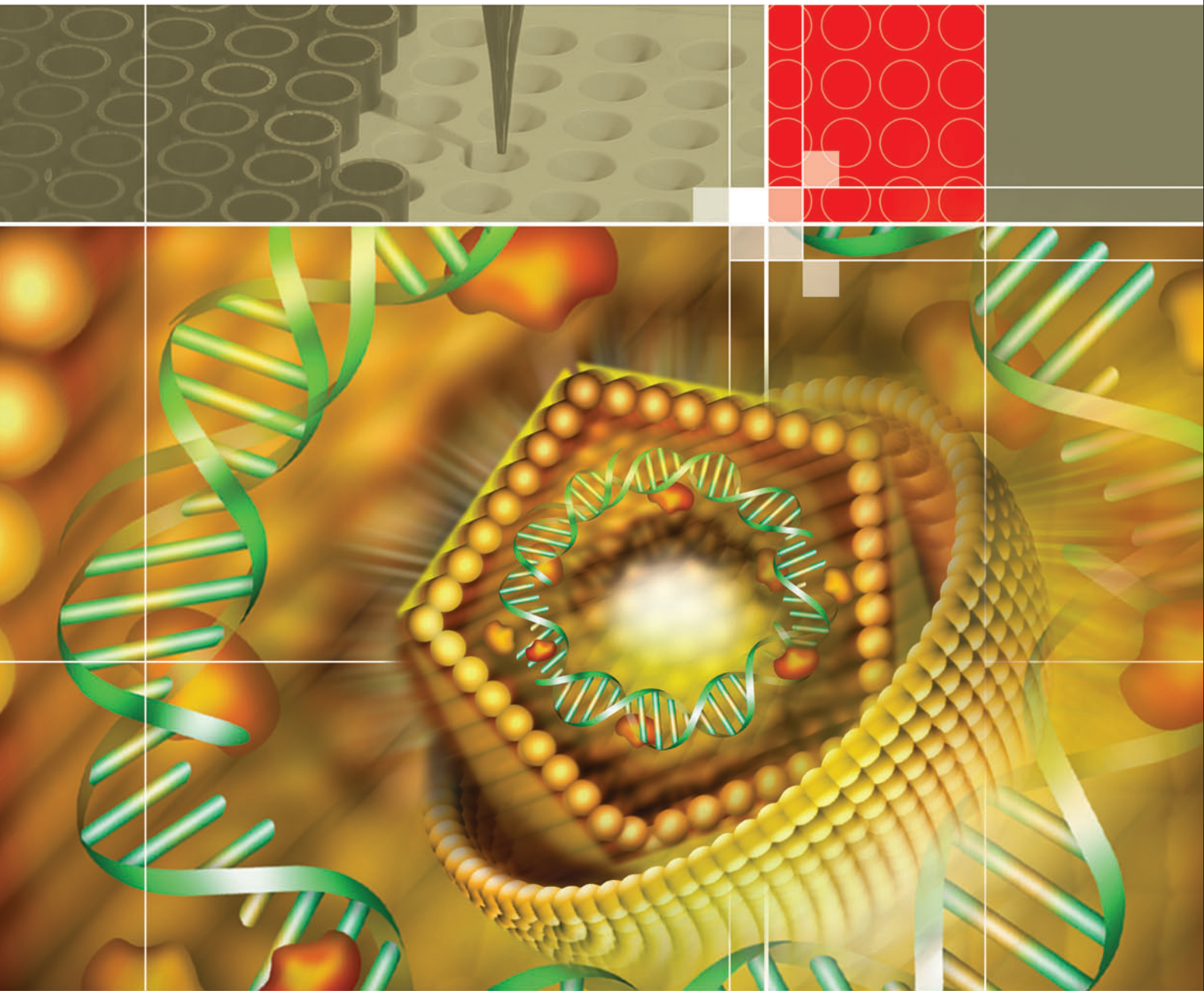

GS HBsAg EIA 3.0

Superior Hepatitis Testing from Bio-Rad

BIO-RAD

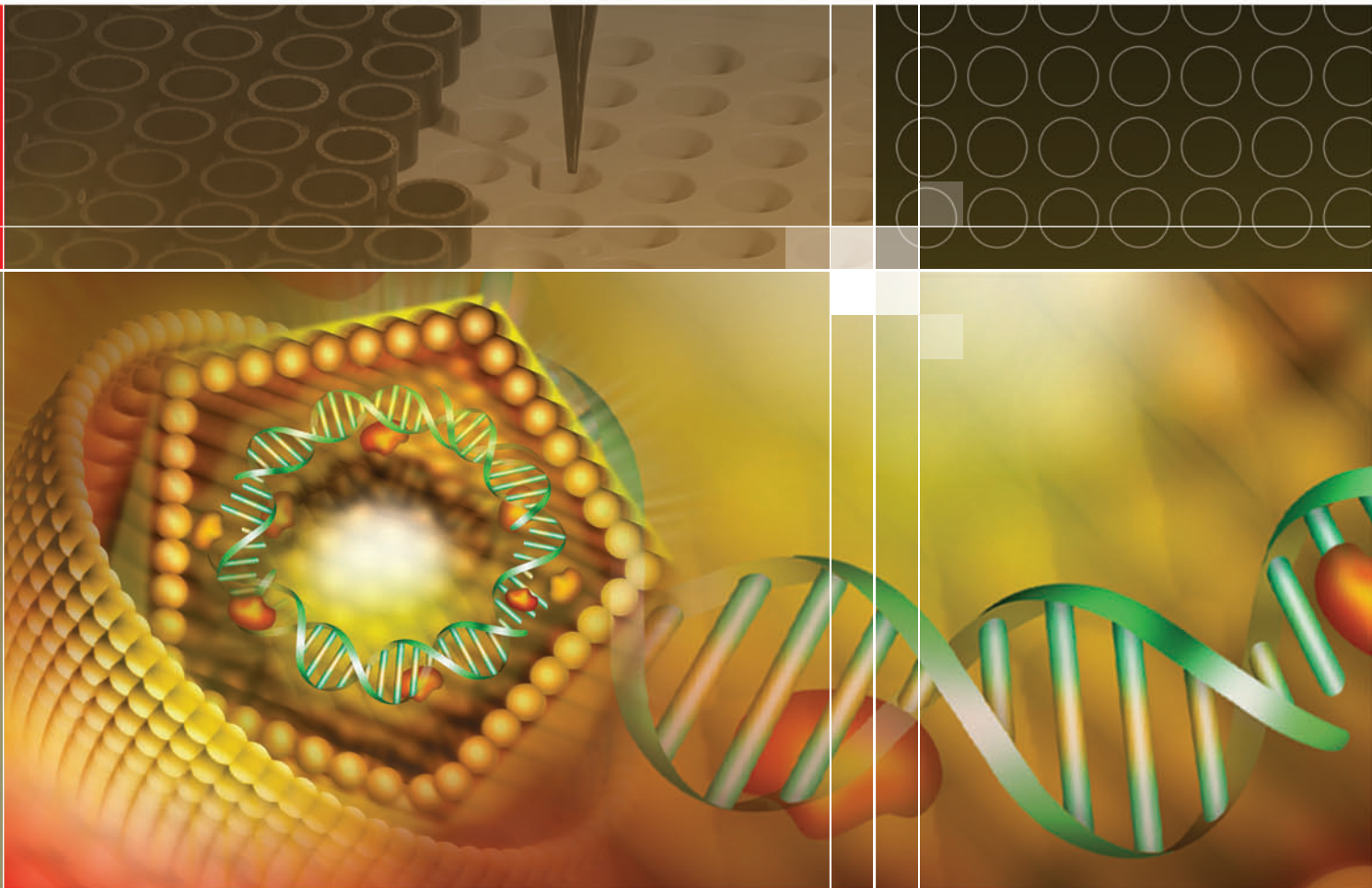

GS HBsAg EIA 3.0

Outstanding sensitivity of less than 0.1 ng/mL of ad/ay HBsAg*

The GS HBsAg EIA 3.0 is an EIA for the detection of Hepatitis B Surface Antigen (HBsAg) in human serum, plasma, and cadaveric serum samples.

GS HBsAg Confirmatory Assay 3.0

For Confirmatory Neutralization of GS HBsAg EIA 3.0 Reactive Specimens

* Calibrated versus PHA 805 commercial panel

Ease-of-Use

- No sample predilution
- Total incubation time = 2.5 hours
- Common reagents can be used with other Bio-Rad EIA assays
- Manual and automated use

Advanced Safety Features

- Bar-coded reagents
- Color-coded reagent labels
- Microplate strip identification
- Color changes during protocol steps

Earliest Hepatitis B Infection Detection Demonstrated on 21 Commercial Seroconversion Panels*

GS HBsAg EIA 3.0 Reactivity on Commercial HBsAg Seroconversion Panels Versus Competitors				
Earliest Detection				
Panel I.D.#	GS HBsAg EIA 3.0		HBsAg EIA Competitor #1 (overnight)	HBsAg EIA Competitor #2
	Static	Shaker		
902	7	7	9	9
903	3	2	5	5
904	2	1	3	3
905	1	1	1	1
906	2	2	2	3
907	6	5	6	6
908	6	4	7	7
909	4	3	4	5
910	3	3	3	4
911	20	20	21	22
912	7	6	8	8
914	2	2	5	5
915	1	1	9	11
917	2	2	3	3
918	2	2	3	3
919	5	4	6	7
920	3	3	3	3
924	3	2	3	4
931	5	3	NT	7
932	10	8	NT	10
40565L	2	1	3	NT

* The sensitivity of GS HBsAg EIA 3.0, procedures A and B, was assessed using 20 seroconversion panels purchased from Boston Biomedica, Inc. (West Bridgewater, MA), and 1 panel from Serologicals, Inc. (Clarkston, GA). The number of the bleed at which HBsAg was detected (S/C ≥ 1.0) by the GS HBsAg EIA 3.0, procedures A and B, is compared to other licensed HBsAg assays (as recorded in the certificates of analysis).


GS HBsAg 3.0 Static Procedure

- Earlier detection on 12 out of 19 seroconversions for HBsAg 3.0 versus competitor assay number 1
- Earlier detection on 16 out of 20 seroconversions for HBsAg 3.0 versus competitor assay number 2

GS HBsAg 3.0 Shaking Procedure

- Earlier detection on 15 out of 19 seroconversions versus competitor assay number 1
- Earlier detection on 18 out of 20 seroconversions versus competitor assay number 2

GS HBsAg EIA 3.0 Testing Procedure

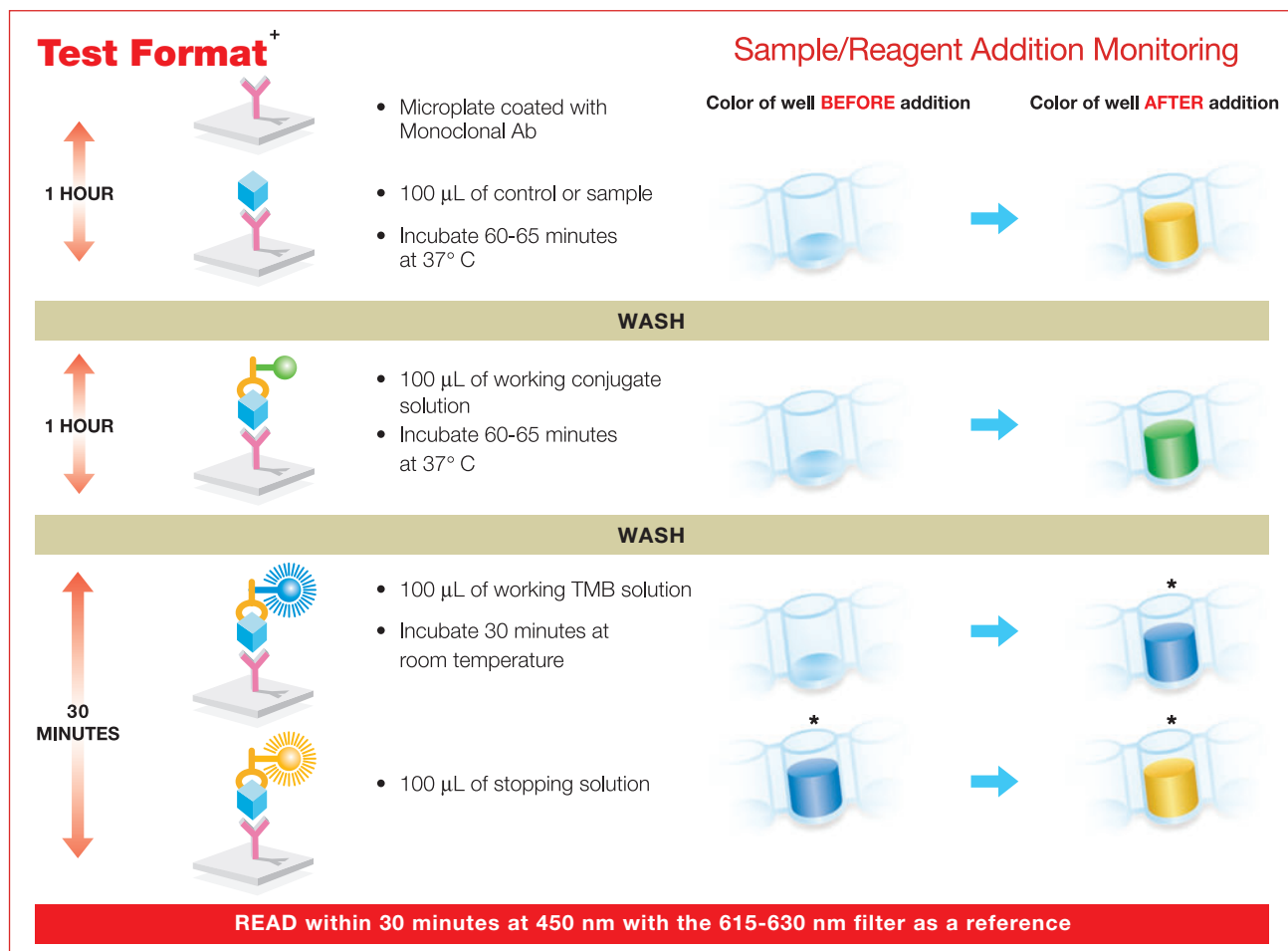

Y Anti-HBsAg Monoclonal Ab HBsAg Y Anti HBsAg Monoclonal Ab Conjugate

+ Static or Shaker.

* Reactive samples after 30 minute incubation.

Ordering Information

Enzyme immunoassay (EIA) for the detection of Hepatitis B Surface Antigen (HBsAg) in Human Serum, Plasma, and Cadaveric Serum Specimens:

Catalog No.	Description	
32591	GS HBsAg EIA 3.0	480 tests
32592	GS HBsAg EIA 3.0	960 tests
25258	GS HBsAg EIA 3.0	4800 tests

For Confirmatory Neutralization of GS HBsAg EIA 3.0 Reactive Specimens:

Catalog No.	Description	
32594	GS HBsAg Confirmatory Assay 3.0.....	25 tests

BIO-RAD

Bio-Rad
Laboratories

For further information, please contact the Bio-Rad office nearest you or visit our website at www.bio-rad.com/diagnostics.

Clinical
Diagnostics Group

Website www.bio-rad.com/diagnostics U.S. 1-800-2BIO-RAD Australia 61-2-9914-2800 Austria 43-1-877-8901 Belgium 32-9-385-5511 Brazil 5521-3237-9400
Canada 1-514-334-4372 China 86-21-64260808 Czech Republic 420-241-430-532 Denmark +45-4452-1000 Finland 358-9-804-22-00 France 33-1-47-95-60-00
Germany +49-(0)89-318-840 Greece 30-210-7774396 Hong Kong 852-2789-3300 Hungary 36-1-455-8800 India 91-124-4029300 Israel 972-3-9636050 Italy +39-02-216091
Japan 81-3-6361-7070 Korea 82-2-3473-4460 Mexico +52(55)5488-7670 The Netherlands +31-318-540666 New Zealand 64-9-415-2280 Norway 47-23-38-41-30
Poland 48-22-3319999 Portugal 351-21-472-7700 Russia 7-495-721-14-04 Singapore 65-6415-3188 South Africa 27-11-442-85-08 Spain 34-91-590-5200
Sweden 46-8-555-127-00 Switzerland 41-61-717-95-55 Thailand 662-651-8311 United Kingdom +44-(0)20-8328-2000